

---

---

**Plastics — Thermogravimetry (TG) of  
polymers —**

Part 3:

**Determination of the activation  
energy using the Ozawa-Friedman plot  
and analysis of the reaction kinetics**

*Plastiques — Thermogravimétrie (TG) des polymères —*

*Partie 3: Détermination de l'énergie d'activation à l'aide du  
graphique d'Ozawa-Friedman et analyse cinétique de la réaction*

STANDARDSISO.COM : Click to view the full PDF of ISO 11358-3:2021



STANDARDSISO.COM : Click to view the full PDF of ISO 11358-3:2021



**COPYRIGHT PROTECTED DOCUMENT**

© ISO 2021

All rights reserved. Unless otherwise specified, or required in the context of its implementation, no part of this publication may be reproduced or utilized otherwise in any form or by any means, electronic or mechanical, including photocopying, or posting on the internet or an intranet, without prior written permission. Permission can be requested from either ISO at the address below or ISO's member body in the country of the requester.

ISO copyright office  
CP 401 • Ch. de Blandonnet 8  
CH-1214 Vernier, Geneva  
Phone: +41 22 749 01 11  
Email: [copyright@iso.org](mailto:copyright@iso.org)  
Website: [www.iso.org](http://www.iso.org)

Published in Switzerland

# Contents

	Page
Foreword .....	iv
Introduction .....	v
<b>1 Scope</b> .....	<b>1</b>
<b>2 Normative references</b> .....	<b>1</b>
<b>3 Terms and definitions</b> .....	<b>1</b>
<b>4 Principle</b> .....	<b>2</b>
<b>5 Apparatus</b> .....	<b>2</b>
<b>6 Test specimens</b> .....	<b>2</b>
<b>7 Mass and temperature calibration</b> .....	<b>2</b>
7.1 Mass calibration .....	2
7.2 Temperature calibration .....	2
<b>8 Procedure</b> .....	<b>2</b>
8.1 General .....	2
8.2 Non-oxidative reactions .....	4
8.3 Oxidative reactions .....	4
<b>9 Expression of results</b> .....	<b>4</b>
9.1 Graphical presentation .....	4
9.2 Determination of activation energy .....	4
9.3 Master curves .....	6
<b>10 Precision</b> .....	<b>9</b>
<b>11 Test report</b> .....	<b>9</b>
<b>Annex A (informative) Mathematical background</b> .....	<b>10</b>
<b>Bibliography</b> .....	<b>14</b>

## Foreword

ISO (the International Organization for Standardization) is a worldwide federation of national standards bodies (ISO member bodies). The work of preparing International Standards is normally carried out through ISO technical committees. Each member body interested in a subject for which a technical committee has been established has the right to be represented on that committee. International organizations, governmental and non-governmental, in liaison with ISO, also take part in the work. ISO collaborates closely with the International Electrotechnical Commission (IEC) on all matters of electrotechnical standardization.

The procedures used to develop this document and those intended for its further maintenance are described in the ISO/IEC Directives, Part 1. In particular, the different approval criteria needed for the different types of ISO documents should be noted. This document was drafted in accordance with the editorial rules of the ISO/IEC Directives, Part 2 (see [www.iso.org/directives](http://www.iso.org/directives)).

Attention is drawn to the possibility that some of the elements of this document may be the subject of patent rights. ISO shall not be held responsible for identifying any or all such patent rights. Details of any patent rights identified during the development of the document will be in the Introduction and/or on the ISO list of patent declarations received (see [www.iso.org/patents](http://www.iso.org/patents)).

Any trade name used in this document is information given for the convenience of users and does not constitute an endorsement.

For an explanation of the voluntary nature of standards, the meaning of ISO specific terms and expressions related to conformity assessment, as well as information about ISO's adherence to the World Trade Organization (WTO) principles in the Technical Barriers to Trade (TBT), see [www.iso.org/iso/foreword.html](http://www.iso.org/iso/foreword.html).

This document was prepared by Technical Committee ISO/TC 61, *Plastics*, Subcommittee SC 5, *Physical-chemical properties*.

This second edition cancels and replaces the first edition (ISO 11358-3:2013), which has been technically revised.

The changes compared to the previous edition are as follows:

- the term "conversion" has been deleted;
- a corresponding reference to ISO 11358-2 has been added in [Clause 3](#);
- details of the gas atmosphere in [8.1](#) and [8.2](#) have been clarified.

A list of all parts in the ISO 11358 series can be found on the ISO website.

Any feedback or questions on this document should be directed to the user's national standards body. A complete listing of these bodies can be found at [www.iso.org/members.html](http://www.iso.org/members.html).

## Introduction

The Ozawa-Friedman plot (logarithm of the rate of mass loss versus the reciprocal of absolute temperature at a given mass loss) is a derivative method that can be applied to data obtained by any mode of temperature change in thermal analysis; e.g. isothermal, constant heating rate, sample-controlled thermal analysis, temperature jump, and repeated temperature scanning.

If controlled rate thermogravimetry (CRTG) is used to study the decomposition of polymers, the Ozawa-Friedman method is typically applied to the analysis of data obtained by CRTG and also to that obtained by the combined use of isothermal thermogravimetry (iso-TG) with conventional linear heating rate thermogravimetry (LHTG), i.e. using a constant heating rate.

STANDARDSISO.COM : Click to view the full PDF of ISO 11358-3:2021

[STANDARDSISO.COM](https://standardsiso.com) : Click to view the full PDF of ISO 11358-3:2021

# Plastics — Thermogravimetry (TG) of polymers —

## Part 3:

# Determination of the activation energy using the Ozawa-Friedman plot and analysis of the reaction kinetics

## 1 Scope

This document specifies an analysis method for determining the activation energy using the Ozawa-Friedman plot. It also specifies the preparation of master plots for verification of the reaction kinetics determined by thermogravimetry.

## 2 Normative references

The following documents are referred to in the text in such a way that some or all of their content constitutes requirements of this document. For dated references, only the edition cited applies. For undated references, the latest edition of the referenced document (including any amendments) applies.

ISO 11358-1, *Plastics — Thermogravimetry (TG) of polymers — Part 1: General principles*

ISO 11358-2, *Plastics — Thermogravimetry (TG) of polymers — Part 2: Determination of activation energy*

## 3 Terms and definitions

For the purposes of this document, the terms and definitions given in ISO 11358-1, ISO 11358-2 and the following apply.

ISO and IEC maintain terminological databases for use in standardization at the following addresses:

- ISO Online browsing platform: available at <https://www.iso.org/obp>
- IEC Electropedia: available at <http://www.electropedia.org/>

### 3.1

#### generalized time

$t_{\text{gen}}$

$$t_{\text{gen}} = \int_0^t \exp\left(-\frac{E_a}{RT}\right) dt$$

where

$E_a$  is the activation energy, expressed in kJ/mol;

$R$  is the gas constant, expressed as 8,314 J/(mol K);

$T$  is the absolute temperature, expressed in Kelvin;

$t$  is time, expressed in minutes.

### 3.2

#### generalized rate of conversion

$dC/dt_{\text{gen}}$

$$\frac{dC}{dt_{\text{gen}}} = \exp\left(\frac{E_a}{RT}\right) \frac{dC}{dt}$$

Note 1 to entry: For the definition of the degree of conversion, refer to ISO 11358-2.

### 3.3

#### master curve

plot that can be used to evaluate the results and investigate the reaction kinetics models

EXAMPLE Conversion versus the generalized time, conversion versus the generalized rate of conversion, generalized time versus the generalized rate of conversion.

## 4 Principle

Test specimens are heated using any temperature profile and the change in mass is measured as a function of temperature and time. At a given conversion, the logarithm of the rate of conversion is plotted versus the reciprocal of the absolute temperature, and the activation energy is calculated from the slope of the straight line fit to the data thus obtained.

At least two of the master curves enable verification of the reaction kinetics analysis.

## 5 Apparatus

The apparatus shall be in accordance with ISO 11358-1.

## 6 Test specimens

Test specimens shall be in the form of powder, pellets, flakes, filaments, or film. The test specimens shall be prepared by cutting the material, as necessary, to a size appropriate for the apparatus (see ISO 11358-1).

## 7 Mass and temperature calibration

### 7.1 Mass calibration

The procedure of mass calibration is given in ISO 11358-1.

### 7.2 Temperature calibration

The procedure of temperature calibration is given in ISO 11358-1.

## 8 Procedure


### 8.1 General

The determination of the rate of conversion  $dC/dt$  is necessary for the analysis in this document. The rate of conversion versus absolute temperature shall be determined.


See ISO 11358-1 and ISO 11358-2.

Perform tests using either an isothermal run, constant rate heating, sample mass controlled rate thermal analysis, temperature jump, repeated temperature scanning, or any combination of the above,

using specimens of similar mass ( $\pm 1\%$ ). For examples of linear heating rate thermogravimetry (LHTG) and sample mass controlled rate thermogravimetry (CRTG), see [Figure 1 a\)](#) and [Figure 1 b\)](#), respectively.

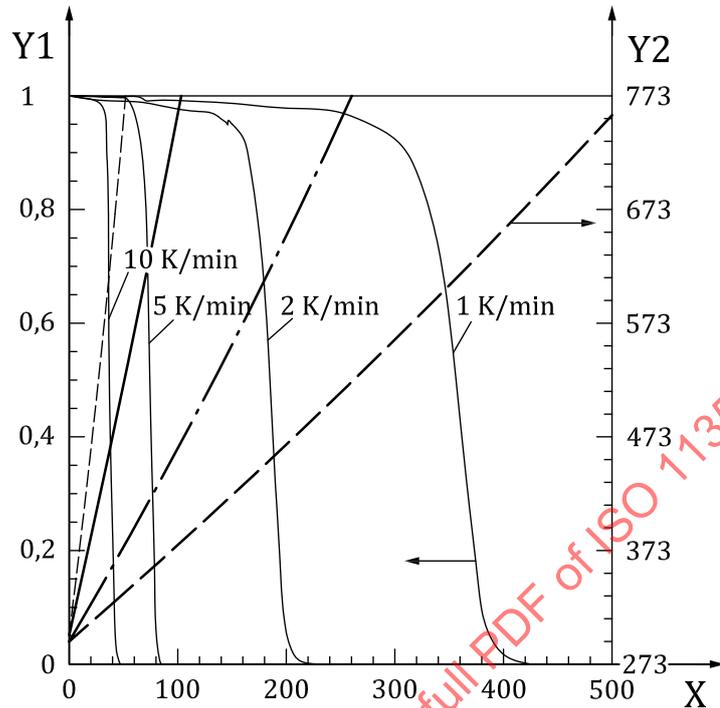

a) Example of linear heating rate thermogravimetry (LHTG) measurements of PMMA

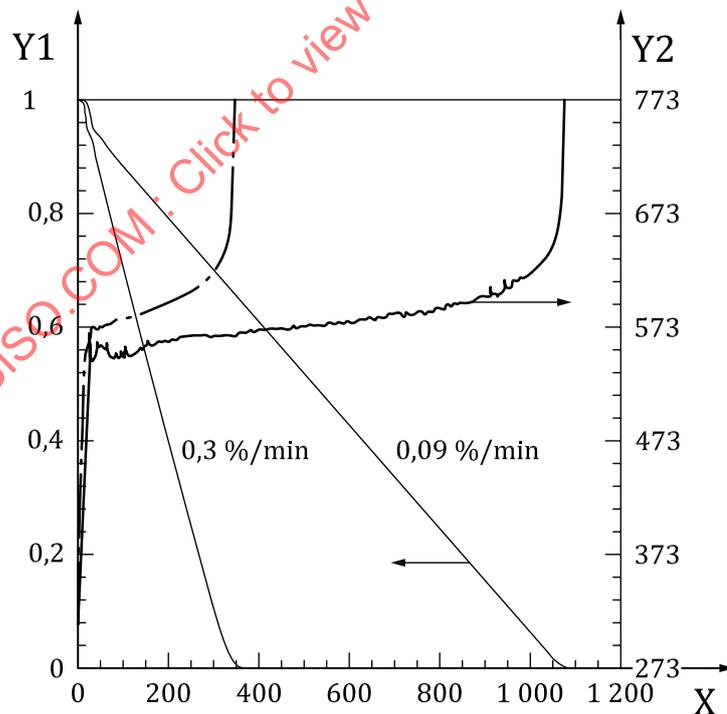

b) Example of sample mass controlled rate thermogravimetry (CRTG) measurements of PMMA

**Key**

X  $t/\text{min}$

Y1  $(1 - C)$

Y2  $T/\text{K}$

NOTE The units of %/min indicate the percentage weight loss per minute of the mass controlled rate.

**Figure 1 — Examples of linear heating rate thermogravimetry (LHTG) measurements and sample mass controlled rate thermogravimetry (CRTG) measurements of PMMA**

Preferably, the specimen mass should be between 1 mg and 10 mg and the temperature scanning rate between  $2 \text{ K min}^{-1}$  and  $10 \text{ K min}^{-1}$ .

Determine the rate of conversion (or rate of change of mass loss fraction with time) at a given conversion (or given mass loss fraction).

## 8.2 Non-oxidative reactions

The inert gas atmosphere shall be in accordance with ISO 11358-2.

## 8.3 Oxidative reactions

The oxidative gas atmosphere shall be in accordance with ISO 11358-2.

## 9 Expression of results

### 9.1 Graphical presentation

For a given conversion (or mass loss fraction), present the thermogravimetry data in the form of a plot of the logarithm of the rate of conversion (or rate of change of mass loss fraction with time) versus the reciprocal of the absolute temperature, i.e. the Ozawa-Friedman plot, see [Figure 2](#).

### 9.2 Determination of activation energy

The conversion (or the mass loss fraction) is not directly related to the quantity of reacted pieces, for example as in the random scission of the main chain of polymers, but the conversion,  $C$ , of a reacting material is expressed as a function of the fraction of a structural quantity, such as a group, a constituent, a broken bond, etc., which is represented by  $\alpha$ <sup>[1]</sup>.

Generally, the rate formulae are as shown in [Formulae \(1\)](#) and [\(2\)](#):

$$C = f(\alpha) \quad (1)$$

and

$$\frac{d\alpha}{dt} = A \exp\left(\frac{-E_a}{RT}\right) g(\alpha) \quad (2)$$

where

$\alpha$  is the reacted fraction at time  $t$ ;

$f(\alpha)$  is an arbitrary function describing the relationship between the conversion  $C$  and the fraction of the reacted structural quantity;

$A$  is the pre-exponential factor;

$g(\alpha)$  is a function describing the reaction mechanism.

By taking the logarithm of both sides of [Formula \(2\)](#), [Formula \(3\)](#) can be obtained as follows:


$$\ln\left(\frac{d\alpha}{dt}\right) = \ln\{Ag(\alpha)\} - \frac{E_a}{RT} \quad (3)$$

When plotted for fixed values of  $\alpha$ , with various conversion rate values determined for different experimental temperature control profiles, the slope of the plot of  $\ln\left(\frac{d\alpha}{dt}\right)$  versus  $\frac{1}{T}$  has the value  $-E_a/R$  and is used to determine the activation energy  $E_a$ , as the first term in the right side of [Formula \(3\)](#) is constant.

As an example, the experimental data presented in [Figures 1a](#) and [1b](#) have been analysed, assuming an  $n^{\text{th}}$  order type reaction and thus  $C = f(\alpha) = \alpha$  (see [Annex A](#)), the results of which are presented in [Table 1](#) and [Figure 2](#).

**Table 1 — Rate of conversion  $dC/dt$  and estimated activation energy  $E_a$  at different conversion levels<sup>a</sup>**

	1 000/T	ln(dC/dt (% min <sup>-1</sup> ))	1 000/T	ln(dC/dt (% min <sup>-1</sup> ))	1 000/T	ln(dC/dt (% min <sup>-1</sup> ))
	(K <sup>-1</sup> )	at C = 0,2	(K <sup>-1</sup> )	at C = 0,5	(K <sup>-1</sup> )	at C = 0,8
CRTG, 0,09 %/min	1,779	-2,425	1,735	-2,411	1,677	-2,448
CRTG, 0,3 %/min	1,731	-1,161	1,693	-1,183	1,641	-1,288
Activation energy (kJmol <sup>-1</sup> )						
Ozawa-Friedman plot	207,3		232,8		261,9	
<sup>a</sup> Determined from the controlled rate thermogravimetry (CRTG) data presented in <a href="#">Figure 1b</a> .						

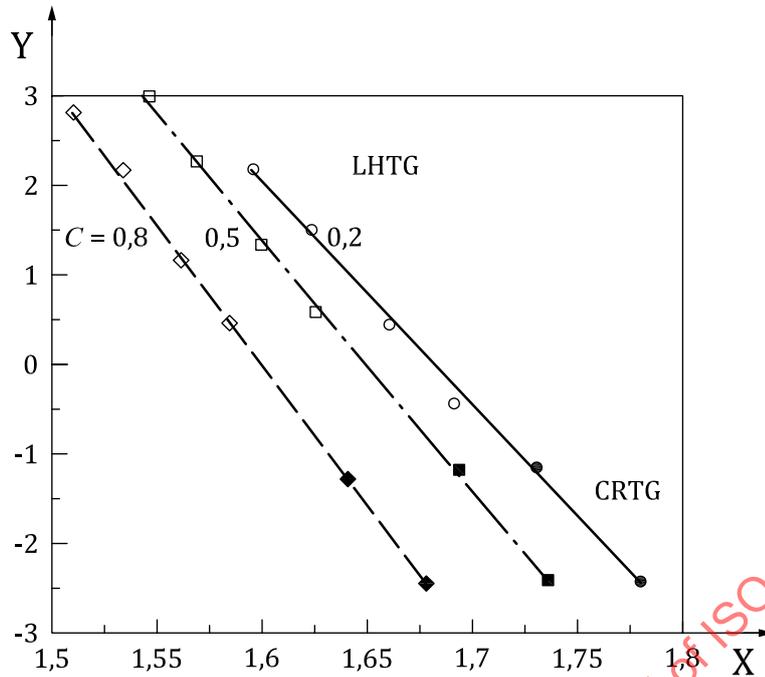

**Key**

X 1 000/T K<sup>-1</sup>  
 Y ln (dC/dt / % min<sup>-1</sup>)

NOTE Determined from [Figure 1 a](#)) (open symbols), and CRTG data from [Table 1](#) (closed symbols) at conversion values of 0,2, 0,5 and 0,8

**Figure 2 — Ozawa-Friedman plot of the LHTG data**

**9.3 Master curves**

Given [Formula \(1\)](#) and [Formula \(2\)](#), the assumed models of the method, [Formula \(4\)](#) is calculated as follows:

$$\frac{dC}{dt} = \frac{df(\alpha)}{d\alpha} \frac{d\alpha}{dt} = \frac{df(\alpha)}{d\alpha} A \exp\left(\frac{-E_a}{RT}\right) g(\alpha) \tag{4}$$

Thus, using the generalized rate of conversion in [3.3](#), [Formula \(5\)](#) is obtained:

$$\frac{dC}{dt_{gen}} = A g(\alpha) \frac{df(\alpha)}{d\alpha} \tag{5}$$

The master curve, the plot of

$$\frac{dC}{dt_{gen}} = \exp\left(\frac{E_a}{RT}\right) \frac{dC}{dt} \tag{6}$$

can be used to demonstrate how consistent the reaction being analysed is with the assumed models of [Formula \(1\)](#) and [Formula \(2\)](#).


Similarly, other relations between C, t<sub>gen</sub> and dC/dt<sub>gen</sub> can be derived. According to [Formula \(6\)](#) and [Formula \(A.7\)](#) in [Annex A](#), the generalized rate of conversion dC/dt<sub>gen</sub> has the following relationship to

the experimental kinetic data via the generalized rate of the reacted fraction  $d\alpha/dt_{\text{gen}}$  for an  $n^{\text{th}}$  order reaction where  $C = f(\alpha) = \alpha$ , as shown in [Table A.1](#):

$$\frac{dC/dt_{\text{gen}}}{(dC/dt_{\text{gen}})_{C=0,5}} = \frac{d\alpha/dt_{\text{gen}}}{(d\alpha/dt_{\text{gen}})_{C=0,5}} = \frac{d\alpha/dt}{(d\alpha/dt)_{C=0,5}} \cdot \frac{\exp(E_a/RT)}{\exp(E_a/RT_{0,5})} \quad (7)$$

where,  $T_{0,5}$  is the reaction temperature at  $C = 0,5$ .

For calculating the experimental values of  $(dC/dt_{\text{gen}})/(dC/dt_{\text{gen}})_{C=0,5}$ , the experimental kinetic data have to be used, for example, see [Figure 3](#).

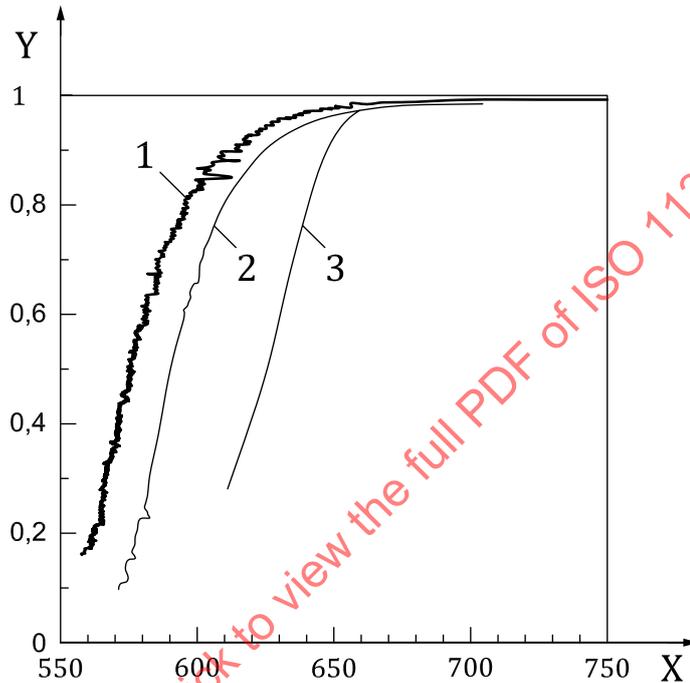

#### Key

- X T/K
- Y Conversion C
- 1 CRTG 0,09 %/min
- 2 CRTG 0,3 %/min
- 3 LHTG 2 K/min

**Figure 3 — Experimental conversion curves of CRTG (0,3 %/min and 0,09 %/min) and LHTG (2 K/min) of PMMA plotted as a function of temperature**

Furthermore, master curves may be used to help identify the reaction mechanism model,  $g(\alpha)$ , and the relation between the conversion and the reaction species,  $f(\alpha)$ , by fitting theoretical models to the experimental data, and revising the model as necessary to achieve a good fit of experimental data with theoretical values.


In the case of  $C = f(\alpha) = \alpha$ , [Formula \(5\)](#) can be written in the form of [Formula \(8\)](#):

$$\frac{dC}{dt_{gen}} = Ag(\alpha) \tag{8}$$

Using the conversion of  $\alpha = 0,5$  as a reference value, [Formula \(9\)](#) can be derived from [Formula \(2\)](#) and [\(7\)](#):

$$\frac{d\alpha / dt_{gen}}{(\frac{d\alpha}{dt_{gen}})_{\alpha=0,5}} = \frac{g(\alpha)}{g(0,5)} \tag{9}$$

[Formula \(9\)](#) states that, at a given reacted fraction  $\alpha$ , the experimentally determined  $(\frac{d\alpha}{dt_{gen}})/(\frac{d\alpha}{dt_{gen}})_{\alpha=0,5}$  and the theoretically calculated  $g(\alpha)/g(0,5)$  are equivalent when an appropriate function  $g(\alpha)$  for describing the reaction process under investigation is applied. [Figure 4](#) compares the theoretical plot of  $g(\alpha)/g(0,5)$  versus  $\alpha$  with the experimental master plots of  $(\frac{d\alpha}{dt_{gen}})/(\frac{d\alpha}{dt_{gen}})_{\alpha=0,5}$  versus  $\alpha$  calculated using the experimental data of [Figure 3](#). In this case an  $n^{\text{th}}$  order reaction with a reaction order  $n = 1,9$  has been assumed and used to generate the theoretical plot and provides a good fit to the experimental master plot data as shown in [Figure 4](#).

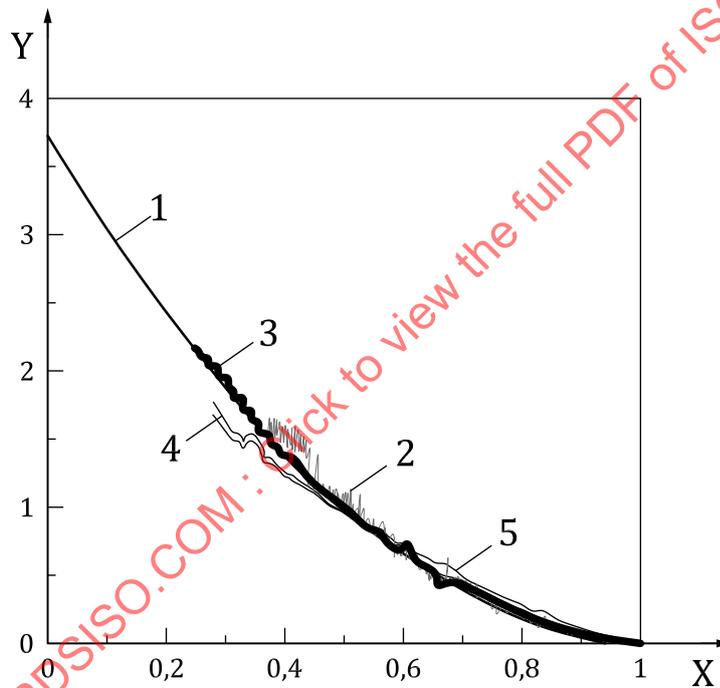

**Key**

- X  $\alpha$
- Y  $(\frac{d\alpha}{dt_{gen}})/(\frac{d\alpha}{dt_{gen}})_{\alpha=0,5}$
- 1  $g(\alpha)/g(0,5)$
- 2 CRTG 0,09 %/min
- 3 CRTG 0,3 %/min
- 4 LHTG 2 °K/min
- 5 LHTG 5 °K/min

NOTE It was assumed that  $C = \alpha$  as for an  $n^{\text{th}}$  order reaction as listed in [Table A.1](#) and thus  $dC/dt_{gen} = d\alpha/dt_{gen}$ .

**Figure 4 — Comparison of the experimental master plots of  $(\frac{d\alpha}{dt_{gen}})/(\frac{d\alpha}{dt_{gen}})_{\alpha=0,5}$  versus  $\alpha$  with the theoretical master curves  $g(\alpha)/g(0,5)$  versus  $\alpha$  for the thermal decomposition of PMMA**

## 10 Precision

The standard deviations of the estimated activation energy  $E_a$  and the generalized rate of conversion  $dC/dt_{\text{gen}}$  are reported for the specific case analysed to be less than 5 %<sup>[4]</sup>.

## 11 Test report

The test report shall include the following information:

- a) a reference to this document, i.e. ISO 11358-3:2021;
- b) the details necessary for complete identification of the material analysed;
- c) the form and dimensions of the test specimen;
- d) the mass of the test specimen;
- e) details of the conditioning of the specimen prior to the test;
- f) the specimen pan size, mass and material of construction;
- g) the atmosphere and gas-flow rate used;
- h) the details of the temperature profile used;
- i) the standard reference material(s) used for temperature calibration;
- j) the activation energy, determined from a linear slope in the plot of the logarithm of the rate of mass loss versus the reciprocal absolute temperature at a given mass loss, and the mass loss (see [Figure 1](#));
- k) the reaction order determined with a fitting function of the kinetics, if required (see [Table 1](#));
- l) any observations regarding equipment, test conditions or test specimen behaviour;
- m) the date of the determination.

## Annex A (informative)

### Mathematical background

In thermal analysis, the property  $P$  of the sample, for example the mass, or the rate of change of the property,  $dP/dt$  (or  $dP/dT$ ) are recorded as a function of time,  $t$ , or temperature,  $T$ . In order to estimate the kinetic parameters by this method, the property measured should be described only by the parameter,  $\alpha$ . Thus, the conversion,  $C$ , which is equal to  $(P-P_0)/(P_\infty-P_0)$ , (where  $P_0$  and  $P_\infty$  are the property values before and after the reaction, respectively), is a function of  $\alpha$  as written in [Formula \(A.1\)](#):

$$C = f(\alpha) \quad (\text{A.1})$$

The parameter  $\alpha$  is assumed to change following ordinary reaction kinetics as shown in [Formula \(A.2\)](#):

$$\frac{d\alpha}{dt} = A \exp\left(-\frac{E_a}{RT}\right) g(\alpha) \quad (\text{A.2})$$

where  $A$ ,  $E_a$  and  $R$  are the pre-exponential factor, the activation energy and the gas constant, respectively.

[Formula \(A.2\)](#) is rearranged as [Formula \(A.3\)](#):

$$\int_0^\alpha \frac{d\alpha}{g(\alpha)} = A \int_0^t \exp\left(-\frac{E_a}{RT}\right) dt \quad (\text{A.3})$$

or [Formula \(A.4\)](#):

$$G(\alpha) = At_{\text{gen}} \quad (\text{A.4})$$

where [Formula \(A.5\)](#) applies

$$G(\alpha) = \int_0^\alpha \frac{d\alpha}{g(\alpha)} \quad (\text{A.5})$$

and  $t_{\text{gen}}$  is the generalized time defined as [Formula \(A.6\)](#):

$$t_{\text{gen}} = \int_0^t \exp\left(-\frac{E_a}{RT}\right) dt \quad (\text{A.6})$$

In principle, from [Formula \(A.1\)](#) and [Formula \(A.5\)](#) and given the functions  $f(\alpha)$  and  $g(\alpha)$ , the relationship between  $P$  and  $At_{\text{gen}}$  can be derived. The relationship between  $C$  and  $At_{\text{gen}}$  is dependent only on the functions  $f(\alpha)$  and  $g(\alpha)$ , i.e. the mechanism of the reaction and the relationship of  $P$  with  $\alpha$  and can be theoretically derived.

Similar theoretical relationships of derivative type can also be derived as shown in [Formula \(A.7\)](#):

$$\frac{dC}{dAt_{\text{gen}}} = \frac{df(\alpha)}{d\alpha} \frac{d\alpha}{dAt_{\text{gen}}} = \frac{df(\alpha)}{d\alpha} \frac{d\alpha}{dt} \frac{dt}{dAt_{\text{gen}}} = g(\alpha) \frac{df(\alpha)}{d\alpha} \quad (\text{A.7})$$

as [Formula \(A.8\)](#) applies:

$$\frac{dAt_{\text{gen}}}{dt} = A \exp\left(-\frac{E_a}{RT}\right) \quad (\text{A.8})$$

Thus, the theoretical relationship between  $dC/dAt_{\text{gen}}$  and  $C$  and between  $dC/dAt_{\text{gen}}$  and  $At_{\text{gen}}$  can be derived and depend only on the functions  $f(\alpha)$  and  $g(\alpha)$ . For some typical cases,  $dC/dAt_{\text{gen}}$  is shown as a function of  $At_{\text{gen}}$  in [Table A.1](#).

STANDARDSISO.COM : Click to view the full PDF of ISO 11358-3:2021