
Determination of the ultimate aerobic biodegradability of plastic materials in an aqueous medium — Method by measuring the oxygen demand in a closed respirometer

Évaluation de la biodégradabilité aérobie ultime des matériaux plastiques en milieu aqueux — Méthode par détermination de la demande en oxygène dans un respiromètre fermé

STANDARDSISO.COM : Click to view the full PDF of ISO 14851:2019



STANDARDSISO.COM : Click to view the full PDF of ISO 14851:2019



COPYRIGHT PROTECTED DOCUMENT

© ISO 2019

All rights reserved. Unless otherwise specified, or required in the context of its implementation, no part of this publication may be reproduced or utilized otherwise in any form or by any means, electronic or mechanical, including photocopying, or posting on the internet or an intranet, without prior written permission. Permission can be requested from either ISO at the address below or ISO's member body in the country of the requester.

ISO copyright office
CP 401 • Ch. de Blandonnet 8
CH-1214 Vernier, Geneva
Phone: +41 22 749 01 11
Fax: +41 22 749 09 47
Email: copyright@iso.org
Website: www.iso.org

Published in Switzerland

Contents

	Page
Foreword.....	iv
Introduction.....	vi
1 Scope.....	1
2 Normative references.....	1
3 Terms and definitions.....	1
4 Principle.....	3
5 Test environment.....	3
6 Reagents.....	3
6.1 Distilled or deionized water.....	4
6.2 Test medium.....	4
6.2.1 Standard test medium.....	4
6.2.2 Optimized test medium.....	4
6.3 Pyrophosphate solution.....	6
6.4 Carbon dioxide absorber.....	6
7 Apparatus.....	6
8 Procedure.....	6
8.1 Test material.....	6
8.2 Reference material.....	7
8.3 Preparation of the inoculum.....	7
8.4 Test.....	8
9 Calculation and expression of results.....	9
9.1 Calculation.....	9
9.2 Expression and interpretation of results.....	10
10 Validity of results.....	10
11 Test report.....	11
Annex A (informative) Theoretical oxygen demand (ThOD).....	12
Annex B (informative) Correction of BOD values for interference by nitrification.....	13
Annex C (informative) Principle of a closed manometric respirometer.....	15
Annex D (informative) Two-phase closed-bottle version of the respirometric test.....	17
Annex E (informative) Example of the determination of a carbon balance.....	20
Annex F (informative) Example of a determination of the amount of water-insoluble polymer remaining at the end of a biodegradation test and the molecular mass of the polymer.....	22
Annex G (informative) Example of the determination of the CO₂ absorbed in the absorbent.....	23
Bibliography.....	25

Foreword

ISO (the International Organization for Standardization) is a worldwide federation of national standards bodies (ISO member bodies). The work of preparing International Standards is normally carried out through ISO technical committees. Each member body interested in a subject for which a technical committee has been established has the right to be represented on that committee. International organizations, governmental and non-governmental, in liaison with ISO, also take part in the work. ISO collaborates closely with the International Electrotechnical Commission (IEC) on all matters of electrotechnical standardization.

The procedures used to develop this document and those intended for its further maintenance are described in the ISO/IEC Directives, Part 1. In particular the different approval criteria needed for the different types of ISO documents should be noted. This document was drafted in accordance with the editorial rules of the ISO/IEC Directives, Part 2 (see www.iso.org/directives).

Attention is drawn to the possibility that some of the elements of this document may be the subject of patent rights. ISO shall not be held responsible for identifying any or all such patent rights. Details of any patent rights identified during the development of the document will be in the Introduction and/or on the ISO list of patent declarations received (see www.iso.org/patents).

ISO (the International Organization for Standardization) is a worldwide federation of national standards bodies (ISO member bodies). The work of preparing International Standards is normally carried out through ISO technical committees. Each member body interested in a subject for which a technical committee has been established has the right to be represented on that committee. International organizations, governmental and non-governmental, in liaison with ISO, also take part in the work. ISO collaborates closely with the International Electrotechnical Commission (IEC) on all matters of electrotechnical standardization.

The procedures used to develop this document and those intended for its further maintenance are described in the ISO/IEC Directives, Part 1. In particular the different approval criteria needed for the different types of ISO documents should be noted. This document was drafted in accordance with the editorial rules of the ISO/IEC Directives, Part 2 (see www.iso.org/directives).

Attention is drawn to the possibility that some of the elements of this document may be the subject of patent rights. ISO shall not be held responsible for identifying any or all such patent rights. Details of any patent rights identified during the development of the document will be in the Introduction and/or on the ISO list of patent declarations received (see www.iso.org/patents).

Any trade name used in this document is information given for the convenience of users and does not constitute an endorsement.

For an explanation on the voluntary nature of standards, the meaning of ISO specific terms and expressions related to conformity assessment, as well as information about ISO's adherence to the World Trade Organization (WTO) principles in the Technical Barriers to Trade (TBT) see the following URL: www.iso.org/iso/foreword.html.

The committee responsible for this document is ISO/TC 61, *Plastics*, Subcommittee SC 14, *Environmental aspects*.

This second edition cancels and replaces the first edition (ISO 14851:1999), which has been technically revised. It also incorporates the Technical Corrigendum ISO 14851:1999/Cor.1:2005. The main changes compared to the previous edition are as follows:

- the footnotes have been renumbered;
- in [Annex C](#), errors in the key to [Figure C.1](#) have been corrected and minor improvements made to the figure itself;
- in scope and [Clause 8](#), soil and compost have been excluded for the inoculums used in this document;

- in 8.4, numbers of test flask for the test material and blank control have been changed from two to three;
- references in this document have been updated for latest active version;
- the Bibliography has been updated.

Any feedback or questions on this document should be directed to the user's national standards body. A complete listing of these bodies can be found at www.iso.org/members.html.

STANDARDSISO.COM : Click to view the full PDF of ISO 14851:2019

Introduction

With the increasing use of plastics, their recovery and disposal have become a major issue. As a first priority, recovery is promoted. Complete recovery of plastics, however, is difficult. For example, plastic litter, which comes mainly from consumers, is difficult to recover completely. Additional examples of plastics which are difficult to recover are fishing tackle, plastic microbeads in personal care products and water-soluble polymers. These plastic materials tend to leak from closed waste-management cycles into the environment. Biodegradable plastics are now emerging as one of the options available to solve such environmental problems. Plastic materials, such as products or packaging, which are sent to composting facilities are expected to be potentially biodegradable. Therefore, it is very important to determine the potential biodegradability of such materials and to obtain an indication of their biodegradability in natural environments.

STANDARDSISO.COM : Click to view the full PDF of ISO 14851:2019

Determination of the ultimate aerobic biodegradability of plastic materials in an aqueous medium — Method by measuring the oxygen demand in a closed respirometer

WARNING — Sewage, activated sludge, soil and compost may contain potentially pathogenic organisms. Therefore, appropriate precautions should be taken when handling them. Toxic test compounds and those whose properties are unknown should be handled with care.

1 Scope

This document specifies a method, by measuring the oxygen demand in a closed respirometer, for the determination of the degree of aerobic biodegradability of plastic materials, including those containing formulation additives. The test material is exposed in an aqueous medium under laboratory conditions to an inoculum from activated sludge.

If an unadapted activated sludge is used as the inoculum, the test simulates the biodegradation processes which occur in a natural aqueous environment; if a mixed or pre-exposed inoculum is used, the method is used to investigate the potential biodegradability of a test material.

The conditions used in this document do not necessarily correspond to the optimum conditions allowing maximum biodegradation to occur, but this document is designed to determine the potential biodegradability of plastic materials or give an indication of their biodegradability in natural environments.

The method enables the assessment of the biodegradability to be improved by calculating a carbon balance (optional, see [Annex E](#)).

The method applies to the following materials.

- Natural and/or synthetic polymers, copolymers or mixtures thereof.
- Plastic materials which contain additives such as plasticizers, colorants or other compounds.
- Water-soluble polymers.
- Materials which, under the test conditions, do not inhibit the microorganisms present in the inoculum. Inhibitory effects can be determined using an inhibition control or by another appropriate method (see, for example, ISO 8192^[2]). If the test material is inhibitory to the inoculum, a lower test concentration, another inoculum or a pre-exposed inoculum can be used.

2 Normative references

The following documents are referred to in the text in such a way that some or all of their content constitutes requirements of this document. For dated references, only the edition cited applies. For undated references, the latest edition of the referenced document (including any amendments) applies.

ISO 8245, *Water quality — Guidelines for the determination of total organic carbon (TOC) and dissolved organic carbon (DOC)*

3 Terms and definitions

For the purposes of this document, the following terms and definitions apply.

ISO and IEC maintain terminological databases for use in standardization at the following addresses:

- ISO Online browsing platform: available at <https://www.iso.org/obp>
- IEC Electropedia: available at <http://www.electropedia.org/>

3.1

ultimate aerobic biodegradation

breakdown of an organic compound by microorganisms in the presence of oxygen into carbon dioxide, water and mineral salts of any other elements present (mineralization) plus new biomass

3.2

activated sludge

biomass produced in the aerobic treatment of waste water by the growth of bacteria and other microorganisms in the presence of dissolved oxygen

3.3

concentration of suspended solids in an activated sludge

amount of solids obtained by filtration or centrifugation of a known volume of *activated sludge* (3.2) and drying at about 105 °C to constant mass

3.4

biochemical oxygen demand

BOD

mass concentration of the dissolved oxygen consumed under specified conditions by the aerobic biological oxidation of a chemical compound or organic matter in water

Note 1 to entry: It is expressed as milligrams of oxygen uptake per milligram or gram of test compound.

3.5

theoretical oxygen demand

ThOD

theoretical maximum amount of oxygen required to oxidize a chemical compound completely, calculated from the molecular formula

Note 1 to entry: It is expressed as milligrams of oxygen uptake per milligram or gram of test compound.

3.6

total organic carbon

TOC

amount of carbon bound in an organic compound

3.7

dissolved organic carbon

DOC

part of the organic carbon in water which cannot be removed by specified phase separation, for example by centrifugation at 40 000 m·s⁻² for 15 min or by membrane filtration using membranes with pores of 0,2 µm to 0,45 µm diameter

3.8

lag phase

time, measured in days, from the start of a test until adaptation and/or selection of the degrading microorganisms is achieved and the degree of biodegradation of a chemical compound or organic matter has increased to about 10 % of the *maximum level of biodegradation* (3.9)

3.9

maximum level of biodegradation

degree of biodegradation, measured in per cent, of a chemical compound or organic matter in a test, above which no further biodegradation takes place during the test

3.10**biodegradation phase**

time, measured in days, from the end of the *lag phase* (3.8) of a test until about 90 % of the *maximum level of biodegradation* (3.9) has been reached

3.11**plateau phase**

time, measured in days, from the end of the *biodegradation phase* (3.10) until the end of a test

3.12**pre-exposure**

pre-incubation of an inoculum in the presence of the chemical compound or organic matter under test, with the aim of enhancing the ability of the inoculum to biodegrade the test material by adaptation and/or selection of the microorganisms

3.13**pre-conditioning**

pre-incubation of an inoculum under the conditions of the subsequent test in the absence of the chemical compound or organic matter under test, with the aim of improving the test by acclimatization of the microorganisms to the test conditions

4 Principle

The biodegradability of a plastic material is determined using aerobic microorganisms in an aqueous system. The test mixture contains an inorganic medium, the organic test material (the sole source of carbon and energy) with a concentration between 100 mg/l and 2 000 mg/l of organic carbon, and activated sludge as the inoculum. The mixture is stirred in closed flasks in a respirometer for a period not exceeding 2 months. The carbon dioxide evolved is absorbed in a suitable absorber in the headspace of the flasks. The consumption of oxygen (BOD) is determined, for example by measuring the amount of oxygen required to maintain a constant volume of gas in the respirometer flasks, or by measuring the change in volume or pressure (or a combination of the two) either automatically or manually. An example of a respirometer is given in [Annex C](#). Alternatively, the two-phase closed-bottle version described in ISO 10708[3] may be used (see [Annex D](#)).

The level of biodegradation is determined by comparing the BOD with the theoretical amount (ThOD) and expressed in per cent. The influence of possible nitrification processes on the BOD has to be considered. The test result is the maximum level of biodegradation determined from the plateau phase of the biodegradation curve. Optionally, a carbon balance may be calculated to give additional information on the biodegradation (see [Annex E](#)). Moreover, also the absorbed carbon dioxide in the adsorber at the end of the test may be determined to give additional information on the biodegradation (see [Annex G](#)).

Unlike ISO 9408[6], which is used for a variety of organic compounds, this document is specially designed for the determination of the biodegradability of plastic materials. The special requirements necessary affect the choice of the inoculum and the test medium, and there is the possibility of improving the evaluation of the biodegradability by calculating a carbon balance.

5 Test environment

Incubation shall take place in the dark or in diffuse light in an enclosure which is free from vapours inhibitory to microorganisms and which is maintained at a constant temperature, preferably between 20 °C and 25 °C, to an accuracy of ± 1 °C, or at any other appropriate temperature depending on the inoculum used and the environment to be assessed.

6 Reagents

Use only reagents of recognized analytical grade.

6.1 Distilled or deionized water

Distilled or deionized water, free of toxic substances (copper in particular) and containing less than 2 mg/l of DOC.

6.2 Test medium

Depending on the purpose of the test, different test media may be used. For example, if simulating a natural environment use the standard test medium (6.2.1). If a test material is used at higher concentrations, use the optimized test medium (6.2.2) with higher buffering capacity and nutrient concentrations.

6.2.1 Standard test medium

6.2.1.1 Solution A

Dissolve the following in water (6.1) and make up to 1 000 ml.

anhydrous potassium dihydrogen phosphate (KH_2PO_4)	8,5 g
anhydrous dipotassium hydrogen phosphate (K_2HPO_4)	21,75 g
disodium hydrogen phosphate dihydrate ($\text{Na}_2\text{HPO}_4 \cdot 2\text{H}_2\text{O}$)	33,4 g
ammonium chloride (NH_4Cl)	0,5 g

The correct composition of the solution can be checked by measuring the pH, which should be 7,4.

6.2.1.2 Solution B

Dissolve 22,5 g of magnesium sulfate heptahydrate ($\text{MgSO}_4 \cdot 7\text{H}_2\text{O}$) in water (6.1) and make up to 1 000 ml.

6.2.1.3 Solution C

Dissolve 36,4 g of calcium chloride dihydrate ($\text{CaCl}_2 \cdot 2\text{H}_2\text{O}$) in water (6.1) and make up to 1 000 ml.

6.2.1.4 Solution D

Dissolve 0,25 g of iron(III) chloride hexahydrate ($\text{FeCl}_3 \cdot 6\text{H}_2\text{O}$) in water (6.1) and make up to 1 000 ml.

Prepare this solution freshly before use to avoid precipitation, or add a drop of concentrated hydrochloric acid (HCl) or a drop of 0,4 g/l aqueous solution of ethylenediaminetetraacetic acid (EDTA).

6.2.1.5 Preparation

To prepare 1 l of test medium, add the following to about 500 ml of water (6.1),

- 10 ml of solution A;
- 1 ml of each of solutions B to D.

Make up to 1 000 ml with water (6.1). Prepare the test medium freshly before use. The solutions A up to C may be stored up to 6 months in the dark at room temperature.

6.2.2 Optimized test medium

This optimized medium is highly buffered and contains more inorganic nutrients. This is necessary to keep the pH constant in the system during the test, even at high concentrations of the test material. The

medium contains about 2 400 mg/l of phosphorus and 50 mg/l of nitrogen and is therefore suitable for concentrations in the test material of up to 2 000 mg/l of organic carbon. If higher test-material concentrations are used, increase the nitrogen content to keep the C:N ratio at about 40:1.

6.2.2.1 Solution A

Dissolve the following in water (6.1) and make up to 1 000 ml.

anhydrous potassium dihydrogen phosphate (KH_2PO_4)	37,5 g
disodium hydrogen phosphate dihydrate ($\text{Na}_2\text{HPO}_4 \cdot 2\text{H}_2\text{O}$)	87,3 g
ammonium chloride (NH_4Cl)	2,0 g

6.2.2.2 Solution B

Dissolve 22,5 g of magnesium sulfate heptahydrate ($\text{MgSO}_4 \cdot 7\text{H}_2\text{O}$) in water (6.1) and make up to 1 000 ml.

6.2.2.3 Solution C

Dissolve 36,4 g of calcium chloride dihydrate ($\text{CaCl}_2 \cdot 2\text{H}_2\text{O}$) in water (6.1) and make up to 1 000 ml.

6.2.2.4 Solution D

Dissolve 0,25 g of iron(III) chloride hexahydrate ($\text{FeCl}_3 \cdot 6\text{H}_2\text{O}$) in water (6.1) and make up to 1 000 ml (see second paragraph of 6.2.1.4).

6.2.2.5 Solution E (trace-element solution, optional)

Dissolve in 10 ml of aqueous HCl solution (25 %, 7,7 mol/l), in the following sequence:

70 mg of ZnCl_2 , 100 mg of $\text{MnCl}_2 \cdot 4\text{H}_2\text{O}$, 6 mg of H_3BO_3 , 190 mg of $\text{CoCl}_2 \cdot 6\text{H}_2\text{O}$, 3 mg of $\text{CuCl}_2 \cdot 2\text{H}_2\text{O}$, 240 mg of $\text{NiCl}_2 \cdot 6\text{H}_2\text{O}$, 36 mg of $\text{Na}_2\text{MoO}_4 \cdot 2\text{H}_2\text{O}$, 33 mg of $\text{Na}_2\text{WO}_4 \cdot 2\text{H}_2\text{O}$ and 26 mg of $\text{Na}_2\text{SeO}_3 \cdot 5\text{H}_2\text{O}$.

Make up to 1 000 ml with water (6.1).

6.2.2.6 Solution F (vitamin solution, optional)

Dissolve in 100 ml of water (6.1) 0,6 mg of biotine, 2,0 mg of niacinamide, 2,0 mg of *p*-aminobenzoate, 1,0 mg of pantothenic acid, 10,0 mg of pyridoxal hydrochloride, 5,0 mg of cyanocobalamine, 2,0 mg of folic acid, 5,0 mg of riboflavin, 5,0 mg of DL-thioctic acid and 1,0 mg of thiamine dichloride or use a solution of 15 mg of yeast extract in 100 ml of water (6.1). Filter the solution for sterilization using membrane filters (see 7.4).

Solutions E and F are optional and are not required if a sufficient concentration of the inoculum is used, for example, activated sludge. It is recommended that 1 ml portions be prepared and kept refrigerated until use.

6.2.2.7 Preparation

To prepare 1 l of test medium, add, to about 800 ml of water (6.1):

- 100 ml of solution A;
- 1 ml of each of solutions B to D and, optionally, E and F.

Make up to 1 000 ml with water (6.1) and measure the pH.

The correct composition of the test medium can be checked by measuring the pH, which should be $7,0 \pm 0,2$.

6.3 Pyrophosphate solution

Dissolve 2,66 g of anhydrous sodium pyrophosphate ($\text{Na}_4\text{P}_2\text{O}_7$) in water (6.1) and make up to 1 000 ml.

6.4 Carbon dioxide absorber

Preferably soda lime pellets or another suitable absorbant.

7 Apparatus

Ensure that all glassware is thoroughly cleaned and, in particular, free from organic or toxic matter.

Required is usual laboratory equipment, plus the following.

7.1 Closed respirometer, including test vessels (glass flasks) fitted with stirrers and all other necessary equipment, and located in a constant-temperature room or in a thermostated apparatus (e.g. water-bath). For an example, see [Annex C](#).

NOTE Any respirometer able to determine with sufficient accuracy the biochemical oxygen demand is suitable, preferably an apparatus which measures and replaces automatically and continuously the oxygen consumed so that no oxygen deficiency and no inhibition of the microbial activity occurs during the degradation process. Instead of an ordinary respirometer, the two-phase closed-bottle version can be used (see [Annex D](#)).

7.2 Analytical equipment for measuring total organic carbon (TOC) and dissolved organic carbon (DOC) according to ISO 8245.

7.3 Analytical equipment for measuring nitrate and nitrite concentrations.

A qualitative test is recommended first to decide if any nitrification has occurred. If there is evidence of nitrate/nitrite in the medium, a quantitative determination using a suitable method (for example ion chromatography) is required.

7.4 Centrifuge, or filtration device with membrane filters (0,45 μm pore size) which neither adsorb nor release organic carbon significantly.

7.5 Analytical balance (usual laboratory equipment).

7.6 pH meter (usual laboratory equipment).

8 Procedure

8.1 Test material

The test material shall be of known mass and contain sufficient carbon to yield a BOD that can be adequately measured by the respirometer used. Calculate from the chemical formula or determine by elemental analysis the ThOD (see [Annex A](#)) and the TOC (using, for example, ISO 8245). Use a test-material concentration of at least 100 mg/l, corresponding to a ThOD of about 170 mg/l or a TOC of about 60 mg/l. Use lower concentrations only if the sensitivity of the respirometer is adequate. The maximum amount of test material is limited by the oxygen supply to the respirometer and the test medium used. When using the optimized test medium (6.2.2), the test-material concentration shall be such that the TOC does not exceed about 2 000 mg/l, i.e. a C:N ratio of about 40:1. If higher concentrations are to be tested, increase the amount of nitrogen in the test medium.

If biodegradation processes in natural environments are to be simulated, the use of the standard test medium and a test-material concentration of 100 mg/l are recommended.

The test material should preferably be used in powder form, but it may also be introduced as films, pieces, fragments or shaped articles. The form and shape of the test material may influence its biodegradability. Similar shapes should preferably be used if different kinds of plastic material are to be compared. If the test material is used in the form of a powder, particles of known, narrow size distribution should be used. A particle-size distribution with the maximum at 250 µm diameter is recommended. Also, the size of the test equipment used may depend on the form of the test material. It should be ascertained that no substantial mechanical aberrations occur due to the test conditions, for example due to the type of stirring mechanism used. Processing of the test material (e.g. the use of powder in the case of composites) should not influence significantly the degradation behaviour of the material. Optionally, record the hydrogen, oxygen, nitrogen, phosphorus and sulfur contents and the molecular mass of a polymeric test material, using for example liquid exclusion chromatography (see, for example, ASTM D 3536–91^[1] or any other applicable standard method). Preferably, plastic materials without additives such as plasticizers should be tested. When the material does contain such additives, information on their biodegradability will be needed to assess the biodegradability of the polymeric material itself.

For details on how to handle poorly water-soluble compounds, see ISO 10634^[7].

8.2 Reference material

Use aniline and/or a well-defined biodegradable polymer (for example microcrystalline cellulose powder, ashless cellulose filters or poly-β-hydroxybutyrate) as a reference material. If possible, the TOC, form and size should be comparable to that of the test material.

As a negative control, a non-biodegradable polymer (e.g. polyethylene) in the same form as the test material can optionally be used.

8.3 Preparation of the inoculum

Activated sludge from a sewage-treatment plant treating predominantly domestic sewage is a suitable source of the inoculum. It is obtained from an active aerobic environment and is available over a wide geographical area in which a broad range of plastic materials needs to be tested. When biodegradation in a specific waste-treatment system is to be determined, collect the inoculum from that environment.

Take a sample of activated sludge collected from a well-operated sewage-treatment plant or a laboratory plant handling predominantly domestic sewage. Mix well, keep the sample under aerobic conditions and use preferably on the day of collection (at least within 72 h).

Before use, determine the concentration of suspended solids (use, for example, ISO 11923^[4]). If necessary, concentrate the sludge by settling so that the volume of sludge added to the test assay is minimal. Add a suitable volume to obtain suspended solids in the range 30 mg/l to 1 000 mg/l in the final mixture.

If the endogenous respiration of the inoculum is too high, stabilize the inoculum by aeration before use. Harmonize the test temperature with the inoculum used.

It may be useful to determine the colony-forming units (cfu) of the inoculum used. The test mixture should preferably contain about 10³ cfu/ml to 10⁶ cfu/ml.

When biodegradation processes in a natural environment are to be simulated or when a carbon balance determination (see [Annex E](#)) is to be carried out, an inoculum concentration of 30 mg/l suspended solids is recommended. As solid matter can interfere with the carbon balance determination, the following procedure for preparing the inoculum is recommended. Take 500 ml of the activated sludge and homogenize for 2 min at medium speed in a blender or in a suitable high-speed mixer. Allow to settle until the supernatant liquid contains no significant amounts of suspended matter, but in any case for at least 30 min. Decant a sufficient volume of the supernatant liquid and add it to the test flasks to

obtain a concentration of a volume fraction of 1 % to a volume fraction of 5 % in the test medium. Avoid carrying over sludge particles.

An inoculum may be pre-conditioned, but normally no pre-exposed inoculum should be used, especially in the case of standard tests simulating biodegradation behaviour in natural environments. Depending on the purpose of the test, a pre-exposed inoculum may also be used, provided this is clearly stated in the test report (e.g. per cent biodegradation = x %, using pre-exposed inocula) and the method of pre-exposure detailed in the test report. Pre-exposed inocula can be obtained from suitable laboratory biodegradation tests (see ISO/TR 15462[8]) conducted under a variety of conditions or from samples collected from locations where relevant environmental conditions exist (e.g. contaminated areas or industrial treatment plants).

8.4 Test

Provide a number of flasks, so that the test includes at least the following:

- a) three test flasks for the test material (symbol F_T);
- b) three test flasks for the blank control (symbol F_B);
- c) two test flasks for checking the inoculum activity using a reference material (symbol F_C).

And, if required:

- d) One flask for checking for possible abiotic degradation or non-biological change in the test material such as by hydrolysis (symbol F_S). The test solution in F_S shall be sterilized, for example by autoclaving or by the addition of a suitable inorganic toxic compound to prevent microbial activity. Use, for example, 5 ml/l of a solution containing 10 g/l of mercury(II) chloride ($HgCl_2$). Add the same amount of the toxic substance during the test if required.
- e) One flask as a negative control (symbol F_N) using a non-biodegradable polymeric substance (e.g. polyethylene) in the same form as the test material.
- f) One flask for checking the possible inhibiting effect of the test material on microbial activity (symbol F_I). Take care that the ratio of carbon in the test and reference material to nitrogen in the medium is at least about C:N = 40:1. Add nitrogen if required.

Add appropriate amounts of the test medium (6.2) and the inoculum (see 8.3) to the test flasks as indicated in Table 1.

Measure the pH in the flasks and adjust to 7 if necessary. Add carbon dioxide absorber (6.4) to the absorber compartments of the respirometer (see Annex C). Add the test material (see 8.1), the reference material and the material for the negative control (see 8.2) to the respective flasks as indicated in Table 1. If a carbon balance is to be run (see Annex E), remove a known sufficient volume of the inoculated test medium from each flask or from additional separate flasks for DOC and biomass determination at the beginning and the end of the incubation period. Consider the removed volume when adjusting the final volume or when calculating the test results.

Place the flasks in a constant-temperature environment (see Clause 5) and allow all vessels to reach the desired temperature. Make any necessary connections, seal the flasks, place them in the respirometer and start the stirrer.

Take the necessary readings on the manometers (if manual) and verify that the recorder of oxygen consumption is functioning properly (automatic respirometer). As an alternative, the two-phase closed-bottle version described in Annex D may be used.

Table 1 — Final distribution of test and reference materials

Flask		Test material	Reference material	Inoculum
F _T	Test	+	-	+
F _T	Test	+	-	+
F _T	Test	+	-	+
F _B	Blank	-	-	+
F _B	Blank	-	-	+
F _B	Blank	-	-	+
F _C	Inoculum check	-	+	+
F _C	Inoculum check	-	+	+
F _S	Abiotic degradation check (optional)	+	-	-
F _I	Inhibition control (optional)	+	+	+
F _N	Negative control (optional)	-	+	+

When a constant level of BOD is attained (plateau phase reached) and no further biodegradation is expected, the test is considered to be completed. The test period should not typically exceed 2 months. However, if significant biodegradation is still observed and the plateau phase has not been reached after this length of time, then the test may be extended, but not longer than 6 months. In the case of long test durations, special attention shall be paid to the technical system (e.g. tightness of the test vessels and connections).

At the end of the test, measure the pH and determine the concentrations of nitrate and nitrite immediately in flasks F_T (see paragraph below), or take suitably preserved samples. Use the values to correct the calculated degree of biodegradation for nitrification (see [Annex B](#)).

Allylthiourea can only inhibit nitrification during short incubation periods, as it is biodegradable. Therefore, addition of allylthiourea to prevent nitrification is not recommended. Experience shows, however, that with low inoculum concentrations (about a volume fraction of 1 %) nitrification will not occur, even during long incubation periods, when no inhibitor is used.

Optionally determine additionally the absorbed carbon dioxide in the adsorber at the end of the test (see [Annex G](#)).

9 Calculation and expression of results

9.1 Calculation

Read the oxygen consumption values for each flask, using the method given by the manufacturer for the appropriate type of respirometer. Calculate the specific biochemical oxygen demand (BOD_S) of the test material as the difference between oxygen consumption in the test flasks F_T and the blanks F_B divided by the concentration of the test material, using [Formula \(1\)](#):

$$\text{BOD}_S = \frac{\text{BOD}_t - \text{BOD}_{Bt}}{\rho_{TC}} \quad (1)$$

where

- BOD_S is the specific BOD, in milligrams per gram of test material;
- BOD_t is the BOD of the flasks F_T containing test material at time *t*, in milligrams per litre;
- BOD_{Bt} is the BOD of the blank F_B at time *t*, in milligrams per litre;
- ρ_{TC} is the concentration of the test material in the reaction mixture of flask F_T, in grams per litre.

Calculate the percentage biodegradation *D_t* as the ratio of the specific biochemical oxygen demand to the theoretical oxygen demand (ThOD, in milligrams per gram of test material), using [Formula \(2\)](#):

$$D_t = \frac{\text{BOD}_S}{\text{ThOD}} \times 100 \quad (2)$$

Calculate in the same way the BOD and percentage biodegradation of the reference material F_C and, if included, the abiotic degradation check F_S, the inhibition control F_I and the negative control F_N.

For calculation of the ThOD, see [Annex A](#). If significant concentrations of nitrite and nitrate are determined, consider the oxygen demand due to nitrification (see [Annex B](#)). If a carbon balance is to be calculated, use information given in [Annex E](#).

9.2 Expression and interpretation of results

Compile a table of the BOD values measured and the percentages of biodegradation for each measurement interval and each test flask. For each vessel, plot a BOD curve and a biodegradation curve in per cent as a function of time. If comparable results are obtained for the triplicate flasks, a mean curve may be plotted.

The maximum level of biodegradation determined as the mean value of the plateau phase of the biodegradation curve or the highest value, such as when the curve decreases or, further on, slowly increases in the plateau phase, characterizes the degree of biodegradation of the test material. If a carbon balance has been determined, the result of this determination characterizes the total degree of biodegradation.

The wettability and the shape of the test material may influence the result obtained, and hence the test procedure may be limited to comparing plastic materials of similar chemical structure.

Information on the toxicity of the test material may be useful in the interpretation of test results showing a low biodegradability.

10 Validity of results

The test is considered valid if:

- a) the degree of biodegradation of the reference material (inoculum check F_C) is > 60 % at the end of the test;
- b) the BOD of the blank F_B at the end of the test does not exceed an upper limiting value obtained by experience (this value depends on the amount of inoculum and is, for example, in the case of 30 mg/l dry matter, about 60 mg/l as interlaboratory tests have shown);
- c) the BOD values of the three blanks F_B and of the three test flasks F_T are within 20 % of the mean at the plateau phase or at the end of the test.

If in flask F_I (inhibition check, if included) the percentage biodegradation is $<25\%$ and no significant degradation of the test material is observed, it can be assumed that the test material is inhibitory.

If in flask F_S (abiotic degradation check, if included) a significant amount ($>10\%$) of BOD is observed, abiotic degradation processes may have taken place.

If flask F_N (negative control) was included, no significant amount of BOD shall be observed.

If these criteria are not fulfilled, repeat the test using another pre-conditioned or pre-exposed inoculum.

11 Test report

The test report shall contain at least the following information:

- a) a reference to this document, i.e. ISO 14851:2019;
- b) all information necessary to identify the test and reference materials, including their TOC, ThOD, chemical composition and formula (if known), shape, form and amount/concentration in the samples tested;
- c) the main test parameters, including test volume, test medium used, incubation temperature and final pH;
- d) the source and amount of the inoculum used, including details of any pre-exposure;
- e) the analytical techniques used, including the principle of the respirometer and the TOC and nitrate/nitrite determinations;
- f) all the test results obtained for the test and reference materials (in tabular and graphical form), including the measured BOD, the percentage biodegradation values, the respective curves of these parameters against time and the nitrate/nitrite concentrations;
- g) the duration of the lag phase, biodegradation phase and maximum level of degradation, as well as the total test duration;

and, optionally, if run or determined:

- h) the results of the abiotic degradation check F_S , the inhibition control F_I and the negative control F_N ;
- i) the results of the carbon balance determination, including for example:
 - 1) the amount of carbon in the test material oxidized to carbon dioxide, estimated from the degree of biodegradation based on the BOD,
 - 2) the increase in DOC in the test medium during the incubation period due to water-soluble substances,
 - 3) the increase in organic carbon in the biomass during the test,
 - 4) the carbon content of the residual polymers at the end of the test,
 - 5) the sum of all the carbon measured, expressed as a percentage of the carbon introduced as the test material;
- j) the colony-forming units (cfu/ml) in the inoculated test mixtures;
- k) any other relevant data (e.g. initial molecular mass of the sample, molecular mass of the residual polymer).

Annex A (informative)

Theoretical oxygen demand (ThOD)

A.1 Calculation of the ThOD

The theoretical oxygen demand (ThOD) of a substance $C_cH_hCl_{cl}N_nS_sP_pNa_{na}O_o$ of relative molecular mass M_r can be calculated, if the elemental composition is known or can be determined by elemental analysis, using [Formula \(A.1\)](#):

$$\text{ThOD} = \frac{16[2c + 0,5(h - cl - 3n) + 3s + 2,5p + 0,5na - o]}{M_r} \quad (\text{A.1})$$

This calculation assumes that carbon is converted to CO_2 , hydrogen to H_2O , phosphorus to P_2O_5 , sulfur to an oxidation state of +6 and halogens eliminated as hydrogen halides. The oxidation of N, P and S has to be checked by analysis. The calculation also assumes that nitrogen is released as ammonium. For the influence of nitrification, see [Annex B](#).

Express the ThOD in milligrams per gram of substance or in milligrams per milligram of substance.

A.2 Example: Poly(β -hydroxybutyric acid) (PHB)

Summary formula: $\text{C}_4\text{H}_6\text{O}_2$, $c = 4$, $h = 6$, $o = 2$; relative molecular mass $M_r = 86$. See [Formula \(A.2\)](#):

$$\text{ThOD} = \frac{16[2 \times 4 + 0,5 \times 6 - 2]}{86} \quad (\text{A.2})$$

ThOD = 1,674 4 mg/mg PHB = 1 674,4 mg/g PHB.

NOTE PHB is a polymer consisting of the β -hydroxybutyrate monomer. For polymerization (ester formation), water is removed, so that the summary formula for PHB is equivalent to that of the monomer minus one H_2O , which is eliminated in the chemical reaction.

A.3 Example: Blend of polyethylene/starch/glycerol

Component	Formula	ThOD mg/g	Amount of component		ThOD mg/flask
			%	mg/flask	
Polyethylene	$(\text{C}_2\text{H}_4)_n$	3 400	50	500	1 700
Starch	$(\text{C}_6\text{H}_{10}\text{O}_5)_n$	1 190	40	400	476
Glycerol	$\text{C}_3\text{H}_8\text{O}_3$	1 200	10	100	120
Total blend			100	1 000	2 296

Annex B (informative)

Correction of BOD values for interference by nitrification

B.1 Influence of nitrification

BOD values can be influenced by nitrification. They need to be corrected if serious errors are to be avoided in the calculation of degrees of biodegradation based on the oxidation of the carbon in a nitrogen-containing test material.

Errors in the case of nitrogen-free substances are normally negligible, because the oxidation of the ammonium in the medium is taken into account by the subtraction of the blank.

Ammonium salts and nitrogen-containing test compounds can be oxidized to nitrite or nitrate during the incubation period of a biodegradation test. Since the reactions are sequential (carried out by different bacterial species), it is possible for the nitrite concentration to increase or decrease. In the latter case, an equivalent concentration of nitrate is formed. The chemical reactions follow [Formulae \(B.1\) to \(B.3\)](#):



Overall:



From these formulae, it can be concluded that:

- for the oxidation of 2 moles (28 g) of ammonia nitrogen (added as NH_4Cl with the inorganic medium) to nitrite, 3 moles (96 g) of oxygen (BOD_{N02}) are needed, resulting in a factor of 3,43 (96/28) mg of oxygen demand per mg of nitrogen;
- for the oxidation of 2 moles (28 g) of ammonia nitrogen to nitrate, 4 moles (128 g) of oxygen (BOD_{N03}) are needed, resulting in a factor of 4,57 (128/28) mg of oxygen demand per mg of nitrogen.

The amount of nitrification can be determined by measuring the nitrate and nitrite concentrations at the end of the test in the medium in flasks F_T . A qualitative test is recommended first to decide if any nitrification has occurred. If there is evidence of nitrate or nitrite, a quantitative determination is required.

The part of the BOD deriving from nitrogen oxidation at the end of the test, BOD_N , is calculated, in milligrams per litre, using [Formula \(B.4\)](#):

$$\text{BOD}_N = (\rho_{\text{N03}} \times 4,57) + (\rho_{\text{N02}} \times 3,43) \quad (\text{B.4})$$

where

ρ_{NO_3} is the measured concentration of nitrate nitrogen in flasks F_T at the end of the test, in milligrams per litre;

ρ_{NO_2} is the measured concentration of nitrite nitrogen in flasks F_T at the end of the test, in milligrams per litre;

4,57 is the factor for the oxygen demand for the formation of nitrate;

3,43 is the factor for the oxygen demand for the formation of nitrite.

The part of the BOD deriving from carbon oxidation at the end of the test, BOD_C , is calculated, in milligrams per litre, using [Formula \(B.5\)](#):

$$\text{BOD}_C = \text{BOD}_G - \text{BOD}_N - \text{BOD}_{Bt} \quad (\text{B.5})$$

where

BOD_G is the measured BOD of flasks F_T at the end of the test, in milligrams per litre;

BOD_{Bt} is the BOD of the blank F_B at the end of the test, in milligrams per litre.

BOD_C corresponds to BOD_t and is used for calculating BODs and D_t [see [Formulae \(1\)](#) and [\(2\)](#) in [9.1](#)].

B.2 Example

Test substance *p*-aminobenzoic acid 2-ethylhexyl ester at a concentration of 100 mg/l in F_T .

ThOD	239 mg/l
Measured BOD_t at the end of the test	199 mg/l
Measured blank BOD_{Bt}	8 mg/l
D_t without correction for nitrification	80 %
Nitrate at the end of the test 15 mg/l	$\rho_{\text{NO}_3} = 3,5 \text{ mg/l}$
Nitrite at the end of the test 1 mg/l	$\rho_{\text{NO}_2} = 0,3 \text{ mg/l}$
BOD_N at the end of the test	17 mg/l
BOD_C	174 mg/l
D_t corrected for nitrification	73 %


Annex C (informative)

Principle of a closed manometric respirometer

The respirometer is set up in a temperature-controlled environment (e.g. a water bath) and contains test vessels each fitted with a magnetic stirrer rod and a container in the headspace for a CO₂ absorber, a coulometric oxygen production unit, a manometer, magnetic stirrers for each vessel, and an external monitoring device and recorder (printer, plotter or computer). The test vessels are filled to about one-third of their volume with the test mixture. See [Figure C.1](#).

Continuous stirring guarantees an equilibrium of oxygen between the aqueous and the gaseous phase. If biodegradation takes place, the microorganisms consume oxygen and produce carbon dioxide which is totally absorbed. The total pressure in the vessels decreases. The pressure drop is detected by the manometer and used to initiate the electrolytic generation of oxygen. When the original pressure is re-established, electrolysis is stopped and the quantity of electricity used, which is proportional to the oxygen consumption, is continuously measured and utilized to indicate the oxygen consumption, in mg/l BOD, on the recorder.

STANDARDSISO.COM : Click to view the full PDF of ISO 14851:2019

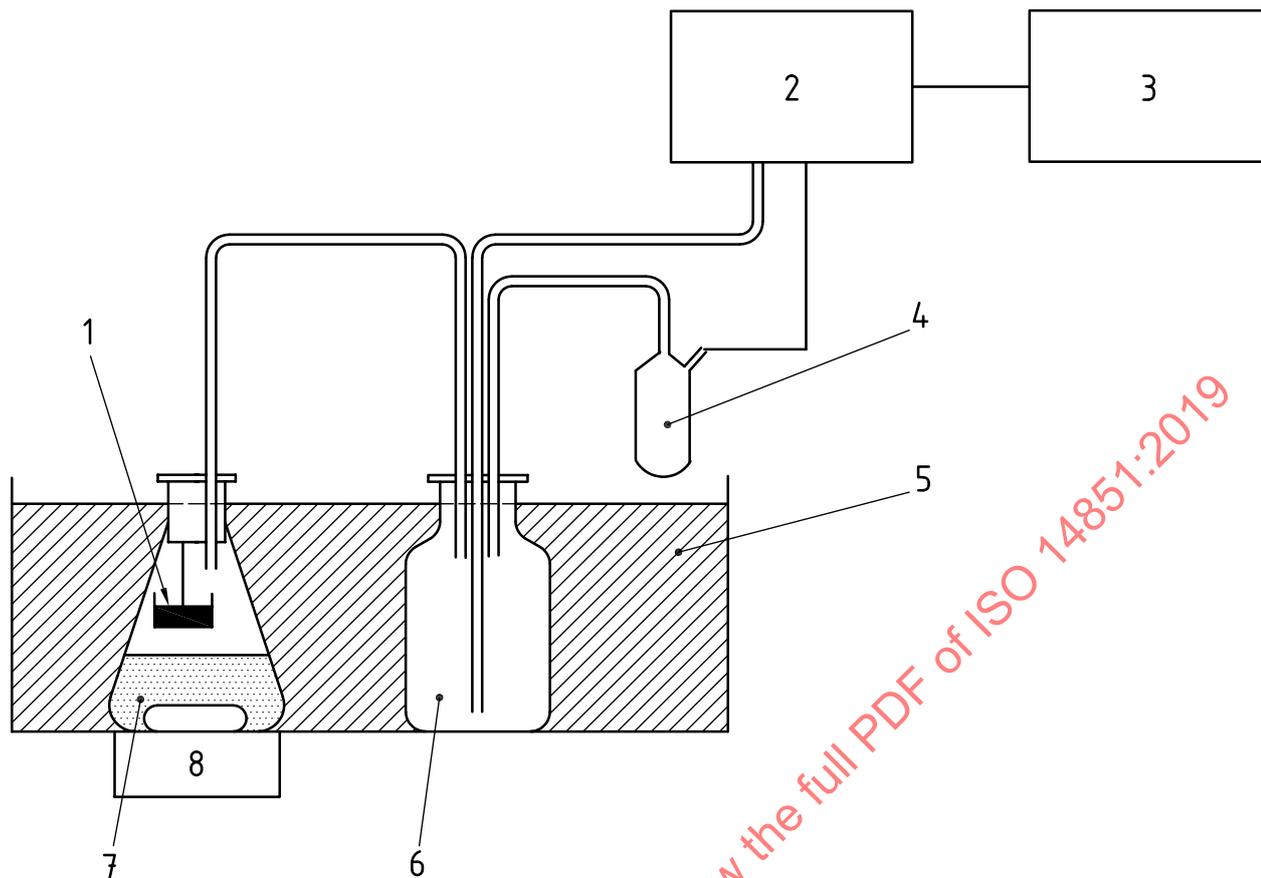

Key

- 1 CO₂ absorber
- 2 monitor
- 3 printer, plotter or computer
- 4 manometer
- 5 water bath
- 6 oxygen-producing unit
- 7 test flask
- 8 stirrer

Figure C.1 — Schematic drawing of a manometric respirometer

Annex D (informative)

Two-phase closed-bottle version of the respirometric test

D.1 Principle

This version can be used as an alternative, for example, if no respirometer is available. The inoculated medium and the test and reference materials are shaken or stirred at 20 °C to 25 °C in closed bottles containing known volumes of aqueous medium and air to ensure steady-state oxygen partitioning between the aqueous and the gas phase. The biodegradation is followed by means of regular measurements of the dissolved-oxygen concentration in the aqueous phase. The total oxygen uptake in the test flasks is calculated from the difference in the measured dissolved-oxygen concentrations in the blank and test flasks divided by the oxygen saturation value under normal conditions and multiplied by the total oxygen content originally present in the aqueous and gas phases. Biodegradability is calculated as the total oxygen uptake divided by the theoretical oxygen demand (ThOD) and expressed as a percentage.

D.2 Special apparatus

D.2.1 Incubation bottles: gas-tight bottles, such as narrow-necked flasks with volumes of 200 ml to 300 ml and with suitable stoppers (e.g. ground-glass stoppers, butyl-rubber stoppers or screw caps), providing shielding from the light (e.g. made of brown glass). Stopper clamps are recommended. Mark each bottle with waterproof markings. If oxygen electrodes with mounted stirrers are not used, provide the bottles with a magnetic stirrer with a PTFE-coated stirrer bar. Either use bottles of standard volume such that the standard deviation from the mean volume for the batch of bottles is less than 1 ml or measure and record the volumes of individual, numbered bottles with an accuracy of 1 ml. Carefully grease the stoppers of the bottles with inert silicone grease to ensure proper closing and easy removal.

D.2.2 Oxygen electrode, preferably with a mounted stirrer, capable of measuring in the range 0 to 10 mg/l to an accuracy of 1 %. The steady-state should be reached within about 1,5 min. Mount the electrode, for example, in an inert stopper which makes a leakproof fit in the ground-glass neck of the incubation bottle or use a technique to measure the oxygen concentration in a circular bypass.

D.2.3 Magnetic stirrer or shaking device.

D.3 Procedure

Set up the incubation bottles as described in 8.4, but with three of each of bottles F_T , F_B and F_C . Place a stirrer bar in each bottle if they are to be stirred rather than shaken. Prepare sufficient test medium, preferably the standard test medium (6.2.1), to perform the complete test. To guarantee a sufficient supply of nutrients, increase the amount of ammonium chloride in solution A (6.2.1.1) by a factor of three to 1,5 g/l. Inoculate the medium in accordance with 8.3, preferably using activated sludge at a concentration of 30 mg/l of suspended solids, mix well and add the mixture to the bottles. Add a volume equal to two-thirds of the volume of the bottle (e.g. 200 ml of liquid to 300 ml bottles). Place the bottles on the shaking device or stir them, and incubate at 20 °C to 25 °C for one week. During this time, the bacteria will use their reserve material and the inoculum will be stabilized. Then aerate the bottles with the help of water-saturated compressed air and an air diffuser for about 15 min. Measure the initial oxygen concentration. Add to the relevant bottles the test or reference material as specified in 8.1 and 8.2. The maximum test material concentration in this test should correspond to 150 mg/l ThOD which corresponds to about 90 mg/l TOC. Stopper all bottles tightly and continue the incubation.

After incubation periods of one week, or more frequently, determine the concentration of dissolved oxygen in each bottle. Maintain the bottles at the incubation temperature kept at a constant value ($\pm 0,5$ °C) during the measurements. Calibrate accurately the oxygen electrode in accordance with the manufacturer's instructions. For oxygen measurements, take each bottle successively and shake vigorously by hand for about 30 s. Place the bottle on a stirrer, without stirring. Remove the stopper and immediately put the oxygen electrode through the neck of the bottle so that the electrode stopper effectively closes the bottle and the electrode tip is well below the liquid surface. Start the stirrer at such a speed that oxygen measurements are possible but no vortex is formed. Use the same stirrer speed throughout a set of measurements and electrode calibrations. Record the oxygen value when it is stable, which should be reached in about 2 min. Use for the calculation of the test result only those bottles in which the oxygen concentration was $>1,5$ mg/l. Alternatively, other appropriate techniques for measuring the dissolved-oxygen concentration in the aqueous phase can be used, such as an oxygen electrode positioned in a closed circular bypass.

Then measure the pH in each bottle and record it. If the pH is below 6,0, adjust to about 7,5 with 0,1 mol/l to 0,5 mol/l sodium hydroxide solution. If the pH is above 8,0, adjust to about 7,5 with 0,1 mol/l to 0,5 mol/l hydrochloric acid solution. Finally, aerate the medium in each bottle with an air-diffuser for 15 min and measure the oxygen concentration again as described previously. Re-stopper the bottles, put them on the shaker and continue the incubation. Correct the measured BOD values at the end of the test for nitrification (see 8.4 and Annex B).

D.4 Calculation of results

Determine the relative oxygen uptake, U_r , measured in the aqueous phase in each bottle, using [Formula \(D.1\)](#):

$$U_r = \frac{C_{Bt} - C_t}{C_s} \quad (D.1)$$

where

C_{Bt} is the mean value, in milligrams per litre, of the dissolved-oxygen concentration in the blank bottles after incubation, at time t ;

C_t is the dissolved-oxygen concentration, in milligrams per litre, in each test bottle after incubation, at time t ;

C_s is the saturation value, in milligrams per litre, for dissolved oxygen.

Use as the saturation value C_s a mean value measured after each aeration or re-aeration of the blank and test bottles. The theoretical value at standard atmospheric pressure (1 013 hPa) and 20 °C is 9,08 mg/l.

Determine, using [Formula \(D.2\)](#), the total oxygen capacity O_c (mg/bottle) of a bottle from the maximum oxygen content in the gas phase and the oxygen content in the aqueous phase at standard pressure and 20 °C:

$$O_c = (0,28 \times V_g) + (0,009 \times V_l) \quad (D.2)$$

where

0,28 is the oxygen content of normal air, in milligrams per millilitre;

V_g is the volume of gas in an incubation bottle, in millilitres;

0,009 is the oxygen content of saturated water, in milligrams per millilitre;

V_l is the volume of liquid in an incubation bottle, in millilitres.

Normally, V_1 will be constant in a test series, except if samples for analysis are taken, but V_g may differ depending on the test bottles used. If the differences between the individual bottles are small, a constant O_c may be used. If the differences are significant (e.g. >2 ml for bottles with a volume of 200 ml), O_c should be calculated for each bottle. If V_1 decreases according to the volume of the sample taken, V_g increases proportionately.

Then calculate the oxygen uptake BOD (mg/bottle) from [Formula \(D.3\)](#):

$$\text{BOD} = U_r \times O_c \quad (\text{D.3})$$

Sum the oxygen uptake Σ BOD (mg/bottle) for all (n) incubation periods using [Formula \(D.4\)](#) to obtain the BOD at the end of the test.

$$\Sigma \text{BOD} = \text{BOD}_1 + \text{BOD}_2 + \dots + \text{BOD}_n \quad (\text{D.4})$$

Finally, calculate the percentage biodegradation as described in [9.1](#).

Use appropriate formulae to determine the abiotic elimination, as well as the biodegradation of the reference material and the inhibition control.

STANDARDSISO.COM : Click to view the full PDF of ISO 14851:2019