

First edition  
2003-04-15

AMENDMENT 2  
2018-06

---

---

**Tobacco — Determination of the  
content of total alkaloids as nicotine  
— Continuous-flow analysis method**

**AMENDMENT 2**

*Tabac — Détermination de la teneur en alcaloïdes totaux exprimés en  
nicotine — Méthode par analyse en flux continu*

AMENDEMENT 2

STANDARDSISO.COM : Click to view the full PDF of ISO 15152:2003/Amd 2:2018



Reference number  
ISO 15152:2003/Amd.2:2018(E)

© ISO 2018

STANDARDSISO.COM : Click to view the full PDF of ISO 15152:2003/Amd 2:2018



**COPYRIGHT PROTECTED DOCUMENT**

© ISO 2018

All rights reserved. Unless otherwise specified, or required in the context of its implementation, no part of this publication may be reproduced or utilized otherwise in any form or by any means, electronic or mechanical, including photocopying, or posting on the internet or an intranet, without prior written permission. Permission can be requested from either ISO at the address below or ISO's member body in the country of the requester.

ISO copyright office  
CP 401 • Ch. de Blandonnet 8  
CH-1214 Vernier, Geneva  
Phone: +41 22 749 01 11  
Fax: +41 22 749 09 47  
Email: [copyright@iso.org](mailto:copyright@iso.org)  
Website: [www.iso.org](http://www.iso.org)

Published in Switzerland

# Tobacco — Determination of the content of total alkaloids as nicotine — Continuous-flow analysis method

## AMENDMENT 2

*Page 1, Clause 4*

Add the following sentence to the end of the second paragraph in Clause 4:

'Alternatively, neutralizing solution C (5.7) shall be pumped through the system and mixed in a 2 litre Büchner flask.'

*Page 4, 5.5*

Delete the second paragraph of 5.5.

*Page 4, 5.6*

Delete the second paragraph of 5.6.

*Page 4*

Insert a new subclause 5.7 as follows and renumber subsequent subclauses:

### **5.7 Alternative neutralizing solution C**

Mix 500 ml of sodium hydroxide (mass fraction of 20 %) with 70 ml of sodium hypochlorite containing a mass fraction of 5 % active chlorine.

*Page 5, 5.9.2, second paragraph*

Delete the expression '(5 times)' from the first sentence.


*Page 6, 7.3, second paragraph*

Delete the expression '(5 times)' from the first sentence.

Replace Annex B with the following:

**Annex B**  
(informative)

**Example of a continuous-flow analyser**

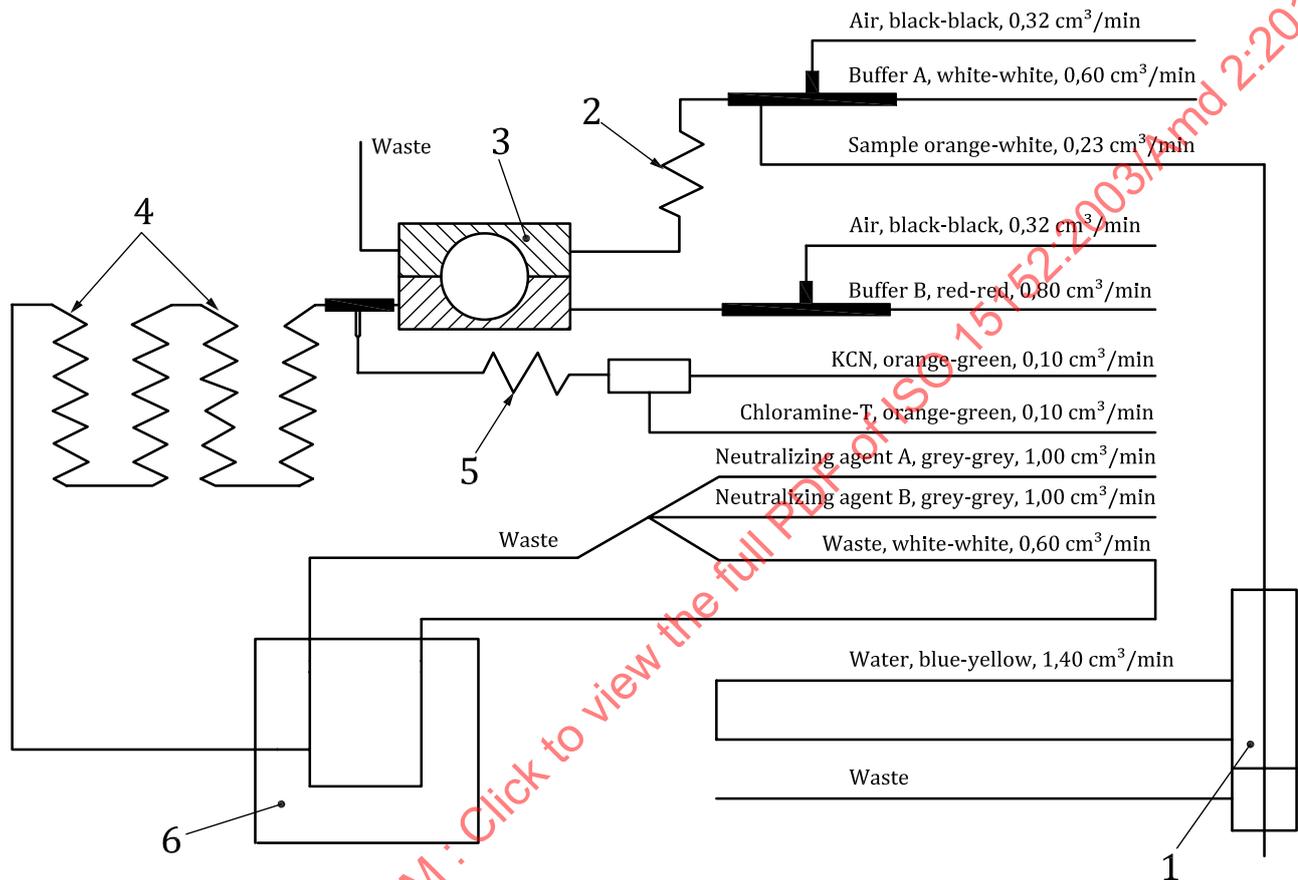

**Key**

- 1 sampler (sample, wash)
- 2 delay coil, 10 turns
- 3 dialyser
- 4 delay coils, 20 turns
- 5 mixing coil
- 6 colorimeter with filter, flow cell and reference cell

NOTE 1 The waste from the colorimeter is only for instruments that require a debubble line.

NOTE 2 Neutralization solution C (see 5.7) can be used instead of a mix of neutralization solutions A and B.

**Figure B.1 — Flow diagram of continuous-flow analyser**