
**Iron ores — Determination of
nickel — Flame atomic absorption
spectrometric method**

*Minerais de fer — Dosage du nickel — Méthode par spectrométrie
d'absorption atomique dans la flamme*

STANDARDSISO.COM : Click to view the full PDF of ISO 15633:2017



STANDARDSISO.COM : Click to view the full PDF of ISO 15633:2017



COPYRIGHT PROTECTED DOCUMENT

© ISO 2017, Published in Switzerland

All rights reserved. Unless otherwise specified, no part of this publication may be reproduced or utilized otherwise in any form or by any means, electronic or mechanical, including photocopying, or posting on the internet or an intranet, without prior written permission. Permission can be requested from either ISO at the address below or ISO's member body in the country of the requester.

ISO copyright office
Ch. de Blandonnet 8 • CP 401
CH-1214 Vernier, Geneva, Switzerland
Tel. +41 22 749 01 11
Fax +41 22 749 09 47
copyright@iso.org
www.iso.org

Contents

	Page
Foreword.....	iv
Introduction.....	v
1 Scope	1
2 Normative references	1
3 Terms and definitions	1
4 Principle	1
5 Reagents	2
6 Apparatus	2
7 Sampling and samples	3
7.1 Laboratory sample.....	3
7.2 Preparation of predried test samples.....	3
8 Procedure	4
8.1 Number of determinations.....	4
8.2 Test portion.....	4
8.3 Blank test and check test.....	4
8.4 Determination.....	4
8.4.1 Decomposition of the test portion.....	4
8.4.2 Removal of iron.....	5
8.4.3 Treatment of the residue.....	5
8.4.4 Preparation of the calibration solutions.....	5
8.4.5 Adjustment of the atomic absorption spectrometer.....	6
8.4.6 Atomic absorption measurements.....	6
9 Expression of results	6
9.1 Calculation of nickel content.....	6
9.2 General treatment of results.....	7
9.2.1 Repeatability and permissible tolerance.....	7
9.2.2 Determination of analytical result.....	7
9.2.3 Between-laboratories precision.....	7
9.2.4 Check for trueness.....	8
9.2.5 Calculation of final result.....	9
9.3 Oxide factor.....	9
10 Test report	9
Annex A (normative) Flowsheet of the procedure for the acceptance of analytical values for test samples	11
Annex B (informative) Derivation of repeatability and permissible tolerance formulae	12
Annex C (informative) Precision data obtained by international analytical trial	13
Bibliography	14

Foreword

ISO (the International Organization for Standardization) is a worldwide federation of national standards bodies (ISO member bodies). The work of preparing International Standards is normally carried out through ISO technical committees. Each member body interested in a subject for which a technical committee has been established has the right to be represented on that committee. International organizations, governmental and non-governmental, in liaison with ISO, also take part in the work. ISO collaborates closely with the International Electrotechnical Commission (IEC) on all matters of electrotechnical standardization.

The procedures used to develop this document and those intended for its further maintenance are described in the ISO/IEC Directives, Part 1. In particular the different approval criteria needed for the different types of ISO documents should be noted. This document was drafted in accordance with the editorial rules of the ISO/IEC Directives, Part 2 (see www.iso.org/directives).

Attention is drawn to the possibility that some of the elements of this document may be the subject of patent rights. ISO shall not be held responsible for identifying any or all such patent rights. Details of any patent rights identified during the development of the document will be in the Introduction and/or on the ISO list of patent declarations received (see www.iso.org/patents).

Any trade name used in this document is information given for the convenience of users and does not constitute an endorsement.

For an explanation on the voluntary nature of standards, the meaning of ISO specific terms and expressions related to conformity assessment, as well as information about ISO's adherence to the World Trade Organization (WTO) principles in the Technical Barriers to Trade (TBT) see the following URL: www.iso.org/iso/foreword.html.

This document was prepared by Technical Committee ISO/TC 102, *Iron ore and direct reduced iron*, Subcommittee SC 02, *Chemical analysis*.

This third edition cancels and replaces the second edition (ISO 15633:2015), of which it constitutes a minor revision with the following changes:

- in [Clause 1](#) “0,1 %” has been replaced by “0,10 %”;
- in [7.1](#) “<160 μm” has been replaced by “–160 μm”;
- in [8.2](#) a new sentence has been included to make reference to ISO 2596;
- in the titles of [9.1](#), [Table B.1](#) and [Figure C.1](#) “sulfur” has been replaced by “nickel”;
- in [9.2.4 Formula \(8\)](#) and the relevant descriptions have been modified to harmonize this subclause across all documents for which ISO/TC 102/SC 2 is responsible.

Introduction

The objective of a proposed revision of ISO 9685:1991 was to extend the lower limit for a flame atomic absorption spectrometric method determination of both chromium and nickel in iron ores down to 0,001 %. However, due to bias, the method for nickel could not be approved for referee purposes.

The 22nd meeting of ISO/TC 102/SC 2 decided to progress this document as a non-referee method.

STANDARDSISO.COM : Click to view the full PDF of ISO 15633:2017

[STANDARDSISO.COM](https://standardsiso.com) : Click to view the full PDF of ISO 15633:2017

Iron ores — Determination of nickel — Flame atomic absorption spectrometric method

WARNING — This document might involve hazardous materials, operations and equipment. This document does not purport to address all of the safety problems associated with its use. It is the responsibility of the user of this document to establish appropriate health and safety practices prior to use.

1 Scope

This document specifies a flame atomic absorption spectrometric method for the determination of the nickel mass fraction of iron ores.

This method is applicable to mass fractions of nickel between 0,001 % and 0,10 % in natural iron ores, iron ore concentrates and agglomerates including sinter products.

This method is not appropriate for referee purposes.

2 Normative references

The following documents are referred to in the text in such a way that some or all of their content constitutes requirements of this document. For dated references, only the edition cited applies. For undated references, the latest edition of the referenced document (including any amendments) applies.

ISO 648, *Laboratory glassware — Single-volume pipettes*

ISO 1042, *Laboratory glassware — One-mark volumetric flasks*

ISO 2596, *Iron ores — Determination of hygroscopic moisture in analytical samples — Gravimetric, Karl Fischer and mass-loss methods*

ISO 3082, *Iron ores — Sampling and sample preparation procedures*

ISO 3696, *Water for analytical laboratory use — Specification and test methods*

ISO 7764, *Iron ores — Preparation of predried test samples for chemical analysis*

ISO 80000-1:2009, *Quantities and units — Part 1: General*

3 Terms and definitions

No terms and definitions are listed in this document.

ISO and IEC maintain terminological databases for use in standardization at the following addresses:

- IEC Electropedia: available at <http://www.electropedia.org/>
- ISO Online browsing platform: available at <http://www.iso.org/obp>

4 Principle

The test portion of iron ore is decomposed by treatment with hydrochloric and nitric acids.

The major portion of iron in the filtrate is removed by extraction with 4-methylpentan-2-one.

The insoluble residue is ignited and silicon dioxide is removed by evaporation with hydrofluoric and sulfuric acids. The residue is fused with a mixture of sodium carbonate and sodium tetraborate, and then dissolved with hydrochloric acid and combined with the main solution.

5 Reagents

During the analysis, use only reagents of recognized analytical grade and only water that conforms to Grade 2 of ISO 3696.

- 5.1 **Sodium carbonate** (Na_2CO_3), anhydrous powder.
- 5.2 **Sodium tetraborate** ($\text{Na}_2\text{B}_4\text{O}_7$), anhydrous powder.
- 5.3 **Lithium tetraborate** ($\text{Li}_2\text{B}_4\text{O}_7$), anhydrous powder.
- 5.4 **Hydrochloric acid**, $\rho = 1,16$ g/ml to 1,19 g/ml.
- 5.5 **Hydrochloric acid**, $\rho = 1,16$ g/ml to 1,19 g/ml, diluted 2 + 1.
- 5.6 **Hydrochloric acid**, $\rho = 1,16$ g/ml to 1,19 g/ml, diluted 1 + 1.
- 5.7 **Hydrochloric acid**, $\rho = 1,16$ g/ml to 1,19 g/ml, diluted 2 + 100.
- 5.8 **Nitric acid**, $\rho = 1,4$ g/ml.
- 5.9 **Nitric acid**, $\rho = 1,4$ g/ml, diluted 1 + 1.
- 5.10 **Hydrofluoric acid**, $\rho = 1,13$ g/ml, 40 % (mass fraction) or $\rho = 1,19$ g/ml, 48 % (mass fraction).
- 5.11 **Sulfuric acid**, $\rho = 1,84$ g/ml.
- 5.12 **Sulfuric acid**, $\rho = 1,84$ g/ml, diluted 1 + 1.
- 5.13 **4-methylpentan-2-one**, (methyl isobutyl ketone, MIBK).
- 5.14 **Nickel standard solution A**, 100 μg Ni/ml.

Dissolve 0,100 0 g of nickel metal [purity > 99,9 % (mass fraction)] (see Note of [5.15](#)) in 30 ml of nitric acid ([5.9](#)). After cooling, transfer quantitatively to a 1 000 ml one-mark volumetric flask, dilute to volume with water and mix.

- 5.15 **Nickel standard solution B**, 10 μg Ni/ml.

Transfer 100,0 ml of nickel standard solution A ([5.14](#)) to a 1 000 ml one-mark volumetric flask. Dilute to volume with water and mix.

NOTE The purity of the metals stated on the certificates does not generally take into account the presence of absorbed gases such as oxygen, carbon monoxide, etc.

6 Apparatus

Ordinary laboratory apparatus including one-mark pipettes and one-mark volumetric flasks complying with the specifications of ISO 648, ISO 1042 and the following.

6.1 Platinum crucible, of minimum capacity 25 ml.

6.2 Muffle furnace.

6.3 Atomic absorption spectrometer, equipped with an air-acetylene burner.

WARNING — Follow the manufacturer's instructions for igniting and extinguishing the air-acetylene flame to avoid possible explosion hazards. Wear tinted safety glasses whenever the burner is in operation.

The atomic absorption spectrometer used in this method shall meet the following criteria.

- a) Minimum sensitivity: the absorbance of the most concentrated calibration solution (see 8.4.4) is at least 0,3.
- b) Graph linearity: the slope of the calibration graph covering the top 20 % of the concentration range (expressed as a change in absorbance) is not less than 0,7 of the value of the slope for the bottom 20 % of the concentration range determined in the same way.
- c) Minimum stability: the standard deviation of the absorbance of the most concentrated calibration solution and that of the zero calibration solution each being calculated from a sufficient number of repetitive measurements are less than 1,5 % and 0,5 %, respectively of the mean value of the absorbance of the most concentrated solution.

The use of a strip-chart recorder and/or digital readout device is recommended to evaluate criteria a), b) and c) and for all subsequent measurements.

NOTE Instrument parameters vary with each instrument. The following parameters were successfully used in several laboratories and they can be used as guidelines. An air-acetylene flame was used.

Hollow-cathode lamp, mA	10
Wavelength, nm	232,0
Air flow rate, l/min	10
Acetylene flow rate, l/min	2,5

In systems where the values shown above for gas flow rates do not apply, the ratio of the gas flow rates can still be a useful guideline.

7 Sampling and samples

7.1 Laboratory sample

For analysis, use a laboratory sample of $-100 \mu\text{m}$ particle size which has been taken and prepared in accordance with ISO 3082. In the case of ores having significant contents of combined water or oxidizable compounds, use a particle size of $-160 \mu\text{m}$.

NOTE A guideline on significant contents of combined water and oxidizable compounds is incorporated in ISO 7764.

Ensure that the sample has not been pulverized in a nickel/chromium pot.

7.2 Preparation of predried test samples

Thoroughly mix the laboratory sample and, taking multiple increments, extract a test sample in such a manner that it is representative of the entire contents of the container. Dry the test sample at $105 \text{ °C} \pm 2 \text{ °C}$ as specified in ISO 7764. This is the predried test sample.

8 Procedure

8.1 Number of determinations

Carry out the analysis at least in duplicate, in accordance with [Annex A](#) and independently on one predried test sample.

NOTE The expression “independently” means that the second and any subsequent result is not affected by the previous result(s). For this particular analytical method, this condition implies that the repetition of the procedure is carried out either by the same operator at a different time or by a different operator including, in either case, appropriate recalibration.

8.2 Test portion

Taking several increments, weigh, to the nearest 0,000 2 g, approximately 1 g of the predried test sample in accordance with [7.2](#). For ores having significant content of combined water or oxidizable compounds, an air-equilibrate test sample shall be prepared in accordance with ISO 2596.

The test portion should be taken and weighed quickly, in order to avoid reabsorption of moisture.

8.3 Blank test and check test

In each run, one blank test and one analysis of a certified reference material of the same type of ore shall be carried out in parallel with the analysis of the ore sample(s) under the same conditions. A predried test sample of the certified reference material shall be prepared as specified in [7.2](#).

The certified reference material should be of the same type as the sample to be analysed and the properties of the two materials should be sufficiently similar to ensure that in either case, no significant changes in the analytical procedure become necessary. Where a certified reference material is not available, a reference material can be used (see [9.2.4](#)).

Where the analysis is carried out on several samples at the same time, the blank value may be represented by one test provided that the procedure is the same and the reagents used are from the same reagent bottles.

Where the analysis is carried out on several samples of the same type of ore at the same time, the analytical value of one certified reference material can be used.

8.4 Determination

8.4.1 Decomposition of the test portion

Transfer the test portion ([8.2](#)) to a 250 ml tall-form beaker. Moisten with a few millilitres of water. Add 25 ml of hydrochloric acid ([5.4](#)), cover with a watch-glass and heat for about 1 h on a hotplate adjusted to give a temperature of 100 °C in a test beaker containing a similar volume and depth of sulfuric acid ([5.11](#)).

If the amount of insoluble residue is high, continue heating on a higher temperature zone of the hotplate without boiling the solution.

After heating the volume of the solution should have decreased to about 10 ml.

Add 5 ml of nitric acid ([5.8](#)) and 0,2 ml of sulfuric acid ([5.12](#)), digest for about 15 min and evaporate the solution almost to dryness.

If the test sample contains significant barium content, the addition of sulfuric acid ([5.12](#)) should be omitted.

Be sure to cover with a watch-glass in order to prevent sublimation of chlorides. During evaporation, shift the glass slightly.

Add 20 ml of hydrochloric acid (5.6) and heat to dissolve the salts. Cool, wash the watch-glass and the walls of the beaker, and filter the solution through a close-texture filter paper containing filter pulp (0,3 g to 0,4 g dry mass) into a 200 ml beaker. Carefully remove all adhering particles with a rubber-tipped glass rod or piece of moistened filter paper and transfer to the filter washing the paper with warm hydrochloric acid (5.7) until it is visibly free from iron, and, again, wash the paper with three or four portions of warm water. Reserve the filtrate and washings in the 200 ml beaker as the main solution. Transfer the filter paper and residue to a platinum crucible (6.1).

8.4.2 Removal of iron

Heat the main solution and evaporate almost to dryness. Dissolve the salts in 15 ml of hydrochloric acid (5.5) and transfer to a 200 ml separating funnel. Rinse the beaker with 20 ml of hydrochloric acid (5.5) and transfer the washings to the funnel.

Add 50 ml of 4-methylpentan-2-one (5.13) to the funnel and shake it for 1 min. Allow the two layers to separate and drain the lower aqueous layer into the original 200 ml beaker. Wash the organic layer by adding 10 ml of hydrochloric acid (5.5) to the funnel and shaking it for 30 s. Allow the two layers to separate and drain the lower aqueous layer into the 200 ml beaker to combine with the first aqueous layer. Discard the organic layer.

Heat the solution gently and expel almost all of the 4-methylpentan-2-one in the solution. Then, add 5 ml of nitric acid (5.8) and evaporate to dryness. Dissolve the salts in 20 ml of hydrochloric acid (5.6).

8.4.3 Treatment of the residue

Dry and burn off the filter paper in the platinum crucible (6.1) at a low temperature (500 °C to 800 °C) and ignite the residue. Cool, moisten with one to two drops of water and add three drops of sulfuric acid (5.11) and 5 ml of hydrofluoric acid (5.10). Evaporate slowly to expel silicon dioxide and continue heating to remove the sulfuric acid. Ignite at 800 °C for several minutes and then cool. Add 1,2 g of a mixture of 0,8 g of sodium carbonate (5.1) and 0,4 g of sodium tetraborate (5.2) to the residue and mix. Heat gently for several minutes, then heat at 1 000 °C for 15 min in a muffle furnace (6.2) or over a pressurized air burner for a time sufficient to produce a clear melt.

Lithium tetraborate (5.3) can be used instead of sodium tetraborate (5.2). In this case, sodium tetraborate should be replaced with lithium tetraborate in the calibration solutions (8.4.4) and the flame-setting solutions (8.4.5).

Allow the crucible to cool, add 10 ml of hydrochloric acid (5.6) and heat gently to dissolve the melt. Combine this solution with the main solution from 8.4.2, heat to remove carbon dioxide and cool.

For samples containing a mass fraction of more than or equal to 0,01 % nickel, transfer to a 100 ml one-mark volumetric flask quantitatively, dilute to volume with water and mix (this is the test solution).

For samples containing a mass fraction of less than 0,01 % nickel, transfer to a 50 ml one-mark volumetric flask quantitatively, dilute to volume with water and mix (this is the test solution).

8.4.4 Preparation of the calibration solutions

Introduce into a series of six 200 ml beakers, 30 ml of hydrochloric acid (5.6), 0,2 ml of sulfuric acid (5.12) and a mixture of 0,8 g of sodium carbonate (5.1) and 0,4 g of sodium tetraborate (5.2). Heat to remove carbon dioxide and cool.

Using pipettes or burettes, add 0 ml, 1,0 ml, 3,0 ml, 5,0 ml, 7,0 ml and 10,0 ml of nickel standard solution A (5.14) for samples containing 0,01 % nickel to 0,10 % nickel or nickel standard solution B (5.15) for samples containing 0,001 % nickel to 0,01 % nickel.

Transfer to six 100 ml one-mark volumetric flasks for samples containing 0,01 % nickel to 0,10 % nickel or 50 ml one-mark volumetric flasks for samples containing 0,001 % nickel to 0,01 % nickel, respectively. Dilute to volume with water and mix. These calibration solutions cover the nickel concentration ranges

0 µg/ml to 10 µg/ml in the case of standard solution A or 0 µg/ml to 2 µg/ml in the case of standard solution B.

8.4.5 Adjustment of the atomic absorption spectrometer

Optimize the response of the instrument as specified in 6.3. Set the wavelength at 232,0 nm to obtain minimum absorbance and adjust the readout to zero absorbance. Light the air-acetylene flame. After 10 min of preheating the burner, aspirate water, and, if necessary, readjust the readout to zero absorbance. Aspirate the calibration solution of highest nickel content (8.4.4) and adjust the fuel flow and burner position to obtain maximum absorbance.

Check that the conditions for zero absorbance have been maintained and evaluate the criteria in 6.3.

Repeat the aspiration of water and the calibration solution of highest nickel content to establish that the absorbance reading is not drifting. Set the reading for water to zero absorbance.

8.4.6 Atomic absorption measurements

Aspirate the calibration solutions (8.4.4) and the test solution (8.4.3) in order of increasing absorption, starting with the zero calibration solution (8.4.4) and the blank test solution. When a stable response has been obtained for each solution, record the readings. Aspirate water between each calibration and test solution.

Correct the absorbance values obtained for the calibration solutions by subtracting the absorbance of the zero calibration solution and prepare a calibration graph by plotting the net absorbance values against micrograms of nickel per millilitre. If the graph is substantially linear, subtract the absorbance obtained for the blank test from the absorbance obtained for the test solution and, using the graph, convert the net absorbance to micrograms of nickel per millilitre.

If any curvature obtained approaches the limit specified in 6.3 b), replot the graph using uncorrected values for all solutions and establish the concentration of the zero calibration solution from the intercept of the graph on the negative side of the concentration axis. Add this value to the nominal concentration values of the calibration solutions and replot the graph to pass through the origin. Determine from the graph the concentration of nickel, in micrograms per millilitre, in the blank test and test solutions, respectively, and correct the concentration of the test solution with the concentration of the blank test.

9 Expression of results

9.1 Calculation of nickel content

The mass fraction of nickel, w_{Ni} , expressed as a percentage by mass is calculated to five decimal places using the [Formulae \(1\)](#) or [\(2\)](#).

a) For < 0,01 % (mass fraction) of nickel content:

$$w_{Ni} \% = \frac{\rho_{Ni} \times 50}{m \times 10^6} \times 100 = \frac{\rho_{Ni} \times 0,005}{m} \quad (1)$$

b) For ≥ 0,01 % (mass fraction) of nickel content:

$$w_{Ni} \% = \frac{\rho_{Ni} \times 100}{m \times 10^6} \times 100 = \frac{\rho_{Ni} \times 0,01}{m} \quad (2)$$

where

w_{Ni} is the mass fraction, in percentage, of nickel in the test portion;

ρ_{Ni} is the concentration, in micrograms per millilitre, of nickel in the final test solution;

m is the mass, in grams, of the test portion.

9.2 General treatment of results

9.2.1 Repeatability and permissible tolerance

The precision of this analytical method is expressed by the regression [Formulae \(3\)](#) to [\(6\)](#).

NOTE Additional information is given in [Annexes B](#) and [C](#).

$$R_d = 0,0043 X^{0,3051} \quad (3)$$

$$P = 0,0215 X^{0,4977} \quad (4)$$

$$\sigma_d = 0,0015 X^{0,3051} \quad (5)$$

$$\sigma_L = 0,0076 X^{0,5242} \quad (6)$$

where

X is the concentration of nickel in the sample, expressed as mass fraction (%), and calculated as follows:

- for the within-laboratory, [Formulae \(3\)](#) and [\(5\)](#): the arithmetic mean of the duplicate values;
- for the between-laboratories, [Formulae \(4\)](#) and [\(6\)](#): the arithmetic mean of the final results ([9.2.5](#)) of the two laboratories;

R_d is the independent duplicate limit;

P is the permissible tolerance between laboratories;

σ_d is the independent duplicate standard deviation;

σ_L is the between-laboratories standard deviation.

9.2.2 Determination of analytical result

Having computed the independent duplicate results according to [Formulae \(1\)](#) or [\(2\)](#), compare them with the independent duplicate limit, R_d , using the procedure given in [Figure A.1](#) and obtain the final laboratory result μ (see [9.2.5](#)).

9.2.3 Between-laboratories precision

Between-laboratories precision is used to determine the agreement between the final results reported by two laboratories. The assumption is that both laboratories followed the procedure described in [9.2.2](#).

Compute the quantity shown in [Formula \(7\)](#):

$$\mu_{1,2} = \frac{\mu_1 + \mu_2}{2} \quad (7)$$

where

- μ_1 is final result reported by laboratory 1;
- μ_2 is the final result reported by laboratory 2;
- $\mu_{1,2}$ is the mean of final results.

Substitute $\mu_{1,2}$ for X in [Formula \(4\)](#) and calculate P .

If $|\mu_1 - \mu_2| \leq P$ the results are in agreement.

9.2.4 Check for trueness

The trueness of the analytical method shall be checked by applying it to a certified reference material (CRM) or a reference material (RM). The procedure is the same as that described above. After confirmation of the precision, the final laboratory result is compared with the reference or certified value, A_c . There are two possibilities.

- a) $|\mu_c - A_c| \leq C$, in which case the difference between the reported result and the reference/certified value is statistically insignificant.
- b) $|\mu_c - A_c| > C$, in which case the difference between the reported result and the reference/certified value is statistically significant.

where

- μ_c is the final result for the certified reference material;
- A_c is the reference/certified value for the CRM/RM;
- C is a value dependent on the type of CRM/RM used.

Certified reference materials used for this purpose should be prepared and certified in accordance with ISO Guide 35.

C should be calculated using [Formula \(8\)](#):

$$C = 2 \sqrt{\frac{s_{Lc}^2 + \frac{s_{Wc}^2}{n_{Wc}}}{N_c} + \sigma_L^2 + \frac{\sigma_d^2}{n}} \quad (8)$$

where

- s_{Lc} is the between-laboratories standard deviation of the certifying laboratories;
- s_{Wc} is the within-laboratory standard deviation of the certifying laboratories;
- n_{Wc} is the average number of replicate determinations in the certifying laboratories;

N_c is the number of certifying laboratories;

n is the number of replicate determinations carried out on the CRM/RM;

σ_L and σ_d are as defined in 9.2.1.

The following procedure should be used when the information on the reference material certificate is incomplete:

- if there are sufficient data to enable the between-laboratories standard deviation to be estimated, delete the expression s_{Wc}^2 / n_{Wc} and regard s_{Lc} as the standard deviation of the laboratory means;
- if the certification has been made by only one laboratory or if the interlaboratory results are missing, use [Formula \(9\)](#):

$$C = 2 \sqrt{2 \sigma_L^2 + \frac{\sigma_d^2}{n}} \quad (9)$$

A CRM certified by only one laboratory should be avoided unless it is known to have an unbiased certified value.

9.2.5 Calculation of final result

The final result is the arithmetic mean of the acceptable analytical values for the test sample, or as otherwise determined by the operations specified in [Annex A](#) calculated to five decimal places and rounded off to the third decimal place according to Rule A in ISO 80000-1:2009, B.3 as follows.

- a) Where the figure in the fourth decimal place is less than five, it is discarded and the figure in the third decimal place is kept unchanged.
- b) Where the figure in the fourth decimal place is five and there is a figure other than 0 in the fifth decimal place or where the figure in the fourth decimal place is greater than five, the figure in the third decimal place is increased by one.
- c) Where the figure in the fourth decimal place is five and the figure 0 is in the fifth decimal place, the five is discarded and the figure in the third decimal place is kept unchanged if it is 0, 2, 4, 6, or 8 and is increased by one if it is 1, 3, 5, 7, or 9.

9.3 Oxide factor

The mass fraction of nickel oxide is calculated using [Formula \(10\)](#).

$$w_{NiO} (\%) = 1,2725 w_{Ni} (\%) \quad (10)$$

10 Test report

The test report shall include the following information:


- a) name and address of the testing laboratory;
- b) date of issue of the test report;
- c) a reference to this document, i.e. ISO 15633;
- d) details necessary for the identification of the sample;
- e) result of the analysis;
- f) reference number of the result;

- g) any characteristics noticed during the determination and any operations not specified in this document which could have had an influence on the result, either for the test sample or for the certified reference material(s).

STANDARDSISO.COM : Click to view the full PDF of ISO 15633:2017

Annex A (normative)

Flowsheet of the procedure for the acceptance of analytical values for test samples

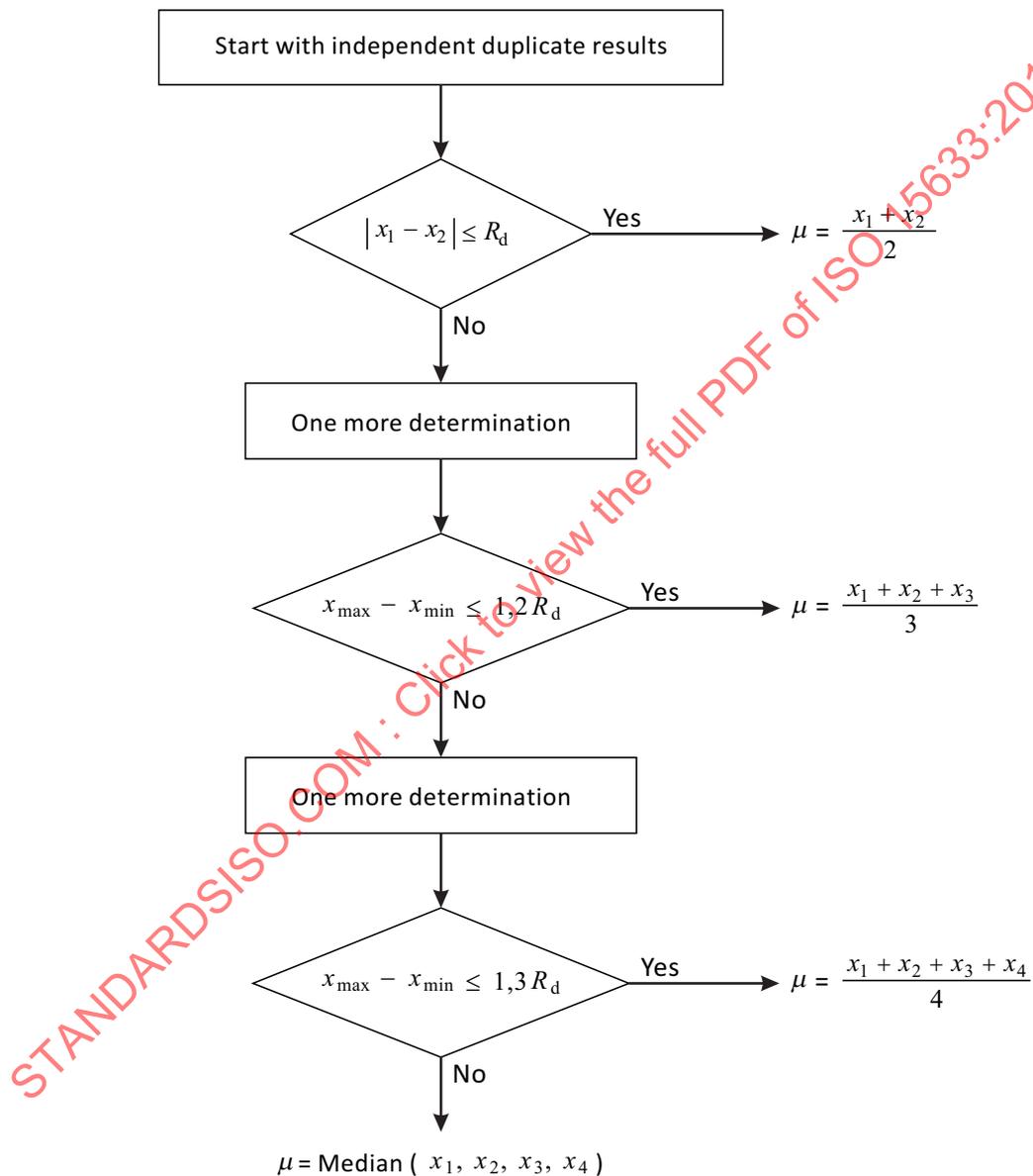

NOTE R_d is as defined in [9.2.1](#).

Figure A.1 — Procedure for the acceptance of analytical values for test samples