

---

---

**Medical electrical equipment —**

Part 2-84:

**Particular requirements for the basic  
safety and essential performance of  
ventilators for the emergency medical  
services environment**

*Appareils électromédicaux —*

*Partie 2-84: Exigences particulières relatives à la sécurité de base  
et aux performances essentielles des ventilateurs utilisés dans  
l'environnement des services médicaux d'urgence*

STANDARDSISO.COM : Click to view the full PDF of ISO 80601-2-84:2023



STANDARDSISO.COM : Click to view the full PDF of ISO 80601-2-84:2023



**COPYRIGHT PROTECTED DOCUMENT**

© ISO 2023

All rights reserved. Unless otherwise specified, or required in the context of its implementation, no part of this publication may be reproduced or utilized otherwise in any form or by any means, electronic or mechanical, including photocopying, or posting on the internet or an intranet, without prior written permission. Permission can be requested from either ISO at the address below or ISO's member body in the country of the requester.

ISO copyright office  
CP 401 • Ch. de Blandonnet 8  
CH-1214 Vernier, Geneva  
Phone: +41 22 749 01 11  
Email: [copyright@iso.org](mailto:copyright@iso.org)  
Website: [www.iso.org](http://www.iso.org)

Published in Switzerland

## Contents

<b>Foreword</b> .....	<b>vi</b>
<b>Introduction</b> .....	<b>viii</b>
<b>201. 1 Scope, object and related standards</b> .....	<b>1</b>
201. 1.1 Scope .....	1
201. 1.2 Object .....	2
201. 1.3 Collateral standards .....	2
201. 1.4 Particular standards .....	3
<b>201. 2 Normative references</b> .....	<b>4</b>
<b>201. 3 Terms and definitions</b> .....	<b>5</b>
<b>201. 4 General requirements</b> .....	<b>23</b>
201. 4.3 <i>Essential performance</i> .....	23
201. 4.3.101 Additional requirements for <i>essential performance</i> .....	23
201. 4.3.102 <i>System recovery</i> .....	23
201. 4.4 Additional requirements for <i>expected service life</i> .....	24
201. 4.5 Alternative <i>risk control</i> measures or test methods for <i>ME equipment</i> or <i>ME system</i> .....	25
201. 4.6 <i>ME equipment</i> or <i>ME system</i> parts that contact the <i>patient</i> .....	25
201. 4.11.101 Additional requirements for pressurized gas input .....	25
201. 4.11.101.1 Overpressure requirement .....	25
201. 4.11.101.2 Compatibility requirement for <i>medical gas pipeline systems</i> .....	25
201. 4.11.101.3 Compatibility requirements for pressure regulators .....	26
<b>201. 5 General requirements for testing of <i>ME equipment</i></b> .....	<b>26</b>
201. 5.101 Additional requirements for general requirements for testing of <i>ME equipment</i> .....	27
201. 5.101.1 <i>EMS ventilator</i> test conditions .....	27
201. 5.101.2 Gas flowrate and leakage specifications .....	27
201. 5.101.3 <i>EMS ventilator</i> testing errors .....	27
<b>201. 6 Classification of <i>ME equipment</i> and <i>ME systems</i></b> .....	<b>27</b>
<b>201. 7 <i>ME equipment</i> identification, marking and documents</b> .....	<b>28</b>
201. 7.1.101 <i>Information to be supplied by the manufacturer</i> .....	28
201. 7.2.3 Consult <i>accompanying documents</i> .....	28
201. 7.2.4.101 Additional requirements for <i>accessories</i> .....	28
201. 7.2.18 External gas source .....	28
201. 7.2.101 Additional requirements for <i>marking</i> on the outside of <i>ME equipment</i> or <i>ME equipment</i> parts .....	29
201. 7.4.2 Control devices .....	30
201. 7.4.3 Units of measurement .....	30
201. 7.4.101 Labelling of units of measurement .....	30
201. 7.9.2.1 General .....	30
201. 7.9.2.1.101 Additional general requirements .....	30
201. 7.9.2.2.101 Additional requirements for warnings and safety notices .....	30
201. 7.9.2.8.101 Additional requirements for start-up <i>procedure</i> .....	31
201. 7.9.2.9.101 Additional requirements for operating instructions .....	32
201. 7.9.2.12 <i>Cleaning, disinfection, and sterilization</i> .....	33
201. 7.9.2.14.101 Additional requirements for <i>accessories</i> , supplementary equipment, used material .....	33
201. 7.9.3.1.101 Additional general requirements .....	33
201. 7.9.3.101 Additional requirements for the <i>technical description</i> .....	34
<b>201. 8 Protection against electrical <i>hazards</i> from <i>ME equipment</i></b> .....	<b>34</b>
<b>201. 9 Protection against <i>mechanical hazards</i> of <i>ME equipment</i> and <i>ME systems</i></b> .....	<b>34</b>

201. 9.4.3.101	Additional requirements for instability from unwanted lateral movement.....	34
201. 9.4.4	Grips and other handling devices .....	35
201. 9.6.2.1.101	Additional requirements for audible acoustic energy.....	35
<b>201. 10</b>	<b>Protection against unwanted and excessive radiation hazards.....</b>	<b>36</b>
<b>201. 11</b>	<b>Protection against excessive temperatures and other hazards.....</b>	<b>36</b>
201. 11.1.2.2	Applied parts not intended to supply heat to a patient .....	37
201. 11.6.6	Cleaning and disinfection of ME equipment or ME system.....	37
201. 11.6.7	Sterilization of ME equipment or ME system.....	38
201. 11.7	Biocompatibility of ME equipment and ME systems .....	38
<b>201. 12</b>	<b>Accuracy of controls and instruments and protection against hazardous outputs .....</b>	<b>38</b>
201. 12.1.101	Volume-control inflation-type.....	38
201. 12.1.102	Pressure-control inflation-type.....	42
201. 12.1.103	Other inflation-types.....	45
201. 12.1.104	Inspiratory volume monitoring equipment.....	45
201. 12.4	Protection against hazardous output.....	45
201. 12.4.101	Oxygen monitor .....	45
201. 12.4.102	Measurement of airway pressure .....	46
201. 12.4.103	Measurement of expired volume and low volume alarm conditions .....	47
201. 12.4.104	Maximum limited pressure protection device.....	48
201. 12.4.105	High airway pressure alarm condition and protection device.....	48
201. 12.4.106	Expiratory end-tidal CO <sub>2</sub> monitoring equipment .....	50
201. 12.4.107	Protection against inadvertent setting of high airway pressure.....	50
<b>201. 13</b>	<b>Hazardous situations and fault conditions for ME equipment .....</b>	<b>50</b>
201. 13.2.101	Additional specific single fault conditions .....	50
201. 13.102	Failure of one gas supply to an EMS ventilator .....	51
201. 13.103	Independence of ventilation control function and related risk control measures .....	51
<b>201. 14</b>	<b>Programmable electrical medical systems (PEMS).....</b>	<b>52</b>
201. 14.1	General.....	52
201. 14.101	Cybersecurity capabilities of EMS ventilators.....	52
<b>201. 15</b>	<b>Construction of ME equipment .....</b>	<b>53</b>
201. 15.3.5.101	Additional requirements for rough handling.....	53
201. 15.4.1	Construction of connectors .....	53
201. 15.101	Mode of operation.....	54
201. 15.102	Delivered oxygen concentration .....	54
201. 15.103	Accessory self-check.....	54
<b>201. 16</b>	<b>ME systems.....</b>	<b>54</b>
201. 16.1.101	Additional general requirements for ME systems .....	54
201. 16.2.101	Additional requirements for power supply.....	55
<b>201. 17</b>	<b>Electromagnetic compatibility of ME equipment and ME systems .....</b>	<b>55</b>
<b>201. 101</b>	<b>Gas connections .....</b>	<b>55</b>
201. 101.1	Protection against reverse gas leakage .....	55
201. 101.2	Connection to a high-pressure inlet .....	56
201. 101.2.1	Low-pressure hose assembly.....	56
201. 101.2.2	Filter .....	56
201. 101.3	VBS connectors .....	56
201. 101.3.1	General.....	56
201. 101.3.2	Other named ports.....	56
201. 101.3.2.1	Patient-connection port.....	56
201. 101.3.2.2	Gas output port and gas return port .....	57
201. 101.3.2.3	Flow-direction-sensitive components .....	57
201. 101.3.2.4	Gas pathway connection port.....	57
201. 101.3.2.5	Gas exhaust port.....	58

201. 101.3.2.6	<i>Gas intake port</i> .....	58
<b>201. 102</b>	<b>Requirements for the VBS and accessories</b> .....	<b>58</b>
201. 102.1	General.....	58
201. 102.2	Labelling.....	58
201. 102.3	Breathing sets.....	59
201. 102.4	Water vapour management.....	59
201. 102.4.1	Humidification system.....	59
201. 102.4.2	<i>Heat and moisture exchanger (HME)</i> .....	59
201. 102.5	<i>Breathing system filters</i> .....	59
201. 102.6	Leakage from complete VBS.....	59
<b>201. 103</b>	<b>Spontaneous breathing during loss of ventilation</b> .....	<b>60</b>
<b>201. 104</b>	<b>Indication of duration of operation</b> .....	<b>60</b>
<b>201. 105</b>	<b>Functional connection</b> .....	<b>60</b>
201. 105.1	General.....	60
201. 105.2	Connection to an electronic health record.....	61
201. 105.3	Connection to a <i>distributed alarm system</i> .....	61
<b>201. 106</b>	<b>Display loops</b> .....	<b>61</b>
201. 106.1	Pressure-volume loops.....	61
201. 106.2	Flow-volume loops.....	61
<b>201. 107</b>	<b>Timed ventilatory pause</b> .....	<b>61</b>
201. 107.1	<i>Expiratory pause</i> .....	62
201. 107.2	<i>Inspiratory pause</i> .....	62
<b>202</b>	<b>Electromagnetic disturbances — Requirements and tests</b> .....	<b>63</b>
<b>206</b>	<b>Usability</b> .....	<b>64</b>
206.101	<i>Primary operating functions</i> .....	64
206.102	Training.....	66
<b>208</b>	<b>General requirements, tests and guidance for alarm systems in medical electrical equipment and medical electrical systems</b> .....	<b>66</b>
<b>212</b>	<b>Requirements for medical electrical equipment and medical electrical systems intended for use in the emergency medical services environment</b> .....	<b>68</b>
<b>Annex C (informative)</b>	<b>Guide to marking and labelling requirements for ME equipment and ME systems</b> .....	<b>71</b>
201.C.1	<i>Marking on the outside of ME equipment, ME systems or their parts</i> .....	71
201.C.2	<i>Accompanying documents, general</i> .....	72
201.C.3	<i>Accompanying documents, instructions for use</i> .....	72
201.C.4	<i>Accompanying documents, technical description</i> .....	74
<b>Annex D (informative)</b>	<b>Symbols on marking</b> .....	<b>76</b>
<b>Annex AA (informative)</b>	<b>Particular guidance and rationale</b> .....	<b>78</b>
AA.1	General guidance.....	78
AA.2	Rationale for particular clauses and subclauses.....	78
<b>Annex BB (informative)</b>	<b>Data interfaces</b> .....	<b>114</b>
BB.1	Background and purpose.....	114
BB.2	Data definition.....	115
<b>Annex CC (informative)</b>	<b>Reference to the IMDRF essential principles and labelling guidances</b> .....	<b>122</b>
<b>Annex DD (informative)</b>	<b>Reference to the essential principles</b> .....	<b>126</b>
<b>Bibliography</b>	.....	<b>129</b>
<b>Terminology — Alphabetized index of defined terms</b>	.....	<b>133</b>

## Foreword

ISO (the International Organization for Standardization) and IEC (the International Electrotechnical Commission) form the specialized system for worldwide standardization. National bodies that are members of ISO or IEC participate in the development of International Standards through technical committees established by the respective organization to deal with particular fields of technical activity. ISO and IEC technical committees collaborate in fields of mutual interest. Other international organizations, governmental and non-governmental, in liaison with ISO and IEC, also take part in the work.

The procedures used to develop this document and those intended for its further maintenance are described in the ISO/IEC Directives, Part 1. In particular, the different approval criteria needed for the different types of document should be noted. This document was drafted in accordance with the editorial rules of the ISO/IEC Directives, Part 2 (see [www.iso.org/directives](http://www.iso.org/directives) or [www.iec.ch/members\\_experts/refdocs](http://www.iec.ch/members_experts/refdocs)).

Attention is drawn to the possibility that some of the elements of this document may be the subject of patent rights. ISO and IEC shall not be held responsible for identifying any or all such patent rights. Details of any patent rights identified during the development of the document will be in the Introduction and/or on the ISO list of patent declarations received (see [www.iso.org/patents](http://www.iso.org/patents)) or the IEC list of patent declarations received (see <https://patents.iec.ch>).

Any trade name used in this document is information given for the convenience of users and does not constitute an endorsement.

For an explanation of the voluntary nature of standards, the meaning of ISO specific terms and expressions related to conformity assessment, as well as information about ISO's adherence to the World Trade Organization (WTO) principles in the Technical Barriers to Trade (TBT) see [www.iso.org/iso/foreword.html](http://www.iso.org/iso/foreword.html). In the IEC, see [www.iec.ch/understanding-standards](http://www.iec.ch/understanding-standards).

This document was prepared jointly by Technical Committee ISO/TC 121, *Anaesthetic and respiratory equipment*, Subcommittee SC 3, *Respiratory devices and related equipment used for patient care*, and Technical Committee IEC/TC 62, *Medical equipment, software, and systems*, Subcommittee SC 62D, *Particular medical equipment, software, and systems*, in collaboration with the European Committee for Standardization (CEN) Technical Committee CEN/TC 215, *Respiratory and anaesthetic equipment*, in accordance with the Agreement on technical cooperation between ISO and CEN (Vienna Agreement).

This second edition cancels and replaces the first edition (ISO 80601-2-84:2020), which has been technically revised. The main changes are as follows:

- alignment with IEC 60601-1:2005+AMD1:2012+AMD2:2020, IEC 60601-1-2:2014+AMD1:2020, IEC 60601-1-6:2010+AMD1:2013+AMD2:2020, IEC 60601-1-8:2006+AMD1:2012+AMD2:2020, and IEC 60601-1-12:2014+AMD1:2020;
- added requirements for a *responsible organization* log
- added requirements for the display legibility for *operators* wearing personal protective equipment;
- added requirements for display during calibration of gas monitors;
- clarified *maximum limited pressure* requirements;
- added requirements for *ventilator system recovery*;
- added requirements and definitions for *cybersecurity*; and

— harmonization with ISO 20417, where appropriate.

A list of all parts in the ISO 80601 series and the IEC 80601 series can be found on the ISO and IEC websites.

Any feedback or questions on this document should be directed to the user's national standards body. A complete listing of these bodies can be found at [www.iso.org/members.html](http://www.iso.org/members.html) and [www.iec.ch/national-committees](http://www.iec.ch/national-committees).

STANDARDSISO.COM : Click to view the full PDF of ISO 80601-2-84:2023

## Introduction

In referring to the structure of this document, the term

- “clause” means one of the five numbered divisions within the table of contents, inclusive of all subdivisions (e.g. Clause 201 includes subclauses 201.7, 201.8, etc.);
- “subclause” means a numbered subdivision of a clause (e.g. 201.7, 201.8 and 201.9 are all subclauses of Clause 201).

References to clauses within this document are preceded by the term “Clause” followed by the clause number. References to subclauses within this particular document are by number only.

In this document, the conjunctive “or” is used as an “inclusive or” so a statement is true if any combination of the conditions is true.

In this document, the following verbal forms are used:

- “shall” indicates a requirement;
- “should” indicates a recommendation;
- “may” indicates a permission;
- “can” is used to describe a possibility or capability; and
- “must” is used to indicate an external constraint.

Annex C contains a guide to the *marking* and labelling requirements in this document.

Annex D contains a summary of the *symbols* referenced in this document.

Requirements in this document have been decomposed so that each requirement is uniquely delineated. This is done to support automated requirements tracking.

# Medical electrical equipment —

Part 2-84:

## Particular requirements for the basic safety and essential performance of ventilators for the emergency medical services environment

### 201.1 Scope, object and related standards

IEC 60601-1:2005+AMD1:2012+AMD2:2020, Clause 1 applies, except as follows:

#### 201.1.1 Scope

*Replacement:*

NOTE 1 There is guidance or rationale for this subclause contained in Clause AA.2.

This document applies to the *basic safety* and *essential performance* of an *EMS ventilator* in combination with its *accessories*, hereafter also referred to as *ME equipment*:

- intended for *patients* who need differing levels of support from *artificial ventilation* including *ventilator-dependent patients*;
- intended to be operated by a *healthcare professional operator*;
- intended for use in the *EMS environment*; and
- intended for *invasive* or *non-invasive ventilation*.

NOTE 2 An *EMS ventilator* can also be used for transport within a *professional healthcare facility*.

An *EMS ventilator* is not considered to use a *physiologic closed loop-control system* unless it uses a *physiological patient variable* to adjust the *artificial ventilation* therapy settings.

This document is also applicable to those *accessories* intended by their *manufacturer* to be connected to the *ventilator breathing system*, or to an *EMS ventilator*, where the characteristics of those *accessories* can affect the *basic safety* or *essential performance* of the *EMS ventilator*.

NOTE 3 If a clause or subclause is specifically intended to be applicable to *ME equipment* only, or to *ME systems* only, the title and content of that clause or subclause will say so. If that is not the case, the clause or subclause applies both to *ME equipment* and to *ME systems*, as relevant.

*Hazards* inherent in the intended physiological function of *ME equipment* or *ME systems* within the scope of this document are not covered by specific requirements in this document except in IEC 60601-1:2005+AMD1:2012+AMD2:2020, 7.2.13 and 8.4.1.

NOTE 4 Additional information can be found in IEC 60601-1:2005+AMD1:2012+AMD2:2020, 4.2.

This document does not specify the requirements for the following:

NOTE 5 See ISO/TR 21954 for guidance on the selection of the appropriate *ventilator* for a given *patient*.

- *ventilators* or *accessories* intended for *ventilator-dependent patients* in critical care applications, which are given in ISO 80601-2-12.

- *ventilators* or *accessories* intended for *ventilator-dependent patients* in the *home healthcare environment*, which are given in ISO 80601-2-72.
- *ventilators* or *accessories* intended for anaesthetic applications, which are given in ISO 80601-2-13.
- *ventilators* or *accessories* intended for ventilatory support equipment (intended only to augment the *ventilation* of spontaneously breathing *patients*), which are given in ISO 80601-2-79 and ISO 80601-2-80.
- obstructive sleep apnoea therapy *ME equipment*, which are given in ISO 80601-2-70.
- user-powered resuscitators, which are given in ISO 10651-4.
- gas-powered emergency resuscitators, which are given in ISO 10651-5.
- *continuous positive airway pressure (CPAP) ME equipment*.
- high-frequency jet *ventilators* (HFJVs), which are given in ISO 80601-2-87.
- high-frequency oscillatory *ventilators* (HFOVs)<sup>[44]</sup>, which are given in ISO 80601-2-87.

NOTE 6 An *EMS ventilator* can incorporate high-frequency jet or high-frequency oscillatory *ventilation-modes*.

- respiratory high-flow therapy equipment, which are given in ISO 80601-2-90.

NOTE 7 An *EMS ventilator* can incorporate high-flow therapy operational mode, but such a mode is only for spontaneously breathing *patients*.

- oxygen therapy constant flow *ME equipment*.
- cuirass or “iron-lung” *ventilators*.

### 201.1.2 Object

*Replacement:*

The object of this particular document is to establish *basic safety* and *essential performance* requirements for an *EMS ventilator*, as defined in 201.3.228, and its *accessories*.

*Accessories* are included because the combination of the *EMS ventilator* and the *accessories* needs to have acceptable *risk*. *Accessories* can have a significant impact on the *basic safety* or *essential performance* of an *EMS ventilator*.

NOTE 1 This document has been prepared to address the relevant *essential principles*<sup>[38]</sup> and labelling<sup>[39]</sup> guidances of the International Medical Devices Regulators Forum (IMDRF) as indicated in Annex CC.

NOTE 2 This document has been prepared to address the relevant *essential principles of safety and performance* of ISO 16142-1:2016 as indicated in Annex DD.

NOTE 3 This document has been prepared to address the relevant general safety and performance requirements of European regulation (EU) 2017/745<sup>[40]</sup>.

### 201.1.3 Collateral standards

*Amendment (add at the end of the subclause):*

This document refers to those applicable collateral standards that are listed in Clause 2 of the general standard and in 201.2 of this document.

NOTE The general standard is IEC 60601-1:2005+AMD1:2012+AMD2:2020.

IEC 60601-1-2:2014+AMD1:2020, IEC 60601-1-6:2010+AMD1:2013+AMD2:2020, IEC 60601-1-8:2016+AMD1:2012+AMD2:2020 and IEC 60601-1-12:2014+AMD1:2020 apply as modified in Clauses 202, 206, 208 and 212 respectively.

IEC 60601-1-3, IEC 60601-1-9, and IEC 60601-1-11 do not apply.

All other published collateral standards in the IEC 60601-1 series apply as published.

### **201.1.4 Particular standards**

*Replacement:*

In the IEC 60601 series, particular standards may modify, replace or delete requirements contained in the general standard, including the collateral standards, as appropriate for the particular *ME equipment* under consideration, and may add other *basic safety* or *essential performance* requirements.

A requirement of a particular standard takes priority over IEC 60601-1:2005+AMD1:2012+AMD2:2020 or the collateral standards.

For brevity, IEC 60601-1:2005+AMD1:2012+AMD2:2020 is referred to in this document as the general standard. Collateral standards are referred to by their document number.

The numbering of clauses and subclauses of this particular document corresponds to those of the general standard with the prefix "201" (e.g. 201.1 in this document addresses the content of Clause 1 of the general standard) or applicable collateral standard with the prefix "2xx" where xx is the final digits of the collateral standard document number (e.g. 202.4 in this document addresses the content of Clause 4 of the IEC 60601-1-2 collateral standard, 208.4 in this document addresses the content of Clause 4 of the IEC 60601-1-8 collateral standard, etc.). The changes to the text of the general standard are specified by the use of the following words:

"Replacement" means that the clause or subclause of IEC 60601-1:2005+AMD1:2012+AMD2:2020 or the applicable collateral standard is replaced completely by the text of this particular document.

"Addition" means that the text of this document is additional to the requirements of IEC 60601-1:2005+AMD1:2012+AMD2:2020 or the applicable collateral standard.

"Amendment" means that the clause or subclause of IEC 60601-1:2005+AMD1:2012+AMD2:2020 or the applicable collateral standard is amended as indicated by the text of this document.

Subclauses, figures or tables that are additional to those of the general standard are numbered starting from 201.101. However, due to the fact that definitions in the general standard are numbered 3.1 through 3.154, additional definitions in this document are numbered beginning from 201.3.201. Additional annexes are lettered AA, BB, etc., and additional items aa), bb), etc.

Subclauses or figures that are additional to those of a collateral standard are numbered starting from 20x, where "x" is the number of the collateral standard (e.g. 202 for IEC 60601-1-2, 208 for IEC 60601-1-8, etc.).

The term "this document" is used to make reference to the general standard, any applicable collateral standards and this particular document taken together.

Where there is no corresponding clause or subclause in this document, the clause or subclause of IEC 60601-1:2005+AMD1:2012+AMD2:2020 or the applicable collateral standard, although possibly not relevant, applies without modification; where it is intended that any part of

IEC 60601-1:2005+AMD1:2012+AMD2:2020 or the applicable collateral standard, although possibly relevant, is not to be applied, a statement to that effect is given in this particular document.

## **201.2 Normative references**

The following documents are referred to in the text in such a way that some or all of their content constitutes requirements of this document. For dated references, only the edition cited applies. For undated references, the latest edition of the referenced document (including any amendments) applies.

ISO 32:1977, *Gas cylinders for medical use — Marking for identification of content*

ISO 3744:2010, *Acoustics — Determination of sound power levels and sound energy levels of noise sources using sound pressure — Engineering methods for an essentially free field over a reflecting plane*

ISO 4871:1996, *Acoustics — Declaration and verification of noise emission values of machinery and equipment*

ISO 5356-1:2015, *Anaesthetic and respiratory equipment — Conical connectors — Part 1: Cones and sockets*

ISO 5359:2014+AMD1:2017, *Anaesthetic and respiratory equipment — Low-pressure hose assemblies for use with medical gases*

ISO 5367:2014, *Anaesthetic and respiratory equipment — Breathing sets and connectors*

ISO 7396-1:2016+AMD1:2017, *Medical gas pipeline systems — Part 1: Pipeline systems for compressed medical gases and vacuum*

ISO 9360-1:2000, *Anaesthetic and respiratory equipment — Heat and moisture exchangers (HMEs) for humidifying respired gases in humans — Part 1: HMEs for use with minimum tidal volumes of 250 ml*

ISO 9360-2:2001, *Anaesthetic and respiratory equipment — Heat and moisture exchangers (HMEs) for humidifying respired gases in humans — Part 2: HMEs for use with tracheostomized patients having minimum tidal volumes of 250 ml*

ISO 10524-1:2018, *Pressure regulators for use with medical gases — Part 1: Pressure regulators and pressure regulators with flow-metering devices*

ISO 10524-3:2019, *Pressure regulators for use with medical gases — Part 3: Pressure regulators integrated with cylinder valves (VIPRS)*

ISO 14937:2009, *Sterilization of health care products — General requirements for characterization of a sterilizing agent and the development, validation and routine control of a sterilization process for medical devices*

ISO 17664-1:2021, *Processing of health care products — Information to be provided by the medical device manufacturer for the processing of medical devices — Part 1: Critical and semi-critical medical devices*

ISO 17664-2:2021, *Processing of health care products — Information to be provided by the medical device manufacturer for the processing of medical devices — Part 2: Non-critical medical devices*

ISO 18562-1:2023, *Biocompatibility evaluation of breathing gas pathways in healthcare applications — Part 1: Evaluation and testing within a risk management process*

ISO 19054:2005+AMD1:2016, *Rail systems for supporting medical equipment*

ISO 19223:2019, *Lung ventilators and related equipment — Vocabulary and semantics*

ISO 20417:2021, *Medical devices — Information to be supplied by the manufacturer*

ISO 23328-1:2003, *Breathing system filters for anaesthetic and respiratory use: — Part 1: Salt test method to assess filtration performance*

ISO 23328-2:2002, *Breathing system filters for anaesthetic and respiratory use: — Part 2: Non-filtration aspects*

ISO 80601-2-55:2018, *Medical electrical equipment — Part 2-55: Particular requirements for the basic safety and essential performance of respiratory gas monitors*

ISO 80601-2-74:2021, *Medical electrical equipment — Part 2-74: Particular requirements for the basic safety and essential performance of respiratory humidifying equipment*

IEC 60601-1:2005+AMD1:2012+AMD2:2020, *Medical electrical equipment — Part 1: General requirements for basic safety and essential performance*

IEC 60601-1-12:2014+AMD1:2020, *Medical Electrical Equipment — Part 1-12: General requirements for basic safety and essential performance — Collateral Standard: Requirements for medical electrical equipment and medical electrical systems used in the emergency medical services environment*

IEC 61672-1:2013, *Electroacoustics — Sound level meters — Part 1: Specifications*

IEC 62366-1:2015+AMD1:2020, *Medical devices — Part 1: Application of usability engineering to medical devices*

IEC 62570:2014, *Standard practice for marking medical devices and other items for safety in the magnetic resonance environment*

IEC Guide 115:2021, *Application of uncertainty of measurement to conformance assessment activities in the electrotechnical sector*

### 201.3 Terms and definitions

For the purposes of this document, the terms and definitions given in IEC 60601-1:2005+AMD1:2012+AMD2:2020, and the following apply.

ISO and IEC maintain terminology databases for use in standardization at the following addresses:

- ISO Online browsing platform: available at <https://www.iso.org/obp>
- IEC Electropedia: available at <https://www.electropedia.org/>

#### 201.3.201 accompanying information

information accompanying or *marked* on a medical device or *accessory* for the user or those accountable for the installation, use, *processing*, maintenance, decommissioning and disposal of the medical device or *accessory*, particularly regarding safe use

Note 1 to entry: The *accompanying information* shall be regarded as part of the medical device or *accessory*.

Note 2 to entry: The *accompanying information* can consist of the label, *marking*, *instructions for use*, *technical description*, installation manual, quick reference guide, etc.

Note 3 to entry: *Accompanying information* is not necessarily a written or printed document but could involve auditory, visual, or tactile materials and multiple media types (e.g., CD/DVD-ROM, USB stick, website).

[SOURCE: ISO 20417:2021, 3.2, modified — deleted note 4.]

**201.3.202**

**acknowledged**

state of an *alarm system* initiated by *operator* action, where the auditory *alarm signal* associated with a currently active *alarm condition* is inactivated until the *alarm condition* no longer exists or until a predetermined time interval has elapsed

Note 1 to entry: *Acknowledged* only affects *alarm signals* that are active at the time of the *operator* action.

[SOURCE: IEC 60601-1-8:2006+AMD1:2012+AMD2:2020, 3.37]

**201.3.203**

**airway pressure**

$P_{aw}$

pressure at the *patient-connection port* or at the distal *outlet* of the equipment where there is no *patient-connection port*

Note 1 to entry: The *airway pressure* can be derived from pressure measurements made anywhere within the equipment.

[SOURCE: ISO 4135:2022, 3.1.4.41.1]

**201.3.204**

**alarm condition delay**

time from the occurrence of a triggering event either in the *patient*, for *physiological alarm conditions*, or in the equipment, for *technical alarm conditions*, to when the *alarm system* determines that an *alarm condition* exists

[SOURCE: IEC 60601-1-8:2006+AMD1:2012+AMD2:2020, 3.2]

**201.3.205**

**alarm limit**

threshold used by an *alarm system* to determine an *alarm condition*

[SOURCE: IEC 60601-1-8:2006+AMD1:2012+AMD2:2020, 3.3]

**201.3.206**

**alarm off**

state of indefinite duration in which an *alarm system* or part of an *alarm system* does not generate *alarm signals*

[SOURCE: IEC 60601-1-8:2006+AMD1:2012+AMD2:2020, 3.4]

**201.3.207**

**alarm paused**

state of limited duration in which the *alarm system* or part of the *alarm system* does not generate *alarm signals*

[SOURCE: IEC 60601-1-8:2006+AMD1:2012+AMD2:2020, 3.5]

**201.3.208**

**alarm settings**

*alarm system* configuration, including but not limited to:

- *alarm limits*;
- the characteristics of any *alarm signal* inactivation states; and
- the values of variables or parameters that determine the function of the *alarm system*

Note 1 to entry: Some algorithmically-determined *alarm settings* can require time to be determined or re-determined.

[SOURCE: IEC 60601-1-8:2006+AMD1:2012+AMD2:2020, 3.8]

### 201.3.209

#### **artificial ventilation**

intermittent elevation of the pressure in the *patient's* airway relative to that in the *lungs* by external means with the intention of augmenting, or totally controlling, the *ventilation* of a *patient*

EXAMPLE Means used to provide *artificial ventilation* are: manual resuscitation; mouth-to-mouth resuscitation; automatic *ventilation*; mechanical *ventilation*.

Note 1 to entry: Common classifications of areas of application of *artificial ventilation* are: emergency; transport; home-care; anaesthesia; critical care; rehabilitation.

Note 2 to entry: Classifications used to denote means used for *artificial ventilation* include: positive-pressure; negative-pressure; gas-powered; *operator*-powered; electrically-powered.

Note 3 to entry: Negative-pressure *ventilation* elevates the relative pressure in the airway by intermittently lowering the pressure in the *lungs*.

[SOURCE: ISO 19223:2019, 3.1.10]

### 201.3.210

#### **assured inflation-type rate**

number of assured *inflation-type* initiations in a specified period of time, expressed as breaths per minute

[SOURCE: ISO 19223:2019, 3.5.2.1, modified — deleted notes.]

### 201.3.211

#### **attack**

attempt to destroy, expose, alter, disable, steal or gain unauthorized access to or make unauthorized use of an asset

[SOURCE: IEC 81001-5-1:2021, 3.5]

### 201.3.212

#### **audio off**

state of indefinite duration in which the *alarm system* or part of the *alarm system* does not generate an auditory *alarm signal*

[SOURCE: IEC 60601-1-8:2006+AMD1:2012+AMD2:2020, 3.12]

### 201.3.213

#### **audio paused**

state of limited duration in which the *alarm system* or part of the *alarm system* does not generate an auditory *alarm signal*

[SOURCE: IEC 60601-1-8:2006+AMD1:2012+AMD2:2020, 3.13]

### 201.3.214

#### **BAP**

quantity by which the baseline *airway pressure* is set to be positively offset from the ambient pressure

[SOURCE: ISO 19223:2019, 3.10.2, modified — deleted notes.]

### 201.3.215

#### **biocompatibility**

ability to be in contact with a living system without producing an unacceptable adverse effect

Note 1 to entry: A medical device or *accessory* may produce some level of adverse effect, but that level may be determined to be acceptable when considering the benefits provided by the medical device or *accessory*.

[SOURCE: ISO 18562-1:2023, 3.6]

**201.3.216**

**body temperature and pressure saturated**

**BTPS**

ambient atmospheric pressure, at a temperature of 37 °C, and a relative humidity of 100 %

[SOURCE: ISO 4135:2022, 3.1.1.7]

**201.3.217**

**breathing system**

pathways through which gas flows to or from the *patient* at respiratory pressures and continuously or intermittently in fluid communication with the *patient's* respiratory tract during any form of *artificial ventilation* or respiratory therapy

[SOURCE: ISO 4135:2022, 3.6.1.1, modified — deleted notes.]

**201.3.218**

**breathing system filter**

**BSF**

device intended to reduce transmission of particulates, including microorganisms, in *breathing systems*

[SOURCE: ISO 23328-2:2002, 3.2]

**201.3.219**

**cleaning**

removal of contaminants to the extent necessary for further *processing* or for *intended use*

Note 1 to entry: *Cleaning* consists of the removal of adherent soil (e.g. blood, protein substances, and other debris) from the surfaces, crevices, serrations, joints, and lumens of a medical device by a manual or automated process that prepares the items for safe handling or further *processing*.

[SOURCE: ISO 17664-2:2021, 3.1, modified — replaced 'and/or' with 'or'.]

**201.3.220**

**communicator**

**COM**

function of the *alarm system* that generates *alarm signals* to notify an *operator* (e.g. to the presence of an *alarm condition*)

Note 1 to entry: A *communicator* can receive an *operator* response.

Note 2 to entry: An *operator* response is not limited to direct *operator* action.

[SOURCE: IEC 60601-1-8:2006+AMD1:2012+AMD2:2020, 3.46, modified — deleted note 3.]

**201.3.221**

**connector**

fitting to join two or more components

EXAMPLE *Connectors* for a *low-pressure hose assembly* are any of a range of mating components intended to maintain gas specificity by the allocation of a set of different diameters to the mating *connectors* for each particular gas.

[SOURCE: ISO 4135:2022, 3.1.4.5]

**201.3.222****CPAP****continuous positive airway pressure**

*ventilation-mode* or sleep-apnoea breathing-therapy mode in which the *patient* breathes continuously at a set *airway-pressure* level, above ambient pressure

[SOURCE: ISO 19223:2019, 3.11.15, modified — deleted notes.]

**201.3.223****cybersecurity**

state where information and systems are protected from unauthorized activities, such as access, use, disclosure, disruption, modification, or destruction to a degree that the related *risks* to violation of confidentiality, integrity, and availability are maintained at an acceptable level throughout the life cycle

[SOURCE: IEC 81001-5-1:2021, 3.30]

**201.3.224** **$\Delta$  inspiratory pressure**

differential *airway pressure* relative to baseline *airway pressure* during an *inflation phase*

Note 1 to entry: In addition to its direct reference, this term or an appropriate *symbol* may be used, in context or by qualification, to designate this concept as a set quantity or a measured quantity.

Note 2 to entry: There is currently no agreed convention as to whether an inspiratory pressure is always to be expressed as an absolute quantity relative to ambient pressure or an absolute quantity for one group of *inflation-types* and relative for another. This has unacceptable *patient-safety* implications that need to be addressed in a vocabulary of *lung ventilation*. The *symbol*,  $\Delta$ , is currently sometimes used as a prefix to make this distinction, and that convention has been adopted as a requirement in this document. Without a prefix, or any other indication, respiratory pressures are always to be considered to be relative to ambient pressure. The addition of a  $\Delta$  prefix, is used to indicate a pressure that is relative to the set *BAP* level. In *ventilation-modes* where there is a second, higher, baseline *airway pressure*, then the prefix for a pressure relative to that higher-pressure level becomes  $\Delta H$ . These prefixes are applicable to relevant terms, *symbols* and displayed values but not to *inflation-types*.

Note 3 to entry: The sum of the set *BAP* level and the  $\Delta$  *inspiratory pressure* equals the inspiratory pressure. This applies to both settings and measurements of this parameter.

[SOURCE: ISO 19223:2019, 3.6.7, modified — deleted note 4.]

**201.3.225****disinfection**

*process* to inactivate viable microorganisms to a level previously specified as being appropriate for a defined purpose

[SOURCE: ISO 17664-2:2021, 3.5]

**201.3.226****distributed alarm system****DAS**

*alarm system* that involves more than one item of equipment of a *ME system* intended for delivery of *alarm conditions* with technical confirmation

Note 1 to entry: The parts of a *distributed alarm system* can be widely separated in distance.

Note 2 to entry: A *distributed alarm system* is intended to notify *operators* of the existence of an *alarm condition*.

Note 3 to entry: For the purposes of this document, technical confirmation means that each element of a *distributed alarm system* confirms or guarantees the successful delivery of the *alarm condition* to the next element

or appropriate *technical alarm conditions* are created as described in IEC 60601-1-8:2006+AMD1:2012+AMD2:2020, 6.11.2.2.1.

[SOURCE: IEC 60601-1-8:2006+AMD1:2012+AMD2:2020, 3.17]

### 201.3.227

#### EMS environment

##### emergency medical services environment

actual conditions and settings, in which *operators* interact with the *ME equipment* or *ME system*, in and around the scene of an emergency outside of a *professional healthcare facility* where a *patient* can be given medical care, basic or advanced life support as well as during professional transport to a *professional healthcare facility* or between *professional healthcare facilities*

EXAMPLE 1 Responding to and providing life support at the scene of an emergency to a *patient* reported as experiencing injury or illness in a pre-hospital setting, and transporting the *patient*, while continuing such life support care, to an appropriate *professional healthcare facility* for further care.

EXAMPLE 2 Providing monitoring, treatment or diagnosis during transport between *professional healthcare facilities*.

Note 1 to entry: For the purposes of this standard, use of equipment intended for the *EMS environment* and temporarily used in the *home healthcare environment* by emergency medical personnel is considered use in the *EMS environment*.

Note 2 to entry: For the purposes of this standard, the *operators* of equipment intended for the *EMS environment* are presumed to be professional medical personnel or personnel with relevant specialized training.

Note 3 to entry: *Professional healthcare facilities* include hospitals, physician offices, freestanding surgical centres, dental offices, freestanding birthing centres, limited care facilities, first aid rooms or rescue rooms and multiple treatment facilities.

Note 4 to entry: Emergency medical services are known by various names in different countries and regions.

Note 5 to entry: For the purposes of this standard, transport includes road, rotary and fixed-wing ambulances.

[SOURCE: IEC 60601-1-12:2014+AMD1:2020, 3.1]

### 201.3.228

#### EMS ventilator

ventilator for emergency medical services environment  
*ventilator* intended for use in the *EMS environment*

### 201.3.229

#### end-expiratory flow

expiratory flow at the point of initiation of an *inflation* or an inspiration

[SOURCE: ISO 19223:2019, 3.7.6, modified — deleted notes.]

### 201.3.230

#### essential function

function or capability that is required to maintain *basic safety*, *essential performance*, a minimum of clinical functionality as specified by the *manufacturer*, and operational availability for the medical device

Note 1 to entry: *Essential functions* include, but are not limited to, the safety instrumented function (*basic safety* and *essential performance*), the control function and the availability of urgently needed functions and such allowing the *operator* to view and manipulate the medical device safely with the most urgently needed

performance (operational availability). The loss of *essential function* is commonly termed loss of protection, loss of control and loss of view respectively.

Note 2 to entry: The term is derived from IEC 62443-4-2:2019, 3.1.20, and has been refined for the purpose and scope of this document.

[SOURCE: IEC/TR 60601-4-5:2021, 3.10]

### 201.3.231

#### **essential principles**

#### **essential principles of safety and performance**

fundamental high-level requirements that when complied with ensure a medical device is safe and performs as intended

[SOURCE: ISO 16142-1:2016, 3.3]

### 201.3.232

#### **exhaust port**

port of the medical equipment or device from which gas is discharged to the atmosphere during *normal use*, either directly or via an anaesthetic gas scavenging system

[SOURCE: ISO 19223:2019, 3.14.2]

### 201.3.233

#### **expiratory pause**

portion of the *expiratory phase* from the end of expiratory flow to the start of inspiratory flow

[SOURCE: ISO 19223:2019, 3.4.4, modified — deleted notes.]

### 201.3.234

#### **expiratory phase**

interval from the start of expiratory flow to the start of inspiratory flow within a *respiratory cycle*

[SOURCE: ISO 19223:2019, 3.4.2, modified — deleted notes.]

### 201.3.235

#### **false positive alarm condition**

presence of an *alarm condition* when no valid triggering event has occurred in the *patient*, the equipment or the *alarm system*

Note 1 to entry: A *false positive alarm condition* can be caused by spurious information produced by the *patient*, the *patient*-equipment interface, other equipment or the *alarm system* itself.

[SOURCE: IEC 60601-1-8:2006+AMD1:2012+AMD2:2020, 3.21]

### 201.3.236

#### **firecall**

method established to provide emergency access to a secure medical device

Note 1 to entry: In an emergency situation, unprivileged users can gain access to key systems to correct the problem. When a *firecall* is used, there is usually a review *process* to ensure that the access was used properly to correct a problem. These methods generally either provide a one-time use user identifier (ID) or one-time password or other suitable measures.

Note 2 to entry: Also referred to as "break glass" feature.

[SOURCE: IEC/TR 60601-4-5:2021, 3.11]

**201.3.237**

**flow-direction-sensitive component**

component or *accessory* through which gas flow has to be in one direction only for proper functioning or *patient* safety

[SOURCE: ISO 4135:2022, 3.1.4.15]

**201.3.238**

**fresh gas**

respirable gas delivered to a *ventilator breathing system*

[SOURCE: ISO 4135:2022, 3.1.1.16, modified — Added 'ventilator' and deleted notes.]

**201.3.239**

**gas intake port**

port through which gas is drawn for use by the *patient*

Note 1 to entry: Gas is drawn at a sub-ambient pressure at a *gas intake port*, in apposition to an *inlet*, at which gas is provided by a medical gas supply system.

[SOURCE: ISO 4135:2022, 3.1.4.21]

**201.3.240**

**gas output port**

port of the *ventilator* through which gas is delivered at respiratory pressures to an *operator*-detachable part of the *ventilator breathing system*

[SOURCE: ISO 19223:2019, 3.14.3]

**201.3.241**

**gas pathway**

interior surfaces, over which gases or liquids that can be inspired, in a medical device bounded by the ports through which gases or liquids enter and leave the medical device including the *patient* interface or the interior surfaces of *enclosures* that are in contact with gases or liquids that can be inspired

Note 1 to entry: *Patient* contact surfaces such as the outer surfaces of a tracheal tube or the cushion of a *mask* are evaluated according to the ISO 10993 series.

EXAMPLE 1 The *ventilator breathing system*, *inlet* filter, gas mixer, blower and internal piping.

EXAMPLE 2 Enclosed chamber of an incubator including the mattress or the inner surface of an oxygen hood.

EXAMPLE 3 The inner surfaces of breathing tubes, tracheal tubes or *masks* and mouthpieces.

[SOURCE: ISO 18562-1:2023, 3.11]

**201.3.242**

**gas return port**

port of the *ventilator* through which gas is returned at respiratory pressures through an *operator*-detachable part of the *ventilator breathing system*, from the *patient-connection port*

[SOURCE: ISO 19223:2019, 3.14.4]

**201.3.243**

**healthcare professional**

<adj> appropriately trained, knowledgeable, and skilled, providing systematic preventive, curative, promotional or rehabilitative healthcare services

[SOURCE: ISO 80601-2-12:2023, 201.3.247]

#### 201.3.244

##### heat and moisture exchanger

##### HME

device intended to retain a portion of the *patient's* expired moisture and heat, and return it to the respiratory tract during inspiration

[SOURCE: ISO 9360-1:2000, 3.1]

#### 201.3.245

##### high-pressure inlet

*inlet* to which gas is supplied at a pressure exceeding 100 kPa above ambient

Note 1 to entry: The phrases 'low-pressure' and 'high-pressure' are used differently in various contexts, including *breathing system* pressures (typically less than 10 kPa), terminal *outlet* pressures (less than 600 kPa), manifold pressures (typically up to 3 000 kPa) and cylinder pressures (typically less than 30 000 kPa).

[SOURCE: ISO 4135:2022, 3.1.4.24]

#### 201.3.246

##### home healthcare environment

dwelling place in which a *patient* lives or other places where *patients* are present, excluding *professional healthcare facility* environments where *operators* with medical training are continually available when *patients* are present

EXAMPLE In a car, bus, train, boat or plane, in a wheelchair or walking outdoors.

Note 1 to entry: *Professional healthcare facilities* include hospitals, physician offices, freestanding surgical centres, dental offices, freestanding birthing centres, limited care facilities, first aid rooms or rescue rooms, multiple treatment facilities and emergency medical services.

Note 2 to entry: For the purpose of this collateral standard, nursing homes are considered *home healthcare environments*.

Note 3 to entry: Other places where a *patient* is present include the outdoor environment, while working and in vehicles.

[SOURCE: IEC 60601-1-11:2015+AMD1:2020, 3.1]

#### 201.3.247

##### humidifier

device that adds water in the form of droplets or vapour, or both, to the inspired gas

Note 1 to entry: This term includes vaporising, bubble-through and ultrasonic *humidifiers* and active *heat and moisture exchangers*.

[SOURCE: ISO 4135:2022, 3.7.2.1]

#### 201.3.248

##### I:E ratio

ratio of the *inspiratory time* to the expiratory time in a *respiratory cycle*

[SOURCE: ISO 19223:2019, 3.4.19, modified — deleted notes.]

**201.3.249**

**immunity**

ability of *ME equipment* or an *ME system* to perform without degradation in the presence of an electromagnetic disturbance

[SOURCE: IEC 60601-1-2:2014+AMD1:2020, 3.8]

**201.3.250**

**inflation**

*ventilator* action intended to increase the volume of gas in the *lungs* by the application of an elevated-pressure waveform to the *patient-connection port* until a specified termination criterion is met

[SOURCE: ISO 19223:2019, 3.3.1, modified — deleted notes.]

**201.3.251**

**inflation phase**

interval from the start of the rise in *airway pressure* resulting from the initiation of an *inflation* to the start of the expiratory flow resulting from its termination

[SOURCE: ISO 19223:2019, 3.4.10, modified — deleted notes.]

**201.3.252**

**inflation-type**

*inflation* characterized by its temporal delivery pattern following initiation, and its termination criteria

[SOURCE: ISO 19223:2019, 3.3.2, modified — deleted notes.]

**201.3.253**

**information supplied by the manufacturer**

information related to the identification and use of a medical device or *accessory*, in whatever form provided, intended to ensure the safe and effective use of the medical device or *accessory*

Note 1 to entry: For the purposes of this document, e-documentation is included in *information supplied by the manufacturer*.

Note 2 to entry: For the purposes of this document, shipping documents and promotional material are excluded from *information supplied by the manufacturer*. However, some authorities having jurisdiction can consider such supplemental information as *information supplied by the manufacturer*.

Note 3 to entry: The primary purpose of *information supplied by the manufacturer* is to identify the medical device and its *manufacturer*, and provide essential information about its safety, performance, and appropriate use to the user or other relevant persons.

[SOURCE: ISO 20417:2021, 3.10, modified — deleted note 4.]

**201.3.254**

**inlet**

opening through which gas or other material is pushed by an elevated upstream pressure

[SOURCE: ISO 4135:2022, 3.1.4.26, modified — deleted note.]

**201.3.255**

**inspiratory pause**

interval from the end of inspiratory flow to the start of expiratory flow within a *respiratory cycle*

[SOURCE: ISO 19223:2019, 3.4.12, modified — deleted notes.]

**201.3.256****inspiratory time** $t_i$ duration of an *inflation phase* or inspiratory phase

[SOURCE: ISO 19223:2019, 3.4.8, modified — deleted notes.]

**201.3.257****inspiratory volume** $V_{\text{insp}}$ volume of gas delivered through the *patient-connection port* during an inspiratory phase or *inflation phase*

[SOURCE: ISO 19223:2019, 3.8.3, modified — deleted notes.]

**201.3.258****instructions for use****IFU**portion of the *accompanying information* that is essential for the safe and effective use of a medical device or *accessory* directed to the user of the medical device

Note 1 to entry: For the purposes of this document, a user can be either a lay user or professional user with relevant specialized training.

Note 2 to entry: For the purposes of this document, instructions for the professional *processing* between uses of a medical device or *accessory* can be included in the *instructions for use*.

Note 3 to entry: For the purposes of this document, information indicated on a graphical user interface (GUI) is considered as appearing on the item.

Note 4 to entry: The *instructions for use*, or portions thereof, can be located on the display of a medical device or *accessory*.

Note 5 to entry: Medical devices or *accessories* that can be used safely and effectively without *instructions for use* are exempted from having *instructions for use* by some authorities having jurisdiction.

[SOURCE: ISO 20417:2021, 3.11, modified — deleted note 6.]

**201.3.259****intelligent alarm system***alarm system* that makes logical decisions based on monitored information without *operator* intervention

EXAMPLE 1 An *alarm system* that changes priority based on the rate of change of a monitored variable.

EXAMPLE 2 An *alarm system* that suppresses an *alarm condition* when a related *alarm condition* of higher priority has recently generated an *alarm signal*.

[SOURCE: IEC 60601-1-8:2006+AMD1:2012+AMD2:2020, 3.24]

**201.3.260****latching alarm signal***alarm signal* that continues to be generated after its triggering event no longer exists until stopped by deliberate *operator* action

[SOURCE: IEC 60601-1-8:2006+AMD1:2012+AMD2:2020, 3.26]

**201.3.261**

**low-pressure hose assembly**

assembly consisting of a flexible hose with permanently attached gas-specific *inlet* and *outlet connectors* and designed to conduct a medical gas at pressures less than 1 400 kPa

Note 1 to entry: The phrases 'low-pressure' and 'high-pressure' are used differently in various contexts, including *breathing system* pressures (typically less than 10 kPa), terminal *outlet* pressures (less than 1 400 kPa), manifold pressures (typically up to 3 000 kPa) and cylinder pressures (typically less than 30 000 kPa).

[SOURCE: ISO 4135:2022, 3.2.3.1]

**201.3.262**

**lung**

each of the pair of compliant organs within the ribcage (thorax), bounded by the terminal bronchiole and the visceral pleura, which during *ventilation* provide gas/blood interfaces that enable oxygen from the gas to pass into the blood and carbon dioxide to be removed

[SOURCE: ISO 19223:2019, 3.1.16, modified — deleted notes.]

**201.3.263**

**marking**

information, in text or graphical format, durably affixed, printed, etched (or equivalent) to a medical device or *accessory*

Note 1 to entry: For the purposes of this document, the term *marked* is used to designate the corresponding act.

Note 2 to entry: For the purposes of this document, *marking* is different from 'direct *marking*' as commonly described in unique device identification (UDI) standards and regulations. A UDI 'direct *marking*' is a type of *marking*.

[SOURCE: ISO 20417:2021, 3.16, modified — deleted note 3.]

**201.3.264**

**mask**

device which provides a non-invasive interface between the *patient's airway* and a *patient-connection port* or other connection to a source of respirable gas

[SOURCE: ISO 4135:2022, 3.8.6.4]

**201.3.265**

**maximum limited pressure**

$P_{LIM,max}$

highest *airway pressure* that can occur during *normal use* or under *single fault condition*

[SOURCE: ISO 4135:2022, 3.1.4.41.3]

**201.3.266**

**maximum working pressure**

$P_{W,max}$

highest *airway pressure* that can be generated by the *ventilator* during *intended use* and *normal condition*

[SOURCE: ISO 19223:2019, 3.13.4, modified — deleted notes.]

**201.3.267**

**medical gas pipeline system**

combination of a supply system, a monitoring and *alarm system* and a pipeline distribution system with terminal units for provision of medical gases or vacuum

[SOURCE: ISO 4135:2022, 3.2.1.1]

### 201.3.268

#### **monitoring equipment**

equipment or part that measures and indicates the value of a variable to the *operator*

Note 1 to entry: *Monitoring equipment* includes devices that are not electrical in operation, such as a pressure gauge.

Note 2 to entry: The value can be displayed continually or intermittently.

Note 3 to entry: The *monitoring equipment* can be primarily intended for detection of an *alarm condition* or for external communication.

[SOURCE: ISO 4135:2022, 3.11.1.3, modified —replaced "user" with "operator".]

### 201.3.269

#### **operator group**

subset of *operators* who are differentiated from other *operators* by factors that are likely to influence their interactions with the *ME equipment*

Note 1 to entry: Attributes of *operator groups* can include age, culture, expertise.

[SOURCE: IEC 60601-1-6:2010+AMD1:2013+AMD2:2020, 3.3]

### 201.3.270

#### **operator interface**

means by which the *operator* and the *ME equipment* interact

Note 1 to entry: The *accompanying documents* are considered part of the *ME equipment* and its *operator interface*.

Note 2 to entry: *Operator interface* includes all the elements of the *ME equipment* with which the *operator* interacts including the physical aspects of the *ME equipment* as well as visual, auditory, tactile displays and is not limited to a software interface.

Note 3 to entry: For the purposes of this standard, the *manufacturer* may treat the combination of *ME equipment* and other equipment as a single *operator interface*.

Note 4 to entry: See IEC 62366-1:2015+AMD1:2020, 3.26.

[SOURCE: IEC 60601-1-6:2010+AMD2:2020, 3.1]

### 201.3.271

#### **operator's position**

intended position of the *operator* with respect to the *alarm signal* generating part of the *alarm system*

Note 1 to entry: A *distributed alarm system* can have multiple *operator's positions*.

[SOURCE: IEC 60601-1-8:2006+AMD1:2012+AMD2:2020, 3.30]

### 201.3.272

#### **outlet**

opening through which gas leaves a device or component

[SOURCE: ISO 4135:2022, 3.1.4.40]

**201.3.273**

**patient-connection port**

port of a *breathing system* intended for connection to an airway device

Note 1 to entry: The *patient-connection port* is the end of the *breathing system* proximal to the *patient*.

Note 2 to entry: The *patient-connection port* is typically a *connector* suitable for connection to an airway device such as a tracheal tube, tracheostomy tube, face *mask* or supralaryngeal airway.

Note 3 to entry: Current product standards typically specify that the *patient-connection port* is required to be in the form of specific standardized *connectors*, for example, a *connector* conforming to ISO 5356-1.

[SOURCE: ISO 4135:2022, 3.1.4.41, modified — deleted note 4.]

**201.3.274**

**PEEP**

**positive end-expiratory pressure**

<actual and measured value> respiratory pressure at the end of an *expiratory phase*

[SOURCE: ISO 19223:2019, 3.10.4, modified — deleted notes.]

**201.3.275**

**physiologic closed-loop control system**

part of *ME equipment* or *ME system* used to adjust a physiologic variable (y) relative to a command variable (c) using a feedback variable (f)

[SOURCE: IEC 60601-1-10:2007+AMD2:2020, 3.19]

**201.3.276**

**physiological alarm condition**

*alarm condition* arising from a monitored *patient*-related variable

EXAMPLE 1 High exhaled anaesthetic agent concentration.

EXAMPLE 2 Low exhaled *tidal volume*.

EXAMPLE 3 Low oxygen saturation measured by pulse oximetry.

EXAMPLE 4 High arterial pressure.

EXAMPLE 5 High heart rate.

[SOURCE: IEC 60601-1-8:2006+AMD1:2012+AMD2:2020, 3.31]

**201.3.277**

**pressure-control**

*inflation-type* that acts to generate a constant inspiratory pressure at a set level, after a set rise time

[SOURCE: ISO 19223:2019, 3.3.4, modified — deleted notes.]

**201.3.278**

**primary operating function**

function that involves user interaction that is related to the safety of the medical device

Note 1 to entry: Often a *primary operating function* is interacted with by a series of tasks that can be broken down into a series of user interactions.

Note 2 to entry: The concept of safety includes loss or degradation of performance resulting in an unacceptable *risk* to the *patient*, including use error that prevents the user from effectively using the medical device to achieve its intended medical purpose. In IEC 60601-1, this is referred to as *essential performance*.

[SOURCE: IEC 62366-1:2015+AMD1:2020, 3.11]

### 201.3.279

#### processing

<preparation of medical device, *accessory*> activity to prepare a new or used medical device and *accessory* for its *intended use*

[SOURCE: ISO 20417:2021, 3.20]

### 201.3.280

#### professional healthcare facility

facility that is continually staffed by suitably trained *healthcare professional operators*

EXAMPLE Hospitals, physician offices, freestanding surgical centres, dental offices, freestanding birthing centres, limited care facilities, first aid rooms or rescue rooms, multiple treatment facilities and emergency medical services

[SOURCE: ISO 80601-2-12:2023, 201.3.285, modified — deleted note.]

### 201.3.281

#### protection device

part or function of *ME equipment* that, without intervention by the *operator*, protects the *patient* from hazardous output due to incorrect delivery of energy or substances

[SOURCE: ISO 80601-2-12:2023, 201.3.286]

### 201.3.282

#### respiratory cycle

complete sequence of respiratory events that leads to an increase, followed by a corresponding decrease, of gas volume in the *lung* regardless of how it is generated

[SOURCE: ISO 19223:2019, 3.4.16, modified — deleted notes.]

### 201.3.283

#### security level

level corresponding to the required set of countermeasures and inherent *cybersecurity* properties of devices and systems for a zone or conduit based on assessment of *risk* for the zone or conduit

[SOURCE: IEC/TR 60601-4-5:2021, 3.23]

### 201.3.284

#### set rate

number of assured *inflations* that are set to occur in a specified period of time, expressed as breaths per minute

[SOURCE: ISO 19223:2019, 3.5.1.1, modified — deleted notes and examples.]

### 201.3.285

#### single use

<medical device, *accessory*> intended by the *manufacturer* to be used on an individual *patient* or specimen during a single *procedure* and then disposed of

Note 1 to entry: A *single use* medical device or *accessory* is not intended by its *manufacturer* to be further *processed* and used again.

[SOURCE: ISO 20417:2021, 3.26]

**201.3.286**

**software item**

any identifiable part of a computer program, i.e., source code, object code, control code, control data, or a collection of these items

Note 1 to entry: Three terms identify the software decomposition. The top level is the software system. the lowest level that is not further decomposed is the software unit. All levels of composition, including the top and bottom levels, can be called *software items*. A *software system*, then, is composed of one or more *software items*, and each *software item* is composed of one or more software units or decomposable *software items*. The responsibility is left to the *manufacturer* to provide the definition and granularity of the *software items* and software units.

[SOURCE: IEC 62304:2006+AMD1:2015, 3.25, modified — deleted note 2.]

**201.3.287**

**spontaneous breath rate**

total number of spontaneous breaths initiated in a specified period of time, expressed as breaths per minute

[SOURCE: ISO 19223:2019, 3.5.1.3, modified — deleted notes.]

**201.3.288**

**standard temperature and pressure, dry**

**STPD**

pressure of 101,325 kPa at a temperature of 20 °C, dry

Note 1 to entry: See also rationale for 201.5.101.2.

[SOURCE: ISO 4135:2022, 3.1.1.8, modified — added note.]

**201.3.289**

**sterile**

free from viable microorganisms

[SOURCE: ISO 20417:2021, 3.28]

**201.3.290**

**sterilization**

*process* used to render product free from viable microorganisms

Note 1 to entry: In a *sterilization process*, the nature of microbial inactivation is exponential and thus the survival of a microorganism on an individual item can be expressed in terms of probability. While this probability can be reduced to a very low number, it can never be reduced to zero.

[SOURCE: ISO 17664-1:2021, 3.17]

**201.3.291**

**symbol**

graphical representation appearing on the label or associated documentation of a medical device that communicates characteristic information without the need for the supplier or receiver of the information to have knowledge of the language of a particular nation or people

Note 1 to entry: The *symbol* can be an abstract pictorial or a graphical representation, or one that uses familiar objects, including alphanumeric characters (with sufficient justification).

[SOURCE: ISO 20417:2021, 3.29]

**201.3.292****system recovery**

method for fault handling via an automatic restart of *PESS* for parts of the *ME equipment* or for the complete *ME equipment*

Note 1 to entry: There is guidance or rationale for this definition contained in Clause AA.2.

[SOURCE: ISO 80601-2-12:2023, 201.3.298]

**201.3.293****technical alarm condition**

*alarm condition* arising from a monitored equipment-related or *alarm system*-related variable

EXAMPLE 1 An electrical, mechanical or other failure.

EXAMPLE 2 A failure of a sensor or component (unsafe voltage, high impedance, signal impedance, artefact, noisy signal, disconnection, calibration error, tubing obstruction, etc.).

EXAMPLE 3 An algorithm that cannot classify or resolve the available data.

[SOURCE: IEC 60601-1-8:2006+AMD1:2012+AMD2:2020, 3.36]

**201.3.294****technical description**

portion of the *accompanying information* directed to the *responsible organization* and *service personnel* that is essential for preparation for the first use and safe use, maintenance or repair as well as *processing*, transport or storage for the *expected service life* of a medical device

Note 1 to entry: The *technical description* may be included in the *instructions for use*.

[SOURCE: ISO 20417:2021, 3.30, modified — replaced 'expected lifetime' with 'expected service life' and deleted note 2.]

**201.3.295****tidal volume**

$V_T$

volume of gas that enters and leaves the *lung* during a breath

[SOURCE: ISO 19223:2019, 3.8.1, modified — deleted notes.]

**201.3.296****total respiratory rate**

number of *respiratory cycles* in a specified period of time, expressed as breaths per minute

[SOURCE: ISO 19223:2019, 3.5.1.2, modified — deleted notes.]

**201.3.297****use environment**

actual conditions and setting in which *operators* interact with the equipment

Note 1 to entry: The conditions of use or attributes of the *use environment* can include hygienic requirements, frequency of use, location, lighting, noise, temperature, mobility, and degree of internationalization. Social attributes such as team versus individual, chaotic versus calm, stress level and length of shift can also play a role.

[SOURCE: IEC 62366-1:2015+AMD1:2020, 3.20, modified — replaced "users" with "operators" and "medical device" with "equipment".]

**201.3.298**

**use scenario**

specific sequence of tasks performed by a specific *operator* in a specific *use environment* and any resulting response of the equipment

[SOURCE: IEC 62366-1:2015+AMD1:2020, 3.22, modified — replaced "user" with "*operator*" and "medical device" with "equipment".]

**201.3.299**

**validation**

confirmation, through the provision of *objective evidence*, that specified requirements have been fulfilled

Note 1 to entry: The *objective evidence* needed for a *validation* is the result of a test or other form of determination such as performing alternative calculations or reviewing documents.

Note 2 to entry: The word "*validated*" is used to designate the corresponding status.

Note 3 to entry: The use conditions for *validation* can be real or simulated.

[SOURCE: ISO 9000:2015, 3.8.13]

**201.3.300**

**ventilation**

cyclical movement of a respirable gas into and out of the *lungs*

Note 1 to entry: This might be by external or spontaneous means, or by a combination of both.

[SOURCE: ISO 19223:2019, 3.1.9, modified — deleted note 2.]

**201.3.301**

**ventilation-mode**

specified manner in which a *ventilator* performs its ventilatory function when connected to a *patient*

[SOURCE: ISO 19223:2019, 3.11.2, modified — deleted notes.]

**201.3.302**

**ventilator breathing system**

**VBS**

pathways through which gas flows to or from the *patient* at respiratory pressures, bounded by the port through which respirable gas enters, the *patient-connection port* and the *gas exhaust port*

Note 1 to entry: These pathways typically extend within and outside the body of the *ventilator*, with those outside being *operator-detachable*.

Note 2 to entry: The port of entry of a respirable gas into the *ventilator breathing system* can be inside the body of the *ventilator* and should not be confused with an external connection port into which respirable gas enters before being reduced to respirable pressures.

[SOURCE: ISO 19223:2019, 3.1.18, modified — deleted notes 3 and 4.]

**201.3.303**

**ventilator-dependent**

<*patient*> dependent upon *artificial ventilation* in order to prevent serious deterioration of health or death

Note 1 to entry: A *ventilator-dependent patient* cannot breathe well enough to maintain life-sustaining levels of oxygen and carbon dioxide in the blood.

EXAMPLE *Patients* with Duchenne muscular dystrophy or other degenerative disease resulting in their unsupported respiratory effort being insufficient to sustain life.

[SOURCE: ISO 4135:2022, 3.1.5.19]

**201.3.304  
volume-control**

*inflation-type* that generates inspiratory flow to a selected flow-waveform, for a set *inspiratory-time*, or until the set volume has been delivered

[SOURCE: ISO 19223:2019, 3.3.3, modified — deleted notes.]

**201.4 General requirements**

IEC 60601-1:2005+AMD1:2012+AMD2:2020, Clause 4 applies, except as follows:

**201.4.3 Essential performance**

*Addition:*

**201.4.3.101 Additional requirements for essential performance**

NOTE There is guidance or rationale for this subclause contained in Clause AA.2.

Additional *essential performance* requirements are found in the subclauses listed in Table 201.101.

**Table 201.101 — Distributed essential performance requirements**

Requirement	Subclause
Delivery of <i>ventilation</i> at the <i>patient-connection port</i> within the <i>alarm limits</i> set by the <i>operator</i> or generation of an <i>alarm condition</i>	
expired volume, if equipped	201.12.4.103
gas supply failure	201.13.102
high <i>airway pressure</i>	201.12.4.105
<i>internal electrical power source</i> depleted	212.8.2
inspiratory oxygen concentration, if equipped	201.12.4.101
<i>system recovery</i>	201.4.3.102
Subclauses 202.4.3.1 and 202.8.1.101 indicates methods of evaluating delivery of <i>ventilation</i> as acceptance criteria following specific tests required by this document.	

**201.4.3.102 System recovery**

NOTE 1 There is guidance or rationale for this subclause contained in Clause AA.2.

- a) Following a malfunction, the *EMS ventilator* shall attempt to perform a *system recovery* to restore *essential performance* of the *EMS ventilator*.

NOTE 2 An *EMS ventilator* or a subassembly function can become disturbed by a malfunction that could jeopardize the *essential performance* of the *EMS ventilator*. Without *system recovery*, the *patient* would have to be disconnected and connected to an alternative means of *ventilation*. *Ventilation* would thus be interrupted until an alternative means of *ventilation* is connected. *System recovery*, as specified in this subclause, attempts

to automatically re-establish *ventilation* in a shorter period. If necessary, the *EMS ventilator* can be replaced later when it is less consequential for the *patient's* therapy.

b) *System recovery* may result in the temporary:

- 1) cessation in the *ventilation* of the *patient*; or
- 2) reduction in the function of *EMS ventilator* subassemblies without impacting the *ventilation* of the *patient*.

EXAMPLE The temporary blanking of the display.

c) During a *system recovery* with a cessation of *ventilation*:

- 1) the *EMS ventilator* shall allow spontaneous *patient* breathing in accordance with 201.103; and
- 2) the *EMS ventilator* shall be equipped with an *alarm system* to indicate *system recovery* with a cessation of *ventilation*.

i) The *system recovery* with a cessation of *ventilation alarm condition* shall be *high priority*.

d) During a *system recovery* without a cessation of *ventilation*, the *EMS ventilator* shall be equipped with an *alarm system* to indicate *system recovery* without a cessation of *ventilation*.

- 1) The *system recovery* without a cessation of *ventilation alarm condition* shall be at least *low priority*.
- 2) A *low priority system recovery* without a cessation of *ventilation alarm condition* need not have an auditory *alarm signal*.

e) Following a *system recovery* without *operator* intervention, the *EMS ventilator* shall attempt to operate with the same system configuration settings, *ventilation* settings and *alarm settings* as before the *system recovery*.

- 1) If the system configuration settings, *ventilation* settings or *alarm settings* are different after the *system recovery*, the *EMS ventilator* shall be equipped with an *alarm system* to indicate any change in settings.
- 2) The change in settings *alarm condition* shall be at least *medium priority*.

f) The duration of a *system recovery* with a cessation of *ventilation* should be as short as practicable to avoid an unacceptable *risk* to the *patient*.

g) The maximum duration of a *system recovery* of *ventilation* shall be disclosed in the *instructions for use*.

Check conformity by inspection of the *instructions for use* and functional testing.

#### 201.4.4 Additional requirements for *expected service life*

*Amendment (add as a second paragraph):*

The *manufacturer* shall:

- a) state the probability of component failure that results in the *EMS ventilator* needing to be taken out of service during the *expected service life*, assuming that the preventative inspection, preventive maintenance and calibration are performed according to the *accompanying documents*; and

b) summarize the methodology used to determine this probability.

*Replacement (replace the compliance check):*

Conformity is checked by inspection of the design documentation and the *risk management file*.

#### **201.4.5 Alternative risk control measures or test methods for ME equipment or ME system**

*Amendment (add prior to the compliance check):*

Subsequent revisions of dated references (new editions or amendments) may be used in substitution of a referenced document provided the *manufacturer* can demonstrate the *hazard* or *hazardous situation* addressed in the dated normative reference is adequately resolved in the subsequent revision.

#### **201.4.6 ME equipment or ME system parts that contact the patient**

NOTE There is guidance or rationale for this subclause contained in Clause AA.2.

*Amendment (add at end of subclause):*

aa) The *VBS* or its parts or *accessories* that can come into contact with the *patient* shall be subject to the requirements for *applied parts* in accordance with this subclause.

*Addition:*

#### **201.4.11.101 Additional requirements for pressurized gas input**

NOTE There is guidance or rationale for this subclause contained in Clause AA.2.

##### **201.4.11.101.1 Overpressure requirement**

a) An *EMS ventilator* shall:

- 1) operate and meet the requirements of this document throughout its *rated* range of input pressure; and
- 2) not cause an unacceptable *risk* under the *single fault condition* of an input pressure of 1 000 kPa.

b) An *EMS ventilator* with a maximum *rated* input pressure in excess of 600 kPa shall not cause an unacceptable *risk* under the *single fault condition* of twice the maximum *rated* input pressure.

NOTE 1 Internal pressure regulators can be required to accommodate the *single fault condition* of maximum input pressure as well as the *rated* range of input pressure.

NOTE 2 Under the *single fault condition* of overpressure, it is desirable for gas to continue to flow to the *VBS*. Under this condition, the flowrate from the *EMS ventilator* is likely to be outside of its specification.

Check conformity by functional testing in *normal use* and under *normal condition* with the most adverse operating settings, by functional testing in *single fault condition* and inspection of the *risk management file* for *risks* associated with overpressure conditions.

##### **201.4.11.101.2 Compatibility requirement for medical gas pipeline systems**

If the *EMS ventilator* is intended to be connected to a *medical gas pipeline system* conforming with ISO 7396-1:2016+AMD1:2017, then:

a) the *rated* range of input pressure shall cover the range specified in ISO 7396-1:2016+AMD1:2017; and

NOTE Taking account of requirements for over-pressure and under-pressure, this corresponds to a range 280 kPa to 600 kPa.

- b) under *normal condition*,
- 1) the maximum input flow required by the *EMS ventilator* for each gas shall not exceed 60 l/min averaged over 10 s, measured at the *gas intake port* with an input pressure of 280 kPa; and
  - 2) any transient input flow shall not exceed 200 l/min averaged for 3 s.
- or:
- 3) the *accompanying documents* shall disclose:
    - i) the maximum input flow required by the *EMS ventilator* for each gas at a pressure of 280 kPa averaged over 10 s, measured at the *gas intake port*;
    - ii) the maximum transient input flow averaged for 3 s required by the *EMS ventilator* for each gas at a pressure of 280 kPa, measured at the *gas intake port*; and
    - iii) a warning to the effect that this *EMS ventilator* is a high-flow device and should only be connected to a pipeline installation designed using a diversity factor that allows for the indicated high flow at a specified number of terminal *outlets*, in order to avoid exceeding the pipeline design flow, thereby minimising the *risk* that the *EMS ventilator* interferes with the operation of adjacent equipment.

Check conformity by functional testing in *normal use* and under *normal condition* with the most adverse operating settings and by inspection of the *accompanying documents*.

EXAMPLE Highest driving gas consumption under worst-case settings for *set rate* and *tidal volume* within the *rated* range for *inlet* pressure.

### 201.4.11.101.3 Compatibility requirements for pressure regulators

If the *EMS ventilator* is intended to be connected to an oxygen gas cylinder via a pressure regulator then:

- a) the *rated* range of input pressure shall cover the range specified in:
  - 1) ISO 10524-1:2018; or
  - 2) ISO 10524-3:2019; and
- b) Under *normal condition*, the 10 s average input flow required by the *EMS ventilator* for oxygen shall not exceed 40 l/min at a pressure of 360 kPa, measured at the *gas intake port*.

Check conformity by functional testing in *normal use* and under *normal condition* with the most adverse operating settings and by inspection of the *accompanying documents*.

EXAMPLE Highest driving gas consumption, highest *fresh gas* delivery and, if provided, the highest *rated* gas consumption at any gas power supply output.

## 201.5 General requirements for testing of *ME equipment*

IEC 60601-1:2005+AMD1:2012+AMD2:2020, Clause 5 applies, except as follows:

*Addition:*

### 201.5.101 Additional requirements for general requirements for testing of *ME equipment*

NOTE There is guidance or rationale for this subclause contained in Clause AA.2.

#### 201.5.101.1 *EMS ventilator test conditions*

a) For testing, the *EMS ventilator*:

- 1) shall be connected to gas supplies as specified for *normal use*;
- 2) except that industrial grade oxygen and air may be substituted for the equivalent medical gas, as appropriate, unless otherwise stated.

NOTE There is guidance or rationale for this list item contained in Clause AA.2.

b) When using substitute gases, care should be taken to ensure that the test gases are oil-free and appropriately dry.

#### 201.5.101.2 Gas flowrate and leakage specifications

NOTE There is guidance or rationale for this subclause contained in Clause AA.2.

All requirements for gas flowrate, volume and leakage in this document

- a) are expressed at *STPD*,
- b) except for those associated with the *VBS*, which are expressed at *BTPS*.

*Correct all test measurements to STPD or BTPS, as appropriate.*

#### 201.5.101.3 *EMS ventilator testing errors*

NOTE 1 There is guidance or rationale for this subclause contained in Clause AA.2.

a) For the purposes of this document, acceptance criteria for declared tolerances of testing shall use:

- 1) *procedure 1* (calculation of uncertainty of measurement) from IEC Guide 115:2021, 4.4.2; or
- 2) *procedure 2* (accuracy method) from IEC Guide 115:2021, 4.4.3.

NOTE 2 This is a change from previous revisions of this document, which required tolerances to be adjusted by subtracting measurement uncertainty from disclosed tolerance values to determine acceptance criteria.

b) Test equipment and methods shall be selected and controlled to ensure that the uncertainty (with coverage factor  $k = 2$ , for confidence of  $\sim 95\%$ ) is no more than 30% of the disclosed tolerance for the parameter being tested.

EXAMPLE If the *manufacturer* wishes to claim a tolerance for delivered volume of  $\pm(10 \text{ ml} + 10\% \text{ of set volume})$  then the uncertainty of the measurement cannot exceed  $\pm(3 \text{ ml} + 3\% \text{ of set volume})$ .

c) The *manufacturer* shall disclose the measurement uncertainty for each disclosed tolerance in the *technical description*.

Check conformity by inspection of the *technical description*.

### 201.6 Classification of *ME equipment* and *ME systems*

IEC 60601-1:2005+AMD1:2012+AMD2:2020, Clause 6 applies.

## 201.7 ME equipment identification, marking and documents

IEC 60601-1:2005+AMD1:2012+AMD2:2020, Clause 7 applies, except as follows:

*Addition:*

### 201.7.1.101 Information to be supplied by the manufacturer

- a) The *information supplied by the manufacturer* of an *EMS ventilator* and its *accessories* shall conform with ISO 20417:2021.
- b) In applying ISO 20417:2021, the terms in this document and those in IEC 60601-1:2005+AMD1:2012+AMD2:2020 shall be used as follows.
  - 1) The term "*accompanying information*" shall assume the same meaning as *accompanying documents*.
  - 2) The term "*medical device*" shall assume the same meaning as *ME equipment*.
  - 3) The term "*user*" shall assume the same meaning as *operator*.
  - 4) The term "*patient*" shall include animals.

Check conformity by application of ISO 20417:2021.

### 201.7.2.3 Consult accompanying documents

*Replacement:*

NOTE There is guidance or rationale for this subclause contained in Clause AA.2.

The *EMS ventilator* shall be *marked* with the *safety sign* for the mandatory action: "follow instructions for use", ISO 7010-M002 (see IEC 60601-1:2005+AMD1:2012+AMD2:2020, Table D.2, *safety sign* 10).

*Addition:*

### 201.7.2.4.101 Additional requirements for accessories

- a) *Accessories* supplied separately shall:
  - 1) fulfil the *marking* requirements of ISO 20417:2021, 6.1.1 c); and
  - 2) be *marked* with an indication of any limitations or adverse effects of the *accessory* on the *basic safety* or *essential performance* of the *EMS ventilator*, if applicable. See also 201.7.2.101.
- b) If *marking* the *accessory* is not practicable, this information may be placed in the *instructions for use*.

Check conformity by inspection, inspection of the *risk management file* for any limitations or adverse effects of the *accessory*, if applicable, and when necessary, inspection of the *instructions for use*.

### 201.7.2.18 External gas source

*Amendment (add before the first dash):*

- aa) the gas name or chemical *symbol* in accordance with ISO 5359:2014+AMD1:2017,
- bb) the *rated* range of gas pressure,
- cc) for oxygen gas inputs, the *rated* range of oxygen concentration,

dd) gas-specific colour coding in accordance with ISO 32:1977, if colour coding is used,

EXAMPLE Colour coding to match the colour of the flexible supply hose or a gas cylinder intended to be attached to the *connector* of the *inlet*.

NOTE In some countries, other colour coding is used.

*Addition:*

### 201.7.2.101 Additional requirements for *marking* on the outside of *ME equipment* or *ME equipment parts*

NOTE 1 There is guidance or rationale for this subclause contained in Clause AA.2.

a) The *EMS ventilator* shall be provided with a *clearly legible* permanently attached *start-up procedure* that is intended to be performed prior to use.

b) The *start-up procedure* may be provided electronically.

EXAMPLE A *start-up procedure* provided on the graphical *operator interface* or provided orally.

c) If applicable, *operator-accessible ME equipment*, parts or *accessories* shall have *clearly legible markings* of the following

1) for an *EMS ventilator* intended to be used in the magnetic resonance (MR) environment, in accordance with IEC 62570:2014,

- i) *symbol 7.3.1-1* (Table 201.D.2.101, *symbol 3*) or *symbol 7.3.1-2* (Table 201.D.2.101, *symbol 4*) of IEC 62570:2014 for an 'MR Safe' *EMS ventilator*, or
- ii) *symbol 7.3.2* (Table 201.D.2.101, *symbol 5*) of IEC 62570:2014 for an 'MR Conditional' *EMS ventilator*,

2) for an *EMS ventilator* not intended to be used in the magnetic resonance (MR) environment, in accordance with IEC 62570:2014, *symbol 7.3.3* (Table 201.D.2.101, *symbol 6*) for an 'MR Unsafe' *EMS ventilator*

3) for *flow-direction-sensitive components* that are *operator-removable* without the use of a *tool*, an arrow indicating the direction:

- i) of the flow; or
- ii) of inspiratory flow if exposed to inspiratory and expiratory flow.

d) If applicable, *operator-accessible ME equipment*, parts or *accessories* shall have *clearly legible markings* of the following on or adjacent to the following items:

1) a warning not to obstruct the *gas intake port*.

EXAMPLE 1 WARNING: Gas Intake — Do not obstruct

i) A *symbol* or *safety sign* evaluated in accordance with IEC 62366-1 as *information for safety* may be used.

2) an arrow indicating the intended direction of gas flow:

- i) for the *gas output port*; and
  - I) *Symbol 0795* of ISO 7000 (Table 201.D.2.101, *symbol 2*) may be used for *gas output port*.
- ii) for the *gas return port*.
  - I) *Symbol 0794* of ISO 7000 (Table 201.D.2.101, *symbol 1*) may be used for *gas return port*.

- 3) for a sampling gas *inlet*, either with the text “Gas sample” or the *Symbol* ISO 7000-0794, (Table 201.D.2.101, *symbol* 1).

EXAMPLE 2 Input to a diverting respiratory gas monitor integrated into the *EMS ventilator*.

#### 201.7.4.2 Control devices

*Amendment (add between the first and second paragraph):*

NOTE There is guidance or rationale for this subclause contained in Clause AA.2.

- aa) An *EMS ventilator* with an oxygen concentration control with an error exceeding 20 % [see 201.12.1.101 b) 3) and 201.12.1.102 b) 3) ] shall not be *marked* with quantitative numerical concentration values.

EXAMPLE “No air mix” or “High” for the setting for the maximum achievable oxygen concentration or “Air mix” or “Low” for the setting with a reduced oxygen concentration.

Check conformity by inspection.

#### 201.7.4.3 Units of measurement

NOTE There is guidance or rationale for this subclause contained in Clause AA.2.

*Amendment (add to the bottom as a new row in Table 1):*

All gas volume, flowrate and leakage specifications:

- aa) shall be expressed at *STPD*; except for those associated with the *VBS* which  
bb) shall be expressed at *BTPS*.

*Addition:*

##### 201.7.4.101 Labelling of units of measurement

Quantitative numeric indications of parameters and settings shall be labelled with the units of measurement, either:

- a) continuously; or  
b) on *operator* action.

#### 201.7.9.2.1 General

*Amendment (add as fifth bullet following the first paragraph):*

- the intended position of the *operator*;

##### 201.7.9.2.1.101 Additional general requirements

The *instructions for use* shall disclose the following:

- a) the intended range of *inspiratory volume*.

Check conformity by inspection.

##### 201.7.9.2.2.101 Additional requirements for warnings and safety notices

NOTE There is guidance or rationale for this subclause contained in Clause AA.2.

a) The *instructions for use* shall include:

- 1) a warning statement to the effect of “Warning: This ventilator is intended to be continuously attended by an operator. Failure to be in close proximity to this ventilator can contribute to patient death or serious injury.”
- 2) a warning statement to the effect of “Warning: Do not cover the ventilator or place in a position that affects proper operation”, including applicable examples.

EXAMPLE 1      **WARNING:** Do not position next to a blanket that blocks the flow of cooling air, thereby causing the ventilator to overheat, thereby interfering with patient ventilation.

EXAMPLE 2      **WARNING:** Do not block the gas intake port thereby interfering with patient ventilation.

- 3) a warning statement to the effect of “Warning: Always have immediate access to an alternative means of ventilation, which is ready for use, in order to reduce the possibility of patient death or additional deterioration of health.”

EXAMPLE 3      **WARNING:** Failure to have an alternative means of ventilation such as a self-inflating, manually powered resuscitator (as specified in ISO 10651-4) with mask can result in patient death if the ventilator fails.

- 4) a warning statement to the effect of “WARNING: Do not add any attachments or accessories to the ventilator that are not listed as intended for use in combination with the ventilator in the instructions for use of the ventilator or accessory, as the ventilator might not function correctly leading to the risk of patient death or additional serious deterioration of health.”

b) If applicable, the *instructions for use* shall include the following:

- 1) a warning statement to the effect of “WARNING: When using nebulisation or humidification, the breathing system filters can require more frequent replacement to prevent increased resistance and blockage.
- 2) a warning statement to the effect of “Warning: The ventilator shall not be used in a hyperbaric chamber. Such use might cause the ventilator to not function correctly causing patient death or additional serious deterioration of health.”
- 3) a warning statement to the effect of “Warning: The ventilator accuracy can be affected by the gas added to the ventilator breathing system by use of a pneumatic nebuliser.
- 4) a warning statement to the effect of “Warning: Do not use the ventilator in explosive environments. Such use might cause an explosion.”

Check conformity by inspection.

#### **201.7.9.2.8.101 Additional requirements for start-up procedure**

NOTE 1 There is guidance or rationale for this subclause contained in Clause AA.2.

NOTE 2 A start-up *procedure* includes a pre-use functional test that is used to determine whether the *EMS ventilator* is ready for use.

- a) The *instructions for use* or the *operator interface* shall include a pre-use checklist.
- b) The *instructions for use* shall disclose a method by which all of the *alarm signals* can be functionally tested by the clinical *operator* to determine if they are operating correctly.

- c) Portions of this test method may:
- 1) be automatically performed by the *EMS ventilator*; or
  - 2) require *operator* action.

EXAMPLE Combination of the power-on self-test routines and *operator* action.

Check conformity by inspection of the *instructions for use*.

### 201.7.9.2.9.101 Additional requirements for operating instructions

NOTE 1 There is guidance or rationale for this subclause contained in Clause AA.2.

The *instructions for use* shall disclose:

- a) a listing of the following pressures:
- 1) the range of *maximum limited pressure* ( $P_{LIM\ max}$ );
  - 2) if provided, the *rated* range to which the *maximum working pressure* ( $P_{W\ max}$ ) can be set, if adjustable;
  - 3) the means by which the *maximum working pressure* is accomplished.

EXAMPLE 1 Pressure cycling, pressure limiting, pressure generation.

- b) the *rated* range of the following characteristics of the assembled *operator*-detachable parts of the *VBS*, over which the accuracies of set and monitored volumes and pressures are maintained:
- 1) inspiratory *gas pathway* resistance,
  - 2) expiratory *gas pathway* resistance, and
  - 3) *VBS* compliance.
    - i) These specifications may be presented in ranges.
    - ii) The accuracies of set and monitored volumes may be presented as a function of these characteristics.

NOTE 2 Compliance and resistance can be non-linear. These characteristics might need to be specified over a range (e.g. at 15 l/min, 30 l/min, 60 l/min and the maximum flowrate or the maximum pressure).

- c) the conditions under which the *EMS ventilator* maintains the accuracy of controlled and displayed variables as disclosed in the *instructions for use*.

EXAMPLE 2 Interval of calibration of a flow sensor.

- d) if end-tidal gas concentration measurement is provided, a summary of effects of respiratory rate on the end-tidal gas concentration measurement accuracy.
- e) if a connection for a *distributed alarm system* is provided, a description of how to connect and test the connection of a *distributed alarm system*.
- f) a cross reference between the *manufacturer*-specific naming of the *EMS ventilator's ventilation-modes* and the *ventilation-mode* systematic coding scheme in Annex E of ISO 19223:2019.

g) whether or not the *EMS ventilator* is intended for use with closed suctioning.

Check conformity by inspection.

#### **201.7.9.2.12 Cleaning, disinfection, and sterilization**

*Amendment: (add after normal use)*

or in *single fault condition*

*Amendment: (add after bulleted list)*

aa) The *instructions for use* shall identify which portions of the *gas pathways* through the *EMS ventilator* can become contaminated with body fluids or by microbial material conveyed by the expired breathing gases during both *normal condition* and *single fault condition*.

*Addition:*

#### **201.7.9.2.14.101 Additional requirements for accessories, supplementary equipment, used material**

NOTE 1 There is guidance or rationale for this subclause contained in Clause AA.2.

- a) The *instructions for use* of the *EMS ventilator* shall identify:
- 1) at least one set of *accessories*; and
  - 2) if applicable, the *ME equipment* necessary for the *EMS ventilator's intended use*.
- b) If applicable, the *instructions for use* of the *EMS ventilator* shall disclose:
- 1) any restrictions on the positioning of components within the *ventilator breathing system*; and  
 EXAMPLE Where such components are *flow-direction-sensitive components*.
  - 2) any reasonably foreseeable adverse effect of any recommended *accessory* on the *essential performance* or *basic safety* of the *EMS ventilator* (additional requirements are found in 201.16).

Check conformity by inspection of the *instructions for use* and inspection of the *risk management file* for any adverse effect of any recommended *accessory*.

#### **201.7.9.3.1.101 Additional general requirements**

NOTE There is guidance or rationale for this subclause contained in Clause AA.2.

- a) The *technical description* shall disclose:
- 1) a summary description of the filtering or smoothing techniques for all measured or computed variables that are displayed or used for control necessary for the *operator* to form a mental model of the operation of the *EMS ventilator*;
  - 2) a pneumatic diagram of the *EMS ventilator*, including a diagram for *operator-detachable parts* of the *ventilator breathing system* either supplied or recommended in the *instructions for use*;
  - 3) a summary description of the means of initiating and terminating the *inflation phase* in each *ventilation-mode* of the *EMS ventilator*; and
  - 4) the outside dimensions and mass of the *EMS ventilator*.

- b) If applicable, the *technical description* shall disclose:
- 1) the essential technical characteristics of each recommended *breathing system filter*; and  
EXAMPLE    Deadspace and resistance.
  - 2) a description of the method used to calculate end-tidal gas readings.

Check conformity by inspection.

### **201.7.9.3.101    Additional requirements for the *technical description***

The *technical description* shall disclose:

- a) a description of a *procedure* for checking the functioning of the *alarm system* for each of the *alarm conditions* specified in this document, if not performed automatically during start-up; and
- b) which checks are performed automatically.

Check conformity by inspection of the *technical description*.

## **201.8    Protection against electrical *hazards* from *ME equipment***

IEC 60601-1:2005+AMD1:2012+AMD2:2020, Clause 8 standard applies.

## **201.9    Protection against *mechanical hazards* of *ME equipment* and *ME systems***

IEC 60601-1:2005+AMD1:2012+AMD2:2020, Clause 9 applies, except as follows:

*Addition:*

### **201.9.4.3.101    Additional requirements for instability from unwanted lateral movement**

NOTE 1    There is guidance or rationale for this subclause contained in Clause AA.2.

- a) An *EMS ventilator* shall include a means by which the *EMS ventilator* can be attached to prevent unwanted movement during transport while in use.
- b) The means shall secure the *EMS ventilator* to withstand accelerations or decelerations of
  - 1) 10 g longitudinal (forward, backward), 10 g transverse (left, right) and 10 g vertical (up, down) for at least 150 ms in each orientation for a road ambulance.
  - 2) 20 g longitudinal (forward, backward), 20 g transverse (left, right), 20 g vertical (down) and 4 g vertical (up) for at least 3 s in each orientation for an air ambulance.
- c) Terminal units and electrical sockets shall not be used as part of the fixation system.
- d) If rail systems are used, they shall conform with ISO 19054:2005+AMD1:2016.

NOTE 2    Rail systems consist of e.g. rail supports, rails, rail clamps, equipment mount holders, equipment mounts, equipment pin holders and equipment pins.

- e) The means to secure a *portable EMS ventilator* to withstand forces caused by accelerations or decelerations shall be operable without the use of a *tool*.

f) The *technical description* shall disclose the technical details of this means for attachment.

Check conformity by inspection and functional testing.

EXAMPLE These tests can be performed by using static loading, acceleration sleds or centrifuges.

#### 201.9.4.4 Grips and other handling devices

Replace list item b) with:

b) A portable EMS ventilator shall be designed such that it can be carried by no more than one hand.

1) It shall be possible for the operator to:

- i) observe the display;
- ii) recognize the *alarm signals*; and
- iii) actuate the controls while holding the EMS ventilator with no more than one hand.

2) A carrying case or handle may be used to fulfill this requirement.

Amendment (add to the end of the compliance check):

Check conformity by functional testing while carrying with no more than one hand.

#### 201.9.6.2.1.101 Additional requirements for audible acoustic energy

a) The A-weighted sound pressure level emitted by the EMS ventilator shall be:

- 1) measured in accordance with ISO 4871:1996 and ISO 3744:2010 using engineering method grade 2; and
- 2) disclosed in the *instructions for use*.

b) The A-weighted sound power level shall be:

- 1) calculated in accordance with 8.2.5 and 8.6 of ISO 3744:2010; and
- 2) disclosed in the *instructions for use*.

c) Check conformity with the following test:

- 1) Place the EMS ventilator on the sound-reflecting plane and attach the least favourable VBS from those indicated in the *instructions for use*.

NOTE 1 The least favourable VBS configuration can vary by *ventilation-mode*, *inflation-type* and flow pattern, as applicable.

- 2) If a *humidifier* is provided with the EMS ventilator, include the *humidifier* filled to the least favourable level in the test.
- 3) Configure the test lung with the compliance and resistance components whose values are indicated in Table 201.102.
  - i) Acoustically isolate the test lung by a suitable means so that any noise caused by the test lung does not interfere with the sound measurement of the EMS ventilator.
  - ii) Connect the *patient-connection port* to the test lung.

- 4) Set the EMS ventilator to the least favourable ventilation-mode, inflation-type and flow pattern, as applicable, that generates ventilation as indicated in Table 201.102.

NOTE 2 The least favourable ventilation-mode, inflation-type and flow pattern can vary by VBS configuration.

- 5) Using a microphone of the sound level meter, conforming with the requirements of type 1 instruments specified in IEC 61672-1:2013, measure the maximum time-weighted sound pressure level using frequency-weighting A and the time-weighting F of the sound level meter (i.e.  $L_{AFmax}$ ) at 10 positions in a hemisphere with a radius from the geometric centre of the EMS ventilator in a free field over a reflecting plane as specified in 8.1.1 of ISO 3744:2010. Average the values in accordance with 8.2.2 of ISO 3744:2010.
- 6) Calculate the A-weighted sound pressure level averaged over the measurement surface in accordance with 8.2.2 of ISO 3744:2010.
- 7) Calculate the A-weighted sound power level in accordance with 8.6 of ISO 3744:2010.
- 8) Confirm that the criteria for background noise specified in 4.2 of ISO 3744:2010 are fulfilled.
- 9) Ensure that the average measured sound pressure level is less than or equal to that disclosed in the instructions for use.
- 10) Ensure that the sound power level is less than that disclosed in the instructions for use.

Table 201.102 — Test conditions for acoustic tests

Adjustable parameter	Test condition		
	For an EMS ventilator intended to provide tidal volume		
	$V_T \geq 300$ ml	$300 \text{ ml} \geq V_T \geq 50$ ml	$V_T \leq 50$ ml
Tidal volume, $V_T^a$	500 ml	150 ml	30 ml
Set rate	10 breaths/min	20 breaths/min	30 breaths/min
I:E ratio	1:2	1:2	1:2
BAP	10 hPa	10 hPa	10 hPa
Linear resistance, $R^{b[43][53][56]}$	5 hPa·(l/s) <sup>-1</sup> ± 10 %	20 hPa·(l/s) <sup>-1</sup> ± 10 %	50 hPa·(l/s) <sup>-1</sup> ± 10 %
Isothermal Compliance, $C^b$	50 ml/hPa ± 5 %	20 ml/hPa ± 5 %	1 ml/hPa ± 5 %

<sup>a</sup>  $V_T$  is measured by means of a pressure sensor at the test lung, where  $V_T = C (P_{max} - P_{min})$ , and  
 $V_T$  is the volume delivered to the test lung  
 $C$  is the Isothermal Compliance of the test lung  
 $P_{max}$  is the maximum pressure measured in the test lung  
 $P_{min}$  is the minimum pressure measured in the test lung.

<sup>b</sup> The accuracy for  $C$  and  $R$  applies over the ranges of the measured parameters.

NOTE There is guidance or rationale for this table contained in Clause AA.2.

### 201.10 Protection against unwanted and excessive radiation hazards

IEC 60601-1:2005+AMD1:2012+AMD2:2020, Clause 10 applies.

### 201.11 Protection against excessive temperatures and other hazards

IEC 60601-1:2005+AMD1:2012+AMD2:2020, Clause 11 applies, except as follows:

**201.11.1.2.2 Applied parts not intended to supply heat to a patient**

NOTE There is guidance or rationale for this subclause contained in Clause AA.2.

*Amendment (add between the existing paragraphs):*

In *normal use* and *single fault conditions* and over the *rated* flowrate range and at the maximum *rated* operating temperature, the temperature of the gas delivered by the *EMS ventilator* at the *patient-connection port*, both with and without each *humidifier* specified for use in the *instructions for use*, when averaged over 120 s shall not exceed:

aa) 70 °C; and

bb) an energy equivalent to 43 °C and 100 % relative humidity (a specific enthalpy not to exceed 197 kJ/m<sup>3</sup> dry air).

Table 201.103 contains examples of combinations of temperature and relative humidity with such a specific enthalpy.

**Table 201.103 — Examples of permissible combinations of temperature and relative humidity**

Temperature °C	Relative humidity %
43	100
44	95
45	90
48	76
50	69
55	52
60	40
65	30
70	23

**201.11.6.6 Cleaning and disinfection of ME equipment or ME system**

NOTE 1 There is guidance or rationale for this subclause contained in Clause AA.2.

*Amendment (add additional requirement as new first paragraph):*

aa) *Gas pathways* through the *EMS ventilator* and its *accessories* not intended for *single use* that can become contaminated with body fluids or by microbial material conveyed by the expired gases during *normal condition* or *single fault condition* shall be designed to allow for:

- 1) *cleaning and disinfection*; or
- 2) *cleaning and sterilization*.

NOTE 2 Additional *sterilization* requirements are found in 11.6.7 of IEC 60601-1:2005+AMD1:2012+AMD2:2020.

bb) *Processing* instructions for the *gas pathways* of the *EMS ventilator* and its *accessories* that can become contaminated shall:

- 1) conform with ISO 17664-1:2021;

- 2) conform with ISO 14937:2009, if applicable; and
- 3) be disclosed in the *instructions for use*.

*Amendment (add additional requirement as the last requirement and replace the compliance test):*

- cc) *EMS ventilator enclosures* shall be designed to allow for surface *cleaning* and *disinfection* to reduce to acceptable levels the *risk* of cross infection of the next *patient*.

NOTE 3 ISO 14159 provides guidance for the design of *enclosures*.

- dd) *Processing* instructions for the *EMS ventilator enclosure* shall:

- 1) conform with ISO 17664-2:2021; and
- 2) be disclosed in the *instructions for use*.

Check conformity by inspection of the *risk management file* to determine the *gas pathways* that can become contaminated in *normal condition* and *single fault condition* and conformity with ISO 17664-1:2021 and ISO 17664-2:2021, as appropriate. When conformity with this document could be affected by the *cleaning* or the *disinfection* of the *EMS ventilator* or its parts or *accessories*, clean and disinfect them for the number of *processing* cycles and in accordance with the methods indicated in the *instructions for use*, including any cooling or drying period. After these *procedures*, ensure that *basic safety* and *essential performance* are maintained. Confirm that the *manufacturer* has evaluated the effects of multiple *processing* cycles and the effectiveness of those cycles.

#### **201.11.6.7 Sterilization of ME equipment or ME system**

*Amendment (add note before compliance test):*

NOTE Additional *sterilization* requirements are found in 11.6.6 of IEC 60601-1:2005+AMD1:2012+AMD2:2020.

#### **201.11.7 Biocompatibility of ME equipment and ME systems**

*Amendment (add after existing text prior to the compliance statement):*

- aa) The *manufacturer* of an *EMS ventilator*, *VBS*, its parts and *accessories* shall address in the *risk management process* the *risks* associated with the leaching or leaking of substances into the *gas pathway*.
- bb) The *gas pathways* shall be evaluated for *biocompatibility* in accordance with to ISO 18562-1:2023.

### **201.12 Accuracy of controls and instruments and protection against hazardous outputs**

IEC 60601-1:2005+AMD1:2012+AMD2:2020, Clause 12 applies, except as follows:

*Addition:*

#### **201.12.1.101 Volume-control inflation-type**

NOTE 1 There is guidance or rationale for this subclause contained in Clause AA.2.

- a) If provided with a *volume-control inflation-type*, then with a *volume-control inflation-type* selected and the *EMS ventilator* operating in *normal condition*, the accuracy as determined for the test settings and conditions specified in this document shall be disclosed in the *instructions for use*, as the maximum bias error and maximum linearity error.

EXAMPLE  $\pm (5 + 10 \% \text{ of the set volume}) \text{ ml}$

- b) This disclosure shall include at least:
- 1) the maximum error of the *inspiratory volume* in relation to the *set tidal volume*;
  - 2) the maximum error of the *PEEP* in relation to the set value of *BAP*; and
  - 3) the maximum error of the inspiratory oxygen ( $\text{FiO}_2$ ) concentration at the *patient-connection port* in relation to the set value.

i) The disclosed accuracies shall include the effects of the *rated* input oxygen concentration.

- c) All of the errors may be reported separately for the following ranges of intended *tidal volume*:

- 1)  $V_T \geq 300 \text{ ml}$
- 2)  $300 \text{ ml} \geq V_T \geq 50 \text{ ml}$
- 3)  $V_T \leq 50 \text{ ml}$

- d) The accuracy of the performance of the *EMS ventilator* shall either be:

- 1) determined for each *VBS* configuration indicated in the *instructions for use*; or
- 2) determined for the worst-case *VBS* configurations indicated in the *instructions for use*.

NOTE 2 The worst-case *VBS* configuration can be different for each error or *nominal inspiratory volume*.

- e) If worst-case *VBS* configurations are used, the rationale for their selection shall be documented in the design documentation.
- f) Check conformity by inspection of the design documentation for the rationale, if applicable, and with the following tests:


NOTE 3 In some cases, the following tests can be carried out simultaneously.

- 1) *Tidal volume* and end-expiratory pressure errors
  - i) Set up the *EMS ventilator* as shown in Figure 201.101.
  - ii) If applicable, determine or input the *VBS* compliance required for compliance correction as indicated in the *instructions for use* and activate this correction. If a *humidifier* is used, fill the *humidifier* to the maximum water level prior to determining the *VBS* compliance.
  - iii) Use the test parameters and settings of the first applicable row (selected by intended *tidal volume*) of Table 201.104.
  - iv) Wait for steady-state conditions to be achieved.
  - v) Determine the *tidal volume*, for example via integration of the flow signal provided by a calibrated flow sensor located at the *patient-connection port* or by the product of the test lung compliance and the measured change of test lung pressure, compensated for temperature effects due to fast compression of the gas, if necessary.

NOTE 4 Additional information on the construction of an isothermal test lung is found in Reference [48].

- vi) Compare the result with the volume setting for the test and the resulting difference with the tolerance indicated in the *instructions for use*.

- vii) If the *EMS ventilator* is equipped with *tidal volume monitoring equipment*, determine the accuracy of the *tidal volume monitoring equipment* by comparing its reading to the *tidal volume* determined in v). Refer to 201.12.1.104.
- viii) Determine the *PEEP* as the average of the *airway pressure* measurements over the last 50 ms of the *expiratory phase*.
- ix) Compare the result with the *BAP* setting for the test and the resulting difference with the tolerance indicated in the *instructions for use*.
- x) Repeat iii) to ix) for 30 consecutive *inflations*.
- xi) Repeat iv) to x) for each applicable row (selected by intended *tidal volume*) of Table 201.104.
- xii) If a *humidifier* is included in the *VBS*, repeat the *tidal volume* tests with the minimum *humidifier* water level without re-determining the *VBS* compliance.
- xiii) Unless it can be demonstrated that the worst-case flow pattern (e.g. constant flow, decelerating flow) has been selected for the tests, repeat ii) to xii) for each flow pattern available on the *EMS ventilator*.
- xiv) If the *EMS ventilator* permits operation without compliance correction, repeat ii) to xiii) without compliance correction.

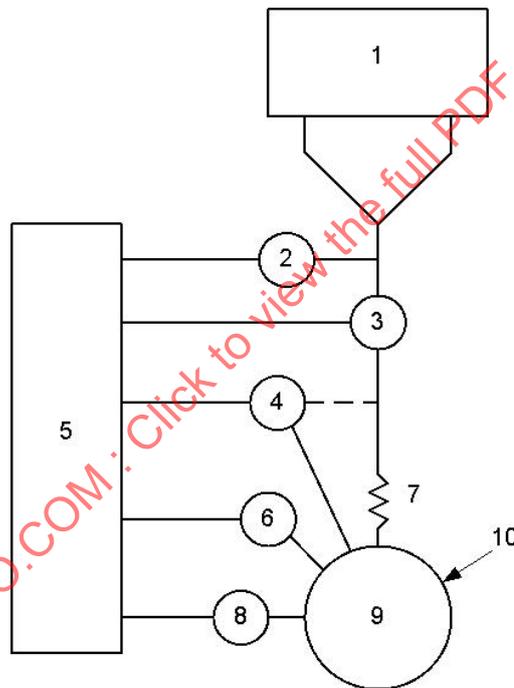

**Key**

- 1 *EMS ventilator* under test
- 2 pressure sensor
- 3 flow sensor, with a 10 % to 90 % rise time of no greater than 10 ms
- 4 oxygen sensor
- 5 data acquisition system, with minimum sample rate of 200 samples/s
- 6 temperature sensor
- 7 test lung resistance ( $R_{lung}$ )
- 8 pressure sensor, with a 10 % to 90 % rise time of no greater than 10 ms
- 9 test lung compliance ( $C_{lung}$ )
- 10 test lung

The oxygen sensor may be placed in the *VBS*.

**Figure 201.101 — Volume-control and pressure-control inflation-type accuracy, typical test setup**

2) Oxygen (FiO<sub>2</sub>) error

- i) The accuracy of the inspiratory oxygen concentration of the gas delivered is assessed by placing the sensor of an oxygen concentration-measuring device at the *patient-connection port* or inside the test lung. If the sensor is located at the *patient-connection port*, the value of the concentration is the flow-weighted average concentration during the *inflation phase*.
- ii) Evaluate the measured oxygen concentration with the acceptance criteria derived from the oxygen setting for the test and the disclosed tolerance.

**Table 201.104 — Volume-control inflation-type testing**

Test number	Test lung parameters		EMS ventilator settings				
	Compliance ml/hPa ±10 %	Linear resistance [43][53][56] hPa·(l/s) <sup>-1</sup> ±10 %	Tidal volume ml	Set rate <sup>a</sup> breaths/min	Inspiratory time s	FiO <sub>2</sub> <sup>b</sup> %	BAP hPa (cmH <sub>2</sub> O)
1	— <sup>c</sup>	— <sup>c</sup>	Maximum settable tidal volume	— <sup>c</sup>	— <sup>c</sup>	— <sup>c</sup>	— <sup>c</sup>
2	50	5	500	20	1	30	5
3	50	20	500	12	1	85	10
4	20	5	500	20	1	85	5
5	20	20	500	20	1	30	10
6	20	20	300	20	1	30	5
7	20	50	300	12	1	85	10
8	10	50	300	20	1	30	10
9	10	10	200	20	1	85	5
10	3	10	50	30	0,6	30	5
11	3	20	50	30	0,6	30	10
12	3	50	50	20	0,6	60	5
13	3	20	30	30	0,6	30	5
14	3	50	30	20	0,6	85	10
15	1	20	30	30	0,6	85	5
16	1	100	30	30	0,6	30	10
17	1	200	20	50	0,4	30	5
18	1	200	15	50	0,4	60	10
19	1	50	10	60	0,4	60	5
20	0,5	50	5	60	0,4	60	10
21	0,5	200	5	30	0,4	30	5
22	0,5	200	5	60	0,4	30	10
23	— <sup>c</sup>	— <sup>c</sup>	Minimum settable tidal volume	— <sup>c</sup>	— <sup>c</sup>	— <sup>c</sup>	— <sup>c</sup>

<sup>a</sup> If the *end-expiratory flow* does not reach zero, reduce the *set rate* until it does.  
<sup>b</sup> In the case that *EMS ventilator* does not provide the settings as specified in column FiO<sub>2</sub>, select the closest available FiO<sub>2</sub> *EMS ventilator* setting.  
<sup>c</sup> The *instructions for use* shall indicate this value.  
 NOTE There is guidance or rationale for this table contained in Clause AA.2.

- iii) Compare the resulting difference with the tolerance indicated in the *instructions for use*. If the sensor is located inside the test lung, care should be taken to allow sufficient time for the gas mixture in the test lung to reach a stable concentration. If the sensor is located in the tubing system, care should be taken to ensure that the measuring device has a sufficiently fast response to permit determination of the flow-weighted average concentration during the *inflation phase*. If the oxygen concentration measuring device has pressure dependencies, compensate for these dependencies.
- iv) Compare each result with the acceptance criteria derived from the oxygen setting for the test in f) 1), above and the tolerance indicated in the *instructions for use*.

**201.12.1.102 Pressure-control inflation-type**

NOTE 1 There is guidance or rationale for this subclause contained in Clause AA.2.

- a) If provided with a *pressure-control inflation-type*, with the *pressure-control inflation-type* selected and the *EMS ventilator* operating in *normal condition*, the accuracy as determined for the test settings and conditions specified in this document shall be disclosed in the *instructions for use*, as the maximum bias error and maximum linearity error.

EXAMPLE ± (3,0 hPa +5 % of the set pressure); (± (3,0 cmH<sub>2</sub>O +5 % of the set pressure) )

- b) This disclosure shall include at least:
  - 1) the maximum error of the *airway pressure* ( $P_{aw}$ ) at the end of the *inflation phase* in relation to the set value;
  - 2) the maximum error of *PEEP* in relation to the set value of *BAP*; and
  - 3) the maximum error of the inspiratory oxygen ( $F_{iO_2}$ ) concentration at the *patient-connection port* in relation to the set value.

- i) The disclosed accuracies shall include the effects of the *rated* input oxygen concentration.

- c) All of the errors may be reported separately for the following ranges of intended *tidal volume*:
  - 1)  $V_T \geq 300$  ml
  - 2)  $300 \text{ ml} \geq V_T \geq 50$  ml
  - 3)  $V_T \leq 50$  ml

- d) The accuracy of the performance of the *EMS ventilator* shall either be:
  - 1) determined for each *VBS* configuration indicated in the *instructions for use*; or
  - 2) determined for the worst-case *VBS* configuration indicated in the *instructions for use*.

NOTE 2 The worst-case *VBS* configuration can be different for each error or each *nominal inspiratory volume* range.

- e) If worst-case *VBS* configurations are used, the rationale for their selection shall be documented in the *risk management file*.
- f) Check conformity by inspection of the *risk management file* for the rationale, if applicable, and with the following tests:

NOTE 3 In some cases, the following tests can be carried out simultaneously.

- 1) End-inspiratory and end-expiratory pressure errors
  - i) Set up the *EMS ventilator* as shown in Figure 201.101.

**Table 201.105 — Pressure-control inflation-type testing**

Test number	Intended tidal volume <sup>a</sup> ml	Test lung parameters		EMS ventilator settings				
		Compliance ml/hPa ±10 %	Linear Resistance [43][53][56] hPa·(l/s) <sup>-1</sup> ±10 %	Set rate <sup>b</sup> breaths/ min	Inspiratory time <sup>c</sup> s	Δ inspiratory pressure <sup>d</sup> hPa	FiO <sub>2</sub> <sup>e</sup> %	BAP hPa (cmH <sub>2</sub> O)
1	— <sup>f</sup>	— <sup>f</sup>	— <sup>f</sup>	— <sup>f</sup>	— <sup>f</sup>	Maximum settable pressure	— <sup>f</sup>	— <sup>f</sup>
2	500	50	5	20	1	10	30	5
3	500	50	20	12	1	15	85	10
4	500	20	5	20	1	25	85	5
5	500	20	20	20	1	25	30	10
6	300	20	20	20	1	15	30	5
7	300	20	50	12	1	25	85	10
8	300	10	50	20	1	30	85	5
9	200	10	10	20	1	25	30	10
10	50	3	10	30	0,6	15	30	5
11	50	3	20	30	0,6	15	30	10
12	50	3	50	20	0,6	25	60	5
13	30	3	20	30	0,6	10	30	5
14	30	3	50	20	0,6	15	85	10
15	30	1	20	30	0,6	30	85	5
16	30	1	100	30	0,6	30	30	10
17	20	1	200	50	0,4	20	30	5
18	15	1	200	50	0,4	15	60	10
19	10	1	50	60	0,4	10	60	5
20	5	0,5	50	60	0,4	15	60	10
21	5	0,5	50	30	0,4	10	30	5
22	5	0,5	200	60	0,4	15	30	10
23	— <sup>f</sup>	— <sup>f</sup>	— <sup>f</sup>	— <sup>f</sup>	— <sup>f</sup>	Minimum settable pressure	— <sup>f</sup>	— <sup>f</sup>

<sup>a</sup> The volume in this column is intended to be used for the selection of the test conditions and parameters based on the intended inspiratory volume of the *EMS ventilator*.

<sup>b</sup> If the *end-expiratory flow* does not reach zero, reduce the *set rate* until it does.

<sup>c</sup> The rise time of the *EMS ventilator* should be set to a value that ensures the set pressure can be reached within the *inspiratory time*.

<sup>d</sup> For the purposes of this test, the set pressure is relative to set *BAP*.

<sup>e</sup> If the *EMS ventilator* does not provide the settings as specified in column *FiO<sub>2</sub>*, select the closest available *FiO<sub>2</sub>* *EMS ventilator* setting.

<sup>f</sup> The *instructions for use* shall indicate this value.

NOTE There is guidance or rationale for this table contained in Clause AA.2.

- ii) If applicable, determine or input the *VBS* compliance required for compliance correction as indicated in the *instructions for use* and activate this correction. If a *humidifier* is used, fill the *humidifier* to the maximum water level prior to determining the *VBS* compliance.
- iii) Use the test parameters and settings of the first applicable row (selected by typical intended *tidal volume*) of Table 201.105.
- iv) Wait until steady-state conditions are achieved.
- v) Determine the *airway pressure* at the end of the *inflation phase* as the average over the preceding 50 ms.
- vi) Compare the result with the pressure setting for the test and the resulting difference with the tolerance indicated in the *instructions for use*.
- vii) Determine the *tidal volume*, for example via integration, of the flow signal provided by a calibrated flow sensor located at the *patient-connection port* or by the product of the test lung compliance and the measured change of test lung pressure, compensated for temperature effects due to fast compression of the gas, if necessary.

NOTE 4 Additional information on the construction of an isothermal test lung is found in Reference [48].

- viii) If the *EMS ventilator* is equipped with *tidal volume monitoring equipment*, determine the accuracy of the *tidal volume monitoring equipment* by comparing its reading to the *tidal volume* determined in vii). Refer to 201.12.1.105.
- ix) Determine the *PEEP* as the average of the *airway pressure* measurements over the last 50 ms of the *expiratory phase*.
- x) Compare the result with the acceptance criteria derived from the *BAP* setting for the test, the measurement uncertainty and the resulting difference with the tolerance indicated in the *instructions for use*.
- xi) Repeat ii) to x) for 30 consecutive *inflations*.
- xii) Repeat ii) to xi) for each applicable row (selected by intended *tidal volume*) of Table 201.105.
- xiii) If a *humidifier* is included in the *VBS*, repeat the *airway pressure* tests with the minimum *humidifier* water level without re-determining the *VBS* compliance.
- xiv) Unless it can be demonstrated that the worst-case flow pattern (e.g. constant flow, decelerating flow) has been selected for the tests, repeat ii) to xiii) for each flow pattern available on the *EMS ventilator*.
- xv) If the *EMS ventilator* permits operation without compliance correction, repeat ii) to xiv) without compliance correction.

## 2) Oxygen (FiO<sub>2</sub>) error

- i) The accuracy of the inspiratory oxygen concentration of the gas delivered is assessed by placing the sensor of an oxygen concentration-measuring device at the *patient-connection port* or inside the test lung. If the sensor is located at the *patient-connection port*, the value of the concentration is the flow-weighted average concentration as a function of flow during the *inflation phase*.
- ii) Evaluate the measured oxygen concentration with the acceptance criteria derived from the oxygen setting for the test and the disclosed tolerance.
- iii) Compare the resulting difference with the tolerance indicated in the *instructions for use*. If the sensor is located inside the test lung, care should be taken to allow sufficient time for the gas mixture in the test lung to reach a stable concentration. If the sensor is located in the tubing system, care should be taken to ensure that the measuring device has a sufficiently fast response to permit determination of the flow-weighted average concentration during the *inflation phase*. If the oxygen concentration-measuring device has pressure dependencies, compensate for these dependencies.

- iv) Compare each result with the acceptance criteria derived from the oxygen setting for the test in i), above and the tolerance indicated in the *instructions for use*.

### 201.12.1.103 Other inflation-types

- a) If other *inflation-types* are provided, then with each other *inflation-type* selected and the *EMS ventilator* operating in *normal condition*,

- 1) the performance at the *patient-connection port*; and
- 2) their acceptance criteria;

as determined by the *manufacturer*, shall be disclosed in the *instructions for use*.

- 3) The disclosed performance and acceptance criteria shall include the effects of the range of the *rated* input oxygen concentration.

- b) All of the acceptance criteria may be separately reported for the following ranges of intended *inspiratory volume*:

- 1)  $V_{\text{insp}} \geq 300$  ml;
- 2)  $300 \text{ ml} \geq V_{\text{insp}} \geq 50$  ml; and
- 3)  $V_{\text{insp}} \leq 50$  ml.

- c) The acceptance criteria of the performance of the *EMS ventilator* shall either be:

- 1) determined for each *VBS* configuration indicated in the *instructions for use*; or
- 2) determined for the worst-case *VBS* configuration indicated in the *instructions for use*.

NOTE 1 The worst-case *VBS* configuration can be different for each error or each *nominal inspiratory volume* range.

- d) If worst-case *VBS* configurations are used, the rationale for their selection shall be documented in the *risk management file*.

Check conformity by inspection of the *instructions for use*, inspection of the *risk management file* for the rationale, if applicable, and with the tests specified by the *manufacturer*.

### 201.12.1.104 Inspiratory volume monitoring equipment

If the *EMS ventilator* is equipped with *inspiratory volume monitoring equipment*, the accuracy of the *inspiratory volume monitoring equipment* shall be disclosed in the *instructions for use*.

Check conformity with the following.

Confirm that the *inspiratory volume monitoring equipment* accuracy as measured in 201.12.1.101 f) 7) and 201.12.1.102 f) 7) is within the accuracy disclosed in the *instructions for use*.

## 201.12.4 Protection against hazardous output

*Addition:*

### 201.12.4.101 Oxygen monitor

NOTE 1 There is guidance or rationale for this subclause contained in Clause AA.2.

- a) The *EMS ventilator*, intended to provide a *inspiratory volume* less than 50 ml [see 201.12.1.101 b) 3), 201.12.1.102 b) 3) and 201.15.102], shall either
- 1) be equipped with *oxygen monitoring equipment* for measurement of the *inspiratory oxygen concentration* that is integral to the *EMS ventilator*; or
  - 2) the *instructions for use* shall contain a statement to the effect that the *EMS ventilator* is to be equipped with *oxygen monitoring equipment* for measurement of the *inspiratory oxygen concentration* before being put into service.
- b) Integrated *oxygen monitoring equipment* shall conform with the following subclauses of ISO 80601-2-55:2018:
- 1) 201.7.4.3;
  - 2) 201.7.9.2.9.101 k);
  - 3) 201.12.1.101;
    - i) None of the display of the *EMS ventilator* shall be blanked when fulfilling the requirement of ISO 80601-2-55:2018, 201.12.1.101.4.1 c).
  - 4) 201.12.1.102;
  - 5) 201.12.1.103; and
  - 6) 208.6.1.2.
- c) Where the *oxygen monitoring equipment* is not an integral part of the *EMS ventilator*, the *instructions for use* shall include information on where to connect the *oxygen monitoring equipment*.
- d) If equipped and the intended *inspiratory volume* is <50 ml, the *oxygen monitoring equipment* shall, in addition, be equipped with an *alarm system* that includes a high oxygen level *alarm condition*.
- e) The high oxygen level *alarm condition*:
- 1) shall be at least *medium priority*, unless
  - 2) an *intelligent alarm system*, based on additional information, determines that the high oxygen level *alarm condition* is suppressed or its priority is changed.

NOTE 2 A low oxygen level *alarm condition* is required by ISO 80601-2-55.

Check conformity by inspection of the *instructions for use* or application of the tests of ISO 80601-2-55:2018.

#### 201.12.4.102 Measurement of *airway pressure*

NOTE There is guidance or rationale for this subclause contained in Clause AA.2.

- a) The *EMS ventilator* shall be equipped with *monitoring equipment* to measure the *airway pressure*.
- b) The site of actual measurement:
- 1) may be anywhere in the *EMS ventilator breathing system*; but
  - 2) the indicated value shall be referenced to the *patient-connection port*.

- c) Under steady-state conditions, the indicated *airway pressure* shall be accurate to within  $\pm$  (2 hPa +8 % of the actual reading); ( $\pm$  (2 cmH<sub>2</sub>O +8 % of the actual reading)).

Check conformity by functional testing.

### 201.12.4.103 Measurement of expired volume and low volume *alarm conditions*

NOTE 1 There is guidance or rationale for this subclause contained in Clause AA.2.

- a) An *EMS ventilator* should either:

- 1) be equipped with *monitoring equipment* for indicating the volume expired through the *patient-connection port*; or
- 2) if not so equipped, the *instructions for use* should include a statement to the effect that the *EMS ventilator* is to be equipped with *monitoring equipment* that conforms with this document before being put into service.

- b) If equipped,

- 1) unless the expired volume *monitoring equipment* is an integral part of the *EMS ventilator*, information on where to connect the expired volume *monitoring equipment* shall be disclosed in the *instructions for use*.
- 2) the accuracy of measurement of expired volumes greater than 100 ml shall be within  $\pm$ 20 % of the actual volume expired through the *patient-connection port*.
- 3) the accuracy of expired volume *monitoring equipment* shall be disclosed in the *instructions for use*.
- 4) the disclosed accuracies shall include the effects of the range of the *rated* input oxygen concentration.

NOTE 2 The use of Oxygen 93 % can reduce the accuracy of measurement due to influence of other gases such as argon on flow measurement.

- 5) the expired volume *monitoring equipment* shall be equipped with an *alarm system* to indicate when
  - i) the low-expired volume *alarm limit* is reached; and
  - ii) the high-expired volume *alarm limit* is reached.
- 6) the low-expired volume and the high-expired volume *alarm conditions*
  - i) shall be at least *medium priority*, unless
  - ii) an *intelligent alarm system*, based on additional information, determines that
    - II) the low-expired volume *alarm condition* or
    - III) the high-expired volume *alarm condition*
 is suppressed or its priority is changed.

- c) The expired volume *monitoring equipment* may be equipped:

- 1) with an *alarm system* that starts with *low priority alarm conditions* to indicate when the expired volume reaches either *alarm limit*; and
- 2) if this state continues, escalates to *medium priority alarm conditions*.

d) The expired volume *alarm limits* may be:

- 1) pre-adjusted;
- 2) *responsible organization*-configurable;

NOTE 3 There is guidance or rationale for this list element contained in Clause AA.2.

- 3) *operator*-adjustable;
- 4) *EMS ventilator*-adjustable; or
- 5) a combination of *operator*-adjustable and *EMS ventilator*-adjustable.

e) If the *alarm limits* are adjustable by the *EMS ventilator*, a summary description of the algorithm that determines the *alarm limit* values shall be disclosed in the *instructions for use*.

NOTE 4 Depending on the type of *ventilation-mode* being used, there can be more than one active *alarm limit*.

Check conformity by functional testing using the test conditions described in Table 201.104 and Table 201.105 selecting the appropriate rows based on intended *tidal volume* and inspection of the *instructions for use*. Select and set up the worst-case *VBS* configuration indicated in the *instructions for use*.

EXAMPLE Minimum or maximum *VBS* compliance.

For testing with a *humidifier*, repeat the tests at the minimum and maximum water levels (2 sets of tests for a *humidifier*).

#### **201.12.4.104 Maximum limited pressure protection device**

NOTE 1 There is guidance or rationale for this subclause contained in Clause AA.2.

A *protection device* shall be provided to prevent the *airway pressure* from exceeding for more than 200 ms the greater of:

- a) 85 hPa; or
- b) 120 % of the *maximum working pressure*.

NOTE 2 This requirement applies in both *normal condition* and *single fault condition*. See the definition of *maximum limited pressure*.

Check conformity by functional testing.

#### **201.12.4.105 High airway pressure alarm condition and protection device**

NOTE 1 There is guidance or rationale for this subclause contained in Clause AA.2.

- a) The *EMS ventilator* shall be equipped with *monitoring equipment* with an *alarm system* to indicate when the high-pressure limit for *airway pressure* is reached during the entire *respiratory cycle*.
- b) The high *airway pressure alarm condition*:
  - 1) shall be *high priority*; unless
  - 2) an *intelligent alarm system*, based on additional information, determines that:

- i) the high *airway pressure alarm condition* is suppressed; or
  - ii) its priority is changed.
- c) The high *airway pressure alarm limit* may be:
- 1) independently adjustable;
  - 2) connected to an adjustable pressure limitation; or
  - 3) related to the set pressure of the *EMS ventilator*.
- d) If the high *airway pressure alarm limit* is independently adjustable, it shall not be possible to set:
- 1) the *alarm limit* to a value less than that of the adjustable pressure limitation; and
  - 2) the adjustable pressure limitation to a value more than that of the *alarm limit*.
- e) Means shall be provided to require the *operator* to perform a deliberate sequence of actions to confirm the setting of the high *airway pressure alarm limit* to values exceeding the lower of:
- 1) 30 hPa (30 cmH<sub>2</sub>O) more than the *operator-set pressure*; or
- NOTE 2 An example of the maximum set pressure is the sum of the set *BAP* and the set  $\Delta$  *inspiratory pressure*.
- NOTE 3 An example for a bi-level positive *airway pressure ventilation-mode*, is the sum of the set *BAP<sub>H</sub>* and the set  $\Delta_H$  *inspiratory pressure*. See ISO 19223:2019, Figure C.33.
- NOTE 4 The *operator-set pressure* does not apply when by design the *ventilator* adjusts the *airway pressure* on a breath-by-breath basis.
- 2) 60 hPa (60 cmH<sub>2</sub>O).
- f) *Patient-generated transient pressure increases* should not cause the high *airway pressure alarm condition*.
- EXAMPLE Transient pressure increase caused by the *patient* coughing.
- g) The high *airway pressure alarm condition delay* shall not exceed 200 ms and the *EMS ventilator* shall:
- 1) act to attempt to cause the pressure to start to decline within that duration; and;
  - 2) act to prevent the pressure from continuing to rise.
- h) In *normal condition*, whenever the high *airway pressure alarm condition* occurs, the *EMS ventilator* shall reduce the *airway pressure* to the set *BAP* level within the lesser of:
- 1) two *respiratory cycles*; or
  - 2) 15 s.
- i) During *single fault condition*, the *EMS ventilator* shall reduce the *airway pressure* to the set *BAP* level or below within no more than 30 s.

Check conformity by functional testing.

#### 201.12.4.106 Expiratory end-tidal CO<sub>2</sub> monitoring equipment

NOTE There is guidance or rationale for this subclause contained in Clause AA.2.

- a) An *EMS ventilator* should either:
- 1) be equipped with *CO<sub>2</sub> monitoring equipment* for the measurement of the expiratory carbon dioxide concentration that is integral to the *EMS ventilator*; or
  - 2) the *instructions for use* should contain a statement to the effect that the *EMS ventilator* is to be equipped with *CO<sub>2</sub> monitoring equipment* for measurement of the expiratory carbon dioxide concentration (e.g. in the expiratory limb or at the *patient-connection port*) before being put into service.
- b) Integrated *CO<sub>2</sub> monitoring equipment* shall conform with the following subclauses of ISO 80601-2-55:2018:
- 1) 201.7.4.3;
  - 2) 201.7.9.2.9.101 k);
  - 3) 201.12.1.101;
    - i) None of the display of the *EMS ventilator* shall be blanked when fulfilling the requirement of ISO 80601-2-55:2018, 201.12.1.101.4.1 c).
  - 4) 201.12.1.102;
  - 5) 201.12.1.103; and
  - 6) for expired CO<sub>2</sub> concentration, 208.6.1.2.
- c) Where the *CO<sub>2</sub> monitoring equipment* is not an integral part of the *EMS ventilator*, the *instructions for use* should include information on where to connect the *CO<sub>2</sub> monitoring equipment*.

Check conformity by inspection.

#### 201.12.4.107 Protection against inadvertent setting of high airway pressure

Means shall be provided to require the *operator* to perform a deliberate sequence of actions to confirm any *airway pressure* settings exceeding 60 hPa (60 cmH<sub>2</sub>O). See also 201.12.4.105 e) 2) and 206.101 b) 6).

Check conformity by functional testing.

### 201.13 Hazardous situations and fault conditions for ME equipment

IEC 60601-1:2005+AMD1:2012+AMD2:2020, Clause 13 applies, except as follows:

*Addition:*

#### 201.13.2.101 Additional specific single fault conditions

NOTE There is guidance or rationale for this subclause contained in Clause AA.2.

An *EMS ventilator* shall be so designed and constructed that the following *single fault conditions* shall not cause an unacceptable *risk*:

- a) disconnection or blockage of the gas delivery to the *patient-connection port* from the *EMS ventilator*;
- b) disconnection or blockage of the gas flow pathway from the *patient-connection port* to the *exhaust port*;
- c) failure to install, removal or failure of an *operator-detachable breathing system filter*; and
- d) interruption of a *functional connection* between parts of the *EMS ventilator* or *ME system*.

EXAMPLE 1 Loss of communication between the *EMS ventilator* and its remote (wired or wireless) control or monitoring module.

EXAMPLE 2 Loss of communication between the *EMS ventilator* and its *distributed alarm system*.

EXAMPLE 3 Loss of communication between the *EMS ventilator* and the means for generating remote *alarm signals*.

Check conformity by functional testing and inspection of *risk management file*.

#### **201.13.102 Failure of one gas supply to an EMS ventilator**

- a) For an *EMS ventilator* with multiple high-pressure gas supplies, following the failure of one gas supply, the *EMS ventilator* shall:
  - 1) continue ventilation using the other gas supply;
  - 2) maintain *basic safety*; and
  - 3) maintain *essential performance*.
- b) The *EMS ventilator* shall be equipped with an *alarm system* that detects a *technical alarm condition* to indicate this gas-supply failure.
- c) The gas supply failure *technical alarm condition*
  - 1) shall be at least *low priority* if *ventilation* is maintained, unless
  - 2) an *intelligent alarm system*, based on additional information, determines that the gas-supply-failure *technical alarm condition* is suppressed; or
  - 3) shall be at least *medium priority* otherwise, unless
  - 4) an *intelligent alarm system*, based on additional information, determines that the gas-supply-failure *technical alarm condition* is suppressed.

Check conformity by functional testing.

#### **201.13.103 Independence of ventilation control function and related risk control measures**

NOTE There is guidance or rationale for this subclause contained in Clause AA.2.

- a) A *single fault condition* shall not cause the simultaneous failure of:
  - 1) a *ventilation control function*; and
  - 2) the corresponding *protection device*.

- b) A *single fault condition* shall not cause either of the following to fail in such a way that the loss of the *ventilation control function* is not detected:
- 1) a *ventilation-control function* and the corresponding *monitoring equipment*; or
  - 2) a *ventilation-control function* and the corresponding *alarm system*

Check conformity by inspection and functional testing.

## 201.14 Programmable electrical medical systems (PEMS)

IEC 60601-1:2005+AMD1:2012+AMD2:2020, Clause 14 applies, except as follows:

Addition:

### 201.14.1 General

Amendment (add prior to the compliance check):

- aa) The *programmable electronic subsystems (PESS)* of an *EMS ventilator* shall be developed with a *design process* conforming with:
- 1) IEC 62304:2006+AMD1:2015; and
  - 2) If provided with an *electronic functional connection*, IEC 81001-5-1:2021.
- bb) The *ventilation control software items* of the *EMS ventilator PESS* without an independent *risk control* measure external to the *PESS* shall be considered as software safety Class C.

### 201.14.101 Cybersecurity capabilities of EMS ventilators

- a) An *EMS ventilator* should contribute to the safe operation related to *cybersecurity*.
- b) *Risk control* measures as specified in IEC/TR 60601-4-5:2021, Clauses 4 to 7, should be implemented, as appropriate, with the following additions to IEC/TR 60601-4-5:2021, 4.6.3.
- 1) *Essential function*: After disconnection of all data interfaces to *IT-network* connections possibly subject to *attack* (via physical or logical disconnection), if necessary, after restarting the *EMS ventilator* once, all clinical functions except those related to the affected remote functionality should be in place and operate as intended.
  - 2) *Firecall functions*: As an *EMS ventilator* is used with continuous *healthcare professional operator* supervision, the user interface need not require authentication for *operators* and thus need not require a *firecall* function. If authentication for *operators* to *operator accessible* settings is in place also for the user interface, a *firecall* function should be able to overrule that *operator* authentication combined with a log entry protected against modifications by the *operator* (i.e. *responsible organization log*). *IT-network* interfaces including those for remote access typically require authentication but should not include *firecall* functions as the *healthcare professional operator* is always close to the *EMS ventilator* and *patient*.
  - 3) Target *security level SL-T*: For user interfaces and for remote setting of clinical functions, the expected *security level* of SL-T, as specified in IEC/TR 60601-4-5:2021, is 1 or better. However, authentication of *operators* is not required for the user interface, only for access via data interfaces. The expected *security level* of SL-T, as specified in IEC/TR 60601-4-5:2021, is 2 or better for secure updates or restorage of software, remotely or locally.

## 201.15 Construction of *ME* equipment

IEC 60601-1:2005+AMD1:2012+AMD2:2020, Clause 15 applies, except as follows:

*Addition:*

### 201.15.3.5.101 Additional requirements for rough handling

NOTE 1 There is guidance or rationale for this subclause contained in Clause AA.2.

- a) An *EMS ventilator* and its parts, including applicable *accessories* shall have adequate mechanical strength to withstand the mechanical stress caused by *normal use*, pushing and rough handling while operating.
- b) For this test, the *EMS ventilator* and its parts, and applicable *accessories*, shall be mounted using the mounting *accessories* indicated in the *accompanying documents*. If more than one mounting system is described in the *accompanying documents*, multiple tests are required.

NOTE 2 An *EMS ventilator* tested and conforming with a more severe requirement is considered to conform with the corresponding requirement of this subclause.

- c) During the following test, an *EMS ventilator* shall maintain *basic safety* and *essential performance*.
- d) During this testing, volume and pressure *alarm condition alarm limits* shall be set to their least sensitive levels.
- e) During the testing, the error of:
  - 1) the *inspiratory volume* of individual *inflations* shall not deviate by more than 35 % of the *inspiratory volume* measured prior to the test;
  - 2) the *inspiratory volume* averaged over a one min interval shall not deviate by more than 25 % of the *inspiratory volume* measured prior to the test;
    - i) One-minute averaged testing need not be performed if the *inspiratory volume* of individual *inflations* does not deviate by more than 25 %.
  - 3) the *PEEP* of individual *inflations* shall not deviate by more than 5,0 hPa (5,0 cmH<sub>2</sub>O) from the *PEEP* measured prior to the test; and
- f) Check conformity by performing the following tests:
  - 1) Broadband random vibration test in accordance with IEC 60068-2-64:2008, using the following conditions:
    - i) acceleration amplitude:
      - 10 Hz to 100 Hz: 1 (m/s<sup>2</sup>)<sup>2</sup>/Hz;
      - 100 Hz to 200 Hz: -3 dB per octave;
      - 200 Hz to 2 000 Hz: 0,5 (m/s<sup>2</sup>)<sup>2</sup>/Hz;
    - ii) duration: 30 min per perpendicular axis (3 total).

### 201.15.4.1 Construction of *connectors*

*Addition:*

aa) *Operator-detachable VBS connectors* are exempt from this requirement.

### 201.15.101 Mode of operation

An *EMS ventilator* shall be suitable for *continuous operation*.

Check conformity by inspection.

### 201.15.102 Delivered oxygen concentration

NOTE There is guidance or rationale for this subclause contained in Clause AA.2.

- a) An *EMS ventilator* shall provide at least two different settings for the inspiratory oxygen concentration.
- b) An *EMS ventilator* shall be capable of delivering gas to the *patient* containing an oxygen concentration over the range of flowrates indicated in the *instructions for use*:
  - 1) of at least 85 % of the input oxygen concentration in one setting; and
  - 2) of no greater than 65 % of the input oxygen concentration in a second setting.
- c) An *EMS ventilator* intended to provide an *inspiratory volume* of less than 50 ml shall be equipped with an oxygen concentration setting with a delivery error of not more than 20 %.
  - 1) Such an *EMS ventilator* shall have an oxygen concentration setting range from room air to 85 % of the input oxygen concentration.

Check conformity by inspection of the *instructions for use* and functional testing.

### 201.15.103 Accessory self-check

- a) An *EMS ventilator* should be equipped with means that allow the determination of whether or not the *VBS* resistance and compliance characteristics are outside the values necessary to maintain normal operation.

NOTE Additional requirements are found in 201.7.9.2.8.101.

- b) This means may require *operator* action.

Check conformity by functional testing.

## 201.16 ME systems

IEC 60601-1:2005+AMD1:2012+AMD2:2020, Clause 16 applies, except as follows:

*Addition:*

### 201.16.1.101 Additional general requirements for ME systems

*Accessories* connected to the *VBS* shall be considered to:

- a) be part of the *EMS ventilator*; or
- b) form an *ME system* with the *EMS ventilator*.

For *accessories* as part of the *EMS ventilator*, check conformity in combination with the *EMS ventilator*.

The testing need not to be repeated for each combination for a family of *accessories* with a family of *EMS ventilators*. See 201.5 for worst-case combinations.

For *accessories* forming an *ME system* with the *EMS ventilator*, apply the requirements of IEC 60601-1:2005+AMD1:2012+AMD2:2020, 16.1. Check conformity by application of the relevant tests of IEC 60601-1:2005+AMD1:2012+AMD2:2020, 16.1, and by application of the relevant tests.

### 201.16.2.101 Additional requirements for power supply

- a) An *EMS ventilator*, which depends on an external gas source to supply energy for *ventilation* shall comprise an *ME system* with the external gas source.
- b) A *portable* gas source shall indicate, as appropriate:
  - 1) the gas source pressure;
  - 2) remaining gas source volume; or
  - 3) for an oxygen concentrator, the remaining run time prior to exhaustion of its source of electrical power.
- c) The indication shall be:
  - 1) visible from the intended *operator's position*; and
  - 2) *clearly legible* under the conditions of 6.1 of IEC 60601-1-12:2014+AMD1:2020.
    - i) A pressure reducer with manometer mounted to the *portable* gas source, where the manometer is *clearly legible* under the conditions of 6.1 of IEC 60601-1-12:2014+AMD1:2020, may be used.

Check conformity by inspection and the relevant tests of IEC 60601-1-12:2014+AMD1:2020.

## 201.17 Electromagnetic compatibility of *ME equipment* and *ME systems*

IEC 60601-1:2005+AMD1:2012+AMD2:2020, Clause 17 applies.

*Addition:*

### 201.101 Gas connections

#### 201.101.1 Protection against reverse gas leakage

NOTE 1 There is guidance or rationale for this subclause contained in Clause AA.2.

For an *EMS ventilator* with two or more *high-pressure inlets*,

- a) means shall be provided to limit reverse gas leakage flowrate from *gas intake ports* into the supply system of the same gas to a flowrate less than 100 ml/min in *normal condition* or *single fault condition*.
- b) means shall be provided to limit cross leakage from gas supplied through one *high-pressure inlet* into the supply system of a different gas to less than 100 ml/h in *normal condition* or *single fault condition*.

Check conformity by functional testing.

## 201.101.2 Connection to a *high-pressure inlet*

### 201.101.2.1 Low-pressure hose assembly

If an *operator-detachable low-pressure hose assembly* is provided for connection between the *EMS ventilator* and either the *medical gas pipeline system* or a pressure regulator, it shall conform with ISO 5359:2014+AMD1:2017.

Check conformity by application of the tests of ISO 5359:2014+AMD1:2017.

### 201.101.2.2 Filter

NOTE 1 There is guidance or rationale for this subclause contained in Clause AA.2.

Each *high-pressure inlet* shall be provided with a filter having a pore size less than or equal to 100 µm.

NOTE 2 The need for particle filtration in oxygen-enriched environments including the proper choice of filter materials is discussed in ISO 15001:2010.

NOTE 3 Depending on the sensitivity against particles of the components used in the *gas pathways* (e.g. flow sensors) to particles, filtration of smaller particles can be needed.

Check conformity by inspection.

## 201.101.3 VBS connectors

### 201.101.3.1 General

NOTE There is guidance or rationale for this subclause contained in Clause AA.2.

a) *Operator-detachable VBS connectors* through which the main flow of gas to or from the *patient* passes in *normal condition*, excluding the *patient-connection port*:

1) shall be:

- i) a 15 mm *connector* conforming with ISO 5356-1:2015;
- ii) a 22 mm *connector* conforming with ISO 5356-1:2015; or
- iii) for an *EMS ventilator* only intended for *tidal volumes* of  $\leq 300$  ml, a 11,5 mm *connector* conforming with ISO 5356-1:2015 for a neonatal or paediatric use *VBS*; or

2) may be a non-conical *connector* that does not engage with a conical *connector* conforming with ISO 5356-1:2015.

Check conformity by application of the tests of ISO 5356-1:2015 and functional testing.

### 201.101.3.2 Other named ports

#### 201.101.3.2.1 Patient-connection port

The *patient-connection port* shall be one of the following:

- a) a 15 mm conical socket *connector* conforming with ISO 5356-1:2015;
- b) a 22 mm conical cone *connector* conforming with ISO 5356-1:2015; or
- c) a coaxial 15 mm/22 mm conical *connector* conforming with ISO 5356-1:2015.

Check conformity by application of the tests of ISO 5356-1:2015.

**201.101.3.2.2 Gas output port and gas return port**

- a) The *gas output port* and the *gas return port* shall be one of the following
- 1) a 22 mm conical cone *connector* conforming with ISO 5356-1:2015.
  - 2) a 15 mm conical cone *connector* conforming with ISO 5356-1:2015.
  - 3) a coaxial 15 mm/22 mm conical *connector* conforming with ISO 5356-1:2015.
  - 4) a non-conical *connector* that does not engage with a conical *connector* conforming with ISO 5356-1:2015.
- b) Notwithstanding this requirement, an *EMS ventilator* only intended for *inspiratory volumes* of ≤300 ml, may be equipped with a *gas output port* and a *gas return port* using a 11,5 mm conical cone *connector* conforming with ISO 5356-1:2015.

Check conformity by application of the tests of ISO 5356-1:2015.

**201.101.3.2.3 Flow-direction-sensitive components**

NOTE There is guidance or rationale for this subclause contained in Clause AA.2.

Any *flow-direction-sensitive component* of the *VBS* detachable without the use of a *tool* shall be so designed that it cannot be fitted in such a way that it presents an unacceptable *risk* to the *patient*.

Check conformity by inspection of *operator-detachable flow-direction-sensitive components* and inspection of the *risk management file*.

**201.101.3.2.4 Gas pathway connection port**

NOTE 1 There is guidance or rationale for this subclause contained in Clause AA.2.

- a) If provided, each *accessory* port shall:
- 1) be provided with a means to secure the *accessory* in position; and
  - 2) be provided with a means to secure closure after removal of the *accessory*.

NOTE 3 This port connects to the *gas pathway* and is generally used e.g. for measuring gas concentrations or for introduction of liquids.

NOTE 4 For the purposes of this document, the temperature probe port specified in ISO 80601-2-74 is not considered an *accessory* port.

- b) The *gas pathway* connection port of a *EMS ventilator*, *VBS* or *accessory* may conform with ISO 80369-7:2021.
- c) A *gas pathway* connection port that conforms with ISO 80369-7:2021 shall:
- 1) be *marked* with either:
    - i) the *symbol* ISO 7000-1641 (see IEC 60601-1:2005+AMD1:2012+AMD2:2020, Table D.1, *symbol* 11); or
    - ii) if the *marking* is the primary *risk control* measure, the *safety sign* ISO 7010-M002 (see IEC 60601-1:2005+AMD1:2012+AMD2:2020, Table D.2, *safety sign* 10).

- 2) include in its *instructions for use* a warning to the effect that “Warning: As this [insert name of medical device here] uses a Luer connector for other than intravascular or hypodermic access, there is a possibility that an inadvertent connection can occur between this [insert name of medical device here] and another medical device or accessory using a Luer connector, which can result in a hazardous situation causing harm to the patient. Special measures need be taken by the user to mitigate these reasonably foreseeable risks.”

Check conformity by inspection and, if applicable, by application of the tests of ISO 80369-7:2021.

#### 201.101.3.2.5 Gas exhaust port

- a) If a *connector* is provided for the gas *exhaust port*, it shall be a 30 mm *connector* conforming with ISO 5356-1:2015.

NOTE A 30-mm *connector* conforming with ISO 5356-1:2015 is suitable for connection to an anaesthesia gas scavenging system (AGSS) that conforms with ISO 80601-2-13.

- b) An *EMS ventilator* shall be designed so that any provided gas *exhaust port* is not obstructed during *intended use*.

Check conformity by inspection and application of the tests of ISO 5356-1:2015.

#### 201.101.3.2.6 Gas intake port

An *EMS ventilator* that incorporates a *gas intake port* shall be designed such that the *gas intake port* cannot be blocked by a flat surface.

EXAMPLE A *gas intake port* incorporating a protective cage, or with multiple intake orifices arranged in a non-planar geometry.

Check conformity by functional testing.

### 201.102 Requirements for the VBS and accessories

#### 201.102.1 General

NOTE There is guidance or rationale for this subclause contained in Clause AA.2.

All *EMS ventilator breathing systems*, their parts and *accessories* that can affect the *basic safety* or *essential performance* of the *EMS ventilator* shall conform with the requirements of this document, whether they are produced by the *manufacturer* of the *EMS ventilator* or by another entity (“third-party manufacturer” or healthcare provider).

Check conformity by the tests of this document.

#### 201.102.2 Labelling

- a) The *accompanying document* provided with each *VBS*, its parts or *accessories*, conforming with 201.102.1, shall include the *model or type reference* of at least one compatible *EMS ventilator*.
- b) Statements shall be included in the *accompanying documents* of each *EMS ventilator breathing system*, part or *accessory* to the effect that:
  - 1) ventilator breathing systems, their parts and accessories are validated for use with specific ventilators;
  - 2) incompatible parts can result in degraded performance; and

- 3) the responsible organization is responsible for ensuring the compatibility of the ventilator and all of the parts used to connect to the patient before use.

Check conformity by inspection of the *accompanying document*.

### 201.102.3 Breathing sets

Breathing sets, other than heated breathing sets, intended for use in the *VBS* shall conform with ISO 5367:2014, 5.3.4.

Check conformity by application of the tests of ISO 5367:2014.

### 201.102.4 Water vapour management

NOTE There is guidance or rationale for this subclause contained in Clause AA.2.

#### 201.102.4.1 Humidification system

NOTE There is guidance or rationale for this subclause contained in Clause AA.2.

Any *humidifier*, including heated breathing tubes, either incorporated into the *EMS ventilator* or recommended for use with the *EMS ventilator*, shall conform with ISO 80601-2-74:2021.

Check conformity by application of the tests of ISO 80601-2-74:2021.

#### 201.102.4.2 Heat and moisture exchanger (HME)

Any *heat and moisture exchanger*, either incorporated into the *VBS* or recommended for use with the *VBS*, shall conform with:

- a) ISO 9360-1:2000; or
- b) ISO 9360-2:2001.

Check conformity by application of the tests of ISO 9360-1:2000 or ISO 9360-2:2001.

### 201.102.5 Breathing system filters

Any *breathing system filter*, either incorporated into the *EMS ventilator* or recommended for use with the *EMS ventilator*, shall conform with the relevant requirements of:

- a) ISO 23328-1:2003; and
- b) ISO 23328-2:2002.

Check conformity by application of the tests of ISO 23328-1:2003 and ISO 23328-2:2002.

### 201.102.6 Leakage from complete *VBS*

NOTE There is guidance or rationale for this subclause contained in Clause AA.2.

Unintended leakage from the *VBS* should not exceed:

- a) 200 ml/min at 50 hPa (50 cmH<sub>2</sub>O) for an *EMS ventilator* intended to provide a *tidal volume* greater than 300 ml;
- b) 100 ml/min at 40 hPa (40 cmH<sub>2</sub>O) for an *EMS ventilator* intended to provide a *tidal volume* between 300 ml and 50 ml; and
- c) 50 ml/min at 20 hPa (20 cmH<sub>2</sub>O) for an *EMS ventilator* intended to provide a *tidal volume* less than 50 ml.

### 201.103 Spontaneous breathing during loss of ventilation

- a) A *protection device* shall be provided to allow spontaneous breathing when normal *ventilation* is compromised as a result of the electrical or pneumatic supply power being outside the values necessary for normal operation.
- b) Under these conditions, the inspiratory and expiratory pressure drop measured at the *patient-connection port* with all recommended *accessories* in place shall not exceed 6,0 hPa (6,0 cmH<sub>2</sub>O) at a flowrate of:
  - 1) 30 l/min for an *EMS ventilator* intended to provide a *tidal volume*,  $V_T \geq 300$  ml;
  - 2) 15 l/min for an *EMS ventilator* intended to provide a *tidal volume*,  $300 \text{ ml} \geq V_T \geq 50$  ml;
  - 3) 2,5 l/min for an *EMS ventilator* intended to provide a *tidal volume*,  $V_T \leq 50$  ml.

NOTE This requirement is intended to allow the *patient* to breathe spontaneously under compromised conditions.

Check conformity by functional testing and measurement of flow, pressure, and resistance at the *patient-connection port* with that combination of *accessories* indicated in the *instructions for use* which produces the highest pressure drop.

### 201.104 Indication of duration of operation

NOTE There is guidance or rationale for this subclause contained in Clause AA.2.

- a) The *EMS ventilator* should have means to indicate visually the cumulative hours of operation of the *EMS ventilator*, either:
  - 1) automatically; or
  - 2) by *operator* action.
- b) The *EMS ventilator* should also have means to indicate visually:
  - 1) the time since the last preventive maintenance; or
  - 2) the time until the next recommended preventive maintenance.

Check conformity by inspection.

### 201.105 Functional connection

#### 201.105.1 General

*Basic safety* and *essential performance* shall be maintained if connections to a *functional connection* of an *EMS ventilator* are:

- a) disrupted; or
- b) if the equipment connected to those parts fails.

Check conformity by functional testing.

**201.105.2 Connection to an electronic health record**

NOTE There is guidance or rationale for this subclause contained in Clause AA.2.

- a) An *EMS ventilator* should be equipped with a *functional connection* that permits data transmission from the *EMS ventilator* to an electronic health record.
- b) The data transmission should be capable of transmitting the information described in Annex BB.

**201.105.3 Connection to a distributed alarm system**

NOTE There is guidance or rationale for this subclause contained in Clause AA.2.

An *EMS ventilator* should be equipped with a *functional connection* that permits connection to a *distributed alarm system*.

**201.106 Display loops****201.106.1 Pressure-volume loops**

- a) If an *EMS ventilator* is provided with the display of pressure-volume loops the graph shall use:
  - 1) *inspiratory volume* on the vertical axis; and
  - 2) *airway pressure* on the horizontal axis.
- b) Positive values shall be on the top and the right of the display.
- c) Increases in *inspiratory volume* shall be positive values.
- d) The volume shall be reset to the origin at the beginning of each breath.

Check conformity by inspection.

**201.106.2 Flow-volume loops**

- a) If an *EMS ventilator* is provided with the display of flow-volume loops, the graph shall use:
  - 1) flowrate on the vertical axis; and
  - 2) *inspiratory volume* on the horizontal axis.
- b) Positive values shall be on the top and the right of the display.
- c) Gas flow to the *patient* (inspiratory flow) and increases in *inspiratory volume* shall be positive values.
- d) The volume shall be reset to the origin at the beginning of each breath.
- e) The *EMS ventilator* may be provided with an additional optional display configuration for the flow-volume loop where gas flow from the *patient* (expiratory flow) is represented as a positive value.

Check conformity by inspection.

**201.107 Timed ventilatory pause**

NOTE There is guidance or rationale for this subclause contained in Clause AA.2.

### 201.107.1 Expiratory pause

- a) An *EMS ventilator* may be equipped with an *operator-controlled* means to pause the *EMS ventilator* in expiration.
- b) If an *EMS ventilator* is equipped with a means to pause the *EMS ventilator* in expiration,
- 1) more than one *expiratory pause* function may be provided.
  - 2) the maximum allowable duration of an *expiratory pause* shall be 60 s.
  - 3) means may be provided to initiate the *expiratory pause* from a *functional connection*.
- c) If an *EMS ventilator* is equipped with a timed means to pause the *EMS ventilator* in expiration,
- 1) the duration of the *expiratory pause* may be:
    - i) *responsible organization-configurable*; or  
NOTE 1 There is guidance or rationale for this list item contained in Clause AA.2.
    - ii) *operator-adjustable*.
  - 2) during the *expiratory pause*, any apnoea-related ventilatory *alarm condition* that would be caused by this *expiratory pause* shall be:
    - i) *audio paused*; or
    - ii) *alarm paused*for the duration of the *expiratory pause*.
  - 3) in addition to the requirements for *alarm signal* inactivation in 6.8.5 of IEC 60601-1-8:2006+AMD1:2012+AMD2:2020, the *EMS ventilator* shall indicate the presence of the *expiratory pause* with at least:
    - i) an *information signal*; or
    - ii) a *low priority alarm condition*.

NOTE 2 An *expiratory pause* can be equivalent to placing the *EMS ventilator* into standby mode or *CPAP* and automatically resuming *ventilation* after a pre-determined duration.

NOTE 3 The *expiratory pause* can be used to permit defibrillator ECG analysis with a deflated *lung*.

Check conformity by inspection and functional testing.

### 201.107.2 Inspiratory pause

- a) An *EMS ventilator* may be equipped with an *operator-controlled* means to pause automatic *ventilation* at end-inspiration.
- b) If an *EMS ventilator* is equipped with a means to pause the *EMS ventilator* in inspiration,
- 1) more than one *inspiratory pause* function may be provided.
  - 2) the maximum duration of a non-adjustable *inspiratory pause* shall be 10 s.
  - 3) the maximum allowable duration of an adjustable *inspiratory pause* shall be 40 s.
  - 4) means may be provided to initiate the *inspiratory pause* from a *functional connection*.

- c) If a *EMS ventilator* is equipped with a timed means to pause the *EMS ventilator* in inspiration,
- 1) the duration of the *inspiratory pause* may be:
    - i) non-adjustable;
    - ii) *responsible organization*-configurable; or  
NOTE There is guidance or rationale for this list item contained in Clause AA.2.
    - iii) *operator*-adjustable.
  - 2) the high-pressure *alarm condition* and *protection device* of 201.12.4.105 shall remain active during an *inspiratory pause*.
  - 3) during the *inspiratory pause*, any:
    - i) apnoea *alarm condition*; or
    - ii) sustained *airway pressure alarm condition*  
that would be caused by this *inspiratory pause* should, for the duration of the *inspiratory pause*, be:
      - iii) *audio paused*; or
      - iv) *alarm paused*.
  - 4) in addition to the requirements for *alarm signal* inactivation of 6.8.5 of IEC 60601-1-8:2006+AMD1:2012+AMD2:2020, the *EMS ventilator* shall indicate the presence of the *inspiratory pause* with
    - i) at least an *information signal*; or
    - ii) a *low priority alarm condition*.

Check conformity by inspection and functional testing.

## 202 Electromagnetic disturbances — Requirements and tests

IEC 60601-1-2:2014+AMD1:2020 applies except as follows:

### 202.4.3.1 Compliance criteria

NOTE There is guidance or rationale for this subclause contained in Clause AA.2.

*Amendment (replace the second dash of 4.3.1 with):*

- the *EMS ventilator* operated using the worst-case conditions and parameters of Table 201.104 or Table 201.105 selected by intended *tidal volume*, as appropriate. During this testing, set the volume and pressure *alarm limits* to their least sensitive levels.

### 202.5.2.2.1 Requirements applicable to all *ME equipment* and *ME systems*

*Amendment (add note to list element b)):*

NOTE The requirements of this document are not considered deviations or allowances.

*Addition:*

### 202.8.1.101 Additional general requirements

NOTE There is guidance or rationale for this subclause contained in Clause AA.2.

- a) The following degradations, if associated with *basic safety* and *essential performance*, shall not be allowed:
- 1) component failures;
  - 2) changes in programmable parameters or settings;
  - 3) reset to default settings;
  - 4) change of *ventilation-mode*;  
EXAMPLE Change of *inflation-type, ventilation-mode, set rate, I:E ratio*.
  - 5) initiation of an unintended operation;
  - 6) during the testing, the error of the:
    - i) *inspiratory volume* of individual *inflations* deviating by more than 35 % of the *inspiratory volume* measured prior to the test;
    - ii) *inspiratory volume* averaged over a one-minute interval deviating by more than 25 % of the *inspiratory volume* measured prior to the test;
      - l) One-minute averaged testing need not be performed if the volume of individual *inflations* does not deviate by more than 25 %.
    - iii) *PEEP* of individual *inflations* deviating by more than 5 hPa from the *PEEP* measured prior to the test; and
    - iv) delivered  $\text{FiO}_2$  averaged over a one min interval deviating by more than the deviation disclosed by the *manufacturer* in the *instructions for use*.
      - l) One-minute averaged testing need not be performed if the  $\text{FiO}_2$  of individual *inflations* does not deviate by more than the deviation disclosed in the *instructions for use*.
- b) The *EMS ventilator* may exhibit temporary degradation of performance (e.g. deviation from the performance indicated in the *instructions for use*) that does not affect *basic safety* or *essential performance*.

## 206 Usability

IEC 60601-1-6:2010+AMD1:2013+AMD2:2020 applies except as follows:

*Addition:*

### 206.101 *Primary operating functions*

- a) For an *EMS ventilator*, the following shall be considered *primary operating functions*:
- 1) setting the *operator-adjustable controls* from the intended *operator's position*;
    - i) setting *alarm limits*;
    - ii) inactivating *alarm signals*;
    - iii) switching between different *ventilation-modes* and *inflation-types*;
    - iv) setting *ventilation control parameters*;  
EXAMPLE 1 *Tidal volume, set rate, BAP, pressure support*
  - 2) observing and identifying the monitored *ventilation* parameters from the intended *operator's position*;

EXAMPLE 2 *Airway pressure and expired volume*

EXAMPLE 3 *FiO<sub>2</sub> or etCO<sub>2</sub>*

3) observing and identifying:

- i) the *alarm signals*; and
- ii) the *alarm signal* inactivation states;

4) recognizing the presence of an *alarm condition* in the intended *use environments* from the intended *operator positions*;

EXAMPLE 4 *The operator in the ambulance next to the patient.*

EXAMPLE 5 *The operator at the scene of an emergency situation with sirens and flashing emergency lights.*

5) setting the *operator*-adjustable controls while a *portable EMS ventilator* is being carried;

- i) setting *alarm limits*;
- ii) inactivating *alarm signals*;
- iii) switching between different *ventilation-modes* and *inflation-types*;
- iv) setting *ventilation* control parameters;

EXAMPLE 6 *Tidal volume, set rate, BAP, pressure support*

6) observing and identifying the monitored *ventilation* parameters while the *portable EMS ventilator* is being carried;

EXAMPLE 7 *Airway pressure and expired volume*

7) configuring

- i) the *VBS*, including
- ii) connection of the detachable parts of the *VBS* to the *EMS ventilator*;

EXAMPLE 8 *HME, tubing, breathing system filter*

8) connecting or disconnecting the *patient-connection port* of the *VBS* to the *patient* interface;

9) starting the *EMS ventilator* from power off;

10) turning off the *EMS ventilator*;

11) performing a basic pre-use functional check of the *EMS ventilator* including the *alarm system*;

12) connecting the *EMS ventilator*:

- i) to each gas supply indicated in the *instructions for use* including, if applicable
  - I) when attached in the ambulance; or
  - II) within its carrying case;
- ii) to each power source indicated in the *instructions for use*, including, if applicable
  - I) when attached in the ambulance; or
  - II) within its carrying case; and

13) *processing* the *EMS ventilator* between *patient* uses.

EXAMPLE 9 Readyng the *EMS ventilator* for use on the next *patient*. See also 201.11.6.6.

b) The following functions, if available, also shall be considered *primary operating functions*:

- 1) starting *ventilation* from standby;
- 2) activating standby;
- 3) attaching the *portable EMS ventilator* to an ambulance;
- 4) connecting and disconnecting a *distributed alarm system*;
- 5) setting of the adjustable high-pressure *alarm limit* to values exceeding 60 hPa (60 cmH<sub>2</sub>O); and
- 6) setting of the *airway pressure* to values exceeding 60 hPa (60 cmH<sub>2</sub>O).

c) The following actions associated with *ventilation* also shall be considered *primary operating functions*:

NOTE For the purposes of this document, the following functions are considered *primary operating functions* even though they might not be performed on the *EMS ventilator's operator interface*.

- 1) if provided, adding medication to the gas flowing into the *patient*; and

EXAMPLE 10 Nebulisation or injecting fluids into an *accessory port connection* of the *VBS*.

- 2) understanding the need for alternative means of *ventilation* as indicated in the *instructions for use*.

## 206.102 Training

NOTE There is guidance or rationale for this subclause contained in Clause AA.2.

For an *EMS ventilator*, in the application of the requirements of IEC 62366-1:2015, 5.6 and 5.8, training shall be considered necessary for both:

- a) the *operator*; and
- b) the designee of the *responsible organization* (e.g. *service personnel* or *processing personnel*).

## 208 General requirements, tests and guidance for *alarm systems in medical electrical equipment and medical electrical systems*

IEC 60601-1-8:2006+AMD1:2012+AMD2:2020 applies except as follows:

### 208.6.3.2.2.1 4 m (distant) visual alarm signals

NOTE 1 There is guidance or rationale for this subclause contained in Clause AA.2.

*Amendment (Add prior the first paragraph):*

*An EMS ventilator shall be equipped with 4 m (distant) visual alarm signals.*

*Addition:*

### 208.6.8.3.101 Additional requirements for global indefinite *alarm signal* inactivation states

An *EMS ventilator* shall not be equipped with a means to initiate a global *alarm off* while connected to a *patient*.

Check conformity by functional testing.

### 208.6.8.4.101 Additional requirements for termination of *alarm signal* inactivation

NOTE 1 There is guidance or rationale for this subclause contained in Clause AA.2.

The duration of *audio paused* and *alarm paused* for the *alarm conditions* required by this document shall not exceed 120 s without *operator* intervention.

NOTE 2 This permits an *operator* to deliberately extend the duration of *audio paused* by no more than 120 s following each direct action.

Check conformity by functional testing.

### 208.6.12.3 Responsible organization alarm system logging

Replacement:

- a) A *EMS ventilator* shall be provided with a *responsible organization* log with a capacity of at least 1 000 events in total.
- b) Viewing the *responsible organization alarm system* log shall be restricted to the *responsible organization* in accordance with IEC 60601-1-8:2006+AMD1:2012+AMD2:2020, 6.7.
- c) The *responsible organization alarm system* log shall include all of the information contained in the *operator alarm system* log, if an *operator alarm system* log is provided.
- d) The *responsible organization alarm system* log shall contain:
  - 1) the *EMS ventilator* settings; and
  - 2) each change of those settings.

EXAMPLE 1 The name of the *ventilator* preset in use and any changes made to it.
- e) The *responsible organization alarm system* log shall contain:
  - 1) the *alarm settings*; and
  - 2) each change of those settings.

EXAMPLE 2 The name of the *alarm preset* in use and any changes made to it.
- f) Means shall not be provided for the *operator* or *responsible organization* to edit or delete entries in the *responsible organization alarm system* log.
- g) The *responsible organization alarm system* log shall be retained when the *alarm system* is powered down.
- h) The *technical description* shall indicate what happens to the contents of the log after the *alarm system* has experienced an unexpected total loss of power (*supply mains* and *internal electrical power source*).

- i) The *technical description* shall indicate:
- 1) the *responsible organization alarm system* log capacity; and
  - 2) what happens to the contents of the log as it reaches capacity.

EXAMPLE 3 The *alarm system* discards the oldest data when the log becomes full.

- j) The *alarm system* should log *technical alarm conditions* for servicing and maintenance purposes. This log should not be resettable or editable by *operator* action.
- k) The log may be provided either within the equipment or remotely through a communications interface.

Check conformity by inspection and functional testing.

## **212 Requirements for *medical electrical equipment* and *medical electrical systems* intended for use in the *emergency medical services environment***

IEC 60601-1-12:2014+AMD1:2020 applies except as follows:

### **212.4.2.2.1 Continuous operating conditions**

Amendment (replace the first dash with):

- aa) where *permanently installed* or *fixed*, a temperature range of 0 °C to +40 °C; or
- bb) where not *permanently installed* or *fixed*, a temperature range of -10 °C to +40 °C.

### **212.8.1.1 Ingress of water or particulate matter into *ME equipment***

NOTE There is guidance or rationale for this subclause contained in Clause AA.2.

Amendment (in the first sentence of the first paragraph):

Replace 'IP33' with 'IP34'.

### **212.8.2 Additional requirements for interruption of the power supply to *ME equipment* and *ME system***

Amendment (replace the paragraph following note 2 with):

- aa) The *EMS ventilator* shall be equipped with a means to determine the power source that is currently powering the *EMS ventilator*.

Amendment (replace the last paragraph with):

- bb) An *EMS ventilator* shall be equipped with an *alarm system* that detects
- 1) a *technical alarm condition* or presents an *information signal* to indicate a switchover to an *internal electrical power source*.
    - i) The switchover to an *internal electrical power source alarm condition* shall be at least a *low priority*.
  - 2) a *technical alarm condition* to indicate when the *internal electrical power source* nears depletion, prior to the loss of all power.

- i) The *internal electrical power source* nears depletion *alarm condition* shall be at least a *medium priority*.
  - ii) As the *internal electrical power source* nears depletion, at least 5 min prior to depletion, the *internal electrical power source* nears depletion *technical alarm condition* shall escalate to *high priority*.
- 3) a *technical alarm condition* to indicate a complete loss of all power sources.
- i) The duration of the *alarm signals* for this *technical alarm condition* shall be at least 120 s.

### 212.8.3 Additional requirements for *internal electrical power source* for *ME equipment*

In the first paragraph, replace '20 min' with '2 h'.

Amendment (add before the compliance check):

- aa) The *instructions for use* shall describe behaviour of recharging of the *internal electrical power source* while the *EMS ventilator* is connected to each *supply mains*.

### 212.10.1.2 Requirements for mechanical strength for *fixed or permanently installed ME equipment* intended for use in a road ambulance

Amendment (replace the second paragraph with the following):

- aa) While performing the tests of IEC 60601-1-12:2014+AMD1:2020, 10.1.2 a), the *EMS ventilator*:
  - 1) shall be operating and maintain *basic safety* and *essential performance*; and
  - 2) may exhibit the temporary degradation of performance of *inspiratory volume* or *PEEP* for a single *inflation* for each shock impulse.
- bb) While performing the test of IEC 60601-1-12:2014+AMD1:2020, 10.1.2 b), the *EMS ventilator* shall be operating. After this test, the *EMS ventilator* shall maintain *basic safety* and *essential performance*.

### 212.10.1.3 Requirements for mechanical strength for *transportable ME equipment*

NOTE There is guidance or rationale for this subclause contained in Clause AA.2.

Amendment (replace the second paragraph with the following):

- aa) While performing the tests of IEC 60601-1-12:2014+AMD1:2020, 10.1.3 a) and IEC 60601-1-12:2014+AMD1:2020, 10.1.3 c), the *EMS ventilator*:
  - 1) shall be operating and maintain *basic safety* and *essential performance*; and
  - 2) may exhibit the temporary degradation of performance of *inspiratory volume* or *PEEP* for a single *inflation* for each shock impulse.
- bb) While performing the test of IEC 60601-1-12:2014+AMD1:2020, 10.1.3 b), the *EMS ventilator* shall be operating. After this test, the *EMS ventilator* shall maintain *basic safety* and *essential performance*.

Replace the fall height list item in c) with:

- fall height:
  - for mass  $\leq 1$  kg, 1,0 m,
  - for mass  $> 1$  kg and  $\leq 50$  kg, 0,75 m,

- for mass > 50 kg, 0,1 m,
- number of falls: 2 in each specified attitude.

**212.10.1.4 Requirements for mechanical strength for *ME equipment* intended for airborne use**

Amendment (replace the second paragraph with the following):

While performing the tests of IEC 60601-1-12:2014, 10.1.4, the *EMS ventilator* shall be operating and maintaining *basic safety* and *essential performance*.

The annexes of the IEC 60601-1:2005+AMD1:2012+AMD2:2020, apply, except as follows:

STANDARDSISO.COM : Click to view the full PDF of ISO 80601-2-84:2023

## Annex C (informative)

### Guide to *marking* and labelling requirements for *ME equipment* and *ME systems*

IEC 60601-1:2005+AMD1:2012+AMD2:2020, Annex C applies, except as follows:

*Addition:*

#### 201.C.1 *Marking on the outside of ME equipment, ME systems or their parts*

Additional requirements for *marking* on the outside of an *EMS ventilator*, its parts and *accessories* are found in Table 201.C.101.

**Table 201.C.101 — *Marking on the outside of an EMS ventilator, its parts or accessories***

Description of marking	Subclause
Arrow indicating the direction of the flow for <i>flow-direction-sensitive components</i> , if applicable	201.7.2.101 c) 3)
Arrow indicating the intended direction of gas flow for the <i>gas output port</i> and the <i>gas return port</i>	201.7.2.101 d) 2)
For a sampling <i>gas inlet</i> , either with the text "Gas sample" or the <i>Symbol</i> ISO 7000-0794	201.7.2.101 d) 3)
For <i>accessories</i> supplied separately, indication of any limitations or adverse effects of the <i>accessory</i> on the <i>basic safety</i> or <i>essential performance</i> of the <i>EMS ventilator</i> , if applicable	201.7.2.4.101 a) 2)
For <i>accessories</i> supplied separately, the <i>marking</i> requirements of ISO 20417:2021, 6.1.1 c)	201.7.2.4.101 a) 1)
For an <i>EMS ventilator</i> intended to be used in the magnetic resonance (MR) environment, <i>symbol</i> 7.3.1-1 (Table 201.D.2.101, <i>symbol</i> 4) or <i>Symbol</i> 7.3.1-2 (Table 201.D.2.101, <i>symbol</i> 5) of IEC 62570:2014 for an 'MR Safe' ventilator	201.7.2.101 c) 1) i)
For an <i>EMS ventilator</i> intended to be used in the magnetic resonance (MR) environment, <i>symbol</i> 7.3.2 (Table 201.D.2.101, <i>symbol</i> 6) of IEC 62570:2014 for an 'MR Conditional' ventilator	201.7.2.101 c) 1) ii)
For an <i>EMS ventilator</i> not intended to be used in the magnetic resonance (MR) environment, <i>symbol</i> 7.3.3 (Table 201.D.2.101, <i>symbol</i> 6) for an 'MR Unsafe' ventilator	201.7.2.101 c) 2)
For an oxygen concentration control with an error exceeding 20 %, not quantitative numerical concentration values	201.7.4.2 aa)
For oxygen gas inputs, the <i>rated</i> range of oxygen concentration	201.7.2.18 cc)
Gas name or chemical <i>symbol</i> for any gas-specific inputs	201.7.2.18 aa)
Gas volume, flowrate and leakage specifications in <i>STPD</i> except for those associated with the <i>VBS</i> in <i>BTPS</i>	201.7.4.3
Gas-specific colour coding for any gas-specific inputs, if applicable	201.7.2.18 dd)01
Mandatory action <i>safety sign</i> : follow <i>instructions for use</i>	201.7.2.3
Quantitative numeric indications of parameters and settings with the units of measurement	201.7.4.101
<i>Rated</i> range of gas pressure	201.7.2.18 bb)
Requirements of ISO 20417	201.7.1.101
Start-up <i>procedure</i> that is intended to be performed prior to use	201.7.2.101 a)
<i>Symbol</i> or <i>safety sign</i> , as appropriate, for ISO 80369-7 to the <i>gas pathway</i>	201.101.3.2.4 c)
Warning not to obstruct the <i>gas intake port</i> , if applicable	201.7.2.101 d) 1)

### 201.C.2 Accompanying documents, general

Additional requirements for general information to be included in the *accompanying documents* of an *EMS ventilator* or its parts are found in Table 201.C.102.

**Table 201.C.102 — Accompanying documents, general**

Description of requirement	Subclause
<i>EMS ventilator</i> is a high-flow device warning, if applicable	201.4.11.101.2 b) 3) iii)
For each <i>VBS</i> and <i>accessory</i> , the <i>model or type reference</i> of at least one compatible <i>EMS ventilator</i>	201.102.2 a)
For each <i>VBS</i> , part and <i>accessory</i> , a statement to the effect that ventilator breathing systems, their parts and accessories are validated for use with specific ventilators	201.102.2 b) 1)
For each <i>VBS</i> , part and <i>accessory</i> , a statement to the effect that incompatible parts can result in degraded performance	201.102.2 b) 2)
For each <i>VBS</i> , part and <i>accessory</i> , a statement to the effect that the <i>responsible organization</i> is responsible for ensuring the compatibility of the ventilator and all of the parts used to connect to the patient before use	201.102.2 b) 3)
Maximum time-weighted average input flow for each gas, if applicable	201.4.11.101.2 b) 3) i)
Maximum transient input flow for each gas, if applicable	201.4.11.101.2 b) 3) ii)
Requirements of ISO 20417	201.7.1.101

### 201.C.3 Accompanying documents, instructions for use

Additional requirements for information to be included in the *instructions for use* of an *EMS ventilator* or its parts are found in Table 201.C.103.

**Table 201.C.103 — Instructions for use**

Description of requirement	Subclause
Accuracy of expired volume <i>monitoring equipment</i> , if so equipped	201.12.4.103 b) 3)
Accuracy of the <i>inspiratory volume monitoring equipment</i> , if so equipped	201.12.1.104
Any reasonably foreseeable <i>adverse effect</i> of any recommended <i>accessory</i> on the <i>basic safety</i> or <i>essential performance</i> of the <i>EMS ventilator</i> , if applicable	201.7.9.2.14.101 b) 2)
At least one set of <i>accessories</i>	201.7.9.2.14.101 a) 1)
A-weighted sound power level emitted by the <i>EMS ventilator</i>	201.9.6.2.101 b) 2)
A-weighted sound pressure level emitted by the <i>EMS ventilator</i>	201.9.6.2.101 a) 2)
Behaviour of recharging of the <i>internal electrical power source</i> while the <i>EMS ventilator</i> is connected to each <i>supply mains</i>	212.8.3 aa)
Conditions under which the <i>EMS ventilator</i> maintains the accuracy of controlled and displayed variables	201.7.9.2.9.101 c)
Cross reference between the <i>manufacturer-specific naming</i> of the <i>EMS ventilator's ventilation-modes</i> and the <i>ventilation-mode systematic coding scheme</i> of ISO 19223	201.7.9.2.9.101 f)
Description of how to connect and test the connection of a <i>distributed alarm system</i> , if provided	201.7.9.2.9.101 e)
Description of the algorithm that determines the <i>alarm limit</i> values of expired volume <i>monitoring equipment</i> , if so equipped	201.12.4.103 e)

Description of requirement	Subclause
Disclosed accuracies of expired volume <i>monitoring equipment</i> include the effect of <i>rated</i> input oxygen concentration, if so equipped	201.12.4.103 b) 4)
Disclosure of any restrictions on the placing of components within the <i>ventilator breathing system</i> , if applicable	201.7.9.2.14.101 b) 1)
FiO <sub>2</sub> accuracy includes effects of <i>rated</i> input oxygen concentration	201.12.1.101 b) 3) i) 201.12.1.102 b) 3) i)
For <i>accessories</i> supplied separately where <i>marking</i> the <i>accessory</i> is not practicable, the requirements of 201.7.2.4.101	201.7.2.4.101 b)
Information on how to connect CO <sub>2</sub> <i>monitoring equipment</i> , unless such equipment is an integral part of the <i>EMS ventilator</i>	201.12.4.106 c)
Information on how to connect O <sub>2</sub> <i>monitoring equipment</i> , unless such equipment is an integral part of the <i>EMS ventilator</i>	201.12.4.101 c)
Information on how to connect the expired volume <i>monitoring equipment</i> , if not so equipped	201.12.4.103 b) 1)
Intended position of the <i>operator</i>	201.7.9.2.1
Intended range of <i>inspiratory volume</i>	201.7.9.2.1.101 a)
Maximum duration of a <i>system recovery</i> of <i>ventilation</i>	201.4.3.102 g)
Maximum error of the <i>airway pressure</i> at the end of the <i>inflation phase</i> in relation to the set value for a <i>pressure-control inflation</i> in <i>normal condition</i>	201.12.1.102 b) 1)
Maximum error of the inspiratory oxygen concentration (FiO <sub>2</sub> ) at the <i>patient-connection port</i> in relation to the set value for a <i>volume-control inflation</i> in <i>normal condition</i> , if provided	201.12.1.101 b) 3)
Maximum error of the inspiratory oxygen concentration (FiO <sub>2</sub> ) at the <i>patient-connection port</i> in relation to the set value for a <i>pressure-control inflation-type</i> in <i>normal condition</i> , if provided	201.12.1.102 b) 3)
Maximum error of the <i>inspiratory volume</i> in relation to the set value for a <i>volume-control inflation</i> in <i>normal condition</i>	201.12.1.101 b) 1)
Maximum error of the <i>PEEP</i> in relation to the set <i>BAP</i> for a <i>pressure-control inflation-type</i> in <i>normal condition</i>	201.12.1.102 b) 2)
Maximum error of the <i>PEEP</i> in relation to the set <i>BAP</i> for a <i>volume-control inflation</i> in <i>normal condition</i>	201.12.1.101 b) 2)
<i>ME equipment</i> necessary for the <i>EMS ventilator's intended use</i> , if applicable	201.7.9.2.14.101 a) 2)
Means by which the <i>maximum working pressure</i> is ensured	201.7.9.2.9.101 a) 3)
Method by which all of the <i>alarm signals</i> can be functionally tested to determine if they are operating correctly	201.7.9.2.8.101 b)
Performance and acceptance criteria for other <i>inflation-types</i> , including effects of <i>rated</i> oxygen concentration, if so equipped	201.12.1.103
Pre-use checklist	201.7.9.2.8.101 a)
<i>Processing</i> instructions for <i>gas pathways</i> of the <i>EMS ventilator</i> and its <i>accessories</i> that can become contaminated	201.11.6.6 bb) 3)
<i>Processing</i> instructions for the <i>enclosure</i> of the <i>EMS ventilator</i> and its <i>accessories</i>	201.11.6.6 dd) 2)
Range of <i>maximum limited pressure</i>	201.7.9.2.9.101 a) 1)
<i>Rated</i> range of expiratory <i>gas pathway</i> resistance over which the accuracies of set and monitored volumes and pressures are maintained	201.7.9.2.9.101 b) 2)
<i>Rated</i> range of inspiratory <i>gas pathway</i> resistance over which the accuracies of set and monitored volumes and pressures are maintained	201.7.9.2.9.101 b) 1)
<i>Rated</i> range of <i>VBS</i> compliance over which the accuracies of set and monitored volumes	201.7.9.2.9.101 b) 3)

Description of requirement	Subclause
and pressures are maintained	
<i>Rated</i> range to which the <i>maximum working pressure</i> can be set, if adjustable	201.7.9.2.9.101 a) 2)
Requirements of ISO 20417	201.7.1.101
Statement to the effect that the ventilator is to be equipped with O <sub>2</sub> monitoring equipment for the measurement of inspiratory oxygen concentration before being put into service, if not so equipped	201.12.4.101 a) 2)
Statement to the effect that the ventilator is to be equipped with monitoring equipment for indicating expired volume at the patient-connection port before being put into service, if not so equipped	201.12.4.103 a) 2)
Statement to the effect that the ventilator is to be equipped with CO <sub>2</sub> monitoring equipment for the measurement of expiratory carbon dioxide concentration before being put into service, if not so equipped	201.12.4.106 a) 2)
Summary of effects of respiratory rate on end-tidal gas reading accuracy, if provided	201.7.9.2.9.101 d)
<i>Ventilation</i> accuracy expressed as bias error and linearity error	201.12.1.101 a) 201.12.1.102 a)
Warning regarding ISO 80369-7 connections to the <i>gas pathway</i> , if equipped	201.101.3.2.4 c) 2)
Warning statement to the effect that do not add any attachments or accessories to the ventilator that are not listed as intended for use in combination with the ventilator in the instructions for use of the ventilator or accessory, as the ventilator might not function correctly leading to the risk of patient death or additional serious deterioration of health	201.7.9.2.2.101 a) 4)
Warning statement to the effect that do not cover the ventilator or place in a position that affects the proper operation, including applicable examples	201.7.9.2.2.101 a) 2)
Warning statement to the effect that do not use the ventilator in explosive environments. Such use might cause an explosion	201.7.9.2.2.101 b) 4)
Warning statement to the effect that the ventilator accuracy can be affected by the gas added by use of a nebuliser, if applicable	201.7.9.2.2.101 b) 3)
Warning statement to the effect that the ventilator is intended to be continuously attended by an operator. Failure to be in close proximity to this ventilator can contribute to patient death or serious injury	201.7.9.2.2.101 a) 1)
Warning statement to the effect that the <i>ventilator</i> shall not be used in a hyperbaric chamber, if applicable	201.7.9.2.2.101 b) 2)
Warning statement to the effect that when using nebulisation or humidification, the breathing system filters can require more frequent replacement to prevent increased resistance and blockage, if applicable	201.7.9.2.2.101 b) 1)
Warning statement to the effect that, always have immediate access to alternative means of ventilation, which is ready for use, in order to reduce the possibility of patient death or additional deterioration of health	201.7.9.2.2.101 a) 3)
Whether or not the <i>EMS ventilator</i> is intended for use with closed suctioning	201.7.9.2.9.101 g)
Which portions of the <i>gas pathways</i> through the <i>EMS ventilator</i> can become contaminated with body fluids or by contaminants carried by expired gases during both <i>normal condition</i> and <i>single fault condition</i>	201.7.9.2.12 aa)

#### 201.C.4 Accompanying documents, technical description

Additional requirements for information to be included in the *technical description* of an *EMS ventilator* or its parts are found in Table 201.C.104.

Table 201.C.104 — *Technical description*

Description of requirement	Subclause
Description of a method for checking the function of the <i>alarm system</i> for each of the <i>alarm conditions</i> specified in this document	201.7.9.3.101 a)
Description of the method used to calculate end-tidal gas readings, if applicable	201.7.9.3.1.101 b) 2)
Description of what happens to the contents of the log after the <i>alarm system</i> has experienced an unexpected total loss of power ( <i>supply mains</i> and <i>internal electrical power source</i> ) for a finite duration	208.6.12.3 h)
Description of what happens to the contents of the <i>responsible organization alarm system</i> log as it reaches capacity	208.6.12.3 i) 2)
Description of which of the <i>alarm system</i> checks are performed automatically	201.7.9.3.101 b)
Disclosure of the uncertainty for each disclosed tolerance	201.5.101.3 c)
Essential technical characteristics of each recommended <i>breathing system filter</i> , if applicable	201.7.9.3.1.101 b) 1)
Outside dimensions and mass of the <i>EMS ventilator</i>	201.7.9.3.1.101 a) 4)
Pneumatic diagram of the <i>EMS ventilator</i> , including a diagram for <i>operator-detachable</i> parts of the <i>EMS ventilator breathing system</i> either supplied or recommended in the <i>instructions for use</i>	201.7.9.3.1.101 a) 2)
Requirements of ISO 20417	201.7.1.101
<i>Responsible organization alarm system</i> log capacity	208.6.12.3 i) 1)
Summary description of the filtering or smoothing techniques for all measured or computed variables that are displayed or used for control	201.7.9.3.1.101 a) 1)
Summary description of the means of initiating and terminating the <i>inflation phase</i> in each <i>ventilation-mode</i> of the <i>EMS ventilator</i>	201.7.9.3.1.101 a) 3)
Technical details of the means for securing the <i>EMS ventilator</i> to avoid unwanted lateral movement	201.9.4.3.101 f)


## Annex D (informative)

### *Symbols on marking*

IEC 60601-1:2005+AMD1:2012+AMD2:2020, Annex D applies, except as follows:

*Addition:*

**Table 201.D.2.101 — Additional *symbols on marking***

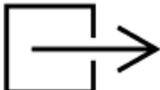No	Symbol	Reference	Title and description
1		IEC 60878:2022 ISO 7000-0794	Input; entrance To identify an entrance, for example exhaust gas entry for measurement (for example of CO <sub>2</sub> value). For electrical (signal) input use <i>symbol</i> IEC 60417-5034.
2		IEC 60878:2022 ISO 7000-0795	Output; exit To identify an exit, for example of a hydraulic pump. For electrical (signal) output use <i>symbol</i> IEC 60417-5035.
3		IEC 60878:2022 <i>Symbol</i> 7.3.1-1 IEC 62570:2014	MR Safe To identify an item which poses no unacceptable <i>risks</i> to the <i>patient</i> , medical staff or other persons within the MR environment. When colour reproduction is not practical, the <i>symbol</i> may be printed in black and white. The use of the coloured icon is strongly encouraged for the added visibility and information provided by the colour.
4		IEC 60878:2022 <i>Symbol</i> 7.3.1-2 IEC 62570:2014	MR Safe Alternative graphical <i>symbol</i> representation. Same meaning as IEC 62570-7.3.1-1.
5		IEC 60878:2022 <i>Symbol</i> 7.3.2 IEC 62570:2014	MR Conditional To identify an item which poses no unacceptable <i>risks</i> within defined conditions to the <i>patient</i> , medical staff or other persons within the MR environment. When color reproduction is not practical, the <i>symbol</i> may be printed in black and white. The use of the colored icon is strongly encouraged for the added visibility and information provided by the color. The MR Conditional <i>symbol</i> may be supplemented by supplementary <i>marking</i> that describes the conditions for which the item has been demonstrated to be MR Conditional.

No	Symbol	Reference	Title and description
6		IEC 60878:2022 Symbol 7.3.3 of IEC 62570:2014	MR Unsafe To identify an item which poses unacceptable <i>risks</i> to the <i>patient</i> , medical staff or other persons within the MR environment. Note - When color reproduction is not practical, the <i>symbol</i> may be printed in black and white. The use of the colored version is strongly encouraged for the added visibility and information provided by the color.

Additional Annexes:

STANDARDSISO.COM : Click to view the full PDF of ISO 80601-2-84:2023

## Annex AA (informative)

### Particular guidance and rationale

#### AA.1 General guidance

This Annex provides rationale for the important requirements of this document and is intended for those who are familiar with the subject of this document but who have not participated in its development. An understanding of the reasons for the main requirements is considered to be essential for its proper application. Furthermore, as clinical practice and technology change, it is believed that rationale for the present requirements will facilitate any revision of this document necessitated by those developments.

#### AA.2 Rationale for particular clauses and subclauses

The numbering of the following rationales corresponds to the numbering of the clauses and subclauses in this document. The numbering is, therefore, not consecutive.

##### — 201.1.1 — Scope

ISO 80601-2-84 *ventilators* are not considered a *physiologic closed-loop control system* due to the fact that parameters monitored during delivery of respiratory gases that are also used to control the delivery of these gases are exclusively physical parameters of the delivered gases. Consequently, these parameters are considered equipment variables as specified in IEC 60601-1-10.

A *ventilator* in a *pressure-control inflation-type* that uses the *breathing system* pressure as a feedback to control *breathing system* pressure is a *closed-loop control system*, but not a *physiologic closed-loop control system*. The *breathing system* pressure is considered both a 'variable' influenced by the *patient* physical conditions and at the same time a 'feedback variable', but it is not a quantity or condition measured from the *patient's* physiology.

The *patient's* physical condition is a disturbance on the closed loop system but the *ventilator* does not adjust the *ventilation* therapy settings based on measurement of these *patient* parameters.

The requirements of this document do not require the *ventilator* to adjust *ventilation* delivery parameters based on the detection in the change of physiological conditions of the *patient*. All automatic adjustments of *ventilator* equipment parameters or generated *alarm conditions* are only based on the measurement of physical variables related to the delivery of breathing gas to the *patient-connection port*. In this sense the *ventilator* ends at the *patient-connection port*, (i.e. has no direct contact to the physiological parameters of the *patient*) and a change in the *patient's* physiological conditions is a disturbance to the *ventilator's* control system that does not act to control the physiological change but continues to control the physical variable(s) to their original objectives.

*Ventilators* create *alarm conditions* when detecting faults in the delivery of breathing gases to the *patient-connection port* but do not adjust therapy setting of the *ventilator*.

The following are examples of medical devices that are considered *physiological closed-loop control system*.

- Insulin infusion pump that adjust the rate of insulin infusion to the *patient* based on the measurement of blood glucose. The physiological feedback mechanism is a blood glucose level monitored by the device.
- External pacemaker that adjusts the pace rate based on the measurement of the cardiac output value. The physiological feedback mechanism is the value detected by the cardiac output monitor.

Unlike a *ventilator*, these devices titrate therapy provided to the *patient* based on the measured physiological parameter. A *ventilator* will not titrate but will either stop *ventilation* or generate an *alarm condition*.

#### — 201.4.3.101 — Additional requirements for *essential performance*

The modern *EMS ventilator* with an active exhalation valve has differing *ventilation-modes* that can consist of multiple *inflation-types*. This is necessary as *patient* response to *ventilation* is unpredictable. *Patient*-initiated breaths or *inflations* where the inspiration is terminated by the *patient* can have characteristics that are different from those that have been set by the *operator*. *Essential performance* as “*ventilation* within the *alarm limits* set by the *operator*” is inclusive of those *inflations* that the *patient* modifies outside of the ventilatory parameters set by the *operator*, but still within the *alarm limits* which are considered safe by the *operator*.

It is expected that the *operator* will set appropriate *alarm limits* which thereby define the *essential performance* for a particular *patient*.

The distributed *essential performance* criteria captured in Table 201.101 have been identified by the committees as the minimum clinical performance necessary to reduce the possibility of exposing the *patient* to unacceptable risk. Conformity criteria for some of the clauses within IEC 60601-1, this document and the other applicable collateral standards include “maintain *essential performance*”. The committees have recognized the difficulty in confirming that all aspects of *essential performance* are maintained when completing longer duration testing.

Footnote a to Table 201.101 indicates methods of evaluating delivery of *ventilation* as acceptance criteria following specific tests required by this document. It is intended to provide criteria that can be used to easily verify that *essential performance* has been maintained. Although the degradation of performance detailed in 202.8.1.101 are associated with *immunity* testing, the same criteria are intended to be used when the conformity criteria from any other clause or subclause require confirmation that *essential performance* is or has been maintained.

Those aspects of *essential performance* that cannot be reasonably linked to the conformity criteria within 202.8.1.101 need to be confirmed via other means. However one need only confirm that the specific requirements indicated in 202.8.1.101 that are likely to have an impact on specific clinical performance are maintained after testing.

#### — 201.4.3.102 *System recovery*

IEC 62304:2006+A1:2015, 5.5.4 d), identifies fault handling as potentially including recovery. In previous standards applicable to *ventilators*, the management of *risk* associated with *single fault conditions* is limited to specifying that a *single fault condition* not causing the simultaneous failure of both a *ventilation-control* function, and the corresponding *protection device*. It is commonplace that a *ventilator* experiencing a *single fault condition* shuts down while providing *alarm signals*. While this is appropriate for some fault modes, for example when a critical hardware component has failed, there are other fault modes where it can be possible for a *system recovery process* to allow *ventilation* to continue.

Many ventilators use one or more microprocessors to control ventilation and to provide monitoring and alarm system functions. Software controlled equipment is inherently susceptible to *single fault conditions* that can be transient in nature. For example, a processing element might stop responding correctly as a result of a software defect, such as a buffer over-run or domain error; memory corruption such as can occur due to ionising radiation; or in input error such as from electromagnetic interference on a *functional connection*.

The risk of harm to the patient is reduced if, in addition to the ventilator providing alarm signals to indicate the loss of ventilation due to a *single fault condition*, it can reset and resume ventilation. This might be implemented by having a watchdog processor that can detect the failure of a ventilation-control software process and restart the ventilation-control processor into a known state.

In some cases, this procedure can also be applicable to *single-fault conditions* that result in non-functional hardware components. An example would be a blower motor that stalls as a result of mechanical shock during operation. If the ventilator control system is able to restart the blower, and hence resume ventilation, this reduces the risk of harm to the patient further below that level that would pertain if the ventilator were to remain inoperable, albeit with alarm signals active.

The manufacturer determines the appropriate system response to a *system recovery process* that fails to restore full functioning. This could occur either because the *single fault condition* is permanent in nature, or because of an ongoing but essentially temporary phenomenon such as continued electromagnetic interference. It is confusing to the operator for the system to continually repeat a *system recovery process* that has failed. In this instance it can be more appropriate to enter a *fail-safe state* with an informative alarm signal.

— **201.4.6 — ME equipment or ME system parts that contact the patient**

Since much of the VBS is likely to be draped over or around the patient, it is likely to come into direct contact with the patient during normal use. Also of concern are electrical hazards should any circuitry be incorporated into the VBS. By ensuring that the gas pathways of the VBS and its parts or accessories are subject to the requirements for applied parts, these issues are addressed by the requirements already in the general standard.

— **201.4.11.101 — Additional requirements for pressurized gas input**

An EMS ventilator designed to be connected to a pressurised gas supply is required to continue to operate reliably throughout its rated range of supply pressures, but these pressures can only be maintained if the EMS ventilator in normal condition does not attempt to draw more flow from the gas source than the gas source is designed to supply. It is also expected that an EMS ventilator should be designed to prevent an unacceptable risk under possible *single fault conditions* of the pressurised gas supply.

Pressurised medical gas supplies, including medical gas pipeline systems and cylinder pressure regulators conforming to current relevant standards, supply medical gases to gas-specific terminal outlets at a pressure that is within an internationally agreed-upon pressure range of 280 kPa to 600 kPa under normal condition. It is expected that EMS ventilators should operate to their declared specification with any supply pressure within this range.

In the case of a pressure regulator failure the gas supply pressure could rise to the pressure regulator's supply pressure — which can be cylinder (tank) pressure. To safeguard against this or similar eventualities, gas-specific medical gas supply systems are required to be provided with a means to limit their output pressure to not more than 1 000 kPa. All gas-powered ME equipment should be designed so as not to present an unacceptable risk if its supply pressure rises up to this value.

*EMS ventilators* with maximum *rated* input pressures exceeding 600 kPa are required to fulfil these conditions at up to twice their maximum *rated* input pressure.

Under the *single fault condition* that the supply pressure of any one gas drops below 280 kPa under steady-state conditions, it is understood that an *EMS ventilator* cannot be expected to continue to operate on this gas. However, it is required that in this case the *EMS ventilator* should detect the unacceptable low pressure, produce an *alarm signal* and also, in the case of two pressurised gas supplies, automatically switch to use the other gas source (oxygen or air) to drive the *EMS ventilator*. This requirement is stated in subclause 201.13.12.

To ensure that the minimum pressure of 280 kPa can be maintained in practice, *medical gas pipeline systems* supplying compressed medical gases through gas-specific terminal *outlets* are designed so that they can maintain this pressure at the input of gas-powered devices whilst supplying steady-state flows up to 60 l/min at a single *outlet* connected directly to the pipeline. Account is taken of the pressure drop in the pipeline supplying the *outlet* and the pressure drop, at 60 l/min, across the terminal unit and the hose assembly connecting the device to the pipeline.

The *medical gas pipeline system* is also required to be capable of supplying sufficient gas that this flow can be drawn from a predetermined number of adjacent terminal units simultaneously. The actual number will have been determined during the design and installation of the *medical gas pipeline system* by the application of a 'diversity factor'; a factor agreed between the supplier and *responsible organization* to be appropriate for each section of the installation according to the designated purpose of each area supplied. Recommended diversity factors are formulated to ensure that the *medical gas pipeline system* is capable of supplying an average flow of 60 l/min to the required proportion of terminal *outlets*. However, if the flow demand from many adjacent *ventilators* exceeds 60 l/min there is an increased possibility that the *ventilator* input pressure could fall below 280 kPa, mainly because of the increased pressure drop across the terminal unit and input hose (also because of the flow-drop characteristic in the case of pressure regulators supplying a single terminal *outlet*).

In addition to steady-state flows of 60 l/min, the switching of the internal pneumatic system and the operation of a *patient-demand* system can result in a *ventilator* requiring transient input flows far in excess of 60 l/min. Because of the compressibility of gas at pipeline pressures and the diameter of piping that is employed in order to minimise pressure drop, such transient demands can generally be accommodated from the gas stored locally within the pipe work of the *medical gas pipeline system*. There can be temporary pressure drops of the input pressure at the *inlet* of the *ventilator* to below 280 kPa, due to transient flows in excess of 200 l/min (over 3 s), but most of these drops will be within the supply hose assemblies specified by the *manufacturer*. *Manufacturers* need to evaluate their own designs to establish whether any consequent transient pressure drop affects the performance of their *ventilator* when used with recommended supply hose configurations and when connected to alternative gas-specific terminal *outlets* such as those fitted to cylinder pressure regulators conforming to ISO 10524-1.

*Ventilators* that can draw greater average or transient flows during *intended use* are permitted, but their *accompanying documents* are required to disclose those flows and warn of the need for a different diversity factor.

The average flow of 60 l/min is greater than the test flow used during the commissioning of *medical gas pipeline systems*. In itself, this should be of no concern because the specific conditions specified for the test do not allow a direct comparison between the two values. The committee responsible for *medical gas pipeline systems* standards, ISO/TC 121/SC 6, in consultation with ISO/TC 121/SC 1 and ISO/TC 121/SC 3, agreed to the 60 l/min average flow value, and also the 200 l/min for up to 3 s transient flows, during the preparation of the first edition of the current series of standards for *medical gas pipeline systems* and were aware of the need to satisfy that specification when finalizing the *medical gas pipeline system* test requirements.

*Manufacturers* should be aware that other standards for supply systems for medical gases permit the fitting of gas-specific terminal *outlets* to spur systems such as pendant supply units. Such subsystems restrict the flow that can be drawn from their terminal *outlets*.

— **201.5.101 — Additional requirements for general requirements for testing of ME equipment**

After due consideration, the committees decided that where this document specifies adjoining ranges for variables as the basis for testing and the declaration of performance, the end value of both ranges should be applicable to both ranges. This means that a *manufacturer* is free to use a round number end value (e.g. 300 ml) in specifications and is not forced to truncate artificially the declared range in order to avoid having to satisfy also the test requirements of the adjacent range. This permits, for example, one *EMS ventilator* to have a declared range *inspiratory volume* of 300 ml to 1000 ml and another 100 ml to 300 ml, with each *EMS ventilator* only being required to be tested for the conditions specified for  $\geq 300$  ml or  $\leq 300$  ml respectively.

— **201.5.101.1 — EMS ventilator test conditions**

— **a) 2)**

Test laboratories do not routinely have access to compressed air and oxygen supplies that are certified to medical pipeline standards. In many jurisdictions, medical grade oxygen is a drug and access is controlled. Enforcing a requirement for type testing to use gases that are approved to medical gas standards would add significant cost, and potentially prevent some test laboratories from being able to perform the *type tests* of this document.

For purposes of testing to this document, there is no significant difference between industrial compressed gases and medical compressed gases.

In practice, the oxygen concentration in industrial grade oxygen is no less than that in medical grade oxygen, although the permitted contaminant profile is different. Similarly, industrial compressed air is typically provided from a compressor that entrains ambient air, with principal components (nitrogen, oxygen, argon and CO<sub>2</sub>) at concentrations that align with the requirements for medical air.

A *ventilator* controls, and measures, flows of gas from air and oxygen *inlets*. The control and measurement technologies are potentially susceptible to gas composition, however the critical parameters for flow measurement and control are density, viscosity, thermal capacity and oxygen concentration. The differences in trace contaminant gases between industrial and medical gases supplies do not result in significant differences in flow measurement and control, and therefore will not change the results of *type test*.

The only likely contaminants in industrial compressed gases that could impact test results would be oil or water in droplet form, as these could potentially damage the *ventilator* parts. Industrial compressed air is usually provided from a compressor that entrains ambient air, which will include water vapour. Dryers and water traps are required to reduce the humidity to prevent condensation in the compressed gas.

Provided the dew point of the industrial gas supplies is below the temperatures of equipment parts during testing, condensation should not occur, and the results of *type tests* remain valid.

— **201.5.101.2 — Gas flowrate and leakage specifications**

Quantities of gas are frequently expressed as the volume that the gas occupies at standardized conditions. Generally, one atmosphere (101,3 kPa) is used as standard pressure. However, several

standard temperatures are used. Whereas 0 °C is used as standard temperature in physics, either 20 °C or 21,2 °C (70 °F) is often used in engineering. In *ventilation*, the gas in the *lungs* has a temperature identical to body temperature (~ 37 °C) irrespective of the temperature of the gas delivered by an *EMS ventilator*. The volume of a given amount of gas increases by about 13,5 % from 0 °C to 37 °C or by 5,8 % from 20 °C to 37 °C.

Gas delivery systems supplying pressurised gas to medical equipment, including *EMS ventilators*, follow engineering conventions and specify gas quantities and flow rates at *STPD* conditions. This practice is followed in this document for all requirements concerning gas input.

However, *EMS ventilators* conforming with this document are likely to be inflating the *patient's lungs* relative to a local atmospheric pressure between 70 kPa and 110 kPa. In addition, the gas in the *lungs* is always saturated with water vapour regardless of the humidity of the gas delivered from an *EMS ventilator*. With a standard temperature of 0 °C, 1 l of gas referenced to *STPD* can expand the *lungs* by 1,8 l at a pressure of 70 kPa. In order to have the values comparable among different *EMS ventilators*, it is essential that the information for all *EMS ventilators* is referenced to the same standard conditions. Because it is the volume of gas and not the number of molecules that expands the *lungs*, *BTPS* is the appropriate set of reference conditions to use.

In *EMS ventilators* a variety of flow transducers are used. Whereas a heated-wire anemometer measures the rate of mass flow of the gas independent of pressure, a pneumotachograph measures the flow of gas at the actual pressure. Therefore the necessary corrections depend on the type of flow transducer. When a pressure correction is required, this can be adequately estimated.

The necessary corrections also depend on the location of the flow transducer in the *VBS*. The humidity of the gas can be zero when the transducer measures the inspired flow inside the *EMS ventilator*. When, however, the flow transducer is located at the Y-piece, the relative humidity can be anything up to 100 %. When a *heat and moisture exchanger* is used for humidification, the output of the flow transducer depends on whether it is located distal or proximal to the *heat and moisture exchanger*. With a blower- or compressor-based *EMS ventilator* that uses ambient air, the humidity of the drawn-in air can be unknown to the *EMS ventilator*. All these effects together will inevitably introduce some errors in the conversion of the measured flow signal to *BTPS* reference conditions. However, these errors are only in the range of several percent. For *inspiratory volumes* greater than 100 ml, the committees came to the conclusion that the permitted inaccuracy of the measurement of  $\pm(4 \text{ ml} + (20 \% \text{ of the measured value})) \text{ ml}$  is sufficiently wide and includes the inaccuracy of the flow transducer and the inaccuracy of the conversion to *BTPS* conditions. However, it remains the responsibility of the *manufacturer* to verify that the accuracy requirements of 201.12.1 and 201.12.4.103 are met.

### — 201.5.101.3 — *EMS ventilator testing errors*

When testing *EMS ventilator* performance several of the test parameters cannot be measured without a significant degree of measurement uncertainty due to limitations of the accuracy that can be achieved, particularly when measuring volumes by the integration of rapidly changing flows.

Because of the relative significance of these uncertainties, it is important that *manufacturers* allow for them when declaring parameter accuracy.

Similarly, it is important for a third-party tester to recognise the significance of the uncertainty in their own measurements when testing to this document.

In practice, this means that, for example, if a *manufacturer* determines that a parameter has an intended tolerance of  $\pm 10 \%$ , but the measurement uncertainty is  $\pm 3 \%$  then the test results are acceptable if, given the uncertainty band for the measured value, the probability of the measured values being within the limit is at least 50 %. In almost all cases, measurement uncertainty has a symmetrical distribution, and the 50 % likelihood criterion is met if the measured value is within

the disclosed limit, in this example, within  $\pm 10\%$  of the setting. If a third-party is testing to this document, they also need to include measurement uncertainty in their testing. The third-party testing organization needs to control measurement uncertainty to the same level as that disclosed for a type test, which in this example is  $\pm 3\%$ .

Note that a third-party testing organization obtaining a measured value outside the limit does not necessarily invalidate the claim – the deviation from the limit is required to be compared to the uncertainty of the measurement to establish the probability of the data representing a true deviation from specification.

See IEC Guide 115 for more information regarding measurement uncertainty.

Furthermore, the *manufacturer* is required to disclose the measurement uncertainty for each declared value in order to provide this information to the *responsible organization* and guidance for a third-party tester as to the needed measurement accuracy when testing to this document.

This is a change from previous revisions of this document, which required tolerances to be adjusted by subtracting measurement uncertainty from disclosed tolerance values to determine acceptance criteria.

#### — 201.7.2.3 — Consult *accompanying documents*

The committees agreed that following the *instructions for use* is a mandatory action for the safe operation of an *EMS ventilator*.

#### — 201.7.2.101 — Additional requirements for *marking on the outside of ME equipment or ME equipment parts*

Although the intended *operator group* for an *EMS ventilator* is primarily paramedical service first responders and emergency care physicians, who are routinely trained in the use of their equipment, their use might not be on a frequent basis. The committees consider that it is important that the *operator* has immediate access to those instructions necessary to place the *EMS ventilator* into use.

This can be in the form of written or illustrated *instructions for use*:

- attached to the *EMS ventilator*;
- in the *instructions for use*;
- provided via a graphical *operator interface* on the *EMS ventilator*; or
- by other means.

One possible example (although not currently implemented on any *EMS ventilator* known to the committees) would be to provide oral prompts in the same fashion as provided by public access AEDs (automated external defibrillators). These use a combination of oral prompts and illuminated controls to guide an *operator* through a life-saving medical *procedure*.

The requirement for legibility in IEC 60601-1:2005+AMD1:2012+AMD2:2020, 7.2 would be deemed to be met if the instructions can be interpreted correctly by the *operator*, whether provided in visual or other format.

#### — 201.7.4.2 — Control devices

See also rationale for 201.15.102.

Many *EMS ventilators* are designed to operate from one single high-pressure oxygen source, and use a venturi injector to entrain air in order to provide one or more 'lower oxygen' settings, in addition to providing *ventilation* with pure oxygen. This design is often used in *EMS ventilators* intended for first responder applications.

*EMS ventilators* according to this design typically do not control oxygen concentration to a specified tolerance.

Because of the physical characteristic of the venturi principle, the oxygen concentration depends significantly on the inspiratory flowrate and the mean *airway pressure*.

*Marking* an *EMS ventilator* using the venturi principle with a quantitative concentration value (such as  $\text{FiO}_2 = 0,40$ , or 40 %) would be potentially misleading owing to the wide tolerance that would apply.

#### — 201.7.4.3 — Units of measurement

Additional information is found in rationale for 201.5.101.2.

#### — 201.7.9.2.2.101 — Additional requirements for warnings and safety notices

##### a) 3)

In case of a severe *ventilator* failure, the *ventilator* stops ventilating the *patient*. There is no back-up mechanism for an alternative means of *ventilation* integrated in the *ventilator*; a *ventilator* is not required to be functional under *single fault condition*. To keep the time for a possible interruption of the *ventilation* of the *patient* as short as possible, an alternative means of *ventilation* has always to be available close to the *ventilator*.

##### a) 4)

Additional attachments or other components or subassemblies increase the resistance and compliance of the *VBS*. The additional flow of a pneumatic nebuliser or the sampling flow of a diverting gas monitor could have a negative impact on the accuracy of the measurement of the flowrates and oxygen concentration.

##### b) 1)

The functionality of *breathing system filters* is affected by a number of aspects of structure, properties and local environment.

At the most basic, a *BSF* is designed to be a filter that removes particles suspended in gas, i.e. a "dry aerosol". The particles primarily targeted in the *VBS* are bacteria or virus particles (although other particles would be subject to retention). The filtering material ("medium") is composed of a matrix of solid material with open passageways to allow gas flow. The passageways in such gas filters are relatively large compared to the bacteria and virus particles that are to be removed. The spatial arrangement of the solid part of the filter medium versus the open spaces in the medium brings the particles in proximity to the surfaces of the medium, where physical forces (electrostatic attraction and Van der Waals forces) attract and bind the particles within the matrix, removing them from the gas flow.

In the practical situation of anaesthesia or respiratory care therapy, environmental factors related to the *patient*, or the therapy can alter the performance of the *BSF* from that which would occur in the simple flow of air with suspended microorganisms through the *BSF*.

One major factor is the presence, phase and amount of moisture present in the airflow.

When there is low humidity in the air (gaseous phase moisture) the gaseous water molecules will generally pass through the filter medium without effect. If there is a sufficiently high relative humidity, some *BSFs* can adsorb or absorb part of this humidity.

If the moisture exists as a liquid aerosol, the water droplets can also be retained by the filter.

The properties of a filter medium that govern the degree to which this interaction with water takes place is its relative affinity for water. A medium which readily attracts water is termed “hydrophilic” and a medium which repels water is termed “hydrophobic”. These properties are, in fact, not discrete, but exist on a continuous scale. Nevertheless, in common parlance filters are grouped into being (relatively) hydrophilic or hydrophobic.

Another example of liquid phase water can be termed “bulk water”. An example of this is the collected condensate that occurs in the expiratory limb of the *VBS*. Depending on the management of the circuit, and the positioning of the *BSF*, this bulk water can actually completely cover and occlude the filter. If a sufficient pressure is applied, the liquid water can be forced through the pores of the filter medium. This requires relatively low pressure for a hydrophilic filter and relatively high pressure for a hydrophobic filter.

The practical consequences of the latter scenario is that if liquid is forced through a hydrophilic *BSF*, gas flow blockage can be relieved, but any microorganisms removed by the filter can be carried past the filter with the liquid stream. In the case of a hydrophobic filter, the pressure in the *VBS* is usually not sufficient to force liquid through the medium, so the microbial retention is not compromised. Airflow occlusion persists, however, until steps are taken to remove the bulk water.

In addition, there can be a temporal aspect to the properties of relative hydrophilicity or hydrophobicity whereby prolonged exposure to water alters these properties during the *expected service life* of the *BSF*. A *BSF* is typically labelled with an *expected service life*, in hours or days, that reflects its ability to perform to its labelled specifications in the clinical environment.

It should be obvious that the potential influence of water on performance differs in anaesthesia and respiratory care applications, although many, if not most, *BSFs* are indicated for use in both applications.

Additional effects on *BSF* functionality can be caused by the introduction of substances other than water or gas into the device. Such substances can originate from the *patient* (e.g. sputum, exudates, blood, vomitus) or substances introduced by the *operator* into the *VBS* (e.g. gross amounts of medications intended to be nebulised for administration through the *VBS*).

The effect of such substances can be an increase in flow resistance of varying degree up to complete occlusion at *ventilator* or physiologic pressures. In the case of nebulised medications, the type of nebuliser, and its operating parameters are variables that affect the likelihood or magnitude of significantly increased *BSF* flow resistance during a prescribed medication regimen. It should be mentioned that accidental introduction of gross amounts of medication from the nebuliser reservoir during *operator* or *patient* manipulation of the *VBS* has been implicated as a source of acute *BSF* blockage.

The cause of increased flow resistance in a *BSF* can be gross blockage of the medium passages, or the effects of surfactant properties of the substances introduced into the *BSF* upon the hydrophobicity of the filter medium. It should be noted that medications indicated for nebulisation can contain surfactant materials that are not identified in the medications’ labelling with respect to their presence or their quantity, and these can change without notice for a given medication. The effect of these substances upon flow resistance differs among individual models and brands of *BSFs*.

The *operator* needs to be aware that these effects can be manifested as increases in the amount of positive *airway pressure* required for an *EMS ventilator*-provided breath, or as an increase in

expiratory flow resistance. These effects can result in a step-wise increase in intrapulmonary pressure that, if not detected, can lead to pneumothorax.

Awareness of the possibility, albeit infrequent or rare, of such significant increases in *BSF* flow resistance, and inclusion in a trouble-shooting scheme for this and other causes of impaired *ventilation* can reduce or eliminate adverse events occurring secondary to *BSF* flow occlusion.

Direct *patient* monitoring, and usage of the appropriate settings for, and prompt attention to, *EMS ventilator alarm conditions* are essential to provide maximum *patient* safety.

Once a *BSF* is recognized to be a source of impaired *ventilation*, simply removing the occluded *BSF* and replacing it with another *BSF* returns *ventilation* to a normal state.

— **201.7.9.2.8.101 — Additional requirements for start-up procedure**

In some designs, adequate checking of the *alarm system* can be performed with a combination of *operator* action and the power-on self-test routines that verify the integrity of the software and the computer controlling the *EMS ventilator*, as well the measuring sensors and the *alarm signal* generation.

— **201.7.9.2.9.101 — Additional requirements for operating instructions**

Some *EMS ventilators* are designed so that they can operate with higher-than-normal tubing circuit compliance and resistance. Thus, knowledge of these *VBS* characteristics is important for the *operator* to be aware of the *EMS ventilator* capability. Also, knowledge of the maximum *VBS* resistance (at *nominal* and maximum flowrates) is important because an occlusion *false positive alarm condition* can be caused by the use of high resistance components in the *VBS*. These characteristics of the *VBS* need to be inclusive of any inhalation or exhalation particle/*BSF*, *humidifier*, *nebuliser*, *water collection vessels* and *connectors* needed for operation.

— **201.7.9.2.14.101 — Additional requirements for accessories, supplementary equipment, used material**

The use of antistatic or electrically conductive materials in the *VBS* is not considered as contributing to any higher degree of safety. On the contrary, the use of such material increases the *risk* of electrical shock to the *patient*.

— **201.7.9.3.1.101 — Additional general requirements**

**a)**

The *manufacturer* is expected to express the description of the *EMS ventilator* in general terms so the reader can understand the important behaviour of the *EMS ventilator* (e.g. mean values and their time specifications, number of breaths and delays, etc.).

— **201.9.4.3.101 — Additional requirements for instability from unwanted lateral movement**

In developing the requirements for this subclause, the committees reviewed the test limits for acceleration from EN 1789<sup>[36]</sup> and from RTCA DO-160G<sup>[37]</sup>.

It was noted that this is a requirement for crash safety — the *EMS ventilator* cannot detach from its attachment in the event of a collision of the vehicle in which it is being used, as this would create a *risk* of traumatic injury both to the *patient* and to attending *operators*.

For road ambulance use, 4.5.9 of EN 1789<sup>[36]</sup> provides requirements for *medical device* restraints to withstand acceleration profiles with peak acceleration in each direction of 10 g. This is applied as an acceleration impulse, involving a peak acceleration in the range of 80 m/s<sup>2</sup> to 120 m/s<sup>2</sup>, (8,2 g to 12,2 g) with a profile with an overall period between 80 ms and 150 ms. It was noted that this is equivalent to being brought to a stop from a speed in the range of 17,3 km/h to 51,8 km/h.

The committees agreed that an appropriate test for the *EMS ventilator* attachment that would ensure conformity to the requirements in EN 1789<sup>[36]</sup> under all allowed test conditions would be to require a test at 10 g for a period of 150 ms. The profile specified in EN 1789 represents an impact at speed of 30 km/h to 32 km/h.

For use in an air vehicle, the most appropriate requirement identified is from 7.2.1 of RTCA DO160<sup>[37]</sup>. The standard includes a table listing acceleration profiles for crash tests applicable to various different types of aircraft. These include different profiles for orientation along and perpendicular to the direction of travel. The committees adopted the worst-case stress for each axis for all types of aircraft, assuming that the mounting could assume a vertical direction, but could not assume forward-to-back or sideways orientation with respect to the direction of travel.

This results in test acceleration values of 20 g vertically downwards; 4 g vertically upwards; and 20 g in any horizontal direction. Testing is required both for multiple impulses (with an impulse duration of not greater than 11 ms) and for sustained acceleration for 3 s.

The committees agreed that an appropriate test for the *EMS ventilator* attachment for an *EMS ventilator* intended for air transport use would be to require a test at 20 g for a period of 3 s horizontal and vertically downwards, and 4 g vertically upwards.

— **201.9.6.2.1.101 — Additional requirements for audible acoustic energy**

— **Table 201.102 — Test conditions for acoustic tests**

After due consideration, the committees decided that where this document specifies adjoining ranges for variables as the basis for testing and the declaration of performance, the end value of both ranges should be applicable to both ranges. This means that a *manufacturer* is free to use a round number end value (e.g. 300 ml) in specifications and is not forced to truncate artificially the declared range in order to avoid having to satisfy also the test requirements of the adjacent range. This permits, for example, one *ventilator* to have a declared *tidal volume* range of 300 ml to 1 000 ml and another 100 ml to 300 ml, with each *ventilator* only being required to be tested for the conditions specified for  $\geq 300$  ml or  $\leq 300$  ml respectively.

— **201.11.1.2.2 — Applied parts not intended to supply heat to a patient**

The human airway has a very significant ability to absorb or deliver heat and moisture. Reference the common practice of sitting in a sauna without *harm* to the respiratory tract<sup>[50]</sup>. Fully saturated gas at 45 °C can be inspired for up to 1 h without damaging the mucosa of the respiratory tract<sup>[58]</sup>. A more recent study reported tolerance of inspired gas temperatures of 46,9 °C to 49,3 °C, 100 % RH (265,6 kJ/kg) for 45 min<sup>[46]</sup>.

Taking into account the enthalpy of inspired gas that has been shown to be tolerated without causing thermal injury to the human airways and the very short exposure times of thermal overshoot from a *ventilator* in clinical practice, the delivered gas energy limit of 197 kJ/m<sup>3</sup> of dry gas when averaged over 120 s can be used.

When considering gas mixtures other than oxygen/air, the following should be observed. Given that most of the energy is contained in the water vapour, the equivalent of air at 43 °C, 100 % RH is the maximum enthalpy that should be allowed. This has a specific volume of 0,9786 m<sup>3</sup>/kg of dry air

and an enthalpy content of 197 kJ/m<sup>3</sup> of dry air. Assuming the volume breathed by the *patient* is the same whatever gas mixture is used, then the safe enthalpy limit is 197 kJ/m<sup>3</sup> of dry gas. This enthalpy per unit volume gives a more relevant measure of the energy delivered to the *patient*.

Sustained delivered gas temperatures above 41 °C, depending on the combination of gas temperature, level of water saturation and *patient* exposure time, can be hazardous. Confirmed by studies conducted by the US Navy Medical Research and Development Command<sup>[58]</sup>, concluded that fully saturated gas at 45 °C can be inspired for up to 1 h without damaging the mucosa of the respiratory tract.

Gas at body temperature and fully saturated (37 °C and 100 % RH) does not transfer thermal energy to or from the *patient* with a normal body temperature of 37 °C. Dry gas at body temperature (37 °C and 0 % RH) draws heat away through evaporation. Gas at 41 °C and fully saturated has the capacity to deliver less than 181,3 kJ/m<sup>3</sup> of dry gas breathed by the *patient*.

The enthalpy content of 197 kJ/m<sup>3</sup> has been in use a long time to limit the energy transfer of humidified breathing air to the respiratory tract of a *patient* with bypassed upper airways and no negative feedback with regard to this limitation was known at the time of the consideration of this document. The committees asked for clinical advice as to whether in addition to the enthalpy limitation, the temperature needed also to be limited. A German group of clinicians after considering the issue and the literature available came to the following conclusion:

The literature shows:<sup>[46][50][58]</sup>

- thermal inhalational traumas with temperatures above 100 °C but with unknown humidity content;
- the very low RH of about 5 % only at 100 °C with an enthalpy content of 197 kJ/m<sup>3</sup>;
- the good experience with the limit of enthalpy content of 197 kJ/m<sup>3</sup> of dry air in humidified breathing gases;
- the physical facts that a blower/turbine type *ventilator* increases the temperature of the gas taken from the environment in the range of 15 °C to 25 °C;
- *ventilators* are used in environments up to 45 °C;
- an additional temperature limitation is intended to limit the temperature under worst-case condition; and
- a sufficient safety margin to protect the *patient* from thermal injuries of its airways.

The clinician group recommended to keep the thermal energy limitation of 197 kJ/m<sup>3</sup> and add a maximum temperature limitation of 70 °C<sup>[50]</sup>, whichever is lower. The committees agreed and confirmed this proposal.

#### — 201.11.6.6 — *Cleaning and disinfection of ME equipment or ME system*

The *essential principles* of ISO 16142-1 require that medical devices are not to be operated or used if their condition could compromise the health and safety of the *patient* on whom they are being used or the employees or third parties interacting with them.

This means that *ventilators*, their *accessories* and parts cannot be used if there is a unacceptable *risk* of the *patient*, *operator* or other person being infected as a result of contact with the *ventilator*, *accessory* or part.

Therefore prior to reuse, *ventilators*, their *accessories* and parts require an appropriate level of *disinfection*, depending on their use, but rarely need to be *sterile*.

Recommendations for hygienic *processing* of *ventilators*, their *accessories* and parts are based on the general hygiene requirements for the *processing* of medical devices and need to take into consideration the special requirements and needs of *patient* care in the clinical environment. The requirements for hygienic *processing* in this document are intended to:

- make the *responsible organization* for *processing* the *ventilator* aware of how to implement these tasks in a responsible manner through appropriate delegation; and
- help all parties involved in the *processing* of *ventilators*, their *accessories* and parts to comply with the *manufacturer's* instructions.

The *cleaning* and *disinfection procedures* of the *manufacturer* are also intended to provide practical support to all those involved in *patient* care in the clinical environment with regard to implementing the hygiene measures required for the *patient's* safety.

The *manufacturer* is encouraged to concentrate on four essential aspects of *cleaning* and *disinfection*:

- a) the external *enclosure* of the *ventilator*;
- b) the *ventilator's* removable breathing circuit, including *accessories* and parts (e.g. *humidifier*, removable flow sensor, *connectors*, water traps, *breathing system filters*);
- c) the internal *gas pathways* that can become contaminated with body fluids or by contaminants carried by expired gases during *normal condition* or *single fault condition*, which generally reside within the *ventilator's enclosure* and are normally replaced or processed between *patients*; and
- d) the internal *gas pathways* that cannot become contaminated with body fluids or by contaminants carried by expired gases during *normal condition* or *single fault condition*, which generally reside within the *ventilator's enclosure* and are not normally removed or processed between *patients*.

Regarding the *cleaning* and *disinfection* or *cleaning* and *sterilization* of a), b) and c), above, *manufacturers* are required to provide detailed, *validated procedures* for ensuring safe and effective *processing* to protect the next *patient*, caregivers, technicians and third parties from pathogenic contamination. *Manufacturers* are required to document these *procedures* in the *accompanying documents* of the *ventilator*. Since the *gas pathways* described in d) above are not contaminated by a *patient*, *manufacturers* are not required to provide *validated procedures* for *processing* them. Item c) and to some extent item d), the *disinfection* of the *ventilator's* internal *gas pathways*, have received renewed attention due the recent outbreaks of contagious diseases like Legionnaires', SARS (severe acute respiratory syndrome) and influenza that affect the respiratory system.

Most modern *ventilators* are designed to permit removal, either for replacement or *processing*, of those portions of the internal *gas pathways* that can become contaminated with body fluids or by contaminants carried by expired gases during *normal condition* or *single fault condition*. *Responsible organizations* need to follow their infection control *procedures* when transferring a *ventilator* from one *patient* to another.

In the event that the *responsible organization* suspects that the internal *gas pathways* of a given *ventilator* might have become contaminated with pathogenic material from the previous *patient*, the committees suggest a three-step *process* in the following order to assess next actions:

- ensure that the breathing circuit and all *accessories* have been removed and dealt with according to applicable *procedures*;

- thoroughly disinfect all outer surfaces of the *ventilator enclosure*, including the outer surfaces of the *gas output port* and *gas return port*; and
- only after executing the previous two steps, swab the inner surfaces of the *gas output port* and *gas return port* and culture the swabs to determine if pathogenic material is present.

If contagious pathogens are detected, the *responsible organization* should follow the *manufacturer's processing procedures* found in the *accompanying documents* to protect the next *patient*, caregivers, technicians and third parties from those detected pathogens. The *responsible organization* should be aware that sensitive and expensive sensors are likely located in the expiratory *gas pathways*. Follow the *accompanying documents* of the affected *ventilators* to ensure that delicate and possibly fragile sensors are not damaged during the *processing procedure*.

Any *ventilator* that has already been used on another *patient* is potentially contaminated with contagious pathogenic microorganisms until proven otherwise. Appropriate handling and *processing procedures* are essential to protect the next person handling the device or the next *patient* on whom the device is used. Hence *ventilators*, their re-usable *accessories* and parts that have been used should, when necessary, undergo a *processing process*, following the *manufacturer's* instructions, prior to reuse on another *patient*.

The following basic considerations need to be addressed by the *manufacturer* when specifying the *processing instructions* of a *ventilator*, its *accessories* or parts:

- protecting the *patient*, the *operator* and the *responsible organization* (including personnel involved in performing the *processing procedure*);
- the limits of the *procedures* used for *processing* (such as the number of *processing cycles*); and
- the necessity to guarantee the proven standardised *procedures* have consistently high and verifiable quality, based on an established quality management system.

The recommended *processing procedures* should be determined by:

- the potential degree and type of contamination of the *ventilator*, *accessories* or parts; and
- the *risk* of infecting another *patient* resulting from their reuse and the type of application of the *ventilator*.

Special consideration of the possible *risk* associated with the contamination of gas-conducting components due to the *patient's* re-breathing under *single fault condition* is required.

On the basis of the above, a *verified* and *validated* documented *processing procedure* needs to be specified in such detail that the outcome is reproducible. An acceptable *residual risk* from the *hazard* of infection for the next *patient* can be assumed if the:

- documented *processing procedure's* effectiveness has been *verified* through appropriate scientific methods by the *manufacturer*; and
- reliability of the documented *processing procedures* has been *verified* in practice through appropriate quality assurance measures by the *responsible organization* carrying out the *processing procedures*.
- When selecting and evaluating the *processing procedures*, the *manufacturer* should consider:
  - the amount and type of pathogenic microorganisms expected to contaminate the *ventilator*, *accessories* or parts;

- the *risk* for the pathogenic microorganisms to be transmitted to the *patient, operator* or other persons; and
- the microorganism's resistance to the recommended *processing procedures*.

The *risks* posed by a processed *ventilator, accessories* or parts are determined by the following factors:

- undesired effects, which can result from:
  - the previous use,
  - the previous *processing*, and
  - transportation and storage;
- the *risks* from subsequent uses, such as the following:
  - residues from the previous use (such as secretions, other body fluids, and drugs);
  - residues from the previous *processing procedures* such as *cleaning* agents, disinfectants and other substances, including their reaction products;
  - changes of physical, chemical or functional properties of the device; and
  - changes in the condition of the material (such as accelerated wear and tear, embrittlement and changed surface conditions, *connectors* and adhesive joints); and
- the *risk* of transmission of any pathogenic microorganisms.

When considering the suitability of the *processing procedures* and the feasibility of the *processing procedure* for the *ventilator, accessories* or parts, the *manufacturer* should consider the following points:

- the *risks* involved in the *processing*;
- the cost effectiveness of the *processing*;
- the practicability of the *processing*;
- the availability of the *cleaning* equipment and the *cleaning* agents specified in the *processing*;
- the efficiency of the *processing*;
- the reproducibility of the *processing*;
- quality management requirements of the *processing*; and
- the environmental impact of the *processing* and the disposal of the *ventilator, accessories* or parts.

The *manufacturer* should verify all *cleaning* agents and *processing* used with regard to their suitability and repeatability with the *ventilator, accessories* or parts, depending on the type of use.

The *responsible organization* should verify that *cleaning* and *disinfection* of the *ventilator*, *accessories* or parts are always carried out in accordance with the *procedures* specified in the *accompanying document*.

The *manufacturer* should specify *validated* automated *cleaning* and *disinfection procedures*. If they are not followed, the effectiveness of the *cleaning* and *disinfection* cannot be guaranteed. Such parameters could include the volume of water used, water pressure, temperature, pH, dosage of *cleaning* agents and disinfectants and residence time.

To ensure the reproducibility of automated *processing procedures*, tests should be carried out on a regular basis.

The *manufacturer* should ensure that the specified *disinfection procedures* are *verified* to be bactericidal, fungicidal and virucidal so that the cleaned and disinfected *ventilator*, *accessories* or parts do not pose an unacceptable *risk* of infection by reproductive pathogenic microorganisms when any of these elements, collectively or individually comes in contact, either directly or indirectly, with the next *patient*, *operator* or another person.

Effective *disinfection* requires that the instructions for the disinfectant, especially with regard to concentration and residence time, are followed.

Following any *processing procedure*, a safety and functional testing of the *ventilator* (as specified by the *manufacturer's* instructions) needs to be carried out. If necessary, safety-relevant functional testing can be carried out directly before use of the *ventilator*.

The extent and type of the tests depends on the *ventilator*, *accessory* or part and these need to be defined in the *accompanying document*.

— **201.12.1.101 — Volume-control inflation-type**

— **b) 3)**

A *ventilator* that is intended for use with a *medical gas pipeline system* containing Oxygen 93 (from an oxygen concentrator) requires a wider *rated* range of input oxygen concentration than a *ventilator* intended only for use with a *medical gas pipeline system* containing oxygen. For a *ventilator* intended for use only with oxygen, the *rated* range of concentration for the oxygen *inlet* would typically be “99,5 % ± 0,5 %” (or alternatively “99,0 % to 100,0 %”). For a *ventilator* intended for use with a *medical gas pipeline system* containing Oxygen 93, the *rated* concentration could be as wide as “90 % to 100 % O<sub>2</sub>” where the *medical gas pipeline system* is backed up with oxygen. The *manufacturer* is required to disclose the effects of variation in the *inlet* oxygen concentration within the *rated* range on delivered oxygen accuracy.

For example, the disclosed tolerance needs to be broadened to account for variation in *inlet* oxygen concentration. A *ventilator* that is indicated for *rated* range of “90 % to 100 % O<sub>2</sub>” and that does not measure the *inlet* concentration, and that has a blending accuracy of ±4,5 %, would have a delivered oxygen concentration tolerance of ±5 % when used with oxygen, but might have a delivered oxygen concentration tolerance of ±10 % when configured for use with Oxygen 93. The variation in *inlet* concentration sums with the blending error term to yield a delivered concentration error.

Alternatively, a *ventilator* with identical blending accuracy of ±4,5 % and that measures the *inlet* concentration with an accuracy of ±1,5 % might have a delivered oxygen concentration tolerance of ±6 % (and a restricted setting range of 21 % to 90 %), when used with Oxygen 93.

The *manufacturer* is not expected to repeat testing with the oxygen *inlet* concentration being varied across the *rated* range. It is expected that the effect of variation in the *inlet* oxygen concentration

that is theoretically derived can be used to adjust the results of testing with a known concentration to determine blending accuracy.

### Table 201.104 — *Volume-control inflation-type testing*

The test lung resistance values used in Tables 201.104 and 201.105 are essentially unchanged from the test cases specified in now withdrawn ASTM F1100:90<sup>[34]</sup>. The only difference has been the addition of five further infant test cases to extend downwards the range of *tidal volume*.

ASTM F1100<sup>[34]</sup> included a note that resistance values are for ‘parabolic resistors’. These are simple devices comprising an orifice plate with a cylindrical hole that are easy to manufacture to good tolerance. However, these have a variation in pressure with flow that follows a quadratic. This results in an effective increase in resistance with flow, much more so than for resistance in the *lung*, and also more than for resistance of a tracheal tube. While linear resistors can be constructed, they typically rely on use of porous media or small-scale laminar-flow structures to achieve a linear pressure-flow profile, and are more difficult to manufacture to close tolerance.

In the first edition of ISO 80601-2-12, the test lung configuration values were retained; however the tables were amended to specify use of linear resistance. At the time, this was thought to be more representative of clinical use, and the increasing availability and use of electronically-controlled test lungs meant that modelling linear resistance was less challenging.

However, it has become clear that in some cases this leads to clinically implausible test cases, particularly those test cases that use R200 resistance (intended to represent a 3,0 mm tracheal tube) with *tidal volumes* of 30 ml or greater. Test lung resistors are calibrated at a test flow of 1 l/s (60 l/min). Both a parabolic and a linear resistor with *nominal* resistance of R200 provide a pressure drop of 200 hPa at 60 l/min test flow. However, the parabolic resistor will have a pressure drop of only 50 hPa at a test flow of 30 l/min, while the linear resistor pressure drop is 100 hPa. The effect is even more pronounced at lower flow rates — the pressure drop at 15 l/min is 12,5 hPa for the parabolic, versus 50 hPa for the linear resistor.

In normal breathing patterns, flow in the *lung* is approximately laminar (owing to the preponderance of the resistance being in the very small-scale passages in the lower levels of the tracheobronchial tree). Typical values for total pulmonary system resistance in mouth-breathing humans vary from approximately 20 hPa·(l/s)<sup>-1</sup> to 25 hPa·(l/s)<sup>-1</sup> in neonates, to less than 5 hPa·(l/s)<sup>-1</sup> in adults<sup>[57]</sup>. These values are significantly increased in restrictive airway disease. However the very high values of resistance found in the *ventilator* test case tables largely represent the resistance of the tracheal tube.


Measurements of pressure drop versus flowrate for seven sizes of tracheal tubes have been published<sup>[49]</sup>. A 7,0 mm tracheal tube has a resistance at 1 l/s that is close to R10 [10 hPa·(l/s)<sup>-1</sup>], however the resistance is lower at lower flow rates, for example at a flowrate of 0,25 l/s the resistance is approximately R4 [4 hPa·(l/s)<sup>-1</sup>].

For smaller sizes of tracheal tubes, the resistance changes significantly as flow increases above a critical value, which represents the onset of turbulence. Figure AA.1 demonstrates this effect.

The critical flow ranges from 0,25 l/s for a 3,5 mm tracheal tube, to 0,63 l/s for a 7,0 mm tracheal tube. The critical flow is evident in Figure AA.1 as a notch in the pressure curve at approximately 12 l/min (0,2 l/s). Below the critical flow, the pressure drop scales approximately as flow to power 1,3, while above the critical flow the pressure drop scales approximately as flow to power 1,8.

Table AA.1 shows the test cases from ISO 80601-2-12:2020, Table AA.1. The first five columns are taken from that table. It should be noted that parabolic resistors are commonly available in a limited number of values, typically Rp5, Rp20, Rp50, and Rp200, and that the table was constructed to use these standard values.

The next two columns show the average flow required (assuming a rectangular flow waveform) and the peak flow required (assuming a worst-case linear decreasing flow pattern, terminating at zero flow). The right hand columns show the worst-case pressure drop for worst-case flow assuming a parabolic resistor, or a linear resistor, respectively.

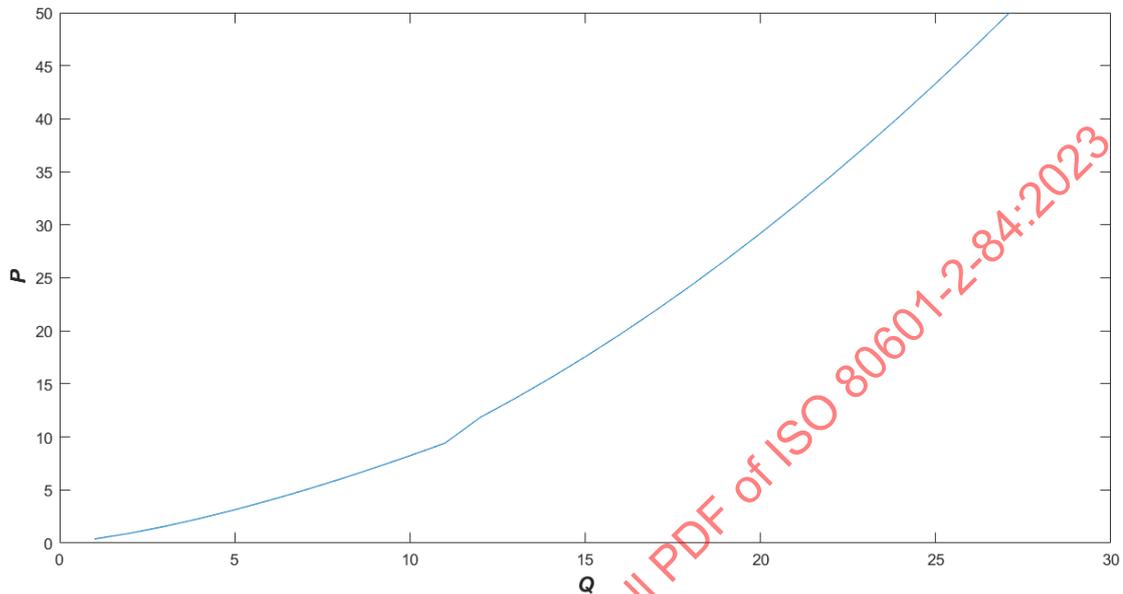

#### Key

- $P$  pressure drop (hPa)  
 $Q$  flowrate (l/min)

**Figure AA.1 — Pressure drop calculation for 3.0 mm tracheal tube, 100 % RH room air at sea level, 37°C, using approach specified by Jarreau<sup>[49]</sup>**

It is clear that changing the resistor specification from parabolic to linear has very significantly altered the pressure drop due to resistance in all test cases with an *inspiratory volume* less than 300 ml. In several cases, such as #11, #13 and #15, the resulting pressure drop is far in excess of the pressure drop appropriate to this *patient* model. However, in other test cases such as #18 through #21, the parabolic resistor is clearly too low a resistance to be appropriate.

Calculated pressure drop across a number of sizes of tracheal tubes at specified flow rates appropriate to that size *patient* are tabulated in Table AA.2, expressed both as actual pressure drop, and as a resistance value [expressed in  $\text{hPa}\cdot(\text{l/s})^{-1}$ ].

The test cases need to account for both the resistance of the tracheal tube, and the resistance of the *patient's* pulmonary system. This suggests that for *tidal volumes* of no greater than 20 ml (test cases #16 to #21) resistance values of  $50 \text{ hPa}\cdot(\text{l/s})^{-1}$  to  $200 \text{ hPa}\cdot(\text{l/s})^{-1}$  remain appropriate. For test cases with *inspiratory volumes* of 30 ml to 50 ml, representing term neonates or infants, resistance values of  $20 \text{ hPa}\cdot(\text{l/s})^{-1}$  to  $100 \text{ hPa}\cdot(\text{l/s})^{-1}$  would be appropriate; for test cases of 50 ml or greater, resistance values should be  $5 \text{ hPa}\cdot(\text{l/s})^{-1}$  to  $50 \text{ hPa}\cdot(\text{l/s})^{-1}$ . Based on this, the resistance values in Table 201.104 and Table 201.105 have been adjusted from the values contained in previous standards.

**Table AA.1 — Flow and pressure drop for linear and parabolic resistors**

Test #	C mL/hPa	R hPa·(l/s) <sup>-1</sup>	V <sub>T</sub> ml	Insp Time s	Mean Flow l/min	Peak Flow l/min	ΔP hPa (parabolic)	Δ hPa (linear)
1	50	5	500	1	30	60	5,0	5,0
2	50	20	500	1	30	60	20,0	20,0
3	20	5	500	1	30	60	5,0	5,0
4	20	20	500	1	30	60	20,0	20,0
5	20	20	300	1	18	36	7,2	12,0
6	20	50	300	1	18	36	18,0	30,0
7	10	50	300	1	18	36	18,0	30,0
8	10	20	200	1	12	24	3,2	8,0
9	3	20	50	0,6	5	10	0,6	3,3
10	3	50	50	0,6	5	10	1,4	8,3
11	3	200	50	0,6	5	10	5,6	33,3
12	3	50	30	0,6	3	6	0,5	5,0
13	3	200	30	0,6	3	6	2,0	20,0
14	1	50	30	0,6	3	6	0,5	5,0
15	1	200	30	0,6	3	6	2,0	20,0
16	1	200	20	0,4	3	6	2,0	20,0
17	1	200	15	0,4	2,25	4,5	1,13	15,0
18	1	50	10	0,4	1,5	3	0,13	2,5
19	0,5	50	5	0,4	0,75	1,5	0,03	1,3
20	0,5	200	5	0,4	0,75	1,5	0,13	5,0
21	0,5	200	5	0,4	0,75	1,5	0,13	5,0

**Table AA.2 — Flow and pressure drop for linear and parabolic resistors**

Tracheal tube diameter mm	Test Case Flow l/min	Pressure Drop hPa	Resistance hPa·(l/s) <sup>-1</sup>
2,5	0,5	0,36	43
2,5	2,0	1,96	59
3,0	1,0	0,4	23
3,0	4,0	2,3	35
4,0	3,0	0,5	9
4,0	10,0	2,5	15
5,0	5,0	0,4	4,5
5,0	24,0	3,6	9,0
7,0	20,0	0,6	1,8
7,0	60,0	3,7	3,7

— 201.12.1.102 — *Pressure-control inflation-type*

b) 3)

See rationale for 201.12.1.102 b) 3).

**Table 201.105 — Pressure-control inflation-type testing**

See rationale for Table 201.104.

— **201.12.4.101 — Oxygen monitor**

After a long discussion, the committees concluded that it is not the length of *intended use* but the *patient* population (i.e. neonatal use) that should be used to make oxygen monitoring mandatory. Monitoring the oxygen level and including level *alarm conditions* is an ineffective *risk control* for an *EMS ventilator* that does not control the oxygen concentration (e.g. venturi entrainment) since the delivered concentration varies widely with flow and pressure.

— **201.12.4.102 — Measurement of airway pressure**

The site in the *VBS* at which pressure is sensed varies from *ventilator* to *ventilator*. Generally, the *manufacturer* chooses one of two strategies:

- measuring the *airway pressure* by direct sampling at the *patient-connection port*: or
- indirectly estimating the pressure at the *patient-connection port* by measuring the pressures at two locations in the *ventilator*, on the inspiratory side of the *VBS* (at the “to *patient*” port) and on the expiratory side of the *VBS* (at the “from *patient*” port), and, after mathematical manipulation, averaging the two values.

Even if the first strategy is chosen, the actual pressure transducer will be located inside of the *ventilator enclosure* with narrow-diameter “plastic” tubing linking the pressure-sampling port at the *patient-connection port* to the sampling nipple on the pressure transducer. And for safety reasons a separate transducer will likely measure the pressure on the inspiratory side at the “to *patient*” port. The displayed *airway pressure*, however, is always expected to estimate accurately the true value that would be measured at the *patient-connection port*. Pressure measurement via the first strategy accurately reflects the true *airway pressure* within the error of the pressure transducer.

If the *manufacturer* selects the second strategy for the prediction of the true *airway pressure*, at least two methodologies can be used to arrive at estimates of this *airway pressure*.

Assuming that during inspiration the gas in the expiratory limb is essentially stagnant, one can conclude that the pressure measured on the expiratory side of the *VBS* reflects the true *airway pressure*. And conversely during exhalation, if one assumes stagnant conditions in the inspiratory limb, the pressure measured on the inspiratory side of the *VBS* can be taken as the *airway pressure*. However, if “bias” or “base” flows during inspiration and exhalation result in significant pressure losses across these individual limbs, these pressure losses need to be estimated. The *airway pressure* on the inspiratory side,  $P_V^I(t)$ , can be approximated by Formula (AA.1):

$$P_V^I(t) = P_I(t) - \dot{V}_I(t) \times R_I \quad (\text{AA.1})$$

where

$P_I(t)$  is measured pressure on the inspiratory side of the *VBS*,

$\dot{V}_I(t)$  is the flow in the inspiratory limb, and

$R_I$  is the resistance of the inspiratory limb.

The *airway pressure* on the expiratory side,  $P_V^E(t)$ , can be approximated by Formula (AA.2):

$$P_Y^E(t) = P_E(t) + \dot{V}_E(t) \times R_E \quad (\text{AA.2})$$

where

$P_E(t)$  is measured pressure on the expiratory side of the *VBS*,

$\dot{V}_E(t)$  is the flow in the expiratory limb, and

$R_E$  is the resistance of the expiratory limb.

Taking the average of the inspiratory and expiratory pressures is shown in Formula (AA.3), which arrives at the best estimate of the *airway pressure*,  $\bar{P}_Y(t)$ .

$$\bar{P}_Y(t) \cong \frac{P_Y^I(t) + P_Y^E(t)}{2} \quad (\text{AA.3})$$

The application of this last method requires a method to estimate  $R_I$  and  $R_E$ . With appropriate algorithms and regular cross checking of the two pressure transducers, the reliability and accuracy of  $\bar{P}_Y(t)$  is assured.

— **201.12.4.103 — Measurement of expired volume and low-volume alarm conditions**

It is desirable to have a fast-responding measurement of volume and narrow *alarm limits*. However, as there is often considerable variation in a *patient's* ventilatory pressures and volumes, narrow *alarm limits* inevitably lead to clinically insignificant *alarm conditions*. As a result, *operators* choose to set wide *alarm limits* to reduce the number of insignificant *alarm conditions* despite the fact that this can compromise *patient* care when there is a prolonged small change in their *ventilation*. Therefore, it is recommended that an *EMS ventilator* be designed to initially use a lower priority *alarm condition*, which escalates to a higher priority if the *alarm limit* violation persists. The initial *alarm condition* priority and the priorities and timing of the escalation should be determined by the severity of the potential *harm* to the *patient* in combination with the length of time that the *operator* has to prevent the *harm* from occurring.

— **d) 2)**

In this document the phrase '*responsible organization-configurable*' is used to describe *ventilator* settings that are prearranged either to a specific value or a range of values by the *responsible organization* for a particular setting. *Responsible organization-configurable* settings are normally protected from change by the normal clinical *operator*. When the *responsible organization-configurable* parameter is constrained to a range, the *operator* typically can then adjust the parameter value within the constrained range.

— **201.12.4.104– Maximum limited pressure protection device**

The value chosen for the *maximum limited pressure*<sup>[47][60]</sup> is a compromise between the need to avoid barotrauma and the need to provide an adequate range of pressure to meet the desire of *operators* to supply both high *PEEP* and high insufflation pressures for specific *patients* who have not responded favourably to optimal *ventilation-management* strategies aimed at *lung* protection and *lung* recruitment. In such cases the *operator* can, as a last resort, elect to ventilate using high *PEEP* and high insufflation pressure.

The predecessor standards ISO 10651-3:1997 specified the *maximum limited pressure* not exceeding 100 hPa or 120 % of the *maximum working pressure*, whichever is greater and EN 794-3:1998+AMD1:2005+AMD2:2009<sup>[35]</sup> specified a *maximum limited pressure* of 100 hPa. These values were based on the clinical need to supply a pressure up to 70 hPa at the *patient-connection port* under emergency and transport conditions for critical *patients*.

The design of pressure limiting in *EMS ventilators* in general is implemented via pneumatic/mechanical pop-off devices. The opening of the pop-off device needs a significant time. The combination of this delay and the relative high resistance of the pop-off device when open, creates a significant increase of pressure (above the opening pressure) during the relief action.

Therefore, an additional pressure margin was needed to provide working pressures up to 100 hPa to the *patient-connection port*.

The technology available to limit the pressure in the *VBS* has changed so that a smaller additional pressure margin is needed to ensure the ability to provide sufficient pressure to the *patient-connection port*.

Considering the *risk* of barotrauma caused by high pressures supplied to the *patient's lungs*, considering the pressures necessary to ventilate most of the *patients* around the world and considering the technology available as of today a survey was initiated to determine:

- the *maximum working pressure* clinically needed to ventilate most *patients*; and

NOTE The consensus was that “to ventilate most *patients*” means that in 90 % to 95 % of the cases *patients* can sufficiently ventilated with this ventilatory pressure.

- the needed pressure difference between the *maximum limited pressure* and the *maximum working pressure* with today's technology.

The following national member bodies participated in this survey: FR, US, GB, DK, NZ, AU, CN, SE, JP, SA, DE and CA.

The conclusion of this survey was that:

- the maximum clinical pressure necessary at the *patient-connection port* is 90 hPa; and
- the maximum needed pressure difference between the *maximum limited pressure* and the *maximum working pressure* at the *patient-connection port* is 15 hPa.

Based on the above, the conclusion of the committees was to limit *maximum limited pressure* to 85 hPa or 120 % of the *maximum working pressure*, whichever is greater.

#### — 201.12.4.105 — High airway pressure alarm condition and protection device

##### e) 1)

The increment of 30 hPa between the *operator-set pressure* and the high *airway pressure alarm limit* is greater than the increment of 20 hPa specified in standards for *lung ventilators* intended for *professional healthcare facility use*. *EMS ventilators* are subject to more stringent environmental conditions, particularly with regards to vibration, and this results in a wider variability in pressure delivery. The permissible increment between the *operator-set pressure* and the high *airway pressure alarm limit* allows for the *operator* to set the *alarm limit* to a value that does not result in false *alarm conditions* due to the influence of vibration, without *risk* of rendering the high *airway pressure alarm condition* ineffective as a control for high *airway pressure*.

##### f)

The high *airway pressure alarm condition* can occur in several scenarios, with different appropriate responses. It should be noted that it is implicit that in *normal condition*, the exhalation resistance of the *VBS* (including breathing circuit, filters, and exhalation valve) is sufficiently low so that the pressure can return to set *BAP* within one *respiratory cycle* of the end of a *ventilator* inspiration.

**Patient cough**

A cough is a transient forced expiration due to involuntary contraction of the abdominal muscles and opening of the glottis and vocal chords, resulting in a very rapid and short duration expiratory flow. This is most commonly triggered by pulmonary irritant receptors, and is a defensive mechanism that ejects foreign matter from the bronchi and trachea.

This action triggers a rapid but transient increase in pressure at the *patient-connection port*. The pressure elevation is generated by the *patient* respiratory muscles. In this scenario, it can be appropriate to transition to exhalation, as this would assist the forced expiration and clearance of foreign material. However, the volume of air ejected has come from the *lungs*, and the normal operation of the *ventilator expiratory phase* control of pressure would be sufficient to restore pressure to a value close to set *BAP*, within the normal exhalation period. There is no clinical reason to open a pressure-relief valve.

**Patient forced expiration**

In some cases, a *patient* can make a forced expiratory effort during the *ventilator inflation phase*, resulting in a pressure increase at the *patient-connection port*. As this is a *patient*-generated pressure, the same rationale applied to cough would apply to this scenario, and transition to exhalation phase is an appropriate and sufficient response. Even in the worst-case scenario, the exhaled volume would be unlikely to exceed twice the normal exhaled volume, and the time required to reduce the pressure to set *BAP* would not significantly exceed one *respiratory cycle*.

**Inappropriate volume for patient compliance**

In a *volume-control* or *volume-targeted ventilation-mode*, the pressure can exceed the high pressure *alarm limit* due to reduced pulmonary compliance. In this scenario it would be preferable to apply a pressure limitation (with a threshold less than the high pressure *alarm limit*). However in the absence of this function, the appropriate action would be to transition to the exhalation phase. As the *inspiratory volume* is not increased, there is no reason why the normal operation of the *ventilator* in the *expiratory phase* would not reduce the *airway pressure* to set *BAP* within one *respiratory cycle*.

**Inspiratory limb occlusion**

In this *single fault condition*, the *patient* is not subject to the high pressure. *Ventilators* that do not incorporate a pressure-monitoring line from the *patient-connection port* might be unable to differentiate this from a high *airway pressure* state. Because the *patient* is not subjected to the high *airway pressure*, an obstruction *alarm condition* is not considered necessary.

**Expiratory limb, breathing system filter or expiratory valve occlusion**

For most software-controlled *ventilators*, this only results in a high *airway pressure alarm condition* if it occurs during the *expiratory phase*, as when the exhalation valve remains closed during *ventilator* inspiration. In some cases (such as a 'pressure limited' *ventilation-mode*), the *patient* pressure can exceed the *alarm limit* if the occlusion prevents the *maximum limited pressure protection device* from operating.

Terminating *outlet flow* — either by terminating the *inflation phase*, or terminating the base or continuous flow — causes the pressure to stop increasing, but might not result in the pressure falling below the *alarm limit*. This constitutes the 'action to cause the pressure to start to decline'. If the pressure fails to return close to set *BAP* within one *respiratory cycle*, this is evidence that the expiratory resistance is excessive, and the *ventilator* should provide an alternative means to relieve the pressure. Typically this would involve opening some form of 'safety valve' in the inspiratory

part of the *VBS*. Assuming that the inspiratory limb of the circuit has comparable resistance to the expiratory limb, it can then take one further *respiratory cycle* to restore pressure at or below set *BAP*. Hence the committees have specified that the pressure is required to be reduced to the set *BAP* or atmospheric pressure within no more than two *respiratory cycles* or 15 s.

### Flow delivery ‘fails open’

In *single fault conditions* such as a mechanical failure of a pressure regulator or a flow control valve, or a failure of software control of a flow source, resulting in uncontrolled flow of gas to the *ventilator outlet*, transition to the *expiratory phase* (which opens the exhalation valve) should be sufficient to control the *risk*. However, the pressure might not reach the limit within 200 ms. If the pressure remains above the high *airway pressure alarm limit* (for example because of a catastrophic flow valve failure) it can then be appropriate to also open some form of pressure-relief valve to divert the *outlet* flow and restore the *patient-connection port* to atmospheric pressure.

### h)

#### Atmospheric pressure

In *normal condition*, termination of *ventilator* inspiration, and the normal action of the expiratory valve, allows pressure to be relieved to *PEEP* within one *respiratory cycle* without loss of control of *PEEP*. *Normal condition* would include any *patient-generated* effect. However, if there is a *single fault condition*, such as an obstruction of the breathing circuit or a flow delivery fault, other means can be required to reduce the pressure, and in this scenario, it is not expected that the *ventilator* is able to maintain control of *PEEP*.

#### Set *BAP* level

As noted above, the worst-case time required to reduce the *airway pressure* to a value less than or equal to set *BAP* is two *respiratory cycles*. To account for the possibility that the mandatory *set rate* has been set to a value significantly less than the *patient's* normal resting breathing rate, the committees have limited the recovery time to a maximum of 15 s. This reflects literature that shows typical breathing rates in adults<sup>[42]</sup>, including geriatric populations<sup>[59]</sup>, that are rarely less than 10 breaths/min, and hence would require no more than 12 s for two *respiratory cycles*.

### — 201.12.4.106 — Expiratory end-tidal CO<sub>2</sub> monitoring equipment

The monitoring of expiratory end-tidal CO<sub>2</sub> is employed clinically as a surrogate for arterial CO<sub>2</sub> tension. It therefore provides an alternative to monitoring expired *tidal volume* in assessing the adequacy of *ventilation* of the *lungs*. However, in the event of an occlusion or leak within the *VBS*, minute *ventilation* can be significantly reduced, while arterial CO<sub>2</sub> rises. This can result in end-tidal CO<sub>2</sub> monitored values that remain within the clinically acceptable range as the *patient tidal volume* is reduced below the level of physiologic dead space. The committees did not believe that safety could be ensured by monitoring of end-tidal CO<sub>2</sub> unless monitoring *inspiratory volume* or expired volume was also in use.

It is not intended that this constrain a *manufacturer* to integrate the CO<sub>2</sub> monitoring equipment into the *ventilator*. This allows for cases such as when the end-tidal CO<sub>2</sub> monitoring equipment is integrated within a vital signs monitor from a separate *manufacturer*.

— **201.13.2.101 — Additional specific *single fault conditions***

**a)**

Disruption of the gas delivery to the *patient* independent of the root cause for the disruption (e.g. by disconnection or blockage of the inspiratory breathing/tubing system) is the most reasonably foreseeable event in the daily practice of *ventilation* around the world that might quickly lead to serious irreversible injury or the death of a *ventilator-dependent patient*.

**b)**

Disruption of the gas flow pathway between the *patient* and the *ventilator* independent of the root cause for the disruption (e.g. by disconnection or blockage of the expiratory breathing / tubing system) is the most reasonably foreseeable event in the daily practice that, depending on the *VBS* in use, might lead to pressure loss. This can consequentially lead to the inability to build up a *VBS* pressure sufficient to ventilate the *patient*, which in further consequence might lead to serious irreversible injury or the death of a *ventilator-dependent patient*.

**c)**

Operation of an *EMS ventilator* without an *operator-detachable breathing system filter* in place is considered reasonably foreseeable when considering those parts of the *VBS* that might become contaminated with body fluids or by contaminants carried by expired gases. If an *EMS ventilator* can operate without the *breathing system filter*, then one has to assume that it has been operated without the *breathing system filter* and therefore those parts of the *VBS* might have been contaminated. Additional information is found in the rationale for 201.11.6.6.

**d)**

Operation of a *ventilator* using an *operator-detachable* remote control or monitoring module is considered as a state-of-the-art option today. Independent of how the communication between the “ventilator module” and the remote control or monitoring module is facilitated (e.g. wired or wireless) this communication needs to be so designed and constructed that a failure or loss of this communication does cause an unacceptable *risk* to the *patient*. Further this communication (e.g. between the “ventilator module”, the remote control or monitoring module, the *distributed alarm system* or a simple remote *alarm signal communicator*) also needs to be designed to be *single fault condition safe*.

— **201.13.103 — Independence of *ventilation control function* and related *risk control measures***

This requirement prevents the use of a monitoring device to control an actuator that would lead to an undetected malfunction of the actuator in case of monitoring failure.

— **201.15.3.5.101 — Additional requirements for rough handling**

An *EMS ventilator* is intended to operate while the *patient* is being transported and is expected to maintain *basic safety* and *essential performance* while it is being moved. Some degradation is permitted, but the *patient* is expected to continue to be adequately and safely ventilated. Rationale for 202.8.1.101 contains additional information regarding appropriate acceptance criteria for *essential performance*.

The committees recognized that during environmental stress (for example, vibration associated with transport in a moving vehicle) *ventilator* performance can degrade. In such a situation, the

duration of altered performance would be at most, few tens of minutes. However, it is imperative for *patient safety* that *essential performance* is maintained.

This document references the vibration test levels defined in IEC 60601-1-12:2014 as survival test. This document provides levels that represent the worst-case exposure in the intended environment.

The committees compared these test levels with those for operating vibration found in standards ISO 10651-3 and EN 794-3<sup>[35]</sup>, currently applicable to *EMS ventilators*; to the vibration levels given in IEC 60601-1-11 for the *home healthcare environment*. The committees also investigated test levels provided for operating vibration in now obsolete FDA guidance (1993 Ventilator Guidance) and appropriate vibration profiles (for example, for lorry transport) in MIL-STD-810G<sup>[54]</sup>. In all cases, the vibration levels are significantly lower than those specified in IEC 60601-1-12. Table AA.3 compares these requirements.

**Table AA.3 — Comparison of vibration requirements**

Document	Frequency coverage Hz	Vibration profile	Vibration amplitude g <sub>rms</sub>
ISO 10651-3	10 to 500	10 Hz to 200 Hz 0,01 g <sup>2</sup> /Hz 200 Hz to 500 Hz 0,003 g <sup>2</sup> /Hz	1,7
EN 794-3	10 to 500	10 Hz to 200 Hz 0,01 g <sup>2</sup> /Hz 200 Hz to 500 Hz 0,003 g <sup>2</sup> /Hz	1,7
IEC 60601-1-11	10 to 2 000	10 Hz to 100 Hz 0,1 g <sup>2</sup> /Hz 100 Hz to 200 Hz -3dB / octave 200 Hz to 2 000 Hz 0,05 g <sup>2</sup> /Hz	10,3
IEC 60601-1-12	10 to 2 000	10 Hz to 100 Hz 0,5 g <sup>2</sup> /Hz 100 Hz to 200 Hz -7dB / octave 200 Hz to 2 000 Hz 0,1 g <sup>2</sup> /Hz	16,0

The committees took regard to statements by *manufacturers* that state-of-the-art homecare *ventilators* are able to pass the testing mandated by ISO 80601-2-72, using vibration profiles defined in IEC 60601-1-11; but that no *manufacturer* was able to confirm that a *ventilator* could operate and provide *essential performance* at the vibration levels stated in IEC 60601-1-12. In particular, at frequencies below 100 Hz, that would be expected to transmit to the *patient*, the test levels in IEC 60601-1-12 are five times greater than in any of the other sources.

Whereas the shock, drop and electromagnetic compatibility for *ME equipment* for the *EMS environment* largely relate to use at an incident scene, operating vibration is primarily a function of motorized transport, and the requirements in a moving ambulance are not deemed to be significantly different from those for a *home healthcare environment ventilator* being used in a car, truck or other motorized vehicle. Accordingly, we have adopted the vibration profile from IEC 60601-1-11 for operational vibration testing.

#### — 201.15.102 — Delivered oxygen concentration

*EMS ventilators* are normally supplied only by one type of supply gas, generally oxygen.

*EMS ventilators* have requirements regarding weight, size, handling, robustness and costs. As a result, an *EMS ventilator* typically is not equipped with accurate oxygen blender functions.


Nevertheless, an *EMS ventilator* is required to provide at least one additional adjustable oxygen concentration with the following intentions.

- High concentrations of oxygen applied for a longer period of time are likely to be toxic for the *patient*. A reduction of the inspiratory oxygen concentration to approximately 60 % reduces the *risk* of toxicity significantly.
- The operating time of *EMS ventilator* supplied via gas cylinders can be substantially extended via a gas-mixing function.

Inspiratory oxygen concentrations of approximately 60 % volume fraction can be easily designed via the injector principle, where high-pressure oxygen is used as driving gas and ambient air is entrained from the atmosphere.

Because of the physical characteristic of the venturi principle, the oxygen concentration depends significantly on the inspiratory flow and the mean *airway pressure*.

Figure AA.2 shows the typical dependency of the resulting inspiratory oxygen concentration from the inspiratory flow at a mean *airway pressure* of e.g. 30 hPa.

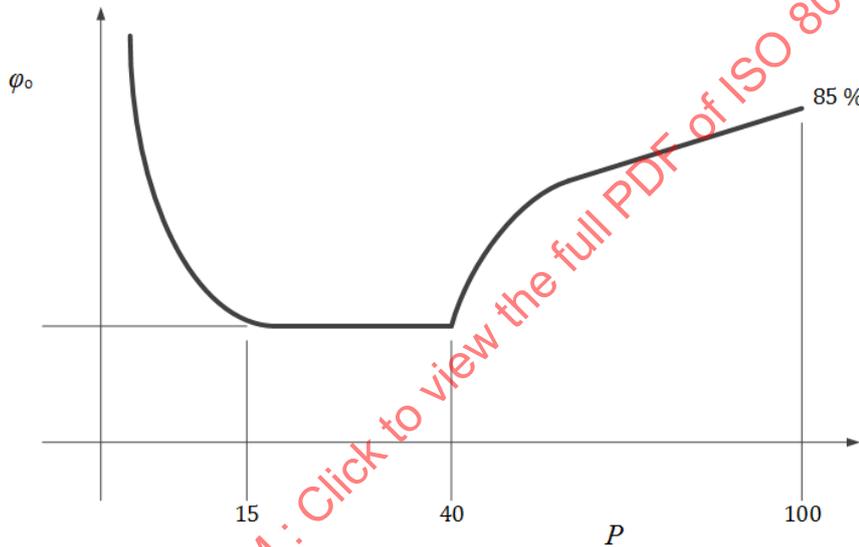


**Key**

- $\varphi_0$  oxygen concentration (volume %)
- $P$  pressure (hPa)

**Figure AA.2 — Entrainment oxygen concentration as a function of flow at 30 hPa pressure**

There is only a limited working range in which the oxygen concentration is on the desired constant level.

At a low flowrate, the venturi is not able to intake sufficient air to achieve the desired oxygen concentration. At high flowrates, the performance limit of the venturi is reached. In both cases, the inspiratory oxygen concentration can rise significantly. Figure AA.2 and Figure AA.3 illustrate this effect.

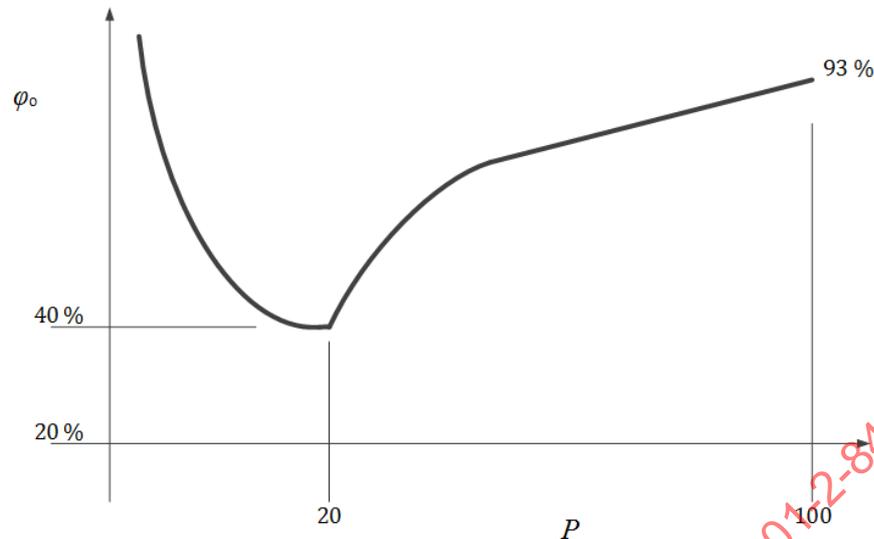**Key**

$\varphi_0$  oxygen concentration (volume %)

$P$  pressure (hPa)

**Figure AA.3 — Entrainment oxygen concentration as a function of flow at 60 hPa pressure**

— **201.101.1 — Protection against reverse gas leakage**

These conditions are necessary to maintain *patient* safety by protecting the *medical gas pipeline system* from contamination via reverse flow.

The basic requirements of this subclause were introduced into standards more than a decade ago because of the *harm* due to reverse gas leakage that was known to have occurred in connection with medical devices that use multiple gas sources.

With devices fitted with multiple *gas intake ports* for the same gas, the *hazardous situation* results from the undetected loss of backup gas supplies due to back leakage into the primary supply. With *gas intake ports* for different gases, the *hazard* is contamination of one gas source by gas from another source. The contamination *hazard* is particularly likely to occur while the medical device is left in a condition where it is connected to the gas supplies but is not drawing flow from the gas supply system.

*Ventilators* are frequently equipped with multiple *gas intake ports* either to achieve a greater flow or to use a local backup supply (e.g. a gas cylinder, in parallel with a *medical gas pipeline systems* supply). With such systems, the backup supply could be depleted prematurely during use or, when connected but not in use, could deplete without detection and not be available when required in an emergency.

With a *ventilator* equipped with more than one different *gas intake port*, even very small leakages from one of the gas systems to the other can cause considerable contamination in a *medical gas pipeline system* over extended periods during which little flow is withdrawn.

More than 10 y of experience has demonstrated that these requirements are effective *risk control* measures.

— **201.101.2.2 — Filter**

The intention of filtration of the gas from the *high-pressure inlet port* is to protect the sensitive components (e.g. flow sensors) of the *ventilator gas pathways* from particles. This gas is provided from a *medical gas pipeline system* or from a gas cylinder.

The standards for high-pressure oxygen compatibility<sup>[12]</sup> and pressure regulators require input filtering that prevents particles greater than 100 µm from entering.

Despite of these requirements in the following cases particles with larger sizes could occur:

- particles collected in *high-pressure inlets* and *port connectors*;
- *high-pressure inlets* of ventilators while disconnected; or
- malfunction of the *medical gas pipeline system*, medical air compressors, oxygen concentrators or filters.

Depending on the design of a specific *ventilator* (e.g. in case that particle-sensitive sensors are used) significantly smaller filter sizes than 100 µm can be required.

— **201.101.3.1 — General**

Non-standard *VBS connectors* can represent an unacceptable *risk* as attempts are made to fit a standard *VBS* to an *EMS ventilator* in an emergency situation. Non-standard *VBS connectors* can cause leaks if used with similar but not compatible *connectors*.

— **201.101.3.2.3 — Flow-direction-sensitive components**

*Flow-direction-sensitive components* commonly used with a *ventilator breathing system* include flow sensors, non-return valves, *breathing system filters* and *heat and moisture exchangers*.

Use errors associated with these components can include the component being placed 'in reverse' (such that the direction of flow is not as intended through it) and can also include the component being placed at an incorrect location within the *ventilator breathing system*.

Each of these use errors can potentially cause *risk of patient harm*. In some cases, this is indirect *harm*. For example, a flow sensor placed in the intended location but incorrect orientation is likely to result in incorrect output from expired volume *monitoring equipment* and can cause incorrect delivered *tidal volume*. In other cases, the *harm* is direct. For example, a non-return valve placed in an incorrect location can cause a complete cessation of inspiratory or expiratory flow, resulting in complete loss of *ventilation*.

At least one case has been reported of post-operative *ventilation* where a non-return valve was used as if it were an adapter to replace an inline nebuliser when it was removed from the circuit. In the reported case, the non-return valve resulted in complete occlusion of the exhalatory pathway, and the *patient* subsequently died.

This standard includes several requirements applicable to *flow-direction-sensitive components*, including *marking* with an arrow for flow direction [201.7.2.101 c) and f)], *instructions for use* [201.7.9.2.14.101 b) 1)], and the constructional requirements of 201.101.3.2.3. However, it is recognised that simply providing a *flow-direction-sensitive component* with mating (rather than identical) *connectors* at *inlet* and *outlet* is not sufficient to prevent all hazardous use errors. *Manufacturers* are expected to consider all foreseeable use errors, including errors of incorrect substitution and incorrect placement as well as incorrect orientation, in their *risk management process*.

#### — 201.101.3.2.4 — *Gas pathway connection port*

The use of Luer taper or Luer-lock *connectors* conforming with ISO 80369-7 are permitted for use in a *VBS*. However, this is only permitted when appropriate information for safety is also provided at the *VBS* connection as a *risk control* measure. Nonetheless, there are several case reports of accidental connection with intravenous or other fluids and parenteral and enteral feeding solutions causing serious morbidity and mortality due to aspiration of these foreign substances into the *lungs*.

The committees expect to change this requirement in the next edition or revision to the R1 *connector* of ISO 80369-2 once ISO 80369-2 is published.

#### — 201.102.1 — **General**

It is the responsibility of the *manufacturer* of a *ventilator breathing system*, its parts or *accessories* to verify that their product conforms with the requirements of this document.

#### — 201.102.4 — **Water vapour management**

Water management refers to the complete *process* by which moisture, in the form of water vapour, is added to the breathing gas delivered to the *patient's lungs* and the *process* by which humidified breathing gas is conducted back to the *ventilator's* expiratory system and exhausted into the room. Intrinsic to this *process* is the necessity to remove bulk water due to condensation of moisture attributable to pressure and temperature changes in the *VBS*. Even if breathing gas reaches the *patient-connection port* without any added moisture, the expired breathing gas directed back to the *ventilator* contains some moisture. Water management in the *VBS* requires attention, whether or not the *VBS* contains an active *humidifier*, with or without heated wires in the inspiratory or the expiratory limbs of the *VBS*, or a passive or an active *HME* at the *patient-connection port*.

Proper management of the *patient's* airway secretions and mucociliary transport system requires that the *ventilator* compensate for the humidity deficit caused by intubation, which bypasses the upper airways where the normal humidification *process* would begin. Excess moisture delivered to the *patient-connection port* can flood the cilia located in the bronchial airways, diminishing their ability to move mucus toward the trachea. On the other hand insufficient humidification of the inspired breathing gas dries the bronchial airways, which leads to thickening of the mucous secretions and likely increased airway resistance or worse. A balanced approach to humidification is needed to maintain healthy cilia. Liquefied mucus can be readily aspirated using a suction catheter.

Optimal humidification of the *patient's* airways results from an understanding of the physics of the techniques chosen to add water vapour to the inspiratory gas stream. Depending on the system selected for delivering humidified breathing gas to the *patient* (for example, active vapour *humidifier* with or without heated wires, conventional *HME* or active *HME*), condensate can accumulate in the inspiratory limb of the *VBS*. If condensation occurs, the *VBS* will need to provide a method by which the liquid can be removed.

In all but the most unusual circumstances, gas leaving the alveoli is saturated at 37 °C. Rainout occurs as the moist gas cools and moves toward the *patient-connection port*, and is conducted back to the *ventilator*. If an *HME* is fitted at the *patient-connection port*, approximately 50 % to 70 % of the water vapour will be trapped in the *HME*. Whatever the configuration of the expiratory limb of the *VBS*, the water vapour content of the exhaled gas will be significant, nearing saturation. Without heated wires, the returning gas cools, causing significant condensation. As in the inspiratory limb, this liquid needs to be removed. The presence of heated wires in the expiratory limb lessens or eliminates condensation before the expired gas enters the *gas return port* of the *ventilator*, but from

this point to the *exhaust port* the gas tends to cool further, so more moisture will condense. The *VBS* needs to include some means to manage this additional condensed water.

— **201.102.4.1 — Humidification system**

An *EMS ventilator* is rarely used with active humidification due to the short-term usage during transport. Normally, an *HME* is used.

— **201.102.6 — Leakage from complete *VBS***

Assuming that the leakage flow can be modelled as if an ideal orifice were producing it, then the leakage flowrate,  $Q_{leak}$ , would follow Formula (AA.4).

$$Q_{leak} = G \times \sqrt{P} \tag{AA.4}$$

where

$G$  is the orifice conductance, and

$P$  is the driving pressure.

Using the leakage limits from ISO 80601-2-12:2020 and Formula (AA.4), the orifice conductance  $G$  can be calculated for each of the *inspiratory volume* ranges. For example, the leakage limit for *inspiratory volumes*  $\geq 300$  ml is 200 ml/min at a pressure of 50 cmH<sub>2</sub>O, which yields a value for  $G$  of 28,28 ml/(min·hPa<sup>1/2</sup>). Conductance values for the other *inspiratory volume* ranges can be similarly calculated. Table AA.4 summarizes these results.

**Table AA.4 — Calculated conductance values by *inspiratory volume* range**


<i>Inspiratory volume</i> range ml	Leakage limit from ISO 80601-2-12:2020 ml/min	Pressure, P hPa (cmH <sub>2</sub> O)	Calculated conductance, G ml/(min·hPa <sup>1/2</sup> )
$V_{insp} \leq 50$	50	20	11,18
$50 \leq V_{insp} \leq 300$	100	40	15,81
$V_{insp} \geq 300$	200	50	28,28

Using these calculated conductance values, it is then possible to find the corresponding *VBS* leakage limit at any pressure. Figure AA.4 demonstrates these relationships.

Using Figure AA.4, one can derive the 60 hPa (60 cmH<sub>2</sub>O) *VBS* leakage flowrate limits. This is the pressure at which most of the *operator-accessible* parts of the *VBS* are specified at for leakage flowrate.

- for  $V_{insp} \leq 50$  ml,  $Q_{leak} = 87$  ml/min
- for  $50 \text{ ml} \leq V_{insp} \leq 300$  ml,  $Q_{leak} = 122$  ml/min
- for  $V_{insp} \geq 300$  ml,  $Q_{leak} = 219$  ml/min

These leakage flowrate limits represent the permissible leakage of the entire *VBS*. It is common to allocate 90 % of the leakage flowrate to the *operator-accessible* parts of the *VBS* and 10 % to the parts of the *VBS* internal to the *ventilator*.

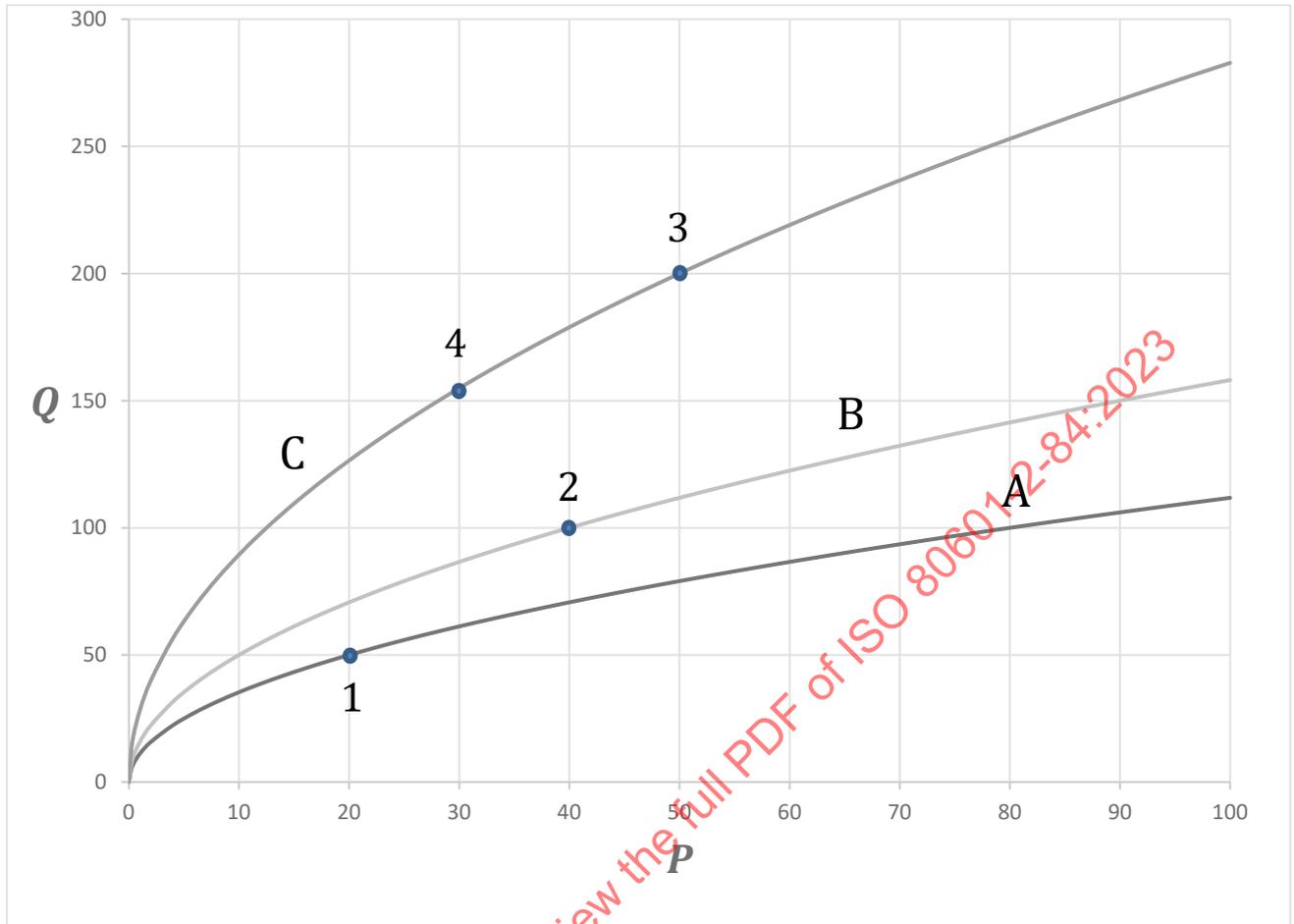**Key**

- $Q$  leakage flowrate in ml/min  
 $P$  pressure in hPa  
 1 leakage limit from ISO 80601-2-12 for  $V_{\text{insp}} \leq 50$  ml  
 2 leakage limit from ISO 80601-2-12 for  $50 \text{ ml} \leq V_{\text{insp}} \leq 300$  ml  
 3 leakage limit from ISO 80601-2-12 for  $V_{\text{insp}} \geq 300$  ml  
 4 leakage limit from ISO 80601-2-13  
 A  $Q_{\text{leak}}$  for  $V_{\text{insp}} \leq 50$  ml  
 B  $Q_{\text{leak}}$  for  $50 \text{ ml} \leq V_{\text{insp}} \leq 300$  ml  
 C  $Q_{\text{leak}}$  for  $V_{\text{insp}} \geq 300$  ml

NOTE This assumes leakage behaves as an orifice according to Formula (AA.4).

**Figure AA.4 — VBS leakage flowrate limits as a function of pressure as specified in ISO 80601-2-12 and ISO 80601-2-13**

— **201.104 — Indication of duration of operation**

*EMS ventilators* require maintenance for continued safe use. A practicable means to ensure that this information is available to the *operator* or the *responsible organization* is to require that the *EMS ventilator* keep track of how long it has been in operation.

— **201.105.2 — Connection to electronic health record**

Electronic documentation of *patient* care interventions is rapidly becoming the standard of care. The primary motivations are to improve the quality of care for an individual *patient* through accurate and complete documentation, and to improve the completeness and accuracy of aggregate data to facilitate continuous quality improvement. In some countries, there is a governmental

directive to provide electronic health records by 2015<sup>[41]</sup>. Electronic data transmission to the electronic health record is essential to meet this requirement. Annex BB contains information that *manufacturers* can find useful as requirements for data interfaces.

The data transmission should be capable of being provided with such a *functional connection* in accordance with AAMI 2700-1<sup>[41]</sup>.

— **201.105.3 — Connection to a *distributed alarm system***

See rationale for 208.6.2.2.1.

— **201.107 — Timed ventilatory pause**

Pausing *artificial ventilation* is necessary for certain clinical *procedures*.

EXAMPLE Permitting defibrillator ECG analysis with a deflated *lung*, measuring respirophasic blood pressure variation, turning the *patient*.

Currently, in order to avoid nuisance *alarm signals* and to avoid cycling the *ventilator* while the *VBS* is disconnected from the *patient*, *operators* usually turn off the *ventilator* and thereby incur the *risk* of undetected prolonged apnoea by subsequently forgetting to turn the *ventilator* back on.

In addition, there are situations where, to permit the minimum disruption of *ventilation*, the initiation of the ventilatory pause needs to come from external equipment. This is particularly important for those *procedures* (e.g. when manual synchronization would be less effective)<sup>[51]</sup>.

As part of the *risk management process*, special attention should be paid to ensuring that the *patient's lungs* remain adequately ventilated when either externally generated or repetitive ventilatory pauses occur.

The *expiratory pause* should be capable of being provided with such a *functional connection* in accordance with AAMI 2700-1<sup>[32]</sup>.

— **201.107.1 — *Expiratory pause***

— c) 1) i)

See the rationale for 201.12.4.103.1 d) 2) in Clause AA.2.

— **201.107.2 — *Inspiratory pause***

— c) 1) ii)

See the rationale for 201.12.4.103.1 d) 2) in Clause AA.2.

— **202.4.3.1 — Compliance criteria**

It is not the intent of the committees to require that the *immunity tests* be performed multiple times (e.g. with *volume-control inflation-type* and *pressure-control inflation-type* at several *inspiratory volumes*), but that the *manufacturer* should determine which *inflation-type* and *inspiratory volume* represents the worst-case condition for a given *immunity test* and use those conditions.

— **202.8.1.101 — Additional general requirements**

The committees recognized that during environmental stress (for example, shock, vibration, electromagnetic disturbances) *ventilator* performance can degrade. In such a situation, the duration