
Ophthalmic instruments — Tonometers

Instruments ophtalmiques — Tonomètres

STANDARDSISO.COM : Click to view the full PDF of ISO 8612:2009



PDF disclaimer

This PDF file may contain embedded typefaces. In accordance with Adobe's licensing policy, this file may be printed or viewed but shall not be edited unless the typefaces which are embedded are licensed to and installed on the computer performing the editing. In downloading this file, parties accept therein the responsibility of not infringing Adobe's licensing policy. The ISO Central Secretariat accepts no liability in this area.

Adobe is a trademark of Adobe Systems Incorporated.

Details of the software products used to create this PDF file can be found in the General Info relative to the file; the PDF-creation parameters were optimized for printing. Every care has been taken to ensure that the file is suitable for use by ISO member bodies. In the unlikely event that a problem relating to it is found, please inform the Central Secretariat at the address given below.

STANDARDSISO.COM : Click to view the full PDF of ISO 8612:2009



COPYRIGHT PROTECTED DOCUMENT

© ISO 2009

All rights reserved. Unless otherwise specified, no part of this publication may be reproduced or utilized in any form or by any means, electronic or mechanical, including photocopying and microfilm, without permission in writing from either ISO at the address below or ISO's member body in the country of the requester.

ISO copyright office
Case postale 56 • CH-1211 Geneva 20
Tel. + 41 22 749 01 11
Fax + 41 22 749 09 47
E-mail copyright@iso.org
Web www.iso.org

Published in Switzerland

Contents

Page

Foreword	iv
1 Scope	1
2 Normative references	1
3 Terms and definitions	1
4 Requirements	2
4.1 General	2
4.2 Design compliance testing (certification)	2
4.3 Verification (instrument compliance)	2
4.4 Construction and function	3
5 Test methods	3
6 Accompanying documents	3
7 Additional information	3
8 Marking	4
Annex A (normative) Reference tonometer and method for determining reference IOP	5
Annex B (normative) Design compliance testing	12

STANDARDSISO.COM : Click to view the full PDF of ISO 8612:2009

Foreword

ISO (the International Organization for Standardization) is a worldwide federation of national standards bodies (ISO member bodies). The work of preparing International Standards is normally carried out through ISO technical committees. Each member body interested in a subject for which a technical committee has been established has the right to be represented on that committee. International organizations, governmental and non-governmental, in liaison with ISO, also take part in the work. ISO collaborates closely with the International Electrotechnical Commission (IEC) on all matters of electrotechnical standardization.

International Standards are drafted in accordance with the rules given in the ISO/IEC Directives, Part 2.

The main task of technical committees is to prepare International Standards. Draft International Standards adopted by the technical committees are circulated to the member bodies for voting. Publication as an International Standard requires approval by at least 75 % of the member bodies casting a vote.

Attention is drawn to the possibility that some of the elements of this document may be the subject of patent rights. ISO shall not be held responsible for identifying any or all such patent rights.

ISO 8612 was prepared by Technical Committee ISO/TC 172, *Optics and photonics*, Subcommittee SC 7, *Ophthalmic optics and instruments*.

This second edition cancels and replaces the first edition (ISO 8612:2001), which has been technically revised.

STANDARDSISO.COM : Click to view the full PDF of ISO 8612:2009

Ophthalmic instruments — Tonometers

1 Scope

This International Standard, together with ISO 15004-1, specifies minimum requirements and the design compliance procedure for tonometers intended for routine clinical use in the estimation of intraocular pressure (IOP).

This International Standard takes precedence over ISO 15004-1, if differences exist.

NOTE The true intraocular pressure is seldom directly measured since it would require invasion of the eye. Since the true IOP cannot be clinically measured, alternative methods are specified for determining a reference IOP (Annex A and Annex B).

2 Normative references

The following referenced documents are indispensable for the application of this document. For dated references, only the edition cited applies. For undated references, the latest edition of the referenced document (including any amendments) applies.

ISO 15004-1, *Ophthalmic instruments — Fundamental requirements and test methods — Part 1: General requirements applicable to all ophthalmic instruments*

IEC 60601-1:2005, *Medical electrical equipment — Part 1: General requirements for basic safety and essential performance*

3 Terms and definitions

For the purposes of this document, the following terms and definitions apply.

3.1

intraocular pressure

IOP

pressure within the eye

NOTE It is expressed in millimetres of mercury (mmHg), where 1 mmHg = 0,133 3 kPa.

3.2

reference tonometer

tonometer as described in Annex A

3.3

test tonometer

verified tonometer used in design compliance testing

3.4

reference IOP

IOP that is measured with a reference tonometer, as specified in Annex A, in accordance with the procedures given in Annex B

**3.5
measured IOP**

IOP reading provided by the test tonometer when used in accordance with the manufacturer's instructions

4 Requirements

4.1 General

4.1.1 The test tonometer shall conform to the general requirements specified in ISO 15004-1.

4.1.2 The test tonometer shall conform to the specific requirements specified in 4.2 to 4.4.

4.2 Design compliance testing (certification)

4.2.1 The manufacturer shall demonstrate, on the basis of design compliance testing as specified in Clause 5, that the test tonometer measurements, when compared to the reference tonometer measurements, meet the requirements as given in Table 1.

A tonometer that meets the requirements of Annex A for a reference tonometer need not undergo design compliance testing.

The requirements are met if not more than 5 % of the paired differences between the reference tonometer reading and the test tonometer reading for each pressure range are greater than the tolerance for that range in Table 1.

NOTE The tolerances given in Table 1 represent 1,96 times the standard deviation allowable for the paired measurement, and so account for not only the allowable error of the tonometer under test but also unavoidable error associated with the reference tonometer.

Table 1 — Requirements for tonometers

IOP range mmHg	Tolerance mmHg	Minimum number of eyes
7 to 16	±5,0	40
> 16 to < 23	±5,0	40
≥ 23	±5,0	40

4.2.2 The manufacturer shall analyse the data, taken in the course of design compliance testing as specified in Clause 5, using the total least squares method for the regression, and make available, as required in Clause 7, the slope, the offset and the standard deviation of the regression line.

4.3 Verification (instrument compliance)

4.3.1 The manufacturer shall develop a method and test apparatus to confirm that the design requirements of 4.2 are met by each manufactured tonometer. Each tonometer shall be verified with this method and apparatus. This method and test apparatus shall be the same as those used to measure and verify the test tonometer in 4.2.

4.3.2 The permissible error of the test apparatus shall be less than or equal to one-half of the permissible tolerance as given in Table 1.

4.4 Construction and function

4.4.1 The surfaces of the tonometer that are intended to come into contact with the cornea shall be:

- a) composed of non-toxic, stable and non-oxidative material which is inert to ocular tissue, tears and appropriate pharmacological agents;
- b) designed either to facilitate disinfection or for single patient use;
- c) free of surface irregularities and imperfections when viewed with unmagnified corrected vision.

4.4.2 The tonometer shall permit the measurement of IOP throughout the range 7 mmHg to 50 mmHg. The scale or display shall either provide a direct measurement of a value whose relationship to IOP is known or give a numerical reading corresponding to the IOP value.

Readings of IOP less than 7 shall be displayed either by their numerical value or by a "low reading" indication. Readings of IOP greater than 50 shall be displayed either by their numerical value or by a "high reading" indication.

5 Test methods

5.1 This International Standard describes type testing (certification) and individual device testing (verification).

5.2 The reference IOP shall be determined as described in Annex A.

5.3 Design compliance testing shall be performed as described in Annex B.

6 Accompanying documents

The tonometer shall be accompanied by documents containing instructions for use together with maintenance procedures and their frequency of application. In particular, this information shall contain:

- a) name and address of the manufacturer;
- b) instructions for effective disinfection of the tonometer where applicable;
- c) contra-indications for the use of the tonometer;
- d) a list of accessories suitable for use with the tonometer;
- e) if appropriate, any additional documents as specified in 7.9 of IEC 60601-1:2005;
- f) a reference to this International Standard, i.e. ISO 8612:2009, if the manufacturer or supplier claims compliance with it.

7 Additional information

Upon request, the manufacturer shall provide information on the operating principles of the certified tonometer.

8 Marking

The tonometer shall be permanently marked with at least the following information:

- a) name of manufacturer or supplier;
- b) name and model of tonometer;
- c) if applicable, marking as required by IEC 60601-1.

STANDARDSISO.COM : Click to view the full PDF of ISO 8612:2009

Annex A (normative)

Reference tonometer and method for determining reference IOP

A.1 Specifications of the reference applanation tonometer

A.1.1 General

The reference tonometer shall be a mechanical-optical applanation tonometer that measures the force required to produce a given area of applanation.

A.1.2 Area of applanation

The area of applanation shall be circular with a diameter of 3,06 mm. The manufacturing tolerance for the diameter of the applanation circle shall be $\pm 0,02$ mm.

A.1.3 Surface of tonometer head

The front (contact) surface of the tonometer head shall have a diameter of at least 6,0 mm and shall be smooth to the touch. When examined by unmagnified corrected vision under direct illumination, the surface shall be free from surface imperfections (free from fissures, cracks and dents) that could damage the eye. Over a central area of not less than 4 mm diameter, the front surface shall be a plane structure with a "peak to valley" deviation from a plane surface of less than 3,0 μm . The outer edge of this tonometer head shall be smoothed. The front surface flatness requirement shall be verified using the method given in A.2.3.6 or an alternate equivalent method.

A.1.4 Measuring force

The measuring force shall be continuously adjustable within a minimum range extending from 0 mN to 49,0 mN, without the use of auxiliary weights. The measured value of the force shall be clearly legible on a linearly divided scale or a digital indication.

The change of force required to move the tonometer head in the opposite direction (reverse span) at the point of transition shall not exceed 0,49 mN.

A.1.5 Display

If lines are used as graduations on the measuring scale, they shall be straight, of equal width, and shall be engraved or otherwise permanently marked. No line shall be wider than 1/4 of the distance between two lines.

If a digital display is used, the increments shall be ≤ 1 mmHg.

One scale unit shall represent either 0,98 mN or 1,96 mN. The main scale graduations shall be numbered with a value. The width of the reference mark shall be no greater than the smallest width of the graduation lines on the measuring scale.

A.1.6 Tolerance for measurement of force

When the tonometer head is adjusted to the verification position, the tolerance for the measured value of the force within the measuring range shall be $\pm 1,5$ % of the nominal value or $\pm 0,49$ mN, whichever is greater, over a temperature range from 15 °C to 30 °C.

A.2 Verification of reference tonometer

A.2.1 General

The reference tonometer shall be verified using the apparatus given in A.2.2 or an alternative equivalent method.

A.2.2 Apparatus

A.2.2.1 Optical limit gauge, consisting of a left line and two right vertical lines divided horizontally by a dashed line, for testing the applanation circle diameter (see Figure A.1). The distance between the right lines corresponds to twice the value of the tolerances of the applanation circle diameter specified in A.1.2.

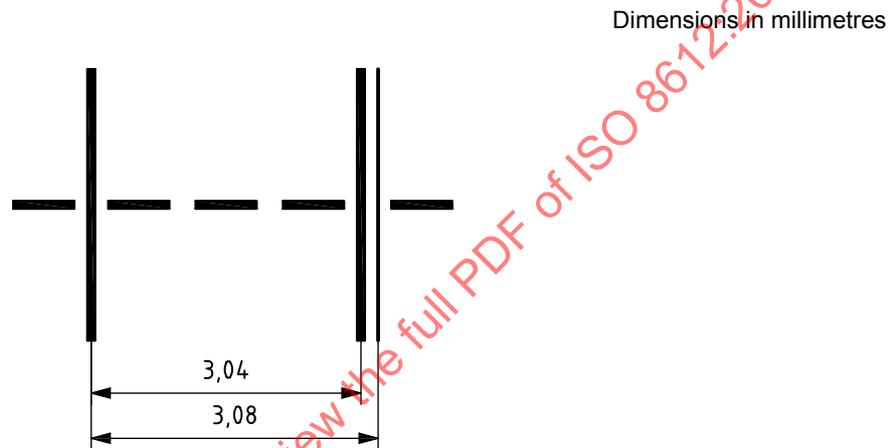


Figure A.1 — Optical limit gauge for verifying diameter of the applanation circle of 3,06 mm

A.2.2.2 Balance, of 0,01 g/scale division sensitivity for a) testing the measuring force, b) testing the reverse span for transitional movement of the tonometer head into the opposite direction, and c) for checking the position of the measuring arm with reference to its freedom of movement at equilibrium of forces.

A.2.2.3 Flatness tester, comprising a low pressure sodium lamp, a glass optical flat (flatness tolerance of $< 1/8$ wave at 589 nm) and a 10× magnifier lens for determining the flatness of the central 4 mm of the reference tonometer head (see Figure A.2).

A.2.3 Verification procedures


A.2.3.1 General

The reference tonometer shall be verified using the method given below or an equivalent method.

A.2.3.2 Diameter of the applanation circle

Substitute the optical limit gauge (A.2.2.1) for the examined eye. Orient the dividing line of the prisms so that it coincides with the dashed line. The action of the doubling prisms is to displace the images of the lines above and below the dividing line by a combined distance equal to the applanation circle diameter, so that they appear to come into coincidence (see Figure A.3). For a tonometer head without doubling prisms, determine the applanation circle diameter using a lined square that is verified by direct comparison with the optical gauge.

The tonometer shall comply with the tolerance requirements if the transposed lower half-line lies within the lateral interval delineated by the upper right line pair (Figure A.3).

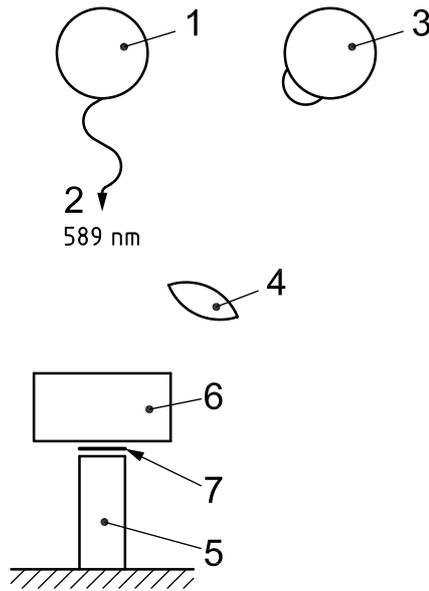

Key

- 1 low pressure sodium lamp, e.g. SOX18 (ANSI L69)
- 2 radiation of wavelength 589 nm
- 3 observer
- 4 magnifier, 10×
- 5 tonometer tip under test
- 6 glass optical flat, 1/8 wave
- 7 interference fringes

Figure A.2 — Optical apparatus for verifying 3 μm flatness of the reference tonometer head

Dimensions in millimetres

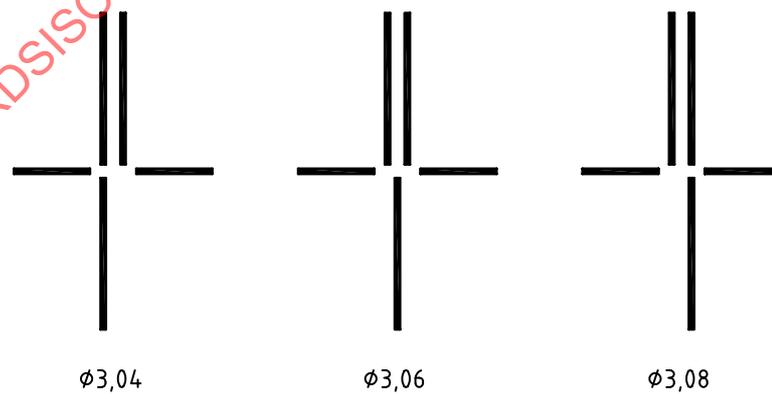


Figure A.3 — Verification of tolerance limits for diameter of the applanation circle, here for diameter values of 3,04 mm, 3,06 mm and 3,08 mm

A.2.3.3 Measuring force, reverse span and central position of the tonometer arm


The measuring force, reverse span and central position of the tonometer are verified by using the tonometer balance system (see Figure A.4) as follows.

NOTE Numbers in parentheses in the following instructions refer to Figure A.4.

Attach the applanation tonometer to the balance system by means of a holding device. Align one of the contact wheels to the tonometer head (1) using a three-coordinate fine adjustment, so that the contact wheel touches the centre of the tonometer head applanation surface and presses against it with enough force such that the balance pointer (3) is aligned to the zero index mark. The tonometer pressure setting for this step shall be in the centre of the measuring range, midway between the front and back overloading stops. Check the balance system for freedom of movement by loading the balance with a 0,5 g weight (2) and ensuring that the pointer swings to the lower balance stop. After ensuring freedom of movement of the balance system, remove the 0,5 g weight. Next, ensure that the balance system responds freely to changes in the tonometer pressure setting as follows:

- adjust the applanation tonometer to a value that corresponds to an IOP of 10 mmHg;
- symmetrically align the tonometer head with the contact wheel without touching it;
- move the tonometer head towards the contact wheel until the balance swings clearly over its zero point;
- adjust the applanation tonometer to a value that corresponds to an IOP of 0 mmHg;
- successively place weights in increments of 1 g on the balance, from 1 g up to the maximum load. Increase the measuring force of the tonometer as the weights are added to the balance, until the pointer again coincides with the zero index mark;
- carry out the test starting with the maximum weight and decreasing the tonometer measuring force;
- adjustment and checking of the applanation tonometer shall be carried out using a low-vibration device.

The tonometer shall comply with the requirements specified in A.1.4 and A.1.6 if the deviations of the various values of the measuring force from nominal values, as displayed by the tonometer scale, do not exceed the tolerance values, and if the change of force required to move the tonometer head in the reverse direction (reverse span) does not exceed the value specified in A.1.4.

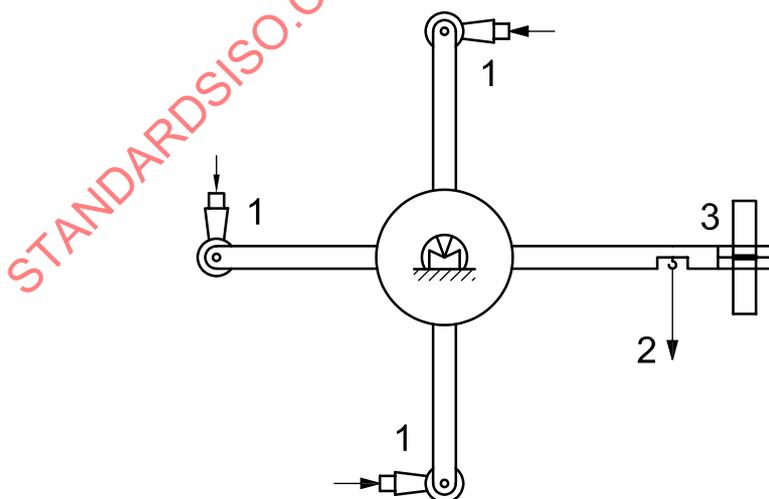

- Key**
- 1 contact wheels and pressure bodies at the different measuring positions
 - 2 balance weight
 - 3 pointer and zero index

Figure A.4 — Arrangement for verifying the measuring force

A.2.3.4 Width of lines on the measuring drum

The tonometer shall comply with the requirements specified in A.1.5, if the measured width is less than or equal to the width of the standard line.

A.2.3.5 Variation of measuring-force scale reading between 15 °C and 30 °C

The test should, in general, be carried out only once for each reference tonometer.

Place the tonometer, mounted on the tonometer balance, in a room in which the temperature can be adjusted.

Monitor the scale readings as a function of temperature.

The tonometer shall comply with the requirements specified in A.1.6 if, throughout the temperature range, the scale readings do not vary by more than the permitted limits.

A.2.3.6 Flatness over central 4 mm diameter of tonometer head

Clean the front and back surfaces of the tonometer head per manufacturer's instructions. The optical flat and the tonometer as shown in Figure A.2 shall be free of any dust or dirt that would affect the contact between the flat and the tonometer head. Place the tonometer head on a clean surface and place the optical flat on top of the tip under the sodium lamp illumination. Allow the sodium lamp to warm up for at least 5 min. Using a lens of not less than 10× magnification observe the fringe pattern appearing at the optical flat-tonometer head tip interface. There shall be ten or less fringes at the interface over the 4 mm central diameter.

A.2.4 Test certificate

The results of verification testing, carried out in accordance with this International Standard, shall be recorded in a test certificate provided by the manufacturer in the countries where no legal verification is available. An example of a test certificate is illustrated in Figure A.5.

A.3 Protocol for using the reference tonometer

The following procedure shall be carried out when using the reference tonometer in the design compliance testing specified in Annex B.

- a) Anaesthetize the selected eye of the subject.
- b) Stain with sodium fluorescein.

NOTE Steps a) and b) can be combined by using an anaesthetic to which sodium fluorescein has already been added.

- c) Set the tonometer drum to a force corresponding to an IOP of 10 mmHg.

If possible, do not touch the eyelid with the fingers in order to open the palpebral aperture. If the palpebral aperture is not wide enough to allow the tonometer cone to make contact, instruct the subject to open his/her eyes wider.

- d) Direct the subject to view a distance fixation point. If distance fixation cannot be maintained and near fixation is used, this fact should be recorded.
- e) Measure the intraocular pressure for the mean of the ocular pulse and remove the tonometer from the eye.

- f) Repeat steps c), d) and e) if the measurement was not valid because:
 - 1) the patient felt a sensation;
 - 2) the eyelid was touched;
 - 3) the fluorescein ring was too broad or too thin;
 - 4) any other circumstances suggest that the measurement may have been inaccurate.
- g) If there is any evidence that the anaesthetic is no longer fully effective, re-administer the anaesthetic.
- h) Repeat steps c), d), e) and f) if a series of multiple measurements is suggested. After the first measurement for a given subject has been taken, all valid subsequent measurements shall be recorded when a series of multiple measurements are desired. If any interruption to the consecutive readings occurs, then recommence the procedure.

STANDARDSISO.COM : Click to view the full PDF of ISO 8612:2009

Tonometer Test Certificate

Tonometer No.: _____ Date: _____

Sender: _____ Tester: _____

Manufacturer: _____ Controller: _____

Diameter of applanation circle, nominal 3,06 mm or _____ $\pm 0,02$ mm

a) Tonometer head No. 1 _____ mm b) Tonometer head No. 2 _____ mm

MEASURING FORCE

Initial state: The tonometer arm shall be adjusted symmetrically in the free space available for movement with reference to the verification position.

Nominal requirement	Error	
	Horizontal	Vertical
9,81 \pm 0,49 mN	_____ mN	_____ mN
19,61 \pm 0,49 mN	_____ mN	_____ mN
29,42 \pm 0,49 mN	_____ mN	_____ mN
39,23 \pm 0,59 mN	_____ mN	_____ mN
49,03 \pm 0,74 mN	_____ mN	_____ mN
58,84 \pm 0,88 mN	_____ mN	_____ mN
68,65 \pm 1,03 mN	_____ mN	_____ mN
78,45 \pm 1,18 mN	_____ mN	_____ mN

Reverse span at the transition of tonometer head movement into the opposite direction:

0,49 mN maximum _____ mN

Position of measuring arm relative to its free play at equilibrium of force:

Middle _____

Tip flatness (10 fringes maximum) _____ fringes

Remarks: _____

Result: _____ Verified _____ Not verified

Signature: _____

Figure A.5 — Example of a test certificate

Annex B (normative)

Design compliance testing

B.1 General

The sponsor of the test shall appoint a monitor to take overall responsibility for the conduct of the test. The monitor shall have no vested or commercial interest in the test tonometer or the outcome of the test.

B.2 Investigators

B.2.1 The monitor shall recruit a minimum of two investigators per centre who will be responsible for the measurements of the IOPs using both the reference tonometer and test tonometer.

B.2.2 The persons conducting the design compliance testing shall have no vested or commercial interest in the test tonometer or the outcome of the test.

NOTE The test can be multi-centred.

B.3 Training

B.3.1 The monitor shall ensure that all investigators are fully trained in the use of both the reference and test tonometers. The investigators shall be made aware that they shall follow the protocol and manufacturer's instructions exactly and at all times.

B.3.2 During the test, the investigators shall not be given any further instruction on the use of the tonometers, nor shall they be given any feedback on the results of the measurements.

B.4 Tonometers

B.4.1 Compliance with the specifications of the test tonometer(s) shall be verified using the methods and apparatus described in 4.3.

B.4.2 The monitor shall ensure that calibration of the reference tonometer is verified at least at the beginning and end of the test.

B.4.3 The test tonometer(s) shall not be subject to non-routine adjustment of the calibration.

B.5 Subjects

B.5.1 Subjects shall be selected according to the following exclusion criteria:

- a) subjects with only one functional eye;
- b) those with one eye having poor or eccentric fixation;
- c) high corneal astigmatism (i.e. those eyes displaying an oval contact image with the Goldmann tonometer);