

First edition  
2015-03-01

Corrected version  
2017-11

---

---

**Health informatics — Personal health  
device communication —**

Part 10442:  
**Device specialization — Strength  
fitness equipment**

*Informatique de santé — Communication entre dispositifs de santé  
personnels —*

*Partie 10442: Spécialisation des dispositifs — Équipement de mise en  
forme musculaire*

STANDARDSISO.COM : Click to view the full PDF of ISO/IEEE 11073-10442:2015



Reference number  
ISO/IEEE 11073-10442:2015(E)

© IEEE 2015

STANDARDSISO.COM : Click to view the full PDF of ISO/IEEE 11073-10442:2015



**COPYRIGHT PROTECTED DOCUMENT**

© IEEE 2015

All rights reserved. Unless otherwise specified, no part of this publication may be reproduced or utilized otherwise in any form or by any means, electronic or mechanical, including photocopying, or posting on the internet or an intranet, without prior written permission. Permission can be requested from either ISO or IEEE at the address below or ISO's member body in the country of the requester.

ISO copyright office  
Ch. de Blandonnet 8 • CP 401  
CH-1214 Vernier, Geneva, Switzerland  
Tel. +41 22 749 01 11  
Fax +41 22 749 09 47  
copyright@iso.org  
www.iso.org

Institute of Electrical and Electronics Engineers, Inc  
3 Park Avenue, New York  
NY 10016-5997, USA

stds.ipr@ieee.org  
www.ieee.org

## Foreword

ISO (the International Organization for Standardization) is a worldwide federation of national standards bodies (ISO member bodies). The work of preparing International Standards is normally carried out through ISO technical committees. Each member body interested in a subject for which a technical committee has been established has the right to be represented on that committee. International organizations, governmental and non-governmental, in liaison with ISO, also take part in the work. ISO collaborates closely with the International Electrotechnical Commission (IEC) on all matters of electrotechnical standardization.

IEEE Standards documents are developed within the IEEE Societies and the Standards Coordinating Committees of the IEEE Standards Association (IEEE-SA) Standards Board. The IEEE develops its standards through a consensus development process, approved by the American National Standards Institute, which brings together volunteers representing varied viewpoints and interests to achieve the final product. Volunteers are not necessarily members of the Institute and serve without compensation. While the IEEE administers the process and establishes rules to promote fairness in the consensus development process, the IEEE does not independently evaluate, test, or verify the accuracy of any of the information or the soundness of any judgments contained in its standards.

The main task of technical committees is to prepare International Standards. Draft International Standards adopted by the technical committees are circulated to the member bodies for voting. Publication as an International Standard requires approval by at least 75 % of the member bodies casting a vote.

Attention is called to the possibility that implementation of this standard may require the use of subject matter covered by patent rights. By publication of this standard, no position is taken with respect to the existence or validity of any patent rights in connection therewith. ISO/IEEE is not responsible for identifying essential patents or patent claims for which a license may be required, for conducting inquiries into the legal validity or scope of patents or patent claims or determining whether any licensing terms or conditions provided in connection with submission of a Letter of Assurance or a Patent Statement and Licensing Declaration Form, if any, or in any licensing agreements are reasonable or non-discriminatory. Users of this standard are expressly advised that determination of the validity of any patent rights, and the risk of infringement of such rights, is entirely their own responsibility. Further information may be obtained from ISO or the IEEE Standards Association.

ISO/IEEE 11073-10442 was prepared by the IEEE 11073 Standards Committee of the IEEE Engineering in Medicine and Biology Society (as IEEE Std 11073-10442-2008). It was adopted by Technical Committee ISO/TC 215, *Health informatics*, in parallel with its approval by the ISO member bodies, under the “fast-track procedure” defined in the Partner Standards Development Organization cooperation agreement between ISO and IEEE. IEEE is responsible for the maintenance of this document with participation and input from ISO member bodies.

This corrected version of ISO/IEEE 11073-10442:2015 incorporates the following corrections:

- replacement of front page footer and deletion of watermark.

STANDARDSISO.COM : Click to view the full PDF of ISO/IEEE 11073-10442:2015

**Health informatics—Personal health device communication**

**Part 10442: Device specialization—  
Strength fitness equipment**

Sponsor

**IEEE 11073™ Standard Committee**

of the

**IEEE Engineering in Medicine and Biology Society**

Approved 26 September 2008

**IEEE-SA Standards Board**

STANDARDSISO.COM : Click to view the full PDF of ISO/IEEE 11073-10442:2015

**Abstract:** Within the context of the ISO/IEEE 11073 family of standards for device communication, this standard establishes a normative definition of the communication between personal strength fitness devices and managers (e.g., cell phones, personal computers, personal health appliances, and set top boxes) in a manner that enables plug-and-play interoperability. It leverages appropriate portions of existing standards including ISO/IEEE 11073 terminology and information models. It specifies the use of specific term codes, formats, and behaviors in telehealth environments restricting optionality in base frameworks in favor of interoperability. This standard defines a common core of communication functionality for personal telehealth strength fitness devices. In this context, strength fitness devices are being used broadly to cover strength fitness devices that measure musculo-skeletal strength-conditioning activities.

**Keywords:** medical device communication, personal health devices, strength fitness equipment

---

The Institute of Electrical and Electronics Engineers, Inc.  
3 Park Avenue, New York, NY 10016-5997, USA

Copyright © 2009 by the Institute of Electrical and Electronics Engineers, Inc.  
All rights reserved. Published 9 January 2009. Printed in the United States of America.

IEEE is a registered trademark in the U.S. Patent & Trademark Office, owned by The Institute of Electrical and Electronics Engineers, Incorporated.

PDF: ISBN 978-0-7381-5822-8 STD95840  
Print: ISBN 978-0-7381-5823-5 STDPD95840

*No part of this publication may be reproduced in any form, in an electronic retrieval system or otherwise, without the prior written permission of the publisher.*

**IEEE Standards** documents are developed within the IEEE Societies and the Standards Coordinating Committees of the IEEE Standards Association (IEEE-SA) Standards Board. The IEEE develops its standards through a consensus development process, approved by the American National Standards Institute, which brings together volunteers representing varied viewpoints and interests to achieve the final product. Volunteers are not necessarily members of the Institute and serve without compensation. While the IEEE administers the process and establishes rules to promote fairness in the consensus development process, the IEEE does not independently evaluate, test, or verify the accuracy of any of the information contained in its standards.

Use of an IEEE Standard is wholly voluntary. The IEEE disclaims liability for any personal injury, property or other damage, of any nature whatsoever, whether special, indirect, consequential, or compensatory, directly or indirectly resulting from the publication, use of, or reliance upon this, or any other IEEE Standard document.

The IEEE does not warrant or represent the accuracy or content of the material contained herein, and expressly disclaims any express or implied warranty, including any implied warranty of merchantability or fitness for a specific purpose, or that the use of the material contained herein is free from patent infringement. IEEE Standards documents are supplied “AS IS.”

The existence of an IEEE Standard does not imply that there are no other ways to produce, test, measure, purchase, market, or provide other goods and services related to the scope of the IEEE Standard. Furthermore, the viewpoint expressed at the time a standard is approved and issued is subject to change brought about through developments in the state of the art and comments received from users of the standard. Every IEEE Standard is subjected to review at least every five years for revision or reaffirmation. When a document is more than five years old and has not been reaffirmed, it is reasonable to conclude that its contents, although still of some value, do not wholly reflect the present state of the art. Users are cautioned to check to determine that they have the latest edition of any IEEE Standard.

In publishing and making this document available, the IEEE is not suggesting or rendering professional or other services for, or on behalf of, any person or entity. Nor is the IEEE undertaking to perform any duty owed by any other person or entity to another. Any person utilizing this, and any other IEEE Standards document, should rely upon the advice of a competent professional in determining the exercise of reasonable care in any given circumstances.

Interpretations: Occasionally questions may arise regarding the meaning of portions of standards as they relate to specific applications. When the need for interpretations is brought to the attention of IEEE, the Institute will initiate action to prepare appropriate responses. Since IEEE Standards represent a consensus of concerned interests, it is important to ensure that any interpretation has also received the concurrence of a balance of interests. For this reason, IEEE and the members of its societies and Standards Coordinating Committees are not able to provide an instant response to interpretation requests except in those cases where the matter has previously received formal consideration. At lectures, symposia, seminars, or educational courses, an individual presenting information on IEEE standards shall make it clear that his or her views should be considered the personal views of that individual rather than the formal position, explanation, or interpretation of the IEEE.

Comments for revision of IEEE Standards are welcome from any interested party, regardless of membership affiliation with IEEE. Suggestions for changes in documents should be in the form of a proposed change of text, together with appropriate supporting comments. Comments on standards and requests for interpretations should be addressed to:

Secretary, IEEE-SA Standards Board  
445 Hoes Lane  
Piscataway, NJ 08854  
USA

Authorization to photocopy portions of any individual standard for internal or personal use is granted by the Institute of Electrical and Electronics Engineers, Inc., provided that the appropriate fee is paid to Copyright Clearance Center. To arrange for payment of licensing fee, please contact Copyright Clearance Center, Customer Service, 222 Rosewood Drive, Danvers, MA 01923 USA; +1 978 750 8400. Permission to photocopy portions of any individual standard for educational classroom use can also be obtained through the Copyright Clearance Center.

## Introduction

This introduction is not part of IEEE Std 11073-10442-2008, Health informatics—Personal health device communication—Part 10442: Device specialization—Strength fitness equipment.

ISO/IEEE 11073 standards enable communication between medical devices and external computer systems. This document uses the optimized framework created in IEEE Std 11073-20601<sup>a</sup> and describes a specific, interoperable communication approach for strength fitness equipment. These standards align with and draw on the existing clinically focused standards to provide easy management of data from either clinical or personal health devices.

## Notice to users

### Laws and regulations

Users of these documents should consult all applicable laws and regulations. Compliance with the provisions of this standard does not imply compliance to any applicable regulatory requirements. Implementers of the standard are responsible for observing or referring to the applicable regulatory requirements. IEEE does not, by the publication of its standards, intend to urge action that is not in compliance with applicable laws, and these documents may not be construed as doing so.

### Copyrights

This document is copyrighted by the IEEE. It is made available for a wide variety of both public and private uses. These include both use, by reference, in laws and regulations, and use in private self-regulation, standardization, and the promotion of engineering practices and methods. By making this document available for use and adoption by public authorities and private users, the IEEE does not waive any rights in copyright to this document.

### Updating of IEEE documents

Users of IEEE standards should be aware that these documents may be superseded at any time by the issuance of new editions or may be amended from time to time through the issuance of amendments, corrigenda, or errata. An official IEEE document at any point in time consists of the current edition of the document together with any amendments, corrigenda, or errata then in effect. In order to determine whether a given document is the current edition and whether it has been amended through the issuance of amendments, corrigenda, or errata, visit the IEEE Standards Association Web site at <http://ieeexplore.ieee.org/xpl/standards.jsp>, or contact the IEEE at the address listed previously.

For more information about the IEEE Standards Association or the IEEE standards development process, visit the IEEE-SA Web site at <http://standards.ieee.org>.

### Errata

Errata, if any, for this and all other standards can be accessed at the following URL: <http://standards.ieee.org/reading/ieee/updates/errata/>. Users are encouraged to check this URL for errata periodically.

### Interpretations

Current interpretations can be accessed at the following URL: <http://standards.ieee.org/reading/ieee/interp/>.

### Patents

Attention is called to the possibility that implementation of this standard may require use of subject matter covered by patent rights. By publication of this standard, no position is taken with respect to the existence or validity of any patent rights in connection therewith. A patent holder or patent applicant has filed a statement of assurance that it will grant licenses under these rights without compensation or under reasonable rates, with reasonable terms and conditions that are demonstrably free of any unfair discrimination to applicants desiring to obtain such licenses. Other Essential Patent Claims may exist for which a statement of assurance has not been received. The IEEE is not responsible for identifying Essential Patent Claims for which a license may be required, for conducting inquiries into the legal validity or scope of Patents Claims, or determining whether any licensing terms or conditions are reasonable or non-discriminatory. Further information may be obtained from the IEEE Standards Association.

<sup>a</sup> For information on references, see Clause 2.

## Participants

At the time this standard was completed, the Personal Health Devices Working Group had the following membership:

**Douglas P. Bogia**, *Chair*

**Eric White**, *Vice Chair*

Karsten Aalders	Jack Harrington	Jayant Parthasarathy
Charles R. Abbruscato	Kai Hassing	Phillip E. Pash
Maher Abuzaid	Hiroshi Hayashi	Thomas Plasa
Manfred Aigner	Torstein Heggebo	Arif Rahman
Murtaza Ali	Ron Hegli	Robert E. Ranslam
Deepak Ayyagari	Rose Higgins	Barry Reinhold
Merat Bagha	Kaoru Hiramatsu	Melvin I. Reynolds
Doug Baird	Allen Hobbs	Jeffrey S. Robbins
David Baker	Alex Holland	Timothy Robertson
Terry Bartlett	Kirsten Howard	Michael B. Robkin
David Bean	Robert Hoy	Bill Saltzstein
Rudy Belliardi	Robert D. Hughes	Stefan Sauer mann
Denis Bettini	Nick Hunn	Naveen Saxena
Ola Björsne	Yutaka Ikeda	Paul S. Schluter
Thomas Blackadar	Philip Isaacson	Lars Schmitt
Marc Blanchet	Ho-In Jeon	Mark Schnell
Douglas P. Bogia	Chris Johnson	Richard A. Schrenker
Terry Bourk	Krishna Jonnalagadda	Aravind Seshagiri
Bernard Burg	Akiyoshi Kabe	Marco Sgroi
Lyle G. Bullock, Jr.	Steve Kahle	Mazen Shihabi
Chris Burns	Tomio Kamioka	Robert Smith
Anthony Butt	Kyung Hee Kang	Motoki Sone
Carole C. Carey	Ulf Karlsson	Emily Sopensky
Randy Carroll	Andy Kaschl	Ryan Spring
Casper Chen	Junzo Kashihara	Nick Steblay
James Cheng	Kohichi Kashiwagi	Lars Steubesand
Silviu Chiricescu	Ralph Kent	John (Ivo) Stivorice
Rick A. Cnossen	Kurt Kermes	Ravi Swami
Moshe Cohen	Ikuo Keshi	Xiaorong Tai
John T. Collins	John Keys	Kunihiro Takiuchi
Cory Condek	Alfred Kloos	Francis Tam
Todd Cooper	Jeongmee Koh	Haruyuyki Tatsumi
Jim DelloStritto	Alexander Kraus	Randy Thomas
Matthew d'Entremont	Falko Kuester	Brad Tipler
Kent Dicks	Nandu Kushalnagar	Bob Tripp
Jakob Ehrensvar	Daniel Lager	Gary Tschautscher
Roger M. Ellingson	Pierre Landau	Masato Tsuchid
Michihiro Enokida	Sungkee Lee	Ken Tubman
Mika Erkkilä	Yonghee Lee	Yoshihiro Uchida
Javier Escayola Calvo	Kathryn A. Lesh	Sunil Unadkat
Leonardo Estevez	Qiong Li	Alpo Värr
Laurent Falconieri	Wei-Jung Lo	Mark Walters
Gear Fisher	Sandra Martinez	Jerry P. Wang
Julie N. Fleischer	Miguel Martínez de Espronceda	Jeff Warner
Joseph W. Forler	Cámara	Toru Watsuji
Eric Freudenthal	Jim McCain	Jeff Webber
Miguel Galarraga	Richard McPartland	Eric White
John Garguilo	Junsei Miyazaki	David L. Whitlinger
Igor Gejdos	Brian Møller	Vernon C. Williams
Chris Gough	Darr Moore	Paul Williamson
Channa Gowda	Joe Morrissey	Jan Wittenber
Niclas Granqvist	Yoshihiko Motohashi	Ariton Xhafa
Jeff Guttmacher	Alex Neefus	Ricky Yang
Christian Habermann	Michael E. Nidd	Done-Sik Yoo
Michael Hagerty	Hiroaki Niwamoto	Thomas Zhao
Rickey L. Hampton	Thomas Norgall	Daidi Zhong
Sten Hanke	Yoshiteru Nozoe	Szymon Zysko
	Mikey Paradis	

The following members of the individual balloting committee voted on this standard. Balloters may have voted for approval, disapproval, or abstention.

Thomas Blackadar  
Douglas P. Bogia  
Lyle G. Bullock, Jr.  
Randy Carroll  
Keith Chow  
Malcolm Clarke  
Rick A. Cnossen

Julie N. Fleischer  
Sergiu Goma  
Randall Groves  
Michael Hagerty  
Kai Hassing  
Werner Hoelzl  
Philip Isaacson

Atsushi Ito  
Piotr Karocki  
Kurt Kermes  
Jayant Parthasarathy  
James E. Smith  
Lars Steubesand  
Eric White

When the IEEE-SA Standards Board approved this standard on 26 September, 2008, it had the following membership:

**Robert M. Grow**, *Chair*  
**Tom A. Prevost**, *Vice Chair*  
**Steve M. Mills**, *Past Chair*  
**Judith Gorman**, *Secretary*

Victor Berman  
Richard DeBlasio  
Andrew Drozd  
Mark Epstein  
Alexander Gelman  
William R. Goldbach  
Arnold M. Greenspan  
Kenneth S. Hanus

James Hughes  
Richard H. Hulet  
Young Kyun Kim  
Joseph L. Koepfinger\*  
John Kulick  
David J. Law  
Glenn Parsons

Ronald C. Petersen  
Chuck Powers  
Narayanan Ramachandran  
Jon Walter Rosdahl  
Anne-Marie Sahazizian  
Malcolm V. Thaden  
Howard L. Wolfman  
Don Wright

\*Member Emeritus

Also included are the following nonvoting IEEE-SA Standards Board liaisons:

Satish K. Aggarwal, *NRC Representative*  
Michael H. Kelley, *NIST Representative*

Don Messina  
*IEEE Standards Program Manager, Document Development*

Kathryn Cush  
*IEEE Standards Program Manager, Technical Program Development*

## Contents

1. Overview .....	1
1.1 Scope .....	1
1.2 Purpose .....	1
1.3 Context .....	2
2. Normative references.....	2
3. Definitions, acronyms, and abbreviations .....	2
3.1 Definitions .....	2
3.2 Acronyms and abbreviations .....	3
4. Introduction to ISO/IEEE 11073 personal health devices .....	3
4.1 General .....	3
4.2 Introduction to IEEE 11073-20601 modeling constructs .....	4
5. Strength fitness device concepts and modalities.....	4
5.1 General concepts.....	4
5.2 Set.....	5
5.3 Repetition .....	5
5.4 Repetition count.....	5
5.5 Resistance.....	5
5.6 Exercise position.....	5
5.7 Exercise laterality .....	5
5.8 Exercise grip.....	5
5.9 Exercise movement.....	5
6. Strength fitness domain information model .....	5
6.1 Overview .....	5
6.2 Class extensions.....	5
6.3 Object instance diagram .....	6
6.4 Types of configuration.....	7
6.5 Medical device system object.....	7
6.6 Numeric objects.....	10
6.7 Real-time sample array objects.....	14
6.8 Enumeration objects .....	14
6.9 PM Store objects.....	20
6.10 Scanner objects .....	20
6.11 Strength fitness information model extensibility rules .....	20
7. Strength fitness service model.....	20
7.1 General .....	20
7.2 Object access services.....	20
7.3 Object access event report services .....	22

8. Strength fitness communication model .....	22
8.1 Overview .....	22
8.2 Communications characteristics .....	22
8.3 Association procedure .....	22
8.4 Configuring procedure .....	24
8.5 Operating procedure .....	24
8.6 Time synchronization .....	24
9. Test associations .....	25
10. Conformance .....	25
10.1 Applicability .....	25
10.2 Conformance specification .....	25
10.3 Levels of conformance .....	25
10.4 Implementation conformance statements .....	26
Annex A (informative) Bibliography .....	30
Annex B (normative) Any additional ASN.1 definitions .....	31
Annex C (normative) Allocation of identifiers .....	32
Annex D (informative) Use cases .....	38
Annex E (informative) Examples .....	40

STANDARDSISO.COM : Click to view the full PDF of ISO/IEEE 11073-10442:2015

## Health informatics—Personal health device communication

# Part 10442: Device specialization— Strength fitness equipment

**IMPORTANT NOTICE:** *This standard is not intended to assure safety, security, health, or environmental protection in all circumstances. Implementers of the standard are responsible for determining appropriate safety, security, environmental, and health practices or regulatory requirements.*

*This IEEE document is made available for use subject to important notices and legal disclaimers. These notices and disclaimers appear in all publications containing this document and may be found under the heading “Important Notice” or “Important Notices and Disclaimers Concerning IEEE Documents.” They can also be obtained on request from IEEE or viewed at <http://standards.ieee.org/IPR/disclaimers.html>.*

## 1. Overview

### 1.1 Scope

Within the context of the ISO/IEEE 11073 family of standards for device communication, this standard establishes a normative definition of the communication between personal strength fitness devices and managers (e.g., cell phones, personal computers, personal health appliances, and set top boxes) in a manner that enables plug-and-play interoperability. It leverages appropriate portions of existing standards, including ISO/IEEE 11073 terminology and information models. It specifies the use of specific term codes, formats, and behaviors in telehealth environments restricting optionality in base frameworks in favor of interoperability. This standard defines a common core of communication functionality for personal telehealth strength fitness devices. In this context, strength fitness devices are being used broadly to cover strength fitness devices that measure musculo-skeletal strength-conditioning activities.

### 1.2 Purpose

This standard addresses a need for an openly defined, independent standard for controlling information exchange to and from personal health devices and managers (e.g., cell phones, personal computers, personal health appliances, and set top boxes). Interoperability is the key to growing the potential market for these devices and to enabling people to be better-informed participants in the management of their health.

## 1.3 Context

See IEEE Std 11073-20601™ for an overview of the environment within which this standard is written.

This document, IEEE Std 11073-10442 defines the device specialization for the strength fitness device, being a specific agent type, and it provides a description of the device concepts, its capabilities, and its implementation according to this standard.

This standard is based on IEEE Std 11073-20601, which in turn draws information from both ISO/IEEE 11073-10201:2004 [B3]<sup>1</sup> and ISO/IEEE 11073-20101:2004 [B4]. The medical device encoding rules (MDER) used within this standard are fully described in IEEE Std 11073-20601.

This standard reproduces relevant portions of the nomenclature found in ISO/IEEE 11073-10101:2004 [B2] and adds new nomenclature codes for the purposes of this standard. Between this standard and IEEE Std 11073-20601, all required nomenclature codes for implementation are documented.

NOTE—In this standard, IEEE Std 11073-104zz is used to refer to the collection of device specialization standards that utilize IEEE Std 11073-20601, where zz can be any number from 01 to 99, inclusive.<sup>2</sup>

## 2. Normative references

The following referenced documents are indispensable for the application of this document (i.e., they must be understood and used, so that each referenced document is cited in text and its relationship to this document is explained). For dated references, only the edition cited applies. For undated references, the latest edition of the referenced document (including any amendments or corrigenda) applies.

IEEE Std 11073-20601™-2008, Health informatics—Personal health device communication—Part 20601: Application profile—Optimized exchange protocol.<sup>3, 4</sup>

See Annex A for all informative material referenced by this standard.

## 3. Definitions, acronyms, and abbreviations

### 3.1 Definitions

For the purposes of this standard, the following terms and definitions apply. The *Authoritative Dictionary of IEEE Standards Terms* [B1] should be referenced for terms not defined in this clause.

**3.1.1 agent:** A node that collects and transmits personal health data to an associated manager.

**3.1.2 class:** In object-oriented modeling, a class describes the attributes, method, and events that objects instantiated from the class utilize.

**3.1.3 compute engine:** *See: manager.*

**3.1.4 device:** A term used to refer to a physical apparatus implementing either an agent or a manager role.

<sup>1</sup> The numbers in brackets correspond to those of the bibliography in Annex A.

<sup>2</sup> Notes in text, tables, and figures are given for information only and do not contain requirements needed to implement the standard.

<sup>3</sup> The IEEE standards or products referred to in this clause are trademarks of the Institute of Electrical and Electronics Engineers, Inc.

<sup>4</sup> IEEE publications are available from the Institute of Electrical and Electronics Engineers, 445 Hoes Lane, Piscataway, NJ 08854, USA (<http://standards.ieee.org>).

IEEE Std 11073-10442-2008  
Health informatics—Personal health device communication  
Part 10442: Device specialization—Strength fitness equipment

**3.1.5 handle:** An unsigned 16-bit number that is locally unique and identifies one of the object instances within an agent.

**3.1.6 manager:** A node receiving data from one or more agent systems. Some examples of managers include a cellular phone, health appliance, set top box, or a computer system.

**3.1.7 obj-handle:** *See:* **handle**.

**3.1.8 object:** In object-oriented modeling, a particular instantiation of a class. The instantiation realizes attributes, methods, and event from the class.

**3.1.9 personal health device:** A device used in personal health applications.

**3.1.10 personal telehealth device:** *See:* **personal health device**.

## 3.2 Acronyms and abbreviations

APDU	application protocol data unit
ASN.1	Abstract Syntax Notation One
DIM	domain information model
EUI-64	extended unique identifier (64 bits)
ICS	implementation conformance statement
MDC	medical device communication
MDER	medical device encoding rules
MDS	medical device system
MOC	medical object class
OID	object identifier
PDU	protocol data unit
PHD	personal health device
RT-SA	real-time sample array
VMO	virtual medical object
VMS	virtual medical system

## 4. Introduction to ISO/IEEE 11073 personal health devices

### 4.1 General

This standard and the remainder of the series of ISO/IEEE 11073 personal health device (PHD) standards fit in the larger context of the ISO/IEEE 11073 series of standards. The full suite of standards enables agents to interconnect and interoperate with managers and with computerized health-care information systems. See IEEE Std 11073-20601 for a description of the guiding principles for this series of ISO/IEEE 11073 personal health device standards.

IEEE Std 11073-20601 supports the modeling and implementation of an extensive set of personal health devices. This standard defines aspects of the strength fitness device. It describes all aspects necessary to implement the application layer services and data exchange protocol between an ISO/IEEE 11073 PHD strength fitness agent and a manager. This standard defines a subset of the objects and functionality contained in IEEE Std 11073-20601 and extends and adds definitions where appropriate. All new definitions are given in Annex B in Abstract Syntax Notation One (ASN.1). Nomenclature codes referenced in this standard, which are not defined in IEEE Std 11073-20601, are normatively defined in Annex C.

## 4.2 Introduction to IEEE 11073-20601 modeling constructs

### 4.2.1 General

The ISO/IEEE 11073 series of standards, and in particular IEEE Std 11073-20601 is based on an object-oriented systems management paradigm. The overall system model is divided into three principal components: the domain information model (DIM), the service model, and the communications model. See IEEE Std 11073-20601 for a detailed description of the modeling constructs.

### 4.2.2 Domain information model

The DIM is a hierarchical model that describes an agent as a set of objects. These objects and their attributes represent the elements that control behavior and report on the status of the agent and data that the agent can communicate to the manager. Communication between the agent and the manager is defined by the application protocol in IEEE Std 11073-20601.

### 4.2.3 Service model

The service model defines the conceptual mechanisms for the data exchange services. Such services are mapped to messages that are exchanged between the agent and the manager. Protocol messages within the ISO/IEEE 11073 series of standards are defined in ASN.1. The messages defined in IEEE Std 11073-20601 can coexist with messages defined in other standard application profiles defined in the ISO/IEEE 11073 series of standards.

### 4.2.4 Communication model

In general, the communication model supports the topology of one or more agents communicating over logical point-to-point connections to a single manager. For each logical point-to-point connection, the dynamic system behavior is defined by a connection state machine as specified in IEEE Std 11073-20601.

### 4.2.5 Implementing the models

An agent implementing this standard shall implement all mandatory elements of the information, service, and communication models as well as all conditional elements where the condition is met. The agent should implement the recommended elements, and it may implement any combination of the optional elements. A manager implementing this standard shall utilize at least one of the mandatory, conditional, recommended, or optional elements. In this context, “utilize” means to use the element as part of the primary function of the manager device. For example, a manager whose primary function is to display data would need to display a piece of data in the element in order to utilize it.

## 5. Strength fitness device concepts and modalities

### 5.1 General concepts

Strength fitness agents measure the extent to which a person can perform a certain motion with a given resistance. The agents in this category are varied, and they range from single exercise machines, such as a leg press, to multi-function agents that can be used to perform a variety of exercises. While these agents may take a variety of forms, the measurements they capture are largely the same.

Agents in this category operate in a store and forward mode. Exercise is measured for a period of time, and sometime after the exercise period has ended, the agent transfers all the measurements to the manager. This period of exercise measurement is called a set.

## 5.2 Set

A set is an episode of a strength training exercise performed on a single agent at the same resistance. A set includes the time that it was started, the muscle group engaged, and optionally, the duration. A set may also have other measurements associated with it, such as the repetition count, details for each repetition, resistance, exercise laterality, exercise grip, exercise position, and exercise movement.

## 5.3 Repetition

A repetition is an instance of a given strength exercise. Each repetition is associated with a set and defines when the repetition was performed. Additionally, the repetition may specify the duration and/or the distance of the repetition.

## 5.4 Repetition count

The repetition count is the number of time an exercise had been performed within an associated set.

## 5.5 Resistance

The resistance is a representation of the amount of force applied to accomplish the exercise.

## 5.6 Exercise position

The position employed to perform the exercise.

## 5.7 Exercise laterality

Exercise laterality identifies which side of the body engaged in the exercise.

## 5.8 Exercise grip

The grip employed to perform the exercise.

## 5.9 Exercise movement

The movement employed in the exercise.

# 6. Strength fitness domain information model

## 6.1 Overview

This clause describes the domain information model for strength fitness equipment.

## 6.2 Class extensions

In this standard, no class extensions are defined with respect to IEEE Std 11073-20601.

### 6.3 Object instance diagram

The object instance diagram of the strength fitness domain information model, defined for the purpose of this standard, is shown in Figure 1.

The objects of the DIM, as shown in Figure 1, are described in 6.4 to 6.11. This includes the medical device system (MDS) object (see 6.5), the numeric objects (see 6.6), the real-time sample array (RT-SA) objects (see 6.7), the enumeration objects (see 6.8), the PM-store objects (see 6.9), and the scanner objects (see 6.10). See 6.11 for rules for extending the information model beyond elements as described in this standard. Each clause that describes an object of the strength fitness contains the following information:

- The nomenclature code used to identify the class of the object. One example of where this code is used is the configuration event, where the object class is reported for each object. This allows the manager to determine whether the class of the object being specified is a numeric, real time sample array, enumeration, scanner, or PM-store class.
- The attributes of the object. Each object has attributes that represent and convey information on the physical agent and its data sources. Each object has a Handle attribute that identifies the object instance within an agent. Attribute values are accessed and modified using methods such as GET and SET. Attribute types are defined using ASN.1. The ASN.1 definitions for new attribute types specific to this standard are in Annex B, and the ASN.1 definitions for existing attribute types referenced in this standard are in IEEE Std 11073-20601.
- The methods available on the object.
- The potential events generated by the object. Data are sent to the manager using events.
- The available services such as getting or setting attributes.

The attributes for each class are defined in tables that specify the name of the attribute, its value, and its qualifier. The qualifiers mean M – Attribute is Mandatory, C – Attribute is Conditional and depends on the condition stated in the Remark or Value column (if IEEE Std 11073-20601 is referenced, then it contains the conditions), R – Attribute is Recommended, NR – Attribute is Not Recommended, and O – Attribute is Optional. Mandatory attributes shall be implemented by an agent. Conditional attributes shall be implemented if the condition applies and may be implemented otherwise. Recommended attributes should be implemented by the agent. Not recommended attributes should not be implemented by the agent. Optional attributes may be implemented by the agent.

The attributes can be either static, meaning that they shall remain unchanged after the configuration is agreed upon, or dynamic, meaning that the attribute may change at some point after configuration.

IEEE Std 11073-10442-2008  
 Health informatics—Personal health device communication  
 Part 10442: Device specialization—Strength fitness equipment

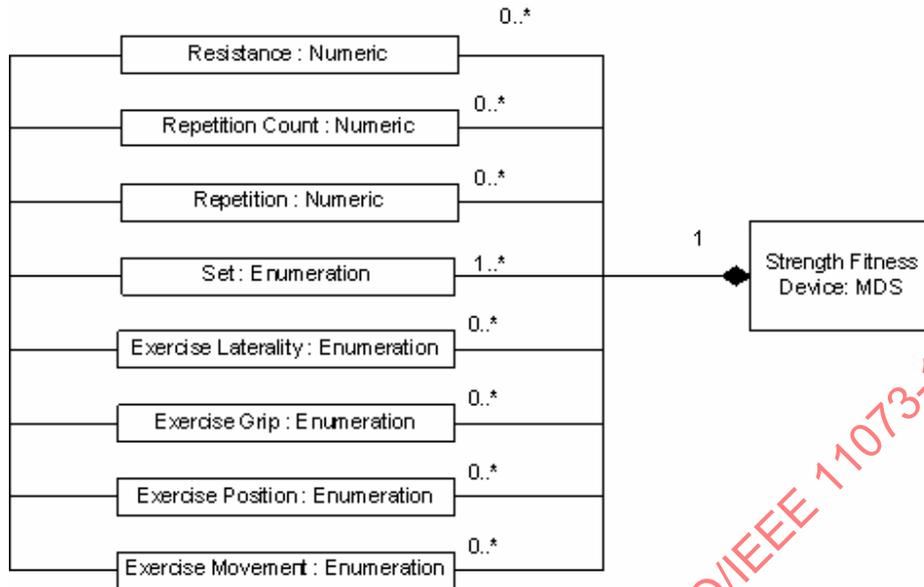


Figure 1—Object containment diagram

## 6.4 Types of configuration

### 6.4.1 General

As specified in IEEE Std 11073-20601, there are two styles of configuration available. Subclauses 6.4.2 and 6.4.3 briefly introduce standard and extended configurations.

### 6.4.2 Standard configuration

This standard does not define any standard configurations since the set of sensors for each agent configuration is likely to vary significantly for each deployment scenario. Therefore, all configurations shall be specified as extended configurations.

### 6.4.3 Extended configuration

In extended configurations, the agent’s configuration is not predefined in a standard. The agent determines which objects, attributes, and values that it wants to use in a configuration and assigns a configuration identifier. When the agent associates with a manager, it negotiates an acceptable configuration. Typically, the manager does not recognize the agent’s configuration on the first connection, so the manager responds that the agent must send its configuration information as a configuration event report. If, however, the manager already understands the configuration, either because it was preloaded in some way or the agent had previously associated with the manager, then the manager responds that the configuration is known and no further configuration information needs to be sent.

## 6.5 Medical device system object

### 6.5.1 MDS object attributes

Table 1 summarizes the attributes of the strength fitness MDS object. The nomenclature code to identify the MDS object class is MDC\_MOC\_VMS\_MDS\_SIMP.

**Table 1—MDS object attributes**

Attribute name	Value	Qual.
Handle	0	M
System-Type-Spec-List	MDC_DEV_SPEC_PROFILE_HF_STRENGTH (OID--Type) 1 (version).	M
System-Model	{“Manufacturer”, “Model”}.	M
System-Id	See IEEE Std 11073-20601.	M
Dev-Configuration-Id	0x4000--0x7FFF.	M
Attribute-Value-Map	See IEEE Std 11073-20601.	C
Production-Specification	See IEEE Std 11073-20601.	O
Mds-Time-Info	See IEEE Std 11073-20601.	C
Date-and-Time	See IEEE Std 11073-20601.	C
Relative-Time	See IEEE Std 11073-20601.	C
HiRes-Relative-Time	See IEEE Std 11073-20601.	C
Date-and-Time-Adjustment	See IEEE Std 11073-20601.	C
Power-Status	onBattery or onMains.	R
Battery-Level	See IEEE Std 11073-20601.	R
Remaining-Battery-Time	See IEEE Std 11073-20601.	R
Reg-Cert-Data-List	See IEEE Std 11073-20601.	O
System-Type-Spec-List	{MDC_DEV_SPEC_PROFILE_HF_STRENGTH, 1}.	M
Confrim-Timeout	See IEEE Std 11073-20601.	O

NOTE—See IEEE Std 11073-20601 for information on whether an attribute is static or dynamic.

In the response to a Get MDS object command, only implemented attributes and their corresponding values are returned.

See IEEE Std 11073-20601 for descriptive explanations of the individual attributes as well as for information on attribute ID and attribute type.

The Dev-Configuration-Id attribute holds a locally unique 16-bit identifier that identifies the device configuration. For a strength fitness agent, this identifier is chosen in the range of extended-config-start to extended-config-end (see IEEE Std 11073-20601) as shown in Table 1.

The agent sends the Dev-Configuration-Id during the Associating state (see 8.3) to identify its configuration for the duration of the association. If the manager already holds the configuration information relating to the Dev-Configuration-Id, it recognizes the Dev-Configuration-Id. Then the Configuration state (see 8.4) is skipped, and the agent and manager enter the Operating state. If the manager does not recognize the Dev-Configuration-Id, the agent and manager enter the Configuring state.

If an agent implements multiple IEEE 11073-104zz specializations System-Type-Spec-List is a list of type/version pairs each referencing the respective device specialization and version of that specialization.

### 6.5.2 MDS object methods

Table 2 defines the methods (actions) of the MDS object. These methods are invoked using the Action service. In Table 2, the Subservice type name column defines the name of the method; the Mode column defines whether the method is invoked as an unconfirmed action (i.e., roiv-cmip-action from IEEE Std 11073-20601) or a confirmed action (i.e., roiv-cmip-confirmed-action); the Subservice type (action-type) column defines the nomenclature code to use in the action-type field of an action request and response (see IEEE Std 11073-20601); the Parameters (action-info-args) column defines the associated ASN.1 data structure (see IEEE Std 11073-20601 and ASN.1 definitions) to use in the action message for the action-info-args field of the request; and the Results (action-info-args) column defines the structure to use in the action-info-args of the response.

**Table 2—MDS object methods**

Service	Subservice type name	Mode	Subservice type (action-type)	Parameters (action-info-args)	Results (action-info-args)
ACTION	Set-Time	Confirmed	MDC_ACT_SET_TIME	SetTimeInvoke	—

*Set-Time*

This method allows the manager to set a real-time clock in the agent with the absolute time. The agent indicates whether the Set-Time command is valid using the mds-time-capab-set-clock bit in the Mds-Time-Info attribute (see IEEE Std 11073-20601). Agents with an internal real-time clock (RTC) shall indicate this capability by also setting the mds-time-capab-real-time-clock bit in the Mds-Time-Info attribute.

Agents following only this device specialization and no others shall send event reports using agent initiated measurement data transmission. Agents following this device specialization as well as others shall send event reports in the appropriate fashion. During the association procedure (see 8.3), DataReqModeCapab shall be set to the appropriate value for the event report style. As a result, the manager shall assume the agent does not support any of the MDS-Data-Request features (see IEEE Std 11073-20601 for additional information). Thus, implementation of the MDS-Data-Request method/action is not required in this standard and is not shown in Table 2.

**6.5.3 MDS object events**

Table 3 defines the events that can be sent by the strength fitness MDS object.

**Table 3—MDS object events**

Service	Subservice type name	Mode	Subservice type (event-type)	Parameters (event-info)	Results (event-reply-info)
EVENT REPORT	MDS-Configuration-Event	Confirmed	MDC_NOTI_CONFIG	ConfigReport	ConfigReportResp
	MDS-Dynamic-Data-Update-Var	Confirmed	MDC_NOTI_SCAN_REPORT_VAR	ScanReportInfoVar	—
	MDS-Dynamic-Data-Update-Fixed	Confirmed	MDC_NOTI_SCAN_REPORT_FIXED	ScanReportInfoFixed	—
	MDS-Dynamic-Data-Update-MP-Var	Confirmed	MDC_NOTI_SCAN_REPORT_MP_VAR	ScanReportInfoMPVar	—
	MDS-Dynamic-Data-Update-MP-Fixed	Confirmed	MDC_NOTI_SCAN_REPORT_MP_FIXED	ScanReportInfoMPFixed	—

— **MDS-Configuration-Event:**

This event is sent by the agent during the configuring state of start up if the manager does not already know the agent's configuration from past associations. The event provides static information about the supported measurement capabilities of the agent.

— **MDS-Dynamic-Data-Update-Var:**

This event provides dynamic measurement data from agent. These data are reported using a generic attribute list variable format. The event is sent as an unsolicited message by the agent (i.e., an agent-initiated measurement data transmission). See 8.5.3 for more information on unsolicited event reporting.

— **MDS-Dynamic-Data-Update-Fixed:**

This event provides dynamic measurement data from the agent. These data are reported in the fixed format defined by the Attribute-Value-Map attribute of the relevant metric objects.

The event is sent as an unsolicited message by the agent (i.e., an agent-initiated measurement data transmission). See 8.5.3 for more information on unsolicited event reporting.

- **MDS-Dynamic-Data-Update-MP-Var:**  
 This is the same as MDS-Dynamic-Data-Update-Var but allows inclusion of data from multiple people.
- **MDS-Dynamic-Data-Update-MP-Fixed:**  
 This is the same as MDS-Dynamic-Data-Update-Fixed but allows inclusion of data from multiple people.

NOTE—IEEE Std 11073-20601 requires that managers support all of the MDS object events listed above.

## 6.5.4 Other MDS services

### 6.5.4.1 GET service

A strength fitness agent shall support the GET service, which is provided by the MDS object to retrieve the values of all implemented MDS object attributes. The GET service can be invoked as soon as the agent receives the Association Response and moves to the Associated state, including the Operating and Configuring substates.

The manager may request the MDS object attributes of the agent, in which the manager shall send the “Remote Operation Invoke | Get” message (see roiv-cmip-get in IEEE Std 11073-20601) with the reserved MDS handle value of 0. The agent shall report its MDS object attributes to the manager using the “Remote Operation Response | Get” message (see rors-cmip-get in IEEE Std 11073-20601). See Table 4 for a summary of the GET service including some message fields.

**Table 4—MDS object GET service**

Service	Subservice type name	Mode	Subservice type	Parameters	Results
GET	<na>	<implied confirmed>	<na>	GetArgumentSimple = (obj-handle = 0), attribute-id-list <optional>	GetResultSimple = (obj-handle = 0), attribute-list

See 8.5.2 for details on the procedure for getting the MDS object attributes.

### 6.5.4.2 SET service

The strength fitness specialization does not require an implementation to support the MDS object SET service.

## 6.6 Numeric objects

### 6.6.1 General

#### 6.6.1.1 Definitions

The numeric objects listed in this subclause represent the numeric observations that can be produced by a strength fitness agent. The nomenclature code to identify a numeric object class is MDC\_MOC\_VMO\_METRIC\_NU. Each object specifies how each of the attributes are to be used, but there are some concepts that apply to each numeric in a general fashion that are described in subclauses 6.6.1.2 through 6.6.1.5.

**6.6.1.2 Dependent attributes**

Sometimes, the interpretation of one attribute value in an object depends on other attribute values in the same object. For example, Unit-Code and Unit-Label-String provide context for the observed values. Whenever a contextual attribute changes, the agent shall report these changes to the manager using an MDS object event (see 6.5.3) prior to reporting any of the dependent values.

**6.6.1.3 Methods**

The objects listed in this clause do not support any methods.

**6.6.1.4 Events**

The objects listed in this clause do not support any events.

**6.6.1.5 Services**

The objects listed in this clause do not support any services.

**6.6.2 Repetition count**

Table 5 defines the attributes for the object that represents the number of repetitions completed in the set. The object may be supported by the agent.

**Table 5—Repetition count attributes**

Attribute name	Value	Qual.
Handle	See IEEE Std 11073-20601.	M
Type	MDC PART PHD HF   MDC HF REP COUNT.	M
Supplemental-Types	See IEEE Std 11073-20601.	NR
Metric-Spec-Small	mss-avail-intermittent   mss-avail-stored-data   mss-updt-a-periodic   mss-msmt-a-periodic   mss-acc-agent-initiated.	M
Metric-Structure-Small	See IEEE Std 11073-20601.	NR
Measurement-Status	See IEEE Std 11073-20601.	O
Metric-Id	See IEEE Std 11073-20601.	O
Metric-Id-List	See IEEE Std 11073-20601.	C
Metric-Id-Partition	See IEEE Std 11073-20601.	O
Unit-Code	See IEEE Std 11073-20601.	O
Attribute-Value-Map	See IEEE Std 11073-20601.	C
Source-Handle-Reference	See IEEE Std 11073-20601.	M
Label-String	See IEEE Std 11073-20601.	O
Unit-LabelString	See IEEE Std 11073-20601.	O
Absolute-Time-Stamp	See IEEE Std 11073-20601.	C
Relative-Time-Stamp	See IEEE Std 11073-20601.	C
HiRes-Time-Stamp	See IEEE Std 11073-20601.	C
Simple-Nu-Observed-Value	See IEEE Std 11073-20601.	C
Compound-Simple-Nu-Observed-Value	See IEEE Std 11073-20601.	NR
Basic-Nu-Observed-Value	See IEEE Std 11073-20601.	C
Compound-Basic-Nu-Observed-Value	See IEEE Std 11073-20601.	NR
Nu-Observed-Value	See IEEE Std 11073-20601.	C
Compound-Nu-Observed-Value	See IEEE Std 11073-20601.	NR
Accuracy	See IEEE Std 11073-20601.	O
Measure-Active-Period	See IEEE Std 11073-20601.	O

NOTE—See IEEE Std 11073-20601 for information on whether an attribute is static or dynamic.

The observed value reported in this object is the number of repetitions completed for a set. Only non-negative whole values shall be used, with zero (0) indicating that no repetitions were completed.

The Source-Handle-Reference attribute shall be valued with the handle of the set object to which this object is associated. Additionally, this object shall use the same timestamp attribute and value as the associated set object.

### 6.6.3 Resistance

Table 6 defines the attributes for the object that represents the resistance used for the set. The object may be supported by the agent.

**Table 6—Resistance attributes**

Attribute name	Value	Qual.
Handle	See IEEE Std 11073-20601.	M
Type	MDC_PART_PHD_HF   MDC_HF_RESISTANCE.	M
Supplemental-Types	See IEEE Std 11073-20601.	NR
Metric-Spec-Small	mss-avail-intermittent   mss-avail-stored-data   mss-updt-aperiodic   mss-msmt-aperiodic   mss-acc-agent-initiated.	M
Metric-Structure-Small	See IEEE Std 11073-20601.	NR
Measurement-Status	See IEEE Std 11073-20601.	O
Metric-Id	See IEEE Std 11073-20601.	O
Metric-Id-List	See IEEE Std 11073-20601.	C
Metric-Id-Partition	See IEEE Std 11073-20601.	O
Unit-Code	Valid values are limited to MDC_DIM_DIMLESS, MDC_DIM_X_G or MDC_DIM_LB.	M
Attribute-Value-Map	See IEEE Std 11073-20601.	C
Source-Handle-Reference	See IEEE Std 11073-20601.	M
Label-String	See IEEE Std 11073-20601.	O
Unit-LabelString	See IEEE Std 11073-20601.	O
Absolute-Time-Stamp	See IEEE Std 11073-20601.	C
Relative-Time-Stamp	See IEEE Std 11073-20601.	C
HiRes-Time-Stamp	See IEEE Std 11073-20601.	C
Simple-Nu-Observed-Value	See IEEE Std 11073-20601.	C
Compound-Simple-Nu-Observed-Value	See IEEE Std 11073-20601.	NR
Basic-Nu-Observed-Value	See IEEE Std 11073-20601.	C
Compound-Basic-Nu-Observed-Value	See IEEE Std 11073-20601.	NR
Nu-Observed-Value	See IEEE Std 11073-20601.	C
Compound-Nu-Observed-Value	See IEEE Std 11073-20601.	NR
Accuracy	See IEEE Std 11073-20601.	O
Measure-Active-Period	See IEEE Std 11073-20601.	O

NOTE—See IEEE Std 11073-20601 for information on whether an attribute is static or dynamic.

IEEE Std 11073-10442-2008  
 Health informatics—Personal health device communication  
 Part 10442: Device specialization—Strength fitness equipment

The observed value reported in this object is the resistance setting for a set. Only non-negative values shall be used, with zero (0) indicating that no resistance was used. When the resistance is measured in weight, either MDC\_DIM\_X\_G or MDC\_DIM\_LB shall be used for the Unit-Code attribute. If resistance is measured in an indexed value, such as setting 1, 2, 3, etc., or some other unit, then the MDC\_DIM\_DIMLESS shall be used with, optionally, the Unit-Label-String attribute providing additional information as needed.

The Source-Handle-Reference attribute shall be valued with the handle of the set object to which this object is associated. Additionally, this object shall use the same timestamp attribute and value as the associated set object.

**6.6.4 Repetition**

Table 7 defines the attributes for the object that represents a repetition in a set. The object may be supported by the agent.

**Table 7—Repetition attributes**

Attribute name	Value	Qual.
Handle	See IEEE Std 11073-20601.	M
Type	MDC_PART_PHD_HF   MDC_HF_REPETITION.	M
Supplemental-Types	See IEEE Std 11073-20601.	NR
Metric-Spec-Small	mss-avail-intermittent   mss-avail-stored-data   mss-updt-aperiodic   mss-msmt-aperiodic   mss-acc-agent-initiated.	M
Metric-Structure-Small	See IEEE Std 11073-20601.	NR
Measurement-Status	See IEEE Std 11073-20601.	O
Metric-Id	See IEEE Std 11073-20601.	O
Metric-Id-List	See IEEE Std 11073-20601.	C
Metric-Id-Partition	See IEEE Std 11073-20601.	O
Unit-Code	Valid values are MDC_DIM_X_M, MDC_DIM_X_INCH.	M
Attribute-Value-Map	See IEEE Std 11073-20601.	C
Source-Handle-Reference	See IEEE Std 11073-20601.	M
Label-String	See IEEE Std 11073-20601.	O
Unit-LabelString	See IEEE Std 11073-20601.	O
Absolute-Time-Stamp	See IEEE Std 11073-20601.	C
Relative-Time-Stamp	See IEEE Std 11073-20601.	C
HiRes-Time-Stamp	See IEEE Std 11073-20601.	C
Simple-Nu-Observed-Value	See IEEE Std 11073-20601.	C
Compound-Simple-Nu-Observed-Value	See IEEE Std 11073-20601.	NR
Basic-Nu-Observed-Value	See IEEE Std 11073-20601.	C
Compound-Basic-Nu-Observed-Value	See IEEE Std 11073-20601.	NR
Nu-Observed-Value	See IEEE Std 11073-20601.	C
Compound-Nu-Observed-Value	See IEEE Std 11073-20601.	NR
Accuracy	See IEEE Std 11073-20601.	O
Measure-Active-Period	See IEEE Std 11073-20601.	O

NOTE—See IEEE Std 11073-20601 for information on whether an attribute is static or dynamic.

The observed value reported in this object is the distance traveled during the repetition, measured as the distance moved from the starting position (i.e., the one-way distance, not the roundtrip cumulative distance). Only non-negative values shall be used. The Measure-Active-Period may be used to specify the amount of time taken to complete the repetition.

The Source-Handle-Reference shall be valued with the handle of the set object to which this object is associated. Additionally, this object shall use the same timestamp attribute as the associated set object.

If the related set object has a Measure-Active-Period defined, the timestamp for this object shall fall within the time range defined by the associated set object's timestamp attribute and shall last for the Measure-Active-Period defined by the related set object's Measure-Active-Period attribute. Additionally, if the repetition object specifies a Measure-Active-Period, the time period defined by the repetition object shall fall completely within the time range defined by the associated set object.

## 6.7 Real-time sample array objects

Real-time sample array objects are not required by this standard.

## 6.8 Enumeration objects

### 6.8.1 General

#### 6.8.1.1 Definitions

The enumeration objects listed in this clause represent the enumeration observations that can be produced by a strength fitness agent. The nomenclature code used to identify an enumeration object class is MDC\_MOC\_VMO\_METRIC\_ENUM. Each object specifies how each of the attributes should be used, but there are some concepts that apply to each enumeration in a general fashion that are described here.

#### 6.8.1.2 Methods

The objects listed in this clause do not support any methods.

#### 6.8.1.3 Events

The objects listed in this clause do not support any events.

#### 6.8.1.4 Services

The objects listed in this clause do not support any services.

### 6.8.2 Set

Table 8 defines the attributes for the object that represents a set. The agent shall include one instance of this object for each measurement period reported.

IEEE Std 11073-10442-2008  
 Health informatics—Personal health device communication  
 Part 10442: Device specialization—Strength fitness equipment

Table 8—Set attributes

Attribute name	Value	Qual.
Handle	See IEEE Std 11073-20601.	M
Type	MDC PART PHD HF   MDC HF SET.	M
Supplemental-Types	See IEEE Std 11073-20601.	NR
Metric-Spec-Small	mss-avail-intermittent   mss-avail-stored-data   mss-updt-aperiodic   mss-msmt-aperiodic   mss-acc-agent-initiated.	M
Metric-Structure-Small	See IEEE Std 11073-20601.	NR
Measurement-Status	See IEEE Std 11073-20601.	O
Metric-Id	See IEEE Std 11073-20601.	O
Metric-Id-List	See IEEE Std 11073-20601.	C
Metric-Id-Partition	See IEEE Std 11073-20601.	O
Unit-Code	See IEEE Std 11073-20601.	O
Attribute-Value-Map	See IEEE Std 11073-20601.	C
Source-Handle-Reference	See IEEE Std 11073-20601.	NR
Label-String	See IEEE Std 11073-20601.	O
Unit-LabelString	See IEEE Std 11073-20601.	O
Absolute-Time-Stamp	See IEEE Std 11073-20601.	C
Relative-Time-Stamp	See IEEE Std 11073-20601.	C
HiRes-Time-Stamp	See IEEE Std 11073-20601.	C
Measure-Active-Period	See IEEE Std 11073-20601.	C
Enum-Observed-Value-Simple-Oid	Valid values are limited to MDC_MUSC *.	M
Enum-Observed-Value-Simple-Bit-Str	See IEEE Std 11073-20601.	NR
Enum-Observed-Value-Basic-Bit-Str	See IEEE Std 11073-20601.	NR
Enum-Observed-Value-Simple-Str	See IEEE Std 11073-20601.	NR
Enum-Observed-Value	See IEEE Std 11073-20601.	NR
Enum-Observed-Value-Partition	MDC_PART_SITES.	M

NOTE—See IEEE Std 11073-20601 for information on whether an attribute is static or dynamic.

The purpose of this object is to define the set at a high level by specifying the muscle group engaged, the date and time of the exercise, and how much time was spent in the activity. This object defines the context for all the objects associated with it.

The muscle group being engaged is defined in the Enum-Observed-Value-Simple-Oid attribute.

The duration of the set, if known, shall be specified in the Measure-Active-Period attribute.

### 6.8.3 Exercise position

Table 9 defines the attributes for the object that represents the user position employed for the exercise defined in the associated set object. This object may be supported by the agent.

**Table 9—Exercise position attributes**

Attribute name	Value	Qual.
Handle	See IEEE Std 11073-20601.	M
Type	MDC_PART_PHD_HF   MDC_HF_EXERCISE_POSITION.	M
Supplemental-Types	See IEEE Std 11073-20601.	NR
Metric-Spec-Small	mss-avail-intermittent   mss-avail-stored-data   mss-updt-aperiodic   mss-msmt-aperiodic   mss-acc-agent-initiated.	M
Metric-Structure-Small	See IEEE Std 11073-20601.	NR
Measurement-Status	See IEEE Std 11073-20601.	O
Metric-Id	See IEEE Std 11073-20601.	O
Metric-Id-List	See IEEE Std 11073-20601.	C
Metric-Id-Partition	See IEEE Std 11073-20601.	O
Unit-Code	See IEEE Std 11073-20601.	O
Attribute-Value-Map	See IEEE Std 11073-20601.	C
Source-Handle-Reference	See IEEE Std 11073-20601.	M
Label-String	See IEEE Std 11073-20601.	O
Unit-LabelString	See IEEE Std 11073-20601.	O
Absolute-Time-Stamp	See IEEE Std 11073-20601.	C
Relative-Time-Stamp	See IEEE Std 11073-20601.	C
HiRes-Time-Stamp	See IEEE Std 11073-20601.	C
Measure-Active-Period	See IEEE Std 11073-20601.	O
Enum-Observed-Value-Simple-Oid	Valid values are limited to MDC_HF_POSITION_*.	M
Enum-Observed-Value-Simple-Bit-Str	See IEEE Std 11073-20601.	NR
Enum-Observed-Value-Basic-Bit-Str	See IEEE Std 11073-20601.	NR
Enum-Observed-Value-Simple-Str	See IEEE Std 11073-20601.	NR
Enum-Observed-Value	See IEEE Std 11073-20601.	NR
Enum-Observed-Value-Partition	See IEEE Std 11073-20601.	NR

NOTE—See IEEE Std 11073-20601 for information on whether an attribute is static or dynamic.

The purpose of this object is to define the position of the exerciser when performing the exercise defined in the associated set. The position is defined in the Enum-Observed-Value-Simple-Oid attribute. This object's timestamp shall be equal to the timestamp of the associated set object. Additionally, the timestamp attribute used for an exercise position object instance shall be the same as that which is used for the set object instance to which it is related. The Source-Handle-Reference attribute shall be equal to the handle of the associated set object.

#### 6.8.4 Exercise laterality

Table 10 defines the attributes for the object that represents laterality of the exercise defined in the set. This object may be supported by the agent.

IEEE Std 11073-10442-2008  
Health informatics—Personal health device communication  
Part 10442: Device specialization—Strength fitness equipment

**Table 10—Exercise laterality attributes**

Attribute name	Value	Qual.
Handle	See IEEE Std 11073-20601.	M
Type	MDC_PART_PHD_HF   MDC_HF_EXERCISE_LATERALITY.	M
Supplemental-Types	See IEEE Std 11073-20601.	NR
Metric-Spec-Small	mss-avail-intermittent   mss-avail-stored-data   mss-updt-a-periodic   mss-msmt-a-periodic   mss-acc-agent-initiated.	M
Metric-Structure-Small	See IEEE Std 11073-20601.	NR
Measurement-Status	See IEEE Std 11073-20601.	O
Metric-Id	See IEEE Std 11073-20601.	O
Metric-Id-List	See IEEE Std 11073-20601.	C
Metric-Id-Partition	See IEEE Std 11073-20601.	O
Unit-Code	See IEEE Std 11073-20601.	O
Attribute-Value-Map	See IEEE Std 11073-20601.	C
Source-Handle-Reference	See IEEE Std 11073-20601.	M
Label-String	See IEEE Std 11073-20601.	O
Unit-LabelString	See IEEE Std 11073-20601.	O
Absolute-Time-Stamp	See IEEE Std 11073-20601.	C
Relative-Time-Stamp	See IEEE Std 11073-20601.	C
HiRes-Time-Stamp	See IEEE Std 11073-20601.	C
Measure-Active-Period	See IEEE Std 11073-20601.	O
Enum-Observed-Value-Simple-Oid	Valid values are limited to MDC_HF_LATERALITY_*.	M
Enum-Observed-Value-Simple-Bit-Str	See IEEE Std 11073-20601.	NR
Enum-Observed-Value-Basic-Bit-Str	See IEEE Std 11073-20601.	NR
Enum-Observed-Value-Simple-Str	See IEEE Std 11073-20601.	NR
Enum-Observed-Value	See IEEE Std 11073-20601.	NR
Enum-Observed-Value-Partition	See IEEE Std 11073-20601.	NR

NOTE—See IEEE Std 11073-20601 for information on whether an attribute is static or dynamic.

The purpose of this object is to define whether the user employed their left, right, or both sides when performing the exercise defined in the associated set. The laterality is defined in the Enum-Observed-Value-Simple-Oid attribute. This object's timestamp shall be equal to the timestamp of the associated set object. Additionally, the timestamp attribute used for an exercise modifier object instance shall be the same as that which is used for the set object instance to which it is related. The Source-Handle-Reference attribute shall be equal to the handle of the associated set object.

6.8.5 Exercise grip

Table 11 defines the attributes for the object that represents the grip employed in the performance of the exercise defined in the associated set. This object may be supported by the agent.

Table 11—Exercise grip attributes

Attribute name	Value	Qual.
Handle	See IEEE Std 11073-20601.	M
Type	MDC_PART_PHD_HF MDC_HF_EXERCISE_GRIP.	M
Supplemental-Types	See IEEE Std 11073-20601.	NR
Metric-Spec-Small	mss-avail-intermittent   mss-avail-stored-data   mss-updt-aperiodic   mss-msmt-aperiodic   mss-acc-agent-initiated.	M
Metric-Structure-Small	See IEEE Std 11073-20601.	NR
Measurement-Status	See IEEE Std 11073-20601.	O
Metric-Id	See IEEE Std 11073-20601.	O
Metric-Id-List	See IEEE Std 11073-20601.	C
Metric-Id-Partition	See IEEE Std 11073-20601.	O
Unit-Code	See IEEE Std 11073-20601.	O
Attribute-Value-Map	See IEEE Std 11073-20601.	C
Source-Handle-Reference	See IEEE Std 11073-20601.	M
Label-String	See IEEE Std 11073-20601.	O
Unit-LabelString	See IEEE Std 11073-20601.	O
Absolute-Time-Stamp	See IEEE Std 11073-20601.	C
Relative-Time-Stamp	See IEEE Std 11073-20601.	C
HiRes-Time-Stamp	See IEEE Std 11073-20601.	C
Measure-Active-Period	See IEEE Std 11073-20601.	O
Enum-Observed-Value-Simple-Oid	Valid values are limited to MDC_HF_GRIP_*.	M
Enum-Observed-Value-Simple-Bit-Str	See IEEE Std 11073-20601.	NR
Enum-Observed-Value-Basic-Bit-Str	See IEEE Std 11073-20601.	NR
Enum-Observed-Value-Simple-Str	See IEEE Std 11073-20601.	NR
Enum-Observed-Value	See IEEE Std 11073-20601.	NR
Enum-Observed-Value-Partition	See IEEE Std 11073-20601.	NR

NOTE—See IEEE Std 11073-20601 for information on whether an attribute is static or dynamic.

The purpose of this object is to define the grip employed by the exerciser when performing the exercise defined in the associated set. The grip is defined in the Enum-Observed-Value-Simple-Oid attribute. This object's timestamp shall be equal to the timestamp of the associated set object. Additionally, the timestamp attribute used for an exercise grip object instance shall be the same as that which is used for the set object

IEEE Std 11073-10442-2008  
 Health informatics—Personal health device communication  
 Part 10442: Device specialization—Strength fitness equipment

instance to which it is related. The Source-Handle-Reference attribute shall be equal to the handle of the associated set object.

**6.8.6 Exercise movement**

Table 12 defines the attributes for the object that represents the movement employed during the exercise defined in the set. This object may be supported by the agent.

**Table 12—Exercise movement attributes**

Attribute name	Value	Qual.
Handle	See IEEE Std 11073-20601.	M
Type	MDC_PART_PHD_HF   MDC_HF_EXERCISE_MOVEMENT.	M
Supplemental-Types	See IEEE Std 11073-20601.	NR
Metric-Spec-Small	mss-avail-intermittent   mss-avail-stored-data   mss-updt-aperiodic   mss-msmt-aperiodic   mss-acc-agent-initiated.	M
Metric-Structure-Small	See IEEE Std 11073-20601.	NR
Measurement-Status	See IEEE Std 11073-20601.	O
Metric-Id	See IEEE Std 11073-20601.	O
Metric-Id-List	See IEEE Std 11073-20601.	C
Metric-Id-Partition	See IEEE Std 11073-20601.	O
Unit-Code	See IEEE Std 11073-20601.	O
Attribute-Value-Map	See IEEE Std 11073-20601.	C
Source-Handle-Reference	See IEEE Std 11073-20601.	O
Label-String	See IEEE Std 11073-20601.	O
Unit-LabelString	See IEEE Std 11073-20601.	O
Absolute-Time-Stamp	See IEEE Std 11073-20601.	C
Relative-Time-Stamp	See IEEE Std 11073-20601.	C
HiRes-Time-Stamp	See IEEE Std 11073-20601.	C
Measure-Active-Period	* See IEEE Std 11073-20601.	O
Enum-Observed-Value-Simple-Oid	Valid values are limited to MDC_HF_MOVEMENT*.	M
Enum-Observed-Value-Simple-Bit-Str	See IEEE Std 11073-20601.	NR
Enum-Observed-Value-Basic-Bit-Str	See IEEE Std 11073-20601.	NR
Enum-Observed-Value-Simple-Str	See IEEE Std 11073-20601.	NR
Enum-Observed-Value	See IEEE Std 11073-20601.	NR
Enum-Observed-Value-Partition	See IEEE Std 11073-20601.	NR

NOTE—See IEEE Std 11073-20601 for information on whether an attribute is static or dynamic.

The purpose of this object is to specify the movement used in the exercise being performed in the associated set. The exercise movement is defined in the Enum-Observed-Value-Simple-Oid attribute.

This object's timestamp shall be equal to the timestamp of the associated set object. Additionally, the timestamp attribute used for an exercise movement object instance shall be the same as that which is used for the set object instance to which it is associated. The Source-Handle-Reference attribute shall be equal to the handle of the associated set object.

### 6.9 PM-store objects

PM-store objects are not required by this standard.

### 6.10 Scanner objects

Scanner objects are not required by this standard.

### 6.11 Strength fitness information model extensibility rules

Beyond the information model as specified in this standard, the device designer may include other vendor-specific metrics and attributes if required. However, any metric or attribute that is implemented should follow the guidelines of this standard as closely as possible.

## 7. Strength fitness service model

### 7.1 General

The service model defines the conceptual mechanisms for data exchange services. These services are mapped to messages that are exchanged between the agent and the manager. Protocol messages within the ISO/IEEE 11073 series of standards are defined in ASN.1. See IEEE Std 11073-20601 for a detailed description of the personal health device service model. Subclauses 7.2 through 7.3 define the specifics of object access and event reporting services for a strength fitness agent according to this standard.

### 7.2 Object access services

The object access services of IEEE Std 11073-20601 are used to access the information objects defined in the domain information model of the strength fitness agent.

The following generic object access services are supported by a strength fitness agent according to this standard:

- GET service: used by the manager to retrieve the values of the agent MDS object attributes. The list of MDS object attributes is given in 6.5.1.
- SET service: used by the manager to set the values of the agent object attributes. There are no settable attributes defined for a strength fitness agent according to this standard.
- Event report service: used by the agent to send configuration reports and measurement data to the manager. The list of event reports for the strength fitness device specialization is given in 6.5.3.
- Action service: used by the manager to invoke actions (or methods) supported by the agent. An example is Set-Time action, which is used to set a real-time clock with the absolute time at the agent.

Table 13 summarizes the object access services supported by the strength fitness equipment object access model.

IEEE Std 11073-10442-2008  
 Health informatics—Personal health device communication  
 Part 10442: Device specialization—Strength fitness equipment

**Table 13—Strength fitness MDS object events**

Service	Subservice type name	Mode	Subservice Type	Parameters	Result	Remarks
GET	<na>	<implied Confirmed>	<na>	GetArgumentSimple = (obj-handle = 0), attribute-id-list <optional>	GetResultSimple = (obj-handle = 0), attribute-list	Allows the manager to retrieve the value of an attribute of an object in the agent.
EVENT REPORT	MDS-Configuration-Event	Confirmed	MDC_NOTI_CONFIG	ConfigReport	ConfigReport Rsp	Configuration report to inform manager of the configuration of the agent.
	MDS-Dynamic-Data-Update-Var	Confirmed	MDC_NOTI_SCAN_REPORT_VAR	ScanReportInfoVar	—	Data report to provide dynamic data to manager for some or all of the agent's objects in variable format.
	MDS-Dynamic-Data-Update-Fixed	Confirmed	MDC_NOTI_SCAN_REPORT_FIXED	ScanReportInfoFixed	—	Data report to provide dynamic data to manager for some or all of the agent's objects in fixed format.
	MDS-Dynamic-Data-Update-MP-Var	Confirmed	MDC_NOTI_SCAN_REPORT_MP_VAR	ScanReportInfoMPVar	—	This is the same as MDS-Dynamic-Data-Update-Var but allows inclusion of data from multiple people.
	MDS-Dynamic-Data-Update-MP-Fixed	Confirmed	MDC_NOTI_SCAN_REPORT_MP_FIXED	ScanReportInfoMPFixed	—	This is the same as MDS-Dynamic-Data-Update-Fixed but allows inclusion of data from multiple people.
ACTION	Set-Time	Confirmed	MDC_ACT_SET_TIME	SetTimeInvoke	—	Manager method to invoke the agent to set time to requested value.

### 7.3 Object access event report services

The event report service (see Table 13) is used by the agent to report its information (e.g., measurements). Event reports in this standard are a property of the MDS object only. The event reports used in this standard are defined in detail in IEEE Std 11073-20601.

The following conditions apply for a strength fitness agent according to this standard:

- Event reports shall be used in confirmed mode.
- Agent-initiated mode shall be supported for measurement data transmission.

An agent designed to operate in an environment where data may be collected from multiple people may use one of the multiple-person event report styles to transmit all the data from each person in a single event. If this functionality is not required, the agent may use the single-person event report styles, which have reduced overhead.

A manager shall support both single-person and multiple-person event reports. A strength fitness agent may support either one or both single-person and multiple-person event reports. The formats for single and multiple person reports are described in IEEE Std 11073-20601.

## 8. Strength fitness communication model

### 8.1 Overview

This clause describes the general communication model and procedures of the strength fitness agent with respect to the general communication model and procedures as defined in IEEE Std 11073-20601. Therefore, the respective parts of IEEE Std 11073-20601 are not reproduced, but rather the specific choices and restrictions with respect to optional elements (e.g., attributes, and procedures) and specific extensions (e.g., nomenclature terms) are specified.

For an illustrative overview of the various message transactions during a typical measurement period, see the sequence diagram for the example use case in Annex D and the corresponding protocol data unit (PDU) examples in Annex E.

### 8.2 Communications characteristics

This standard defines no additional constraints on the communication characteristics defined in IEEE Std 11073-20601.

### 8.3 Association procedure

#### 8.3.1 General

Unless otherwise stated, the association procedure for a strength fitness agent and manager according to this standard shall be pursued as specified in IEEE Std 11073-20601.

### 8.3.2 Agent procedure—association request

In the association request sent by the agent to the manager:

- The version of the association procedure used by the agent shall be set to the *assoc-version1* (i.e., *assoc-version* = 0x80000000).
- The DataProtoList structure element of the data protocol identifier shall be set to *data-PROTO-ID-20601* (i.e., *data-PROTO-ID* = 0x5079).
- The *data-PROTO-INFO* field shall contain a PhdAssociationInformation structure that shall contain the following parameter values:
  - 1) The version of the data exchange protocol shall be set to *protocol-version1* (i.e., *protocol-version* = 0x80000000).
  - 2) At least the MDER shall supported (i.e., *encoding-rules* = 0x8000).
  - 3) The version of the nomenclature used shall be set to *nom-version1* (i.e., *nomenclature-version* = 0x80000000).
  - 4) The field *functional-units* may have the test association bits set but shall not have any other bits set.
  - 5) The field *system-type* shall be set to *sys-type-agent* (i.e., *system-type* = 0x00800000).
  - 6) The *system-id* field shall be set to the value of the System-Id attribute of the MDS object of the agent. The manager may use this field to determine the identity of the agent with which it is associating and, optionally, to implement a simple access restriction policy.
  - 7) The *dev-config-id* field shall be set to the value of the Dev-Configuration-Id attribute of the MDS object of the agent.
  - 8) If the agent supports only the strength fitness specialization, then the field indicating the data request modes (*data-req-mode-capab*) supported by the strength fitness agent shall be set to *data-req-supp-init-agent*.
  - 9) If the agent supports only the strength fitness specialization, then the *data-req-init-manager-count* shall be set to zero, and the *data-req-init-agent-count* shall be set to 1.

### 8.3.3 Manager procedure—association response

In the association response message sent by the manager:

- The *result* field shall be set to an appropriate response from those defined in IEEE Std 11073-20601. For example, if all other conditions of the association protocol are satisfied, *accepted* is returned when the manager recognizes the *dev-config-id* of the agent and *accepted-unknown-config* otherwise.
- In the DataProtoList structure element, the data protocol identifier shall be set to *data-PROTO-ID-20601* (i.e., *data-PROTO-ID* = 0x5079).
- The *data-PROTO-INFO* field shall be filled in with a PhdAssociationInformation structure which shall contain the following parameter values:
  - 1) The version of the data exchange protocol shall be set to *protocol-version1* (i.e., *protocol-version* = 0x80000000).
  - 2) The manager shall respond with a single selected encoding rule that is supported by both agent and manager. The manager shall support at least the MDER.
  - 3) The version of the nomenclature used shall be set to *nom-version1* (i.e., *nomenclature-version* = 0x80000000).
  - 4) The field *functional-units* shall have all bits reset except for those relating to a test association.
  - 5) The field *system-type* shall be set to *sys-type-manager* (i.e., *system-type* = 0x80000000).

- 6) The *system-id* field shall contain the unique system id of the manager device, which shall be a valid extended unique identifier (64 bits) (EUI-64) type identifier.
- 7) The field *dev-config-id* shall be manager-config-response (0).
- 8) The field *data-req-mode-capab* shall be 0.
- 9) The field *data-req-init-\** count shall be 0.

## 8.4 Configuring procedure

The agent enters the Configuring state if it receives an association response of accepted-unknown-config. In this case, the configuration procedure as specified in IEEE Std 11073-20601 shall be followed.

## 8.5 Operating procedure

### 8.5.1 General

The communication of data and status information about the strength fitness agent occurs during the operating state. If not stated otherwise, the operating procedure for a strength fitness agent according to this standard shall be pursued as specified in IEEE Std 11073-20601.

### 8.5.2 GET strength fitness MDS attributes

See Table 4 for a summary of the GET service including some message fields.

If the manager leaves the attribute-id-list field in the roiv-cmip-get service message empty, the agent shall respond with a rors-cmip-get service message, in which the attribute-list contains a list with the values of all implemented attributes of the MDS object.

If the manager requests specific MDS object attributes indicated by the elements in attribute-id-list and the agent supports this capability, the agent shall respond with a rors-cmip-get service message, in which the attribute-list contains a list of the values of the requested attributes of the MDS object that are implemented. If this capability is not implemented, the agent shall respond with a “Remote Operation Error Result” (roer) service message (see IEEE Std 11073-20601), with the error-value field is set to no-such-action (9).

### 8.5.3 Measurement data transmission

See Table 3 for a summary of the event report services available for measurement data transfer.

Measurement data transfer for a strength fitness agent of this standard shall always be initiated by the strength fitness agent (see agent-initiated measurement data transmission in IEEE Std 11073-20601). To limit the amount of data being transported within an application protocol data unit (APDU), the agent shall not include more than 25 temporarily stored measurements in a single event report. If more than 25 pending measurements are available for transmission, they shall be sent using multiple event reports. If multiple measurements are available, up to 25 measurements should be transmitted within a single event report. Alternatively, they may be transmitted using a single event report for each measurement. However, the former strategy is recommended to reduce overall message size and power consumption.

## 8.6 Time synchronization

Time synchronization between a strength fitness agent and a manager may be used to coordinate the clocks used when reporting physiological events. Note that the mechanism for synchronizing an agent to a manager is outside the scope of this standard. If time synchronization is used, then this shall be reported in the Mds-Time-Info attribute of the MDS object.

## 9. Test associations

A cardiovascular fitness and activity monitor standard does not define any test associations.

## 10. Conformance

### 10.1 Applicability

This standard shall be used in conjunction with IEEE Std 11073-20601.

An implementation or a system can conform to the following elements of this standard:

- Domain information model class hierarchy and object definitions (object attributes, notifications, methods and data type definitions)
- Nomenclature code values
- Protocol and service models
- Communication service model (association and configuration)

However, conformance to these elements alone does not provide interoperability of applications or medical devices.

### 10.2 Conformance specification

This standard offers levels of optionality in how the model is applied for an implementation, particularly in the following areas:

- Information model of a specific device
- Use of attributes, value ranges, and access

To support interoperability of applications and systems, an implementation based on this standard shall provide specific details about the way that the definitions of this standard are applied.

These specifications have to be provided in the form of a set of implementation conformance statements. An implementation conformance statement (ICS) is a form of data sheet that discloses details of a specific implementation and specifies which features are provided.

### 10.3 Levels of conformance

#### 10.3.1 General

This standard defines the following levels of conformance for applications.

#### 10.3.2 Conformance level 1

The application uses the object definitions (object hierarchy, data type definitions) and the nomenclature scheme defined in this standard. All mandatory features defined in the object definition tables and in the

ICS tables are implemented. Furthermore, any conditional or optional features that are implemented shall follow the requirements in this standard.

No private (nonstandard) extensions to the object model are used (no additional objects, object attributes, notifications, object methods, or object relations). No extensions to the nomenclature codes are used (no IDs from the private code space).

### 10.3.3 Conformance level 2

Conformance level 2 is the same as conformance level 1, except the application uses private extensions to the object model or the application uses private extensions to the nomenclature (private codes) or the application uses other nomenclature schemes (external nomenclature) in addition to the standard nomenclature.

## 10.4 Implementation conformance statements

### 10.4.1 General format

The ICSs have to be supplied in the form of tables. Subclauses 10.4.2 through 10.4.5 contain templates for these ICS tables. The tables have to be filled out and provided as an overall conformance statement document.

Generally, the column headings of an ICS table contain the following information:

Index	Feature	Reference	Status	Support	Comment
-------	---------	-----------	--------	---------	---------

The table column headings have the following meaning:

- Index is an identifier (e.g., a tag) of a specific feature
- Feature briefly describes the characteristic for which a conformance statement shall be made.
- Reference is a reference to the definition of the feature (may be empty).
- Status specifies the conformance requirement (e.g., mandatory or recommended)—in some cases, this standard does not specify conformance requirements but still wants a definition of the status of a particular feature.
- Support is filled out by the implementer and specifies the characteristics of the feature in the implementation.
- Comment contains additional information provided by the implementer.

More complex expressions or specific lists of items are defined in the specific ICS table.

### 10.4.2 General implementation conformance statement

In a general top-level ICS, the implementer specifies the versions/revisions that are supported by the implementation as well as some high-level system behavior definitions.

Table 14 shows the general ICSs.

IEEE Std 11073-10442-2008  
 Health informatics—Personal health device communication  
 Part 10442: Device specialization—Strength fitness equipment

Table 14—IEEE 11073-10442 general ICSs' table

Index <sup>a</sup>	Feature	Reference	Status	Support	Comment
GEN 11073-10442-1	Implementation Description	—	Identification of the device/application. Description of functionality.		
GEN 11073-10442-2	Standards followed and their revisions	(standard documents)	(set of existing revisions)	(set of supported revision)	
GEN 11073-10442-3	Conformance Adherence - Level 1 -	See 10.3.2	Base conformance declaration that device meets the following IEEE Std 11073-10442 conformance requirements: a) All minimum mandatory (shall) requirements. b) All conditional and optional requirements that were implemented conform to this standard. It is expected that these would be fully described as part of the device specializations conformance statements.	Yes/No (No is not expected as No implies that the implementation is nonconformant)	
GEN 11073-10442-4	Conformance Adherence - Level 2 -	See 10.3.3	In addition to GEN 11073-10442-3, the device implements extensions.	Yes/No	
GEN 11073-10442-5	Object Containment Tree	See 6.3	Provide Object Containment Diagram showing relations between object instances used by the application. A conforming implementation uses only object relations as defined in the DIM.		
GEN 11073-10442-6	Nomenclature document used and revision	(standard documents)	(set of existing revisions)	(set of supported revision)	
GEN 11073-10442-7	Data Structure Encoding	—	—	description of encoding method(s) for ASN.1 data structures	
GEN 11073-10442-8	Use of Private Objects	—	Does the implementation use objects that are not defined in the DIM?	Yes/No (If yes: explain in DIM MOC ICS)	
GEN 11073-10442-9	Use of Private Nomenclature Extensions	—	Does the implementation use private extensions to the nomenclature? Private nomenclature extensions are <i>only</i> allowed if the standard nomenclature does not include the specific terms required by the application.	Yes/No (If yes: explain in the appropriate ICS)	
GEN 11073-10442-10	11073-20601 Conformance		Provide the conformance report required by IEEE Std 11073-10442.		

<sup>a</sup>The prefix GEN11073-10471- is used for the index in the general ICSs table.

**10.4.3 DIM medical object class (MOC) implementation conformance statement**

The DIM MOC ICS template shown in Table 15 is used to define which objects (not base classes) are used by the implementation. The following table is a template only. For each object supported by the implementation, one row shall be filled out.

**Table 15—Template for DIM MOC ICS table**

Index	Feature	Reference	Status	Support	Comment
MOC-[1-n]	Object Name and Object ID	Reference to the clause in the standard where the object is defined.	Implemented	Specify restrictions (e.g., maximum number of supported instances)	

If the implementation uses private objects, these objects should also be specified here. A separate definition can be appended to the conformance statement that can be referenced in the Reference column.

The Support column should indicate specific restrictions about the object implementation.

In addition to the MOC ICS table, an object containment diagram (class instance diagram) should be provided that allows reviewing the class hierarchy used by the implementation.

**10.4.4 MOC attribute implementation conformance statement**

For each supported object as defined in the DIM MOC ICS, a MOC Attribute ICS has to be provided that defines which attributes are used/supported by the implementation, including any inherited attributes. Table 16 is a template only.

**Table 16—Template for MOC attribute ICS table**

Index	Feature	Reference	Status	Support	Comment
ATTR-x-n	Attribute Name and Attribute ID	Reference to the clause in the standard where the object is defined.	M = Mandatory / C = Conditional / O = Optional (as per definition in Attribute Definition Tables)	Implemented? Yes/No	

If the implementation uses private attributes, these attributes should also be specified here. A separate definition can be appended to the conformance statement that can be referenced in the Reference column.

The x in the Index column is the ID of the managed object for which the table is supplied (this is the Index of the managed object as specified in the DIM MOC ICS). There is one separate table for each supported managed object.

The n in the Index column is just a serial number (1..m).

The Support column should also contain attribute value ranges (if applicable), hints about specific restrictions for attribute access, or attribute availability and information, if the attribute value is static or dynamic in the implementation.

NOTE—The attribute definition tables in the standard define a minimum mandatory set of attributes for each object.

**10.4.5 MOC notification implementation conformance statement**

The MOC notification ICS specifies all implemented notifications (typically in the form of the event report service) that are emitted by the agent. Table 17 is a template only. One table has to be provided for each object that supports special object notifications.

**Table 17—Template for MOC notification ICS table**

Index	Feature	Reference	Status	Support	Comment
NOTI-x-n	Notification Name and Notification ID	Reference to the clause in the standard where the event is defined.		Describe how attribute is accessed (e.g., Get, Set, sent in Configuration event report, sent in a data event report). Describe any specific restrictions.	

The x in the Index column is the ID of the managed object for which the table is supplied (this is the Index of the managed object as specified in the DIM MOC ICS). There is one separate table for each managed object that supports specific object notifications (events).

The n in the Index column is just a serial number (1..m).

STANDARDSISO.COM : Click to view the full PDF of ISO/IEEE 11073-10442:2015

## Annex A

(informative)

## Bibliography

[B1] IEEE 100™, *The Authoritative Dictionary of IEEE Standards Terms*, Seventh Edition. New York, Institute of Electrical and Electronic Engineers, Inc.<sup>5,6</sup>

[B2] ISO/IEEE 11073-10101:2004, Health informatics — Point-of-care medical device communication — Part 10101: Nomenclature.<sup>7</sup>

[B3] ISO/IEEE 11073-10201:2004, Health informatics — Point-of-care medical device communication — Part 10201: Domain information model.

[B4] ISO/IEEE 11073-20101:2004, Health informatics — Point-of-care medical device communication — Part 20101: Application profile — Base standard.

[B5] ITU-T Rec. X.680-2002, Information technology — Abstract Syntax Notation One (ASN.1): Specification of basic notation.<sup>8</sup>

<sup>5</sup>IEEE publications are available from the Institute of Electrical and Electronics Engineers, 445 Hoes Lane, Piscataway, NJ 08854, USA (<http://standards.ieee.org>).

<sup>6</sup>The IEEE standards or products referred to in this clause are trademarks of the Institute of Electrical and Electronics Engineers, Inc.

<sup>7</sup>ISO/IEEE publications are available from the ISO Central Secretariat, Case Postale 56, 1 rue de Varembé, CH-1211, Genève 20, Switzerland/Suisse (<http://www.iso.ch/>). ISO/IEEE publications are also available in the United States from the Institute of Electrical and Electronics Engineers, 445 Hoes Lane, Piscataway, NJ 08854, USA (<http://standards.ieee.org>).

<sup>8</sup>ITU publications are available from the International Telecommunications Union, Place des Nations, 1211 Geneva 20, Switzerland (<http://www.itu.in/>).

IEEE Std 11073-10442-2008  
Health informatics—Personal health device communication  
Part 10442: Device specialization—Strength fitness equipment

## Annex B

(normative)

### Any additional ASN.1 definitions

No additional ASN.1 definitions are defined.

STANDARDSISO.COM : Click to view the full PDF of ISO/IEEE 11073-10442:2015

## Annex C

(normative)

### Allocation of identifiers

This annex contains the nomenclature codes used in this document and not found in IEEE Std 11073-20601. For those not contained in this annex, the normative definition is found in IEEE Std 11073-20601.

The format used here follows the one defined in ISO/IEEE 11073-10101.

```

/* Partition codes
#define MDC_PART_SITES 7 /* Body Sites */
#define MDC_PART_PHD_HF 129 /* Health and Fitness */

/*****
* From Dimensions (MDC_PART_DIM)
*****/
#define MDC_DIM_DIMLESS 512 /* dimension-less */
#define MDC_DIM_X_M 1280 /* m */
#define MDC_DIM_X_INCH 1376 /* in */
#define MDC_DIM_X_G 1728 /* g */
#define MDC_DIM_LB 1760 /* lb */

/*****
* From Body Site (MDC_PART_SITES)
*****/

#define MDC_MUSC_SKELETAL 248 /* */
#define MDC_MUSC_HEAD 252 /* */
#define MDC_MUSC_HEAD_EYE 256 /* */
#define MDC_MUSC_HEAD_RECT_SUP 260 /* */
#define MDC_MUSC_HEAD_RECT_INF 264 /* */
#define MDC_MUSC_HEAD_RECT_MED 268 /* */
#define MDC_MUSC_HEAD_RECT_LAT 272 /* */
#define MDC_MUSC_HEAD_OBLIQ_SUP 276 /* */
#define MDC_MUSC_HEAD_OBLIQ_INF 280 /* */
#define MDC_MUSC_HEAD_FACIAL 284 /* */
#define MDC_MUSC_HEAD_OCCIPITOFRONT_VENTER 288 /* */
#define MDC_MUSC_HEAD_ORBIC_OCUL 292 /* */
#define MDC_MUSC_HEAD_ORBIC_OCUL_PARS_ORBIT 296 /* */
#define MDC_MUSC_HEAD_AURIC_POST 300 /* */
#define MDC_MUSC_HEAD_ORBIC_ORIS 304 /* */
#define MDC_MUSC_HEAD_DEPRESSOR_ANGUL_ORIS 308 /* */
#define MDC_MUSC_HEAD_RISOR 312 /* */
#define MDC_MUSC_HEAD_ZYGOMATIC_MAJOR 316 /* */
#define MDC_MUSC_HEAD_ZYGOMATIC_MINOR 320 /* */
#define MDC_MUSC_HEAD_LEVATOR_LAB_SUP 324 /* */
#define MDC_MUSC_HEAD_LEVATOR_LAB_SUP_AL_NASI 328 /* */
#define MDC_MUSC_HEAD_DEPRESSOR_LAB_INF 332 /* */
#define MDC_MUSC_HEAD_LEVATOR_ANGUL_ORIS 336 /* */
#define MDC_MUSC_HEAD_BUCCINATOR 340 /* */
#define MDC_MUSC_HEAD_MENTAL 344 /* */
    
```

IEEE Std 11073-10442-2008  
 Health informatics—Personal health device communication  
 Part 10442: Device specialization—Strength fitness equipment

```

#define MDC_MUSC_HEAD_MASSETER 348 /* */
#define MDC_MUSC_HEAD_TEMPOR 352 /* */
#define MDC_MUSC_HEAD_PTERTYGOID 356 /* */
#define MDC_MUSC_HEAD_PTERTYGOID_LAT 360 /* */
#define MDC_MUSC_HEAD_PTERTYGOID_MED 364 /* */
#define MDC_MUSC_HEAD_LING 368 /* */
#define MDC_MUSC_HEAD_GENIOGLOSS 372 /* */
#define MDC_MUSC_HEAD_LARING 376 /* */
#define MDC_MUSC_HEAD_CRICOTHYROID 380 /* */
#define MDC_MUSC_HEAD_THYROARYTEROID 384 /* */
#define MDC_MUSC_NECK 388 /* */
#define MDC_MUSC_NECK_PLATYSMA 392 /* */
#define MDC_MUSC_NECK_CAPT_LONG 396 /* */
#define MDC_MUSC_NECK_STERNOCLEIDOMASTOID 400 /* */
#define MDC_MUSC_NECK_DIGRASTRIC 404 /* */
#define MDC_MUSC_NECK_DIGRASTRIC_VENTER_ANT 408 /* */
#define MDC_MUSC_NECK_DIGRASTRIC_VENTER_POST 412 /* */
#define MDC_MUSC_NECK_MYLOHYOID 416 /* */
#define MDC_MUSC_TRUNK 420 /* */
#define MDC_MUSC_BACK 424 /* */
#define MDC_MUSC_BACK_UPPER 428 /* */
#define MDC_MUSC_BACK_LOWER 432 /* */
#define MDC_MUSC_BACK_TRAPEZ 436 /* */
#define MDC_MUSC_BACK_LASTISSIM_DORS 440 /* */
#define MDC_MUSC_BACK_RHOMB_MAJOR 444 /* */
#define MDC_MUSC_BACK_RHOMB_MINOR 448 /* */
#define MDC_MUSC_BACK_SCAP_LEVATOR 452 /* */
#define MDC_MUSC_BACK_SERRAT_POST 456 /* */
#define MDC_MUSC_BACK_SPLEN_CAPT 460 /* */
#define MDC_MUSC_BACK_SPLEN_CERVIC 464 /* */
#define MDC_MUSC_BACK_SPLEN 468 /* */
#define MDC_MUSC_BACK_SPINAL_ERECTOR 472 /* */
#define MDC_MUSC_BACK_SPINAL 476 /* */
#define MDC_MUSC_BACK_SPINAL_THORAC 480 /* */
#define MDC_MUSC_BACK_SPINAL_CERVIC 484 /* */
#define MDC_MUSC_BACK_SPINAL_CAPIT 488 /* */
#define MDC_MUSC_BACK_SEMISPINAL 492 /* */
#define MDC_MUSC_BACK_SEMISPINAL_THOR 496 /* */
#define MDC_MUSC_BACK_SEMISPINAL_CERV 500 /* */
#define MDC_MUSC_BACK_SEMISPINAL_CAPIT 504 /* */
#define MDC_MUSC_BACK_MULTIFID 508 /* */
#define MDC_MUSC_BACK_INTERSPINAL 512 /* */
#define MDC_MUSC_BACK_INTERSPINAL_CERVIC 516 /* */
#define MDC_MUSC_BACK_INTERSPINAL_THORAC 520 /* */
#define MDC_MUSC_BACK_INTERSPINAL_LUMBOR 524 /* */
#define MDC_MUSC_THORAX 528 /* */
#define MDC_MUSC_THORAX_PECTORAL_MAJOR 532 /* */
#define MDC_MUSC_THORAX_PECTORAL_MINOR 536 /* */
#define MDC_MUSC_THORAX_SUBCLAV 540 /* */
#define MDC_MUSC_THORAX_SERRAT_ANT 544 /* */
#define MDC_MUSC_THORAX_INTERCOSTAL 548 /* */
#define MDC_MUSC_THORAX_DIAPHRAGM 552 /* */
#define MDC_MUSC ABDOM 556 /* */
#define MDC_MUSC ABDOM ABDOMIN 560 /* */
#define MDC_MUSC ABDOM OBLIQ_EXT 564 /* */

```

```

#define MDC_MUSC_ABDOM_OBLIQ_INT 568 /* */
#define MDC_MUSC_ABDOM_ABDOM_TRANSVERS 572 /* */
#define MDC_MUSC_ABDOM_LUMBOR_QUADRAT 576 /* */
#define MDC_MUSC_ABDOM_PELV 580 /* */
#define MDC_MUSC_ABDOM_PUBORECT 584 /* */
#define MDC_MUSC_ABDOM_COCCYG 588 /* */
#define MDC_MUSC_ABDOM_ANI_SPHINCTER 592 /* */
#define MDC_MUSC_ABDOM_ANI_SPHINCTER_EXT 596 /* */
#define MDC_MUSC_UPEXT 600 /* */
#define MDC_MUSC_UPEXT_DELTOID 604 /* */
#define MDC_MUSC_UPEXT_SUPRASPINAT 608 /* */
#define MDC_MUSC_UPEXT_INFRASPINAT 612 /* */
#define MDC_MUSC_UPEXT_TERES_MINOR 616 /* */
#define MDC_MUSC_UPEXT_TERES_MAJOR 620 /* */
#define MDC_MUSC_UPEXT_SUBSCAP 624 /* */
#define MDC_MUSC_UPEXT_BRACHI_BICEPS 628 /* */
#define MDC_MUSC_UPEXT_BRACHIAL 632 /* */
#define MDC_MUSC_UPEXT_CORACOBRACH 636 /* */
#define MDC_MUSC_UPEXT_BRACH_TRICEPS 640 /* */
#define MDC_MUSC_UPEXT_BRACH_TRICEPS_CAP_LONG 644 /* */
#define MDC_MUSC_UPEXT_BRACH_TRICEPS_CAP_LAT 648 /* */
#define MDC_MUSC_UPEXT_BRACH_TRICEPS_CAP_MED 652 /* */
#define MDC_MUSC_UPEXT_ANCON 656 /* */
#define MDC_MUSC_UPEXT_PRONATOR 660 /* */
#define MDC_MUSC_UPEXT_FLEX_CARPI_RADIAL 664 /* */
#define MDC_MUSC_UPEXT_PALMAR_LONG 668 /* */
#define MDC_MUSC_UPEXT_FLEX_CARPI_ULNAR 672 /* */
#define MDC_MUSC_UPEXT_FLEX_DIGIT_SUPERF 676 /* */
#define MDC_MUSC_UPEXT_FLEX_DIGIT_PROFUND 680 /* */
#define MDC_MUSC_UPEXT_FLEX_POLLIC_LONG 684 /* */
#define MDC_MUSC_UPEXT_PRONATOR_QUADRAT 688 /* */
#define MDC_MUSC_UPEXT_BRACHIORADIAL 692 /* */
#define MDC_MUSC_UPEXT_EXTENS_CARP_RADIAL_LONG 696 /* */
#define MDC_MUSC_UPEXT_EXTENS_CARP_RADIAL_BREV 700 /* */
#define MDC_MUSC_UPEXT_EXTENS_DIGIT 704 /* */
#define MDC_MUSC_UPEXT_EXTENS_DIGIT_MIN 708 /* */
#define MDC_MUSC_UPEXT_EXTENS_CARP_ULNAR 712 /* */
#define MDC_MUSC_UPEXT_SUPINATOR 716 /* */
#define MDC_MUSC_UPEXT_ABDUC_POLLIC_LONG 720 /* */
#define MDC_MUSC_UPEXT_EXTENS_POLLIC_BREV 724 /* */
#define MDC_MUSC_UPEXT_EXTENS_POLLIC_LONG 728 /* */
#define MDC_MUSC_UPEXT_EXTENS_INDIC 732 /* */
#define MDC_MUSC_UPEXT_PALMAR_BREV 736 /* */
#define MDC_MUSC_UPEXT_ABDUC_POLLIC_BREV 740 /* */
#define MDC_MUSC_UPEXT_FLEX_POLLIC_BREV 744 /* */
#define MDC_MUSC_UPEXT_OPPON_POLLIC 748 /* */
#define MDC_MUSC_UPEXT_ADDUC_POLLIC 752 /* */
#define MDC_MUSC_UPEXT_ABDUC_DIGIT_MIN 756 /* */
#define MDC_MUSC_UPEXT_FLEX_DIGIT_BREV_MIN 760 /* */
#define MDC_MUSC_UPEXT_OPPON_DIGIT_MIN 764 /* */
#define MDC_MUSC_UPEXT_LUMBRICAL 768 /* */
#define MDC_MUSC_UPEXT_INTEROSS_DORSAL 772 /* */
#define MDC_MUSC_UPEXT_INTEROSS_PALMAR 776 /* */
#define MDC_MUSC_LOEXT_HIP_THIGH 780 /* */
#define MDC_MUSC_LOEXT_LEG 784 /* */

```